

(19) World Intellectual Property Organization
International Bureau(43) International Publication Date
2 November 2006 (02.11.2006)

PCT

(10) International Publication Number
WO 2006/116127 A2(51) International Patent Classification:
C12Q 1/68 (2006.01)(74) Agents: CASIMIR, David et al.; Medlen & Carroll, 101
Howard Street, Suite 350, San Francisco, CA 94105 (US).(21) International Application Number:
PCT/US2006/015160(81) Designated States (unless otherwise indicated, for every
kind of national protection available): AE, AG, AL, AM,
AT, AU, AZ, BA, BB, BG, BR, BW, BY, BZ, CA, CH, CN,
CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, EG, ES, FI,
GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE,
KG, KM, KN, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV,
LY, MA, MD, MG, MK, MN, MW, MX, MY, NA, NG, NI,
NO, NZ, OM, PG, PH, PL, PT, RO, RU, SC, SD, SE, SG,
SK, SL, SM, SY, TJ, TM, TN, TR, TT, TZ, UA, UG, US,
UZ, VC, VN, YU, ZA, ZM, ZW.

(22) International Filing Date: 21 April 2006 (21.04.2006)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:
60/674,118 21 April 2005 (21.04.2005) US
60/705,631 3 August 2005 (03.08.2005) US
60/732,539 1 November 2005 (01.11.2005) US
60/773,124 13 February 2006 (13.02.2006) US(84) Designated States (unless otherwise indicated, for every
kind of regional protection available): ARIPO (BW, GH,
GM, KE, LS, MW, MZ, NA, SD, SL, SZ, TZ, UG, ZM,
ZW), Eurasian (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM),
European (AT, BE, BG, CH, CY, CZ, DE, DK, EE, ES, FI,
FR, GB, GR, HU, IE, IT, LT, LV, LU, MC, NL, PL, PT,
RO, SE, SI, SK, TR), OAPI (BF, BJ, CF, CG, CI, CM, GA,
GN, GQ, GW, ML, MR, NE, SN, TD, TG).(71) Applicant (for all designated States except US): ISIS
PHARMACEUTICALS, INC. [US/US]; 1896 Ruther-
ford Rd., Carlsbad, California 92008 (US).

(72) Inventors; and

(75) Inventors/Applicants (for US only): SAMPATH, Ran-
garajan [US/US]; 12223 Mannix Road, San Diego, Cal-
ifornia 92129 (US); HALL, Thomas, A. [US/US]; 5239
Wohlford St., Oceanside, California 92056 (US); FSHOO,
Mark, W. [US/US]; 615 Glenmont Dr., Solana Beach, Cal-
ifornia 92075 (US).

Published:

— without international search report and to be republished
upon receipt of that reportFor neo-letter codes and other abbreviations, refer to the "Guid-
ance Notes on Codes and Abbreviations" appearing at the begin-
ning of each regular issue of the PCT Gazette.

(54) Title: COMPOSITIONS FOR USE IN IDENTIFICATION OF BACTERIA

(57) Abstract: The present invention provides compositions, kits and methods for rapid identification and quantification of bacteria
by molecular mass and base composition analysis.

WO 2006/116127 A2

COMPOSITIONS FOR USE IN IDENTIFICATION OF BACTERIA

STATEMENT OF GOVERNMENT SUPPORT

[01] This invention was made with United States Government support under CDC contract RO1 CI000099-01. The United States Government has certain rights in the invention.

FIELD OF THE INVENTION

[02] The present invention provides compositions, kits and methods for rapid identification and quantification of bacteria by molecular mass and base composition analysis.

BACKGROUND OF THE INVENTION

[03] A problem in determining the cause of a natural infectious outbreak or a bioterrorist attack is the sheer variety of organisms that can cause human disease. There are over 1400 organisms infectious to humans; many of these have the potential to emerge suddenly in a natural epidemic or to be used in a malicious attack by bioterrorists (Taylor et al. Philos. Trans. R. Soc. London B. Biol. Sci., 2001, 356, 983-989). This number does not include numerous strain variants, bioengineered versions, or pathogens that infect plants or animals.

[04] Much of the new technology being developed for detection of biological weapons incorporates a polymerase chain reaction (PCR) step based upon the use of highly specific primers and probes designed to selectively detect certain pathogenic organisms. Although this approach is appropriate for the most obvious bioterrorist organisms, like smallpox and anthrax, experience has shown that it is very difficult to predict which of hundreds of possible pathogenic organisms might be employed in a terrorist attack. Likewise, naturally emerging human disease that has caused devastating consequence in public health has come from unexpected families of bacteria, viruses, fungi, or protozoa. Plants and animals also have their natural burden of infectious disease agents and there are equally important biosafety and security concerns for agriculture.

[05] A major conundrum in public health protection, biodefense, and agricultural safety and security is that these disciplines need to be able to rapidly identify and characterize infectious agents, while there is no existing technology with the breadth of function to meet this need. Currently used methods for identification of bacteria rely upon culturing the bacterium to effect isolation from other

organisms and to obtain sufficient quantities of nucleic acid followed by sequencing of the nucleic acid, both processes which are time and labor intensive.

[06] Mass spectrometry provides detailed information about the molecules being analyzed, including high mass accuracy. It is also a process that can be easily automated. DNA chips with specific probes can only determine the presence or absence of specifically anticipated organisms. Because there are hundreds of thousands of species of benign bacteria, some very similar in sequence to threat organisms, even arrays with 10,000 probes lack the breadth needed to identify a particular organism.

[07] The present invention provides oligonucleotide primers and compositions and kits containing the oligonucleotide primers, which define bacterial bioagent identifying amplicons and, upon amplification, produce corresponding amplification products whose molecular masses provide the means to identify bacteria, for example, at and below the species taxonomic level.

SUMMARY OF THE INVENTION

[08] The present invention provides compositions, kits and methods for rapid identification and quantification of bacteria by molecular mass and base composition analysis.

[09] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 456.

[10] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1261.

[11] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 456 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1261.

[12] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 288.

[13] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1269.

[14] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 288 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1269.

[15] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 698.

[16] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1420.

[17] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 698 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1420.

[18] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 217.

[19] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1167

[20] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 217 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1167.

[21] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 399.

[22] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1041.

[23] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 399 and an

oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1041.

[24] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 430.

[25] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1321.

[26] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 430 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1321.

[27] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 174.

[28] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 853.

[29] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 174 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 853.

[30] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 172.

[31] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1360.

[32] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 172 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1360.

[33] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 456 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1261.

[34] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 456 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1261 and further comprising one or more primer pairs wherein each member of said one or more primer pairs is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from the group of primer pairs represented by SEQ ID NOs: 288:1269, 698:1420, 217:1167, 399:1041, 430:1321, 174:853, and 172:1360.

[35] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 681.

[36] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1022.

[37] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 681 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1022.

[38] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 315.

[39] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1379.

[40] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 315 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1379.

[41] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 346.

[42] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 955.

[43] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 346 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 955.

[44] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 504.

[45] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1409.

[46] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 504 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1409.

[47] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 323.

[48] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1068.

[49] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 323 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1068.

[50] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 479.

[51] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 938.

[52] Another embodiment is an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 479 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 938.

[53] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 681 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1022.

[54] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 681 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1022 and further comprising one or more primer pairs wherein each member of said one or more primer pairs is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from the group of primer pairs represented by SEQ ID NOs: 315:1379, 346:955, 504:1409, 323:1068, 479:938.

[55] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 583.

[56] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 923.

[57] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 583 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 923.

[58] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 454.

[59] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1418.

[60] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 454 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1418.

[61] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 250.

[62] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 902.

[63] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 250 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 902.

[64] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 384.

[65] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 878.

[66] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 384 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 878.

[67] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 694.

[68] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1215.

[69] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 694 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1215.

[70] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 194.

[71] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1173.

[72] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 194 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1173.

[73] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 375.

[74] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 890.

[75] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 375 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 890.

[76] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 656.

[77] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1224.

[78] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID

NO: 656 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1224.

[79] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 618.

[80] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1157.

[81] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 618 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1157.

[82] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 302.

[83] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 852.

[84] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 302 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 852.

[85] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 199.

[86] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 889.

[87] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 199 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 889.

[88] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 596.

[89] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1169.

[90] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 596 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1169.

[91] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 150.

[92] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1242.

[93] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 150 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1242.

[94] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 166.

[95] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1069.

[96] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 166 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1069.

[97] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 166.

[98] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1168.

[99] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 166 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1168.

[100] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 583 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 923 and further comprising one or more primer pairs wherein each member of said one or more primer pairs is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from the group of primer pairs represented by SEQ ID NOs: 454:1418, 250:902, 384:878, 694:1215, 194:1173, 375:890, 656:1224, 618:1157, 302:852, 199:889, 596:1169, 150:1242, 166:1069 and 166:1168.

[101] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 437.

[102] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1137.

[103] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 437 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1137.

[104] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 530.

[105] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 891.

[106] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID

NO: 530 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 891.

[107] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 474.

[108] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 869.

[109] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 474 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 869.

[110] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 268.

[111] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1284.

[112] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 268 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1284.

[113] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 418.

[114] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1301.

[115] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 418 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1301.

- [116] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 318.
- [117] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1300.
- [118] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 318 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1300.
- [119] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 440.
- [120] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1076.
- [121] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 440 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1076.
- [122] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 219.
- [123] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1013.
- [124] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 219 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1013.
- [125] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 437 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence

identity with SEQ ID NO: 1137 and further comprising one or more primer pairs wherein each member of said one or more primer pairs is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from the group of primer pairs represented by SEQ ID NOs: 530:891, 474:869, 268:1284, 418:1301, 318:1300, 440:1076 and 219:1013.

[126] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 325.

[127] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1163.

[128] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 325 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1163.

[129] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 278.

[130] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1039.

[131] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 278 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1039.

[132] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 465.

[133] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1037.

[134] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID

NO: 465 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1037.

[135] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 148.

[136] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1172.

[137] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 148 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1172.

[138] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 190.

[139] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1254.

[140] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 190 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1254.

[141] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 266.

[142] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1094.

[143] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 266 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1094.

[144] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 508.

[145] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1297.

[146] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 508 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1297.

[147] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 259.

[148] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1060.

[149] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 259 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1060.

[150] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 325 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1163 and further comprising one or more primer pairs wherein each member of said one or more primer pairs is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from the group of primer pairs represented by SEQ ID NOs: 278:1039: 465:1037, 148:1172, 190:1254, 266:1094, 508:1297 and 259:1060.

[151] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 376.

[152] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1265.

[153] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 376 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1265.

[154] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 267.

[155] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1341.

[156] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 267 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1341.

[157] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 705.

[158] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1056.

[159] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 705 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1056.

[160] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 710.

[161] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1259.

[162] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID

NO: 710 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1259.

[163] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 374.

[164] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1111.

[165] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 374 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1111.

[166] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 545.

[167] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 978.

[168] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 545 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 978.

[169] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 249.

[170] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1095.

[171] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 249 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1095.

- [172] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 195.
- [173] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1376.
- [174] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 195 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1376.
- [175] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 311.
- [176] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1014.
- [177] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 311 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1014.
- [178] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 365.
- [179] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1052.
- [180] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 365 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1052.
- [181] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 527.

[182] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1071.

[183] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 527 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1071.

[184] One embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 490.

[185] Another embodiment is an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1182.

[186] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 490 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1182.

[187] Another embodiment is a kit comprising an oligonucleotide primer pair including an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 376 and an oligonucleotide primer 14 to 35 nucleobases in length having at least 70% sequence identity with SEQ ID NO: 1265 and further comprising one or more primer pairs wherein each member of said one or more primer pairs is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from the group of primer pairs represented by SEQ ID NOs: 267:1341, 705:1056, 710:1259, 374:1111, 545:978, 249:1095, 195:1376, 311:1014, 365:1052, 527:1071 and 490:1182.

[188] In some embodiments, either or both of the primers of a primer pair composition contain at least one modified nucleobase such as 5-propynyluracil or 5-propynylcytosine for example.

[189] In some embodiments, either or both of the primers of the primer pair comprises at least one universal nucleobase such as inosine for example.

[190] In some embodiments, either or both of the primers of the primer pair comprises at least one non-templated T residue on the 5'-end.

[191] In some embodiments, either or both of the primers of the primer pair comprises at least one non-template tag.

[192] In some embodiments, either or both of the primers of the primer pair comprises at least one molecular mass modifying tag.

[193] In some embodiments, the present invention provides primers and compositions comprising pairs of primers, and kits containing the same, and methods for use in identification of bacteria. The primers are designed to produce amplification products of DNA encoding genes that have conserved and variable regions across different subgroups and genotypes of bacteria.

[194] Some embodiments are kits that contain one or more of the primer pair compositions. In some embodiments, each member of the one or more primer pairs of the kit is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from any of the primer pairs listed in Table 2.

[195] Some embodiments of the kits contain at least one calibration polynucleotide for use in quantitation of bacteria in a given sample, and also for use as a positive control for amplification.

[196] Some embodiments of the kits contain at least one anion exchange functional group linked to a magnetic bead.

[197] In some embodiments, the present invention also provides methods for identification of bacteria. Nucleic acid from the bacterium is amplified using the primers described above to obtain an amplification product. The molecular mass of the amplification product is measured. Optionally, the base composition of the amplification product is determined from the molecular mass. The molecular mass or base composition is compared with a plurality of molecular masses or base compositions of known analogous bacterial identifying amplicons, wherein a match between the molecular mass or base composition and a member of the plurality of molecular masses or base compositions identifies the bacterium. In some embodiments, the molecular mass is measured by mass spectrometry in a modality such as electrospray ionization (ESI) time of flight (TOF) mass spectrometry or ESI Fourier transform ion cyclotron resonance (FTICR) mass spectrometry, for example. Other mass spectrometry techniques can also be used to measure the molecular mass of bacterial bioagent identifying amplicons.

[198] In some embodiments, the present invention is also directed to a method for determining the presence or absence of a bacterium in a sample. Nucleic acid from the sample is amplified using the composition described above to obtain an amplification product. The molecular mass of the amplification product is determined. Optionally, the base composition of the amplification product is determined from the molecular mass. The molecular mass or base composition of the amplification product is compared with the known molecular masses or base compositions of one or more known analogous bacterial bioagent identifying amplicons, wherein a match between the molecular mass or base composition of the amplification product and the molecular mass or base composition of one or more known bacterial bioagent identifying amplicons indicates the presence of the bacterium in the sample. In some embodiments, the molecular mass is measured by mass spectrometry.

[199] In some embodiments, the present invention also provides methods for determination of the quantity of an unknown bacterium in a sample. The sample is contacted with the composition described above and a known quantity of a calibration polynucleotide comprising a calibration sequence. Nucleic acid from the unknown bacterium in the sample is concurrently amplified with the composition described above and nucleic acid from the calibration polynucleotide in the sample is concurrently amplified with the composition described above to obtain a first amplification product comprising a bacterial bioagent identifying amplicon and a second amplification product comprising a calibration amplicon. The molecular masses and abundances for the bacterial bioagent identifying amplicon and the calibration amplicon are determined. The bacterial bioagent identifying amplicon is distinguished from the calibration amplicon based on molecular mass and comparison of bacterial bioagent identifying amplicon abundance and calibration amplicon abundance indicates the quantity of bacterium in the sample. In some embodiments, the base composition of the bacterial bioagent identifying amplicon is determined.

[200] In some embodiments, the present invention provides methods for detecting or quantifying bacteria by combining a nucleic acid amplification process with a mass determination process. In some embodiments, such methods identify or otherwise analyze the bacterium by comparing mass information from an amplification product with a calibration or control product. Such methods can be carried out in a highly multiplexed and/or parallel manner allowing for the analysis of as many as 300 samples per 24 hours on a single mass measurement platform. The accuracy of the mass determination methods in some embodiments of the present invention permits allows for the ability to discriminate between different bacteria such as, for example, various genotypes and drug resistant strains of *Staphylococcus aureus*.

BRIEF DESCRIPTION OF THE DRAWINGS

[201] The foregoing summary of the invention, as well as the following detailed description of the invention, is better understood when read in conjunction with the accompanying drawings which are included by way of example and not by way of limitation.

[202] **Figure 1:** process diagram illustrating a representative primer pair selection process.

[203] **Figure 2:** process diagram illustrating an embodiment of the calibration method.

[204] **Figure 3:** common pathogenic bacteria and primer pair coverage. The primer pair number in the upper right hand corner of each polygon indicates that the primer pair can produce a bioagent identifying amplicon for all species within that polygon.

[205] **Figure 4:** a representative 3D diagram of base composition (axes A, G and C) of bioagent identifying amplicons obtained with primer pair number 14 (a precursor of primer pair number 348 which targets 16S rRNA). The diagram indicates that the experimentally determined base compositions of the clinical samples (labeled NHRC samples) closely match the base compositions expected for *Streptococcus pyogenes* and are distinct from the expected base compositions of other organisms.

[206] **Figure 5:** a representative mass spectrum of amplification products indicating the presence of bioagent identifying amplicons of *Streptococcus pyogenes*, *Neisseria meningitidis*, and *Haemophilus influenzae* obtained from amplification of nucleic acid from a clinical sample with primer pair number 349 which targets 23S rRNA. Experimentally determined molecular masses and base compositions for the sense strand of each amplification product are shown.

[207] **Figure 6:** a representative mass spectrum of amplification products representing a bioagent identifying amplicon of *Streptococcus pyogenes*, and a calibration amplicon obtained from amplification of nucleic acid from a clinical sample with primer pair number 356 which targets rplB. The experimentally determined molecular mass and base composition for the sense strand of the *Streptococcus pyogenes* amplification product is shown.

[208] **Figure 7:** a representative mass spectrum of an amplified nucleic acid mixture which contained the Ames strain of *Bacillus anthracis*, a known quantity of combination calibration polynucleotide (SEQ ID NO: 1464), and primer pair number 350 which targets the capC gene on

the virulence plasmid pX02 of *Bacillus anthracis*. Calibration amplicons produced in the amplification reaction are visible in the mass spectrum as indicated and abundance data (peak height) are used to calculate the quantity of the Ames strain of *Bacillus anthracis*.

DEFINITIONS

[209] As used herein, the term "abundance" refers to an amount. The amount may be described in terms of concentration which are common in molecular biology such as "copy number," "pfu or plate-forming unit" which are well known to those with ordinary skill. Concentration may be relative to a known standard or may be absolute.

[210] As used herein, the term "amplifiable nucleic acid" is used in reference to nucleic acids that may be amplified by any amplification method. It is contemplated that "amplifiable nucleic acid" also comprises "sample template."

[211] As used herein the term "amplification" refers to a special case of nucleic acid replication involving template specificity. It is to be contrasted with non-specific template replication (i.e., replication that is template-dependent but not dependent on a specific template). Template specificity is here distinguished from fidelity of replication (i.e., synthesis of the proper polynucleotide sequence) and nucleotide (ribo- or deoxyribo-) specificity. Template specificity is frequently described in terms of "target" specificity. Target sequences are "targets" in the sense that they are sought to be sorted out from other nucleic acid. Amplification techniques have been designed primarily for this sorting out. Template specificity is achieved in most amplification techniques by the choice of enzyme. Amplification enzymes are enzymes that, under conditions they are used, will process only specific sequences of nucleic acid in a heterogeneous mixture of nucleic acid. For example, in the case of Q β replicase, MDV-1 RNA is the specific template for the replicase (D.L. Kacian et al., Proc. Natl. Acad. Sci. USA 69:3038 [1972]). Other nucleic acid will not be replicated by this amplification enzyme. Similarly, in the case of T7 RNA polymerase, this amplification enzyme has a stringent specificity for its own promoters (Chamberlin et al., Nature 228:227 [1970]). In the case of T4 DNA ligase, the enzyme will not ligate the two oligonucleotides or polynucleotides, where there is a mismatch between the oligonucleotide or polynucleotide substrate and the template at the ligation junction (D.Y. Wu and R. B. Wallace, Genomics 4:560 [1989]). Finally, Taq and Pfu polymerases, by virtue of their ability to function at high temperature, are found to display high specificity for the sequences bounded and thus defined by the primers; the high temperature results in thermodynamic conditions that favor primer hybridization with the target sequences and not hybridization with non-target sequences (H.A. Erlich (ed.), PCR Technology, Stockton Press [1989]).

[212] As used herein, the term "amplification reagents" refers to those reagents (deoxyribonucleotide triphosphates, buffer, etc.), needed for amplification, excluding primers, nucleic acid template, and the amplification enzyme. Typically, amplification reagents along with other reaction components are placed and contained in a reaction vessel (test tube, microwell, etc.).

[213] As used herein, the term "analogous" when used in context of comparison of bioagent identifying amplicons indicates that the bioagent identifying amplicons being compared are produced with the same pair of primers. For example, bioagent identifying amplicon "A" and bioagent identifying amplicon "B", produced with the same pair of primers are analogous with respect to each other. Bioagent identifying amplicon "C", produced with a different pair of primers is not analogous to either bioagent identifying amplicon "A" or bioagent identifying amplicon "B".

[214] As used herein, the term "anion exchange functional group" refers to a positively charged functional group capable of binding an anion through an electrostatic interaction. The most well known anion exchange functional groups are the amines, including primary, secondary, tertiary and quaternary amines.

[215] The term "bacteria" or "bacterium" refers to any member of the groups of eubacteria and archaeobacteria.

[216] As used herein, a "base composition" is the exact number of each nucleobase (for example, A, T, C and G) in a segment of nucleic acid. For example, amplification of nucleic acid of *Staphylococcus aureus* strain carrying the lukS-PV gene with primer pair number 2095 (SEQ ID NOs: 456:1261) produces an amplification product 117 nucleobases in length from nucleic acid of the lukS-PV gene that has a base composition of A35 G17 C19 T46 (by convention - with reference to the sense strand of the amplification product). Because the molecular masses of each of the four natural nucleotides and chemical modifications thereof are known (if applicable), a measured molecular mass can be deconvoluted to a list of possible base compositions. Identification of a base composition of a sense strand which is complementary to the corresponding antisense strand in terms of base composition provides a confirmation of the true base composition of an unknown amplification product. For example, the base composition of the antisense strand of the 139 nucleobase amplification product described above is A46 G19 C17 T35.

[217] As used herein, a "base composition probability cloud" is a representation of the diversity in base composition resulting from a variation in sequence that occurs among different isolates of a given

species. The "base composition probability cloud" represents the base composition constraints for each species and is typically visualized using a pseudo four-dimensional plot.

[218] In the context of this invention, a "bioagent" is any organism, cell, or virus, living or dead, or a nucleic acid derived from such an organism, cell or virus. Examples of bioagents include, but are not limited, to cells, (including but not limited to human clinical samples, bacterial cells and other pathogens), viruses, fungi, protists, parasites, and pathogenicity markers (including but not limited to: pathogenicity islands, antibiotic resistance genes, virulence factors, toxin genes and other bioregulating compounds). Samples may be alive or dead or in a vegetative state (for example, vegetative bacteria or spores) and may be encapsulated or bioengineered. In the context of this invention, a "pathogen" is a bioagent which causes a disease or disorder.

[219] As used herein, a "bioagent division" is defined as group of bioagents above the species level and includes but is not limited to, orders, families, classes, clades, genera or other such groupings of bioagents above the species level.

[220] As used herein, the term "bioagent identifying amplicon" refers to a polynucleotide that is amplified from a bioagent in an amplification reaction and which 1) provides sufficient variability to distinguish among bioagents from whose nucleic acid the bioagent identifying amplicon is produced and 2) whose molecular mass is amenable to a rapid and convenient molecular mass determination modality such as mass spectrometry, for example.

[221] As used herein, the term "biological product" refers to any product originating from an organism. Biological products are often products of processes of biotechnology. Examples of biological products include, but are not limited to: cultured cell lines, cellular components, antibodies, proteins and other cell-derived biomolecules, growth media, growth harvest fluids, natural products and biopharmaceutical products.

[222] The terms "biowarfare agent" and "bioweapon" are synonymous and refer to a bacterium, virus, fungus or protozoan that could be deployed as a weapon to cause bodily harm to individuals. Military or terrorist groups may be implicated in deployment of biowarfare agents.

[223] In context of this invention, the term "broad range survey primer pair" refers to a primer pair designed to produce bioagent identifying amplicons across different broad groupings of bioagents. For example, the ribosomal RNA-targeted primer pairs are broad range survey primer pairs which have the capability of producing bacterial bioagent identifying amplicons for essentially all known bacteria. With

respect to broad range primer pairs employed for identification of bacteria, a broad range survey primer pair for bacteria such as 16S rRNA primer pair number 346 (SEQ ID NOs: 202:1110) for example, will produce an bacterial bioagent identifying amplicon for essentially all known bacteria.

[224] The term "calibration amplicon" refers to a nucleic acid segment representing an amplification product obtained by amplification of a calibration sequence with a pair of primers designed to produce a bioagent identifying amplicon.

[225] The term "calibration sequence" refers to a polynucleotide sequence to which a given pair of primers hybridizes for the purpose of producing an internal (i.e: included in the reaction) calibration standard amplification product for use in determining the quantity of a bioagent in a sample. The calibration sequence may be expressly added to an amplification reaction, or may already be present in the sample prior to analysis.

[226] The term "clade primer pair" refers to a primer pair designed to produce bioagent identifying amplicons for species belonging to a clade group. A clade primer pair may also be considered as a "speciating" primer pair which is useful for distinguishing among closely related species.

[227] The term "codon" refers to a set of three adjoined nucleotides (triplet) that codes for an amino acid or a termination signal.

[228] In context of this invention, the term "codon base composition analysis," refers to determination of the base composition of an individual codon by obtaining a bioagent identifying amplicon that includes the codon. The bioagent identifying amplicon will at least include regions of the target nucleic acid sequence to which the primers hybridize for generation of the bioagent identifying amplicon as well as the codon being analyzed, located between the two primer hybridization regions.

[229] As used herein, the terms "complementary" or "complementarity" are used in reference to polynucleotides (i.e., a sequence of nucleotides such as an oligonucleotide or a target nucleic acid) related by the base-pairing rules. For example, for the sequence "5'-A-G-T-3'," is complementary to the sequence "3'-T-C-A-5'." Complementarity may be "partial," in which only some of the nucleic acids' bases are matched according to the base pairing rules. Or, there may be "complete" or "total" complementarity between the nucleic acids. The degree of complementarity between nucleic acid strands has significant effects on the efficiency and strength of hybridization between nucleic acid strands. This is of particular importance in amplification reactions, as well as detection methods that depend upon binding between nucleic acids. Either term may also be used in reference to individual nucleotides, especially within the

context of polynucleotides. For example, a particular nucleotide within an oligonucleotide may be noted for its complementarity, or lack thereof, to a nucleotide within another nucleic acid strand, in contrast or comparison to the complementarity between the rest of the oligonucleotide and the nucleic acid strand.

[230] The term "complement of a nucleic acid sequence" as used herein refers to an oligonucleotide which, when aligned with the nucleic acid sequence such that the 5' end of one sequence is paired with the 3' end of the other, is in "antiparallel association." Certain bases not commonly found in natural nucleic acids may be included in the nucleic acids of the present invention and include, for example, inosine and 7-deazaguanine. Complementarity need not be perfect; stable duplexes may contain mismatched base pairs or unmatched bases. Those skilled in the art of nucleic acid technology can determine duplex stability empirically considering a number of variables including, for example, the length of the oligonucleotide, base composition and sequence of the oligonucleotide, ionic strength and incidence of mismatched base pairs. Where a first oligonucleotide is complementary to a region of a target nucleic acid and a second oligonucleotide has complementary to the same region (or a portion of this region) a "region of overlap" exists along the target nucleic acid. The degree of overlap will vary depending upon the extent of the complementarity.

[231] In context of this invention, the term "division-wide primer pair" refers to a primer pair designed to produce bioagent identifying amplicons within sections of a broader spectrum of bioagents. For example, primer pair number 352 (SEQ ID NOs: 687:1411), a division-wide primer pair, is designed to produce bacterial bioagent identifying amplicons for members of the *Bacillus* group of bacteria which comprises, for example, members of the genera *Streptococci*, *Enterococci*, and *Staphylococci*. Other division-wide primer pairs may be used to produce bacterial bioagent identifying amplicons for other groups of bacterial bioagents.

[232] As used herein, the term "concurrently amplifying" used with respect to more than one amplification reaction refers to the act of simultaneously amplifying more than one nucleic acid in a single reaction mixture.

[233] As used herein, the term "drill-down primer pair" refers to a primer pair designed to produce bioagent identifying amplicons for identification of sub-species characteristics or confirmation of a species assignment. For example, primer pair number 2146 (SEQ ID NOs: 437:1137), a drill-down *Staphylococcus aureus* genotyping primer pair, is designed to produce *Staphylococcus aureus* genotyping amplicons. Other drill-down primer pairs may be used to produce bioagent identifying amplicons for *Staphylococcus aureus* and other bacterial species.

[234] The term "duplex" refers to the state of nucleic acids in which the base portions of the nucleotides on one strand are bound through hydrogen bonding to their complementary bases arrayed on a second strand. The condition of being in a duplex form reflects on the state of the bases of a nucleic acid. By virtue of base pairing, the strands of nucleic acid also generally assume the tertiary structure of a double helix, having a major and a minor groove. The assumption of the helical form is implicit in the act of becoming duplexed.

[235] As used herein, the term "etiology" refers to the causes or origins, of diseases or abnormal physiological conditions.

[236] The term "gene" refers to a DNA sequence that comprises control and coding sequences necessary for the production of an RNA having a non-coding function (e.g., a ribosomal or transfer RNA), a polypeptide or a precursor. The RNA or polypeptide can be encoded by a full length coding sequence or by any portion of the coding sequence so long as the desired activity or function is retained.

[237] The terms "homology," "homologous" and "sequence identity" refer to a degree of identity. There may be partial homology or complete homology. A partially homologous sequence is one that is less than 100% identical to another sequence. Determination of sequence identity is described in the following example: a primer 20 nucleobases in length which is otherwise identical to another 20 nucleobase primer but having two non-identical residues has 18 of 20 identical residues ($18/20 = 0.9$ or 90% sequence identity). In another example, a primer 15 nucleobases in length having all residues identical to a 15 nucleobase segment of a primer 20 nucleobases in length would have $15/20 = 0.75$ or 75% sequence identity with the 20 nucleobase primer. In context of the present invention, sequence identity is meant to be properly determined when the query sequence and the subject sequence are both described and aligned in the 5' to 3' direction. Sequence alignment algorithms such as BLAST, will return results in two different alignment orientations. In the Plus/Plus orientation, both the query sequence and the subject sequence are aligned in the 5' to 3' direction. On the other hand, in the Plus/Minus orientation, the query sequence is in the 5' to 3' direction while the subject sequence is in the 3' to 5' direction. It should be understood that with respect to the primers of the present invention, sequence identity is properly determined when the alignment is designated as Plus/Plus. Sequence identity may also encompass alternate or modified nucleobases that perform in a functionally similar manner to the regular nucleobases adenine, thymine, guanine and cytosine with respect to hybridization and primer extension in amplification reactions. In a non-limiting example, if the 5-propynyl pyrimidines propyne C and/or propyne T replace one or more C or T residues in one primer which is otherwise identical to another primer in sequence and length, the two primers will have 100% sequence identity with each other. In another non-limiting example, Inosine (I) may be used as a replacement for G or T and effectively

hybridize to C, A or U (uracil). Thus, if inosine replaces one or more C, A or U residues in one primer which is otherwise identical to another primer in sequence and length, the two primers will have 100% sequence identity with each other. Other such modified or universal bases may exist which would perform in a functionally similar manner for hybridization and amplification reactions and will be understood to fall within this definition of sequence identity.

[238] As used herein, "housekeeping gene" refers to a gene encoding a protein or RNA involved in basic functions required for survival and reproduction of a bioagent. Housekeeping genes include, but are not limited to genes encoding RNA or proteins involved in translation, replication, recombination and repair, transcription, nucleotide metabolism, amino acid metabolism, lipid metabolism, energy generation, uptake, secretion and the like.

[239] As used herein, the term "hybridization" is used in reference to the pairing of complementary nucleic acids. Hybridization and the strength of hybridization (i.e., the strength of the association between the nucleic acids) is influenced by such factors as the degree of complementarity between the nucleic acids, stringency of the conditions involved, and the T_m of the formed hybrid. "Hybridization" methods involve the annealing of one nucleic acid to another, complementary nucleic acid, i.e., a nucleic acid having a complementary nucleotide sequence. The ability of two polymers of nucleic acid containing complementary sequences to find each other and anneal through base pairing interaction is a well-recognized phenomenon. The initial observations of the "hybridization" process by Marmur and Lane, Proc. Natl. Acad. Sci. USA 46:453 (1960) and Doty et al., Proc. Natl. Acad. Sci. USA 46:461 (1960) have been followed by the refinement of this process into an essential tool of modern biology.

[240] The term "*in silico*" refers to processes taking place via computer calculations. For example, electronic PCR (ePCR) is a process analogous to ordinary PCR except that it is carried out using nucleic acid sequences and primer pair sequences stored on a computer formatted medium.

[241] As used herein, "intelligent primers" are primers that are designed to bind to highly conserved sequence regions of a bioagent identifying amplicon that flank an intervening variable region and, upon amplification, yield amplification products which ideally provide enough variability to distinguish individual bioagents, and which are amenable to molecular mass analysis. By the term "highly conserved," it is meant that the sequence regions exhibit between about 80-100%, or between about 90-100%, or between about 95-100% identity among all, or at least 70%, at least 80%, at least 90%, at least 95%, or at least 99% of species or strains.

[242] The "ligase chain reaction" (LCR; sometimes referred to as "Ligase Amplification Reaction" (LAR) described by Barany, Proc. Natl. Acad. Sci., 88:189 (1991); Barany, PCR Methods and Applic., 1:5 (1991); and Wu and Wallace, Genomics 4:560 (1989) has developed into a well-recognized alternative method for amplifying nucleic acids. In LCR, four oligonucleotides, two adjacent oligonucleotides which uniquely hybridize to one strand of target DNA, and a complementary set of adjacent oligonucleotides, that hybridize to the opposite strand are mixed and DNA ligase is added to the mixture. Provided that there is complete complementarity at the junction, ligase will covalently link each set of hybridized molecules. Importantly, in LCR, two probes are ligated together only when they base-pair with sequences in the target sample, without gaps or mismatches. Repeated cycles of denaturation, hybridization and ligation amplify a short segment of DNA. LCR has also been used in combination with PCR to achieve enhanced detection of single-base changes. However, because the four oligonucleotides used in this assay can pair to form two short ligatable fragments, there is the potential for the generation of target-independent background signal. The use of LCR for mutant screening is limited to the examination of specific nucleic acid positions.

[243] The term "locked nucleic acid" or "LNA" refers to a nucleic acid analogue containing one or more 2'-O-, 4'-C-methylene- β -D-ribofuranosyl nucleotide monomers in an RNA mimicking sugar conformation. LNA oligonucleotides display unprecedented hybridization affinity toward complementary single-stranded RNA and complementary single- or double-stranded DNA. LNA oligonucleotides induce A-type (RNA-like) duplex conformations. The primers of the present invention may contain LNA modifications.

[244] As used herein, the term "mass-modifying tag" refers to any modification to a given nucleotide which results in an increase in mass relative to the analogous non-mass modified nucleotide. Mass-modifying tags can include heavy isotopes of one or more elements included in the nucleotide such as carbon-13 for example. Other possible modifications include addition of substituents such as iodine or bromine at the 5 position of the nucleobase for example.

[245] The term "mass spectrometry" refers to measurement of the mass of atoms or molecules. The molecules are first converted to ions, which are separated using electric or magnetic fields according to the ratio of their mass to electric charge. The measured masses are used to identify the molecules.

[246] The term "microorganism" as used herein means an organism too small to be observed with the unaided eye and includes, but is not limited to bacteria, virus, protozoans, fungi; and ciliates.

[247] The term "multi-drug resistant" or "multiple-drug resistant" refers to a microorganism which is resistant to more than one of the antibiotics or antimicrobial agents used in the treatment of said microorganism.

[248] The term "multiplex PCR" refers to a PCR reaction where more than one primer set is included in the reaction pool allowing 2 or more different DNA targets to be amplified by PCR in a single reaction tube.

[249] The term "non-template tag" refers to a stretch of at least three guanine or cytosine nucleobases of a primer used to produce a bioagent identifying amplicon which are not complementary to the template. A non-template tag is incorporated into a primer for the purpose of increasing the primer-duplex stability of later cycles of amplification by incorporation of extra G-C pairs which each have one additional hydrogen bond relative to an A-T pair.

[250] The term "nucleic acid sequence" as used herein refers to the linear composition of the nucleic acid residues A, T, C or G or any modifications thereof, within an oligonucleotide, nucleotide or polynucleotide, and fragments or portions thereof, and to DNA or RNA of genomic or synthetic origin which may be single or double stranded, and represent the sense or antisense strand

[251] As used herein, the term "nucleobase" is synonymous with other terms in use in the art including "nucleotide," "deoxynucleotide," "nucleotide residue," "deoxynucleotide residue," "nucleotide triphosphate (NTP)," or deoxynucleotide triphosphate (dNTP).

[252] The term "nucleotide analog" as used herein refers to modified or non-naturally occurring nucleotides such as 5-propynyl pyrimidines (i.e., 5-propynyl-dTTP and 5-propynyl-dTCP), 7-deaza purines (i.e., 7-deaza-dATP and 7-deaza-dGTP). Nucleotide analogs include base analogs and comprise modified forms of deoxyribonucleotides as well as ribonucleotides.

[253] The term "oligonucleotide" as used herein is defined as a molecule comprising two or more deoxyribonucleotides or ribonucleotides, preferably at least 5 nucleotides, more preferably at least about 13 to 35 nucleotides. The exact size will depend on many factors, which in turn depend on the ultimate function or use of the oligonucleotide. The oligonucleotide may be generated in any manner, including chemical synthesis, DNA replication, reverse transcription, PCR, or a combination thereof. Because mononucleotides are reacted to make oligonucleotides in a manner such that the 5' phosphate of one mononucleotide pentose ring is attached to the 3' oxygen of its neighbor in one direction via a phosphodiester linkage, an end of an oligonucleotide is referred to as the "5'-end" if its 5' phosphate is not

linked to the 3' oxygen of a mononucleotide pentose ring and as the "3'-end" if its 3' oxygen is not linked to a 5' phosphate of a subsequent mononucleotide pentose ring. As used herein, a nucleic acid sequence, even if internal to a larger oligonucleotide, also may be said to have 5' and 3' ends. A first region along a nucleic acid strand is said to be upstream of another region if the 3' end of the first region is before the 5' end of the second region when moving along a strand of nucleic acid in a 5' to 3' direction. All oligonucleotide primers disclosed herein are understood to be presented in the 5' to 3' direction when reading left to right. When two different, non-overlapping oligonucleotides anneal to different regions of the same linear complementary nucleic acid sequence, and the 3' end of one oligonucleotide points towards the 5' end of the other, the former may be called the "upstream" oligonucleotide and the latter the "downstream" oligonucleotide. Similarly, when two overlapping oligonucleotides are hybridized to the same linear complementary nucleic acid sequence, with the first oligonucleotide positioned such that its 5' end is upstream of the 5' end of the second oligonucleotide, and the 3' end of the first oligonucleotide is upstream of the 3' end of the second oligonucleotide, the first oligonucleotide may be called the "upstream" oligonucleotide and the second oligonucleotide may be called the "downstream" oligonucleotide.

[254] In the context of this invention, a "pathogen" is a bioagent which causes a disease or disorder.

[255] As used herein, the terms "PCR product," "PCR fragment," and "amplification product" refer to the resultant mixture of compounds after two or more cycles of the PCR steps of denaturation, annealing and extension are complete. These terms encompass the case where there has been amplification of one or more segments of one or more target sequences.

[256] The term "peptide nucleic acid" ("PNA") as used herein refers to a molecule comprising bases or base analogs such as would be found in natural nucleic acid, but attached to a peptide backbone rather than the sugar-phosphate backbone typical of nucleic acids. The attachment of the bases to the peptide is such as to allow the bases to base pair with complementary bases of nucleic acid in a manner similar to that of an oligonucleotide. These small molecules, also designated anti gene agents, stop transcript elongation by binding to their complementary strand of nucleic acid (Nielsen, et al. *Anticancer Drug Des.* 8:53-63). The primers of the present invention may comprise PNAs.

[257] The term "polymerase" refers to an enzyme having the ability to synthesize a complementary strand of nucleic acid from a starting template nucleic acid strand and free dNTPs.

[258] As used herein, the term "polymerase chain reaction" ("PCR") refers to the method of K.B. Mullis U.S. Patent Nos. 4,683,195, 4,683,202, and 4,965,188, hereby incorporated by reference, that

describe a method for increasing the concentration of a segment of a target sequence in a mixture of genomic DNA without cloning or purification. This process for amplifying the target sequence consists of introducing a large excess of two oligonucleotide primers to the DNA mixture containing the desired target sequence, followed by a precise sequence of thermal cycling in the presence of a DNA polymerase. The two primers are complementary to their respective strands of the double stranded target sequence. To effect amplification, the mixture is denatured and the primers then annealed to their complementary sequences within the target molecule. Following annealing, the primers are extended with a polymerase so as to form a new pair of complementary strands. The steps of denaturation, primer annealing, and polymerase extension can be repeated many times (i.e., denaturation, annealing and extension constitute one "cycle"; there can be numerous "cycles") to obtain a high concentration of an amplified segment of the desired target sequence. The length of the amplified segment of the desired target sequence is determined by the relative positions of the primers with respect to each other, and therefore, this length is a controllable parameter. By virtue of the repeating aspect of the process, the method is referred to as the "polymerase chain reaction" (hereinafter "PCR"). Because the desired amplified segments of the target sequence become the predominant sequences (in terms of concentration) in the mixture, they are said to be "PCR amplified." With PCR, it is possible to amplify a single copy of a specific target sequence in genomic DNA to a level detectable by several different methodologies (e.g., hybridization with a labeled probe; incorporation of biotinylated primers followed by avidin-enzyme conjugate detection; incorporation of ^{32}P -labeled deoxynucleotide triphosphates, such as dCTP or dATP, into the amplified segment). In addition to genomic DNA, any oligonucleotide or polynucleotide sequence can be amplified with the appropriate set of primer molecules. In particular, the amplified segments created by the PCR process itself are, themselves, efficient templates for subsequent PCR amplifications.

[259] The term "polymerization means" or "polymerization agent" refers to any agent capable of facilitating the addition of nucleoside triphosphates to an oligonucleotide. Preferred polymerization means comprise DNA and RNA polymerases.

[260] As used herein, the terms "pair of primers," or "primer pair" are synonymous. A primer pair is used for amplification of a nucleic acid sequence. A pair of primers comprises a forward primer and a reverse primer. The forward primer hybridizes to a sense strand of a target gene sequence to be amplified and primes synthesis of an antisense strand (complementary to the sense strand) using the target sequence as a template. A reverse primer hybridizes to the antisense strand of a target gene sequence to be amplified and primes synthesis of a sense strand (complementary to the antisense strand) using the target sequence as a template.

[261] The primers are designed to bind to highly conserved sequence regions of a bioagent identifying amplicon that flank an intervening variable region and yield amplification products which ideally provide enough variability to distinguish each individual bioagent, and which are amenable to molecular mass analysis. In some embodiments, the highly conserved sequence regions exhibit between about 80-100%, or between about 90-100%, or between about 95-100% identity, or between about 99-100% identity. The molecular mass of a given amplification product provides a means of identifying the bioagent from which it was obtained, due to the variability of the variable region. Thus design of the primers requires selection of a variable region with appropriate variability to resolve the identity of a given bioagent. Bioagent identifying amplicons are ideally specific to the identity of the bioagent.

[262] Properties of the primers may include any number of properties related to structure including, but not limited to: nucleobase length which may be contiguous (linked together) or non-contiguous (for example, two or more contiguous segments which are joined by a linker or loop moiety), modified or universal nucleobases (used for specific purposes such as for example, increasing hybridization affinity, preventing non-templated adenylation and modifying molecular mass) percent complementarity to a given target sequences.

[263] Properties of the primers also include functional features including, but not limited to, orientation of hybridization (forward or reverse) relative to a nucleic acid template. The coding or sense strand is the strand to which the forward priming primer hybridizes (forward priming orientation) while the reverse priming primer hybridizes to the non-coding or antisense strand (reverse priming orientation). The functional properties of a given primer pair also include the generic template nucleic acid to which the primer pair hybridizes. For example, identification of bioagents can be accomplished at different levels using primers suited to resolution of each individual level of identification. Broad range survey primers are designed with the objective of identifying a bioagent as a member of a particular division (e.g., an order, family, genus or other such grouping of bioagents above the species level of bioagents). In some embodiments, broad range survey intelligent primers are capable of identification of bioagents at the species or sub-species level. Other primers may have the functionality of producing bioagent identifying amplicons for members of a given taxonomic genus, clade, species, sub-species or genotype (including genetic variants which may include presence of virulence genes or antibiotic resistance genes or mutations). Additional functional properties of primer pairs include the functionality of performing amplification either singly (single primer pair per amplification reaction vessel) or in a multiplex fashion (multiple primer pairs and multiple amplification reactions within a single reaction vessel).

[264] As used herein, the terms "purified" or "substantially purified" refer to molecules, either nucleic or amino acid sequences, that are removed from their natural environment, isolated or separated,

and are at least 60% free, preferably 75% free, and most preferably 90% free from other components with which they are naturally associated. An "isolated polynucleotide" or "isolated oligonucleotide" is therefore a substantially purified polynucleotide.

[265] The term "reverse transcriptase" refers to an enzyme having the ability to transcribe DNA from an RNA template. This enzymatic activity is known as reverse transcriptase activity. Reverse transcriptase activity is desirable in order to obtain DNA from RNA viruses which can then be amplified and analyzed by the methods of the present invention.

[266] The term "ribosomal RNA" or "rRNA" refers to the primary ribonucleic acid constituent of ribosomes. Ribosomes are the protein-manufacturing organelles of cells and exist in the cytoplasm. Ribosomal RNAs are transcribed from the DNA genes encoding them.

[267] The term "sample" in the present specification and claims is used in its broadest sense. On the one hand it is meant to include a specimen or culture (e.g., microbiological cultures). On the other hand, it is meant to include both biological and environmental samples. A sample may include a specimen of synthetic origin. Biological samples may be animal, including human, fluid, solid (e.g., stool) or tissue, as well as liquid and solid food and feed products and ingredients such as dairy items, vegetables, meat and meat by-products, and waste. Biological samples may be obtained from all of the various families of domestic animals, as well as feral or wild animals, including, but not limited to, such animals as ungulates, bear, fish, lagomorphs, rodents, etc. Environmental samples include environmental material such as surface matter, soil, water, air and industrial samples, as well as samples obtained from food and dairy processing instruments, apparatus, equipment, utensils, disposable and non-disposable items. These examples are not to be construed as limiting the sample types applicable to the present invention. The term "source of target nucleic acid" refers to any sample that contains nucleic acids (RNA or DNA). Particularly preferred sources of target nucleic acids are biological samples including, but not limited to blood, saliva, cerebral spinal fluid, pleural fluid, milk, lymph, sputum and semen.

[268] As used herein, the term "sample template" refers to nucleic acid originating from a sample that is analyzed for the presence of "target" (defined below). In contrast, "background template" is used in reference to nucleic acid other than sample template that may or may not be present in a sample. Background template is often a contaminant. It may be the result of carryover, or it may be due to the presence of nucleic acid contaminants sought to be purified away from the sample. For example, nucleic acids from organisms other than those to be detected may be present as background in a test sample.

[269] A "segment" is defined herein as a region of nucleic acid within a target sequence.

[270] The "self-sustained sequence replication reaction" (3SR) (Guatelli et al., Proc. Natl. Acad. Sci., 87:1874-1878 [1990], with an erratum at Proc. Natl. Acad. Sci., 87:7797 [1990]) is a transcription-based *in vitro* amplification system (Kwok et al., Proc. Natl. Acad. Sci., 86:1173-1177 [1989]) that can exponentially amplify RNA sequences at a uniform temperature. The amplified RNA can then be utilized for mutation detection (Fahy et al., PCR Meth. Appl., 1:25-33 [1991]). In this method, an oligonucleotide primer is used to add a phage RNA polymerase promoter to the 5' end of the sequence of interest. In a cocktail of enzymes and substrates that includes a second primer, reverse transcriptase, RNase H, RNA polymerase and ribo- and deoxyribonucleoside triphosphates, the target sequence undergoes repeated rounds of transcription, cDNA synthesis and second-strand synthesis to amplify the area of interest. The use of 3SR to detect mutations is kinetically limited to screening small segments of DNA (e.g., 200-300 base pairs).

[271] As used herein, the term "sequence alignment" refers to a listing of multiple DNA or amino acid sequences and aligns them to highlight their similarities. The listings can be made using bioinformatics computer programs.

[272] In context of this invention, the term "speciating primer pair" refers to a primer pair designed to produce a bioagent identifying amplicon with the diagnostic capability of identifying species members of a group of genera or a particular genus of bioagents. Primer pair number 2249 (SEQ ID NOs: 430:1321), for example, is a speciating primer pair used to distinguish *Staphylococcus aureus* from other species of the genus *Staphylococcus*.

[273] As used herein, a "sub-species characteristic" is a genetic characteristic that provides the means to distinguish two members of the same bioagent species. For example, one viral strain could be distinguished from another viral strain of the same species by possessing a genetic change (e.g., for example, a nucleotide deletion, addition or substitution) in one of the viral genes, such as the RNA-dependent RNA polymerase. Sub-species characteristics such as virulence genes and drug-resistance are responsible for the phenotypic differences among the different strains of bacteria.

[274] As used herein, the term "target" is used in a broad sense to indicate the gene or genomic region being amplified by the primers. Because the present invention provides a plurality of amplification products from any given primer pair (depending on the bioagent being analyzed), multiple amplification products from different specific nucleic acid sequences may be obtained. Thus, the term "target" is not used to refer to a single specific nucleic acid sequence. The "target" is sought to be sorted out from other nucleic acid sequences and contains a sequence that has at least partial complementarity with an

oligonucleotide primer. The target nucleic acid may comprise single- or double-stranded DNA or RNA. A "segment" is defined as a region of nucleic acid within the target sequence.

[275] The term "template" refers to a strand of nucleic acid on which a complementary copy is built from nucleoside triphosphates through the activity of a template-dependent nucleic acid polymerase. Within a duplex the template strand is, by convention, depicted and described as the "bottom" strand. Similarly, the non-template strand is often depicted and described as the "top" strand.

[276] As used herein, the term " T_m " is used in reference to the "melting temperature." The melting temperature is the temperature at which a population of double-stranded nucleic acid molecules becomes half dissociated into single strands. Several equations for calculating the T_m of nucleic acids are well known in the art. As indicated by standard references, a simple estimate of the T_m value may be calculated by the equation: $T_m = 81.5 + 0.41(\% \text{ G+C})$, when a nucleic acid is in aqueous solution at 1 M NaCl (see e.g., Anderson and Young, Quantitative Filter Hybridization, in Nucleic Acid Hybridization (1985). Other references (e.g., Allawi, H. T. & SantaLucia, J., Jr. Thermodynamics and NMR of internal G.T mismatches in DNA. Biochemistry 36, 10581-94 (1997) include more sophisticated computations which take structural and environmental, as well as sequence characteristics into account for the calculation of T_m .

[277] The term "triangulation genotyping analysis" refers to a method of genotyping a bioagent by measurement of molecular masses or base compositions of amplification products, corresponding to bioagent identifying amplicons, obtained by amplification of regions of more than one gene. In this sense, the term "triangulation" refers to a method of establishing the accuracy of information by comparing three or more types of independent points of view bearing on the same findings. Triangulation genotyping analysis carried out with a plurality of triangulation genotyping analysis primers yields a plurality of base compositions that then provide a pattern or "barcode" from which a species type can be assigned. The species type may represent a previously known sub-species or strain, or may be a previously unknown strain having a specific and previously unobserved base composition barcode indicating the existence of a previously unknown genotype.

[278] As used herein, the term "triangulation genotyping analysis primer pair" is a primer pair designed to produce bioagent identifying amplicons for determining species types in a triangulation genotyping analysis.

[279] The employment of more than one bioagent identifying amplicon for identification of a bioagent is herein referred to as "triangulation identification." Triangulation identification is pursued by

analyzing a plurality of bioagent identifying amplicons produced with different primer pairs. This process is used to reduce false negative and false positive signals, and enable reconstruction of the origin of hybrid or otherwise engineered bioagents. For example, identification of the three part toxin genes typical of *B. anthracis* (Bowen et al., J. Appl. Microbiol., 1999, 87, 270-278) in the absence of the expected signatures from the *B. anthracis* genome would suggest a genetic engineering event.

[280] In the context of this invention, the term "unknown bioagent" may mean either: (i) a bioagent whose existence is known (such as the well known bacterial species *Staphylococcus aureus* for example) but which is not known to be in a sample to be analyzed, or (ii) a bioagent whose existence is not known (for example, the SARS coronavirus was unknown prior to April 2003). For example, if the method for identification of coronaviruses disclosed in commonly owned U.S. Patent Serial No. 10/829,826 (incorporated herein by reference in its entirety) was to be employed prior to April 2003 to identify the SARS coronavirus in a clinical sample, both meanings of "unknown" bioagent are applicable since the SARS coronavirus was unknown to science prior to April, 2003 and since it was not known what bioagent (in this case a coronavirus) was present in the sample. On the other hand, if the method of U.S. Patent Serial No. 10/829,826 was to be employed subsequent to April 2003 to identify the SARS coronavirus in a clinical sample, only the first meaning (i) of "unknown" bioagent would apply since the SARS coronavirus became known to science subsequent to April 2003 and since it was not known what bioagent was present in the sample.

[281] The term "variable sequence" as used herein refers to differences in nucleic acid sequence between two nucleic acids. For example, the genes of two different bacterial species may vary in sequence by the presence of single base substitutions and/or deletions or insertions of one or more nucleotides. These two forms of the structural gene are said to vary in sequence from one another. In the context of the present invention, "viral nucleic acid" includes, but is not limited to, DNA, RNA, or DNA that has been obtained from viral RNA, such as, for example, by performing a reverse transcription reaction. Viral RNA can either be single-stranded (of positive or negative polarity) or double-stranded.

[282] The term "virus" refers to obligate, ultramicroscopic, parasites that are incapable of autonomous replication (i.e., replication requires the use of the host cell's machinery). Viruses can survive outside of a host cell but cannot replicate.

[283] The term "wild-type" refers to a gene or a gene product that has the characteristics of that gene or gene product when isolated from a naturally occurring source. A wild-type gene is that which is most frequently observed in a population and is thus arbitrarily designated the "normal" or "wild-type" form of the gene. In contrast, the term "modified", "mutant" or "polymorphic" refers to a gene or gene product

that displays modifications in sequence and or functional properties (i.e., altered characteristics) when compared to the wild-type gene or gene product. It is noted that naturally-occurring mutants can be isolated; these are identified by the fact that they have altered characteristics when compared to the wild-type gene or gene product.

[284] As used herein, a “wobble base” is a variation in a codon found at the third nucleotide position of a DNA triplet. Variations in conserved regions of sequence are often found at the third nucleotide position due to redundancy in the amino acid code.

DETAILED DESCRIPTION OF EMBODIMENTS

A. Bioagent Identifying Amplicons

[285] The present invention provides methods for detection and identification of unknown bioagents using bioagent identifying amplicons. Primers are selected to hybridize to conserved sequence regions of nucleic acids derived from a bioagent, and which bracket variable sequence regions to yield a bioagent identifying amplicon, which can be amplified and which is amenable to molecular mass determination. The molecular mass then provides a means to uniquely identify the bioagent without a requirement for prior knowledge of the possible identity of the bioagent. The molecular mass or corresponding base composition signature of the amplification product is then matched against a database of molecular masses or base composition signatures. A match is obtained when an experimentally-determined molecular mass or base composition of an analyzed amplification product is compared with known molecular masses or base compositions of known bioagent identifying amplicons and the experimentally determined molecular mass or base composition is the same as the molecular mass or base composition of one of the known bioagent identifying amplicons. Alternatively, the experimentally-determined molecular mass or base composition may be within experimental error of the molecular mass or base composition of a known bioagent identifying amplicon and still be classified as a match. In some cases, the match may also be classified using a probability of match model such as the models described in U.S. Serial No. 11/073,362, which is commonly owned and incorporated herein by reference in entirety. Furthermore, the method can be applied to rapid parallel multiplex analyses, the results of which can be employed in a triangulation identification strategy. The present method provides rapid throughput and does not require nucleic acid sequencing of the amplified target sequence for bioagent detection and identification.

[286] Despite enormous biological diversity, all forms of life on earth share sets of essential, common features in their genomes. Since genetic data provide the underlying basis for identification of bioagents by the methods of the present invention, it is necessary to select segments of nucleic acids which ideally provide enough variability to distinguish each individual bioagent and whose molecular mass is amenable to molecular mass determination.

[287] Unlike bacterial genomes, which exhibit conservation of numerous genes (i.e. housekeeping genes) across all organisms, viruses do not share a gene that is essential and conserved among all virus families. Therefore, viral identification is achieved within smaller groups of related viruses, such as members of a particular virus family or genus. For example, RNA-dependent RNA polymerase is present in all single-stranded RNA viruses and can be used for broad priming as well as resolution within the virus family.

[288] In some embodiments of the present invention, at least one bacterial nucleic acid segment is amplified in the process of identifying the bacterial bioagent. Thus, the nucleic acid segments that can be amplified by the primers disclosed herein and that provide enough variability to distinguish each individual bioagent and whose molecular masses are amenable to molecular mass determination are herein described as bioagent identifying amplicons.

[289] In some embodiments of the present invention, bioagent identifying amplicons comprise from about 45 to about 150 nucleobases (i.e. from about 45 to about 200 linked nucleosides), although both longer and short regions may be used. One of ordinary skill in the art will appreciate that the invention embodies compounds of 45, 46, 47, 48, 49, 50, 51, 52, 53, 54, 55, 56, 57, 58, 59, 60, 61, 62, 63, 64, 65, 66, 67, 68, 69, 70, 71, 72, 73, 74, 75, 76, 77, 78, 79, 80, 81, 82, 83, 84, 85, 86, 87, 88, 89, 90, 91, 92, 93, 94, 95, 96, 97, 98, 99, 100, 101, 102, 103, 104, 105, 106, 107, 108, 109, 110, 111, 112, 113, 114, 115, 116, 117, 118, 119, 120, 121, 122, 123, 124, 125, 126, 127, 128, 129, 130, 131, 132, 133, 134, 135, 136, 137, 138, 139, 140, 141, 142, 143, 144, 145, 146, 147, 148, 149, and 150 nucleobases in length, or any range therewithin.

[290] It is the combination of the portions of the bioagent nucleic acid segment to which the primers hybridize (hybridization sites) and the variable region between the primer hybridization sites that comprises the bioagent identifying amplicon. Thus, it can be said that a given bioagent identifying amplicon is "defined by" a given pair of primers.

[291] In some embodiments, bioagent identifying amplicons amenable to molecular mass determination which are produced by the primers described herein are either of a length, size or mass compatible with the particular mode of molecular mass determination or compatible with a means of providing a predictable fragmentation pattern in order to obtain predictable fragments of a length compatible with the particular mode of molecular mass determination. Such means of providing a predictable fragmentation pattern of an amplification product include, but are not limited to, cleavage with chemical reagents, restriction enzymes or cleavage primers, for example. Thus, in some

embodiments, bioagent identifying amplicons are larger than 150 nucleobases and are amenable to molecular mass determination following restriction digestion. Methods of using restriction enzymes and cleavage primers are well known to those with ordinary skill in the art.

[292] In some embodiments, amplification products corresponding to bioagent identifying amplicons are obtained using the polymerase chain reaction (PCR) that is a routine method to those with ordinary skill in the molecular biology arts. Other amplification methods may be used such as ligase chain reaction (LCR), low-stringency single primer PCR, and multiple strand displacement amplification (MDA). These methods are also known to those with ordinary skill.

B. Primers and Primer Pairs

[293] In some embodiments, the primers are designed to bind to conserved sequence regions of a bioagent identifying amplicon that flank an intervening variable region and yield amplification products which provide variability sufficient to distinguish each individual bioagent, and which are amenable to molecular mass analysis. In some embodiments, the highly conserved sequence regions exhibit between about 80-100%, or between about 90-100%, or between about 95-100% identity, or between about 99-100% identity. The molecular mass of a given amplification product provides a means of identifying the bioagent from which it was obtained, due to the variability of the variable region. Thus, design of the primers involves selection of a variable region with sufficient variability to resolve the identity of a given bioagent. In some embodiments, bioagent identifying amplicons are specific to the identity of the bioagent.

[294] In some embodiments, identification of bioagents is accomplished at different levels using primers suited to resolution of each individual level of identification. Broad range survey primers are designed with the objective of identifying a bioagent as a member of a particular division (e.g., an order, family, genus or other such grouping of bioagents above the species level of bioagents). In some embodiments, broad range survey intelligent primers are capable of identification of bioagents at the species or sub-species level. Examples of broad range survey primers include, but are not limited to: primer pair numbers: 346 (SEQ ID NOs: 202:1110), 347 (SEQ ID NOs: 560:1278), 348 (SEQ ID NOs: 706:895), and 361 (SEQ ID NOs: 697:1398) which target DNA encoding 16S rRNA, and primer pair numbers 349 (SEQ ID NOs: 401:1156) and 360 (SEQ ID NOs: 409:1434) which target DNA encoding 23S rRNA.

[295] In some embodiments, drill-down primers are designed with the objective of identifying a bioagent at the sub-species level (including strains, subtypes, variants and isolates) based on sub-species characteristics which may, for example, include single nucleotide polymorphisms (SNPs), variable

number tandem repeats (VNTRs), deletions, drug resistance mutations or any other modification of a nucleic acid sequence of a bioagent relative to other members of a species having different sub-species characteristics. Drill-down intelligent primers are not always required for identification at the sub-species level because broad range survey intelligent primers may, in some cases provide sufficient identification resolution to accomplishing this identification objective. Examples of drill-down primers include, but are not limited to: confirmation primer pairs such as primer pair numbers 351 (SEQ ID NOs: 355:1423) and 353 (SEQ ID NOs: 220:1394), which target the pX01 virulence plasmid of *Bacillus anthracis*. Other examples of drill-down primer pairs are found in sets of triangulation genotyping primer pairs such as, for example, the primer pair number 2146 (SEQ ID NOs: 437:1137) which targets the *arcC* gene (encoding carbamate kinase) and is included in an 8 primer pair panel or kit for use in genotyping *Staphylococcus aureus*, or in other panels or kits of primer pairs used for determining drug-resistant bacterial strains, such as, for example, primer pair number 2095 (SEQ ID NOs: 456:1261) which targets the *pv-luk* gene (encoding Panton-Valentine leukocidin) and is included in an 8 primer pair panel or kit for use in identification of drug resistant strains of *Staphylococcus aureus*.

[296] A representative process flow diagram used for primer selection and validation process is outlined in Figure 1. For each group of organisms, candidate target sequences are identified (200) from which nucleotide alignments are created (210) and analyzed (220). Primers are then designed by selecting appropriate priming regions (230) to facilitate the selection of candidate primer pairs (240). The primer pairs are then subjected to *in silico* analysis by electronic PCR (ePCR) (300) wherein bioagent identifying amplicons are obtained from sequence databases such as GenBank or other sequence collections (310) and checked for specificity *in silico* (320). Bioagent identifying amplicons obtained from GenBank sequences (310) can also be analyzed by a probability model which predicts the capability of a given amplicon to identify unknown bioagents such that the base compositions of amplicons with favorable probability scores are then stored in a base composition database (325). Alternatively, base compositions of the bioagent identifying amplicons obtained from the primers and GenBank sequences can be directly entered into the base composition database (330). Candidate primer pairs (240) are validated by testing their ability to hybridize to target nucleic acid by an *in vitro* amplification by a method such as PCR analysis (400) of nucleic acid from a collection of organisms (410). Amplification products thus obtained are analyzed by gel electrophoresis or by mass spectrometry to confirm the sensitivity, specificity and reproducibility of the primers used to obtain the amplification products (420).

[297] Many of the important pathogens, including the organisms of greatest concern as biowarfare agents, have been completely sequenced. This effort has greatly facilitated the design of primers for the detection of unknown bioagents. The combination of broad-range priming with division-wide and drill-down priming has been used very successfully in several applications of the technology, including

environmental surveillance for biowarfare threat agents and clinical sample analysis for medically important pathogens.

[298] Synthesis of primers is well known and routine in the art. The primers may be conveniently and routinely made through the well-known technique of solid phase synthesis. Equipment for such synthesis is sold by several vendors including, for example, Applied Biosystems (Foster City, CA). Any other means for such synthesis known in the art may additionally or alternatively be employed.

[299] In some embodiments primers are employed as compositions for use in methods for identification of bacterial bioagents as follows: a primer pair composition is contacted with nucleic acid (such as, for example, bacterial DNA or DNA reverse transcribed from the rRNA) of an unknown bacterial bioagent. The nucleic acid is then amplified by a nucleic acid amplification technique, such as PCR for example, to obtain an amplification product that represents a bioagent identifying amplicon. The molecular mass of each strand of the double-stranded amplification product is determined by a molecular mass measurement technique such as mass spectrometry for example, wherein the two strands of the double-stranded amplification product are separated during the ionization process. In some embodiments, the mass spectrometry is electrospray Fourier transform ion cyclotron resonance mass spectrometry (ESI-FTICR-MS) or electrospray time of flight mass spectrometry (ESI-TOF-MS). A list of possible base compositions can be generated for the molecular mass value obtained for each strand and the choice of the correct base composition from the list is facilitated by matching the base composition of one strand with a complementary base composition of the other strand. The molecular mass or base composition thus determined is then compared with a database of molecular masses or base compositions of analogous bioagent identifying amplicons for known viral bioagents. A match between the molecular mass or base composition of the amplification product and the molecular mass or base composition of an analogous bioagent identifying amplicon for a known viral bioagent indicates the identity of the unknown bioagent. In some embodiments, the primer pair used is one of the primer pairs of Table 2. In some embodiments, the method is repeated using one or more different primer pairs to resolve possible ambiguities in the identification process or to improve the confidence level for the identification assignment.

[300] In some embodiments, a bioagent identifying amplicon may be produced using only a single primer (either the forward or reverse primer of any given primer pair), provided an appropriate amplification method is chosen, such as, for example, low stringency single primer PCR (LSSP-PCR). Adaptation of this amplification method in order to produce bioagent identifying amplicons can be accomplished by one with ordinary skill in the art without undue experimentation.

[301] In some embodiments, the oligonucleotide primers are broad range survey primers which hybridize to conserved regions of nucleic acid encoding the hexon gene of all (or between 80% and 100%, between 85% and 100%, between 90% and 100% or between 95% and 100%) known bacteria and produce bacterial bioagent identifying amplicons.

[302] In some cases, the molecular mass or base composition of a bacterial bioagent identifying amplicon defined by a broad range survey primer pair does not provide enough resolution to unambiguously identify a bacterial bioagent at or below the species level. These cases benefit from further analysis of one or more bacterial bioagent identifying amplicons generated from at least one additional broad range survey primer pair or from at least one additional division-wide primer pair. The employment of more than one bioagent identifying amplicon for identification of a bioagent is herein referred to as triangulation identification.

[303] In other embodiments, the oligonucleotide primers are division-wide primers which hybridize to nucleic acid encoding genes of species within a genus of bacteria. In other embodiments, the oligonucleotide primers are drill-down primers which enable the identification of sub-species characteristics. Drill down primers provide the functionality of producing bioagent identifying amplicons for drill-down analyses such as strain typing when contacted with nucleic acid under amplification conditions. Identification of such sub-species characteristics is often critical for determining proper clinical treatment of viral infections. In some embodiments, sub-species characteristics are identified using only broad range survey primers and division-wide and drill-down primers are not used.

[304] In some embodiments, the primers used for amplification hybridize to and amplify genomic DNA, and DNA of bacterial plasmids.

[305] In some embodiments, various computer software programs may be used to aid in design of primers for amplification reactions such as *Primer Premier 5* (Premier Biosoft, Palo Alto, CA) or *OLIGO* Primer Analysis Software (Molecular Biology Insights, Cascade, CO). These programs allow the user to input desired hybridization conditions such as melting temperature of a primer-template duplex for example. In some embodiments, an *in silico* PCR search algorithm, such as (ePCR) is used to analyze primer specificity across a plurality of template sequences which can be readily obtained from public sequence databases such as GenBank for example. An existing RNA structure search algorithm (Macke et al., Nucl. Acids Res., 2001, 29, 4724-4735, which is incorporated herein by reference in its entirety) has been modified to include PCR parameters such as hybridization conditions, mismatches, and thermodynamic calculations (SantaLucia, Proc. Natl. Acad. Sci. U.S.A., 1998, 95, 1460-1465, which is incorporated herein by reference in its entirety). This also provides information on primer specificity of

the selected primer pairs. In some embodiments, the hybridization conditions applied to the algorithm can limit the results of primer specificity obtained from the algorithm. In some embodiments, the melting temperature threshold for the primer template duplex is specified to be 35°C or a higher temperature. In some embodiments the number of acceptable mismatches is specified to be seven mismatches or less. In some embodiments, the buffer components and concentrations and primer concentrations may be specified and incorporated into the algorithm, for example, an appropriate primer concentration is about 250 nM and appropriate buffer components are 50 mM sodium or potassium and 1.5 mM Mg²⁺.

[306] One with ordinary skill in the art of design of amplification primers will recognize that a given primer need not hybridize with 100% complementarity in order to effectively prime the synthesis of a complementary nucleic acid strand in an amplification reaction. Moreover, a primer may hybridize over one or more segments such that intervening or adjacent segments are not involved in the hybridization event. (e.g., for example, a loop structure or a hairpin structure). The primers of the present invention may comprise at least 70%, at least 75%, at least 80%, at least 85%, at least 90%, at least 95% or at least 99% sequence identity with any of the primers listed in Table 2. Thus, in some embodiments of the present invention, an extent of variation of 70% to 100%, or any range therewithin, of the sequence identity is possible relative to the specific primer sequences disclosed herein. Determination of sequence identity is described in the following example: a primer 20 nucleobases in length which is identical to another 20 nucleobase primer having two non-identical residues has 18 of 20 identical residues (18/20 = 0.9 or 90% sequence identity). In another example, a primer 15 nucleobases in length having all residues identical to a 15 nucleobase segment of primer 20 nucleobases in length would have 15/20 = 0.75 or 75% sequence identity with the 20 nucleobase primer.

[307] Percent homology, sequence identity or complementarity, can be determined by, for example, the Gap program (Wisconsin Sequence Analysis Package, Version 8 for UNIX, Genetics Computer Group, University Research Park, Madison WI), using default settings, which uses the algorithm of Smith and Waterman (Adv. Appl. Math., 1981, 2, 482-489). In some embodiments, complementarity of primers with respect to the conserved priming regions of viral nucleic acid is between about 70% and about 75% 80%. In other embodiments, homology, sequence identity or complementarity, is between about 75% and about 80%. In yet other embodiments, homology, sequence identity or complementarity, is at least 85%, at least 90%, at least 92%, at least 94%, at least 95%, at least 96%, at least 97%, at least 98%, at least 99% or is 100%.

[308] In some embodiments, the primers described herein comprise at least 70%, at least 75%, at least 80%, at least 85%, at least 90%, at least 92%, at least 94%, at least 95%, at least 96%, at least 98%,

or at least 99%, or 100% (or any range therewithin) sequence identity with the primer sequences specifically disclosed herein.

[309] One with ordinary skill is able to calculate percent sequence identity or percent sequence homology and able to determine, without undue experimentation, the effects of variation of primer sequence identity on the function of the primer in its role in priming synthesis of a complementary strand of nucleic acid for production of an amplification product of a corresponding bioagent identifying amplicon.

[310] In one embodiment, the primers are at least 13 nucleobases in length. In another embodiment, the primers are less than 36 nucleobases in length.

[311] In some embodiments of the present invention, the oligonucleotide primers are 13 to 35 nucleobases in length (13 to 35 linked nucleotide residues). These embodiments comprise oligonucleotide primers 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34 or 35 nucleobases in length, or any range therewithin. The present invention contemplates using both longer and shorter primers. Furthermore, the primers may also be linked to one or more other desired moieties, including, but not limited to, affinity groups, ligands, regions of nucleic acid that are not complementary to the nucleic acid to be amplified, labels, etc. Primers may also form hairpin structures. For example, hairpin primers may be used to amplify short target nucleic acid molecules. The presence of the hairpin may stabilize the amplification complex (see e.g., TAQMAN MicroRNA Assays, Applied Biosystems, Foster City, California).

[312] In some embodiments, any oligonucleotide primer pair may have one or both primers with less than 70% sequence homology with a corresponding member of any of the primer pairs of Table 2 if the primer pair has the capability of producing an amplification product corresponding to a bioagent identifying amplicon. In other embodiments, any oligonucleotide primer pair may have one or both primers with a length greater than 35 nucleobases if the primer pair has the capability of producing an amplification product corresponding to a bioagent identifying amplicon.

[313] In some embodiments, the function of a given primer may be substituted by a combination of two or more primers segments that hybridize adjacent to each other or that are linked by a nucleic acid loop structure or linker which allows a polymerase to extend the two or more primers in an amplification reaction.

[314] In some embodiments, the primer pairs used for obtaining bioagent identifying amplicons are the primer pairs of Table 2. In other embodiments, other combinations of primer pairs are possible by combining certain members of the forward primers with certain members of the reverse primers. An example can be seen in Table 2 for two primer pair combinations of forward primer 16S_EC_789_810_F (SEQ ID NO: 206), with the reverse primers 16S_EC_880_894_R (SEQ ID NO: 796), or 16S_EC_882_899_R or (SEQ ID NO: 818). Arriving at a favorable alternate combination of primers in a primer pair depends upon the properties of the primer pair, most notably the size of the bioagent identifying amplicon that would be produced by the primer pair, which preferably is between about 45 to about 150 nucleobases in length. Alternatively, a bioagent identifying amplicon longer than 150 nucleobases in length could be cleaved into smaller segments by cleavage reagents such as chemical reagents, or restriction enzymes, for example.

[315] In some embodiments, the primers are configured to amplify nucleic acid of a bioagent to produce amplification products that can be measured by mass spectrometry and from whose molecular masses candidate base compositions can be readily calculated.

[316] In some embodiments, any given primer comprises a modification comprising the addition of a non-templated T residue to the 5' end of the primer (i.e., the added T residue does not necessarily hybridize to the nucleic acid being amplified). The addition of a non-templated T residue has an effect of minimizing the addition of non-templated adenosine residues as a result of the non-specific enzyme activity of *Taq* polymerase (Magnuson et al., *Biotechniques*, 1996, 21, 700-709), an occurrence which may lead to ambiguous results arising from molecular mass analysis.

[317] In some embodiments of the present invention, primers may contain one or more universal bases. Because any variation (due to codon wobble in the 3rd position) in the conserved regions among species is likely to occur in the third position of a DNA (or RNA) triplet, oligonucleotide primers can be designed such that the nucleotide corresponding to this position is a base which can bind to more than one nucleotide, referred to herein as a "universal nucleobase." For example, under this "wobble" pairing, inosine (I) binds to U, C or A; guanine (G) binds to U or C, and uridine (U) binds to U or C. Other examples of universal nucleobases include nitroindoles such as 5-nitroindole or 3-nitropyrrole (Loakes et al., *Nucleosides and Nucleotides*, 1995, 14, 1001-1003), the degenerate nucleotides dP or dK (Hill *et al.*), an acyclic nucleoside analog containing 5-nitroindazole (Van Aerschot et al., *Nucleosides and Nucleotides*, 1995, 14, 1053-1056) or the purine analog 1-(2-deoxy- β -D-ribofuranosyl)-imidazole-4-carboxamide (Sala et al., *Nucl. Acids Res.*, 1996, 24, 3302-3306).

[318] In some embodiments, to compensate for the somewhat weaker binding by the wobble base, the oligonucleotide primers are designed such that the first and second positions of each triplet are occupied by nucleotide analogs that bind with greater affinity than the unmodified nucleotide. Examples of these analogs include, but are not limited to, 2,6-diaminopurine which binds to thymine, 5-propynyluracil (also known as propynylated thymine) which binds to adenine and 5-propynylcytosine and phenoxazines, including G-clamp, which binds to G. Propynylated pyrimidines are described in U.S. Patent Nos. 5,645,985, 5,830,653 and 5,484,908, each of which is commonly owned and incorporated herein by reference in its entirety. Propynylated primers are described in U.S. Pre-Grant Publication No. 2003-0170682, which is also commonly owned and incorporated herein by reference in its entirety. Phenoxazines are described in U.S. Patent Nos. 5,502,177, 5,763,588, and 6,005,096, each of which is incorporated herein by reference in its entirety. G-clamps are described in U.S. Patent Nos. 6,007,992 and 6,028,183, each of which is incorporated herein by reference in its entirety.

[319] In some embodiments, primer hybridization is enhanced using primers containing 5-propynyl deoxy-cytidine and deoxy-thymidine nucleotides. These modified primers offer increased affinity and base pairing selectivity.

[320] In some embodiments, non-template primer tags are used to increase the melting temperature (T_m) of a primer-template duplex in order to improve amplification efficiency. A non-template tag is at least three consecutive A or T nucleotide residues on a primer which are not complementary to the template. In any given non-template tag, A can be replaced by C or G and T can also be replaced by C or G. Although Watson-Crick hybridization is not expected to occur for a non-template tag relative to the template, the extra hydrogen bond in a G-C pair relative to an A-T pair confers increased stability of the primer-template duplex and improves amplification efficiency for subsequent cycles of amplification when the primers hybridize to strands synthesized in previous cycles.

[321] In other embodiments, propynylated tags may be used in a manner similar to that of the non-template tag, wherein two or more 5-propynylcytidine or 5-propynyluridine residues replace template matching residues on a primer. In other embodiments, a primer contains a modified internucleoside linkage such as a phosphorothioate linkage, for example.

[322] In some embodiments, the primers contain mass-modifying tags. Reducing the total number of possible base compositions of a nucleic acid of specific molecular weight provides a means of avoiding a persistent source of ambiguity in determination of base composition of amplification products. Addition of mass-modifying tags to certain nucleobases of a given primer will result in simplification of *de novo* determination of base composition of a given bioagent identifying amplicon from its molecular mass.

[323] In some embodiments of the present invention, the mass modified nucleobase comprises one or more of the following: for example, 7-deaza-2'-deoxyadenosine-5'-triphosphate, 5-iodo-2'-deoxyuridine-5'-triphosphate, 5-bromo-2'-deoxyuridine-5'-triphosphate, 5-bromo-2'-deoxycytidine-5'-triphosphate, 5-iodo-2'-deoxycytidine-5'-triphosphate, 5-hydroxy-2'-deoxyuridine-5'-triphosphate, 4-thiothymidine-5'-triphosphate, 5-aza-2'-deoxyuridine-5'-triphosphate, 5-fluoro-2'-deoxyuridine-5'-triphosphate, O6-methyl-2'-deoxyguanosine-5'-triphosphate, N2-methyl-2'-deoxyguanosine-5'-triphosphate, 8-oxo-2'-deoxyguanosine-5'-triphosphate or thiothymidine-5'-triphosphate. In some embodiments, the mass-modified nucleobase comprises ^{15}N or ^{13}C or both ^{15}N and ^{13}C .

[324] In some embodiments, multiplex amplification is performed where multiple bioagent identifying amplicons are amplified with a plurality of primer pairs. The advantages of multiplexing are that fewer reaction containers (for example, wells of a 96- or 384-well plate) are needed for each molecular mass measurement, providing time, resource and cost savings because additional bioagent identification data can be obtained within a single analysis. Multiplex amplification methods are well known to those with ordinary skill and can be developed without undue experimentation. However, in some embodiments, one useful and non-obvious step in selecting a plurality candidate bioagent identifying amplicons for multiplex amplification is to ensure that each strand of each amplification product will be sufficiently different in molecular mass that mass spectral signals will not overlap and lead to ambiguous analysis results. In some embodiments, a 10 Da difference in mass of two strands of one or more amplification products is sufficient to avoid overlap of mass spectral peaks.

[325] In some embodiments, as an alternative to multiplex amplification, single amplification reactions can be pooled before analysis by mass spectrometry. In these embodiments, as for multiplex amplification embodiments, it is useful to select a plurality of candidate bioagent identifying amplicons to ensure that each strand of each amplification product will be sufficiently different in molecular mass that mass spectral signals will not overlap and lead to ambiguous analysis results.

C Determination of Molecular Mass of Bioagent Identifying Amplicons

[326] In some embodiments, the molecular mass of a given bioagent identifying amplicon is determined by mass spectrometry. Mass spectrometry has several advantages, not the least of which is high bandwidth characterized by the ability to separate (and isolate) many molecular peaks across a broad range of mass to charge ratio (m/z). Thus mass spectrometry is intrinsically a parallel detection scheme without the need for radioactive or fluorescent labels, since every amplification product is identified by its molecular mass. The current state of the art in mass spectrometry is such that less than femtomole quantities of material can be readily analyzed to afford information about the molecular contents of the

sample. An accurate assessment of the molecular mass of the material can be quickly obtained, irrespective of whether the molecular weight of the sample is several hundred, or in excess of one hundred thousand atomic mass units (amu) or Daltons.

[327] In some embodiments, intact molecular ions are generated from amplification products using one of a variety of ionization techniques to convert the sample to gas phase. These ionization methods include, but are not limited to, electrospray ionization (ESI), matrix-assisted laser desorption ionization (MALDI) and fast atom bombardment (FAB). Upon ionization, several peaks are observed from one sample due to the formation of ions with different charges. Averaging the multiple readings of molecular mass obtained from a single mass spectrum affords an estimate of molecular mass of the bioagent identifying amplicon. Electrospray ionization mass spectrometry (ESI-MS) is particularly useful for very high molecular weight polymers such as proteins and nucleic acids having molecular weights greater than 10 kDa, since it yields a distribution of multiply-charged molecules of the sample without causing a significant amount of fragmentation.

[328] The mass detectors used in the methods of the present invention include, but are not limited to, Fourier transform ion cyclotron resonance mass spectrometry (FT-ICR-MS), time of flight (TOF), ion trap, quadrupole, magnetic sector, Q-TOF, and triple quadrupole.

D. Base Compositions of Bioagent Identifying Amplicons

[329] Although the molecular mass of amplification products obtained using intelligent primers provides a means for identification of bioagents, conversion of molecular mass data to a base composition signature is useful for certain analyses. As used herein, "base composition" is the exact number of each nucleobase (A, T, C and G) determined from the molecular mass of a bioagent identifying amplicon. In some embodiments, a base composition provides an index of a specific organism. Base compositions can be calculated from known sequences of known bioagent identifying amplicons and can be experimentally determined by measuring the molecular mass of a given bioagent identifying amplicon, followed by determination of all possible base compositions which are consistent with the measured molecular mass within acceptable experimental error. The following example illustrates determination of base composition from an experimentally obtained molecular mass of a 46-mer amplification product originating at position 1337 of the 16S rRNA of *Bacillus anthracis*. The forward and reverse strands of the amplification product have measured molecular masses of 14208 and 14079 Da, respectively. The possible base compositions derived from the molecular masses of the forward and reverse strands for the *B. anthracis* products are listed in Table 1.

Table 1

Possible Base Compositions for *B. anthracis* 46mer Amplification Product

Calc. Mass Forward Strand	Mass Error Forward Strand	Base Composition of Forward Strand	Calc. Mass Reverse Strand	Mass Error Reverse Strand	Base Composition of Reverse Strand
14208.2935	0.079520	A1 G17 C10 T18	14079.2624	0.080600	A0 G14 C13 T19
14208.3160	0.056980	A1 G20 C15 T10	14079.2849	0.058060	A0 G17 C18 T11
14208.3386	0.034440	A1 G23 C20 T2	14079.3075	0.035520	A0 G20 C23 T3
14208.3074	0.065560	A6 G11 C3 T26	14079.2538	0.089180	A5 G5 C1 T35
14208.3300	0.043020	A6 G14 C8 T18	14079.2764	0.066640	A5 G8 C5 T27
14208.3525	0.020480	A6 G17 C13 T10	14079.2989	0.044100	A5 G11 C11 T19
14208.3751	0.002060	A6 G20 C18 T2	14079.3214	0.021560	A5 G14 C16 T11
14208.3439	0.029060	A11 G8 C1 T26	14079.3440	0.000980	A5 G17 C21 T3
14208.3665	0.006520	A11 G11 C6 T18	14079.3129	0.030140	A10 G5 C4 T27
14208.3890	0.016020	A11 G14 C11 T10	14079.3354	0.007600	A10 G8 C9 T19
14208.4116	0.038560	A11 G17 C16 T2	14079.3579	0.014940	A10 G11 C14 T11
14208.4030	0.029980	A16 G8 C4 T18	14079.3805	0.037480	A10 G14 C19 T3
14208.4255	0.052520	A16 G11 C9 T10	14079.3494	0.006360	A15 G2 C2 T27
14208.4481	0.075060	A16 G14 C14 T2	14079.3719	0.028900	A15 G5 C7 T19
14208.4395	0.066480	A21 G5 C2 T18	14079.3944	0.051440	A15 G8 C12 T3
14208.4620	0.089020	A21 G8 C7 T10	14079.4170	0.073980	A15 G11 C17 T3
-	-	-	14079.4084	0.065400	A20 G2 C5 T19
-	-	-	14079.4309	0.087940	A20 G5 C10 T13

[330] Among the 16 possible base compositions for the forward strand and the 18 possible base compositions for the reverse strand that were calculated, only one pair (shown in **bold**) are complementary base compositions, which indicates the true base composition of the amplification product. It should be recognized that this logic is applicable for determination of base compositions of any bioagent identifying amplicon, regardless of the class of bioagent from which the corresponding amplification product was obtained.

[331] In some embodiments, assignment of previously unobserved base compositions (also known as "true unknown base compositions") to a given phylogeny can be accomplished via the use of pattern classifier model algorithms. Base compositions, like sequences, vary slightly from strain to strain within species, for example. In some embodiments, the pattern classifier model is the mutational probability model. On other embodiments, the pattern classifier is the polytope model. The mutational probability model and polytope model are both commonly owned and described in U.S. Patent application Serial No. 11/073,362 which is incorporated herein by reference in entirety.

[332] In one embodiment, it is possible to manage this diversity by building "base composition probability clouds" around the composition constraints for each species. This permits identification of

organisms in a fashion similar to sequence analysis. A “pseudo four-dimensional plot” can be used to visualize the concept of base composition probability clouds. Optimal primer design requires optimal choice of bioagent identifying amplicons and maximizes the separation between the base composition signatures of individual bioagents. Areas where clouds overlap indicate regions that may result in a misclassification, a problem which is overcome by a triangulation identification process using bioagent identifying amplicons not affected by overlap of base composition probability clouds.

[333] In some embodiments, base composition probability clouds provide the means for screening potential primer pairs in order to avoid potential misclassifications of base compositions. In other embodiments, base composition probability clouds provide the means for predicting the identity of a bioagent whose assigned base composition was not previously observed and/or indexed in a bioagent identifying amplicon base composition database due to evolutionary transitions in its nucleic acid sequence. Thus, in contrast to probe-based techniques, mass spectrometry determination of base composition does not require prior knowledge of the composition or sequence in order to make the measurement.

[334] The present invention provides bioagent classifying information similar to DNA sequencing and phylogenetic analysis at a level sufficient to identify a given bioagent. Furthermore, the process of determination of a previously unknown base composition for a given bioagent (for example, in a case where sequence information is unavailable) has downstream utility by providing additional bioagent indexing information with which to populate base composition databases. The process of future bioagent identification is thus greatly improved as more BCS indexes become available in base composition databases.

E. Triangulation Identification

[335] In some cases, a molecular mass of a single bioagent identifying amplicon alone does not provide enough resolution to unambiguously identify a given bioagent. The employment of more than one bioagent identifying amplicon for identification of a bioagent is herein referred to as “triangulation identification.” Triangulation identification is pursued by determining the molecular masses of a plurality of bioagent identifying amplicons selected within a plurality of housekeeping genes. This process is used to reduce false negative and false positive signals, and enable reconstruction of the origin of hybrid or otherwise engineered bioagents. For example, identification of the three part toxin genes typical of *B. anthracis* (Bowen et al., J. Appl. Microbiol., 1999, 87, 270-278) in the absence of the expected signatures from the *B. anthracis* genome would suggest a genetic engineering event.

[336] In some embodiments, the triangulation identification process can be pursued by characterization of bioagent identifying amplicons in a massively parallel fashion using the polymerase chain reaction (PCR), such as multiplex PCR where multiple primers are employed in the same amplification reaction mixture, or PCR in multi-well plate format wherein a different and unique pair of primers is used in multiple wells containing otherwise identical reaction mixtures. Such multiplex and multi-well PCR methods are well known to those with ordinary skill in the arts of rapid throughput amplification of nucleic acids. In other related embodiments, one PCR reaction per well or container may be carried out, followed by an amplicon pooling step wherein the amplification products of different wells are combined in a single well or container which is then subjected to molecular mass analysis. The combination of pooled amplicons can be chosen such that the expected ranges of molecular masses of individual amplicons are not overlapping and thus will not complicate identification of signals.

F. Codon Base Composition Analysis

[337] In some embodiments of the present invention, one or more nucleotide substitutions within a codon of a gene of an infectious organism confer drug resistance upon an organism which can be determined by codon base composition analysis. The organism can be a bacterium, virus, fungus or protozoan.

[338] In some embodiments, the amplification product containing the codon being analyzed is of a length of about 35 to about 200 nucleobases. The primers employed in obtaining the amplification product can hybridize to upstream and downstream sequences directly adjacent to the codon, or can hybridize to upstream and downstream sequences one or more sequence positions away from the codon. The primers may have between about 70% to 100% sequence complementarity with the sequence of the gene containing the codon being analyzed.

[339] In some embodiments, the codon base composition analysis is undertaken

[340] In some embodiments, the codon analysis is undertaken for the purpose of investigating genetic disease in an individual. In other embodiments, the codon analysis is undertaken for the purpose of investigating a drug resistance mutation or any other deleterious mutation in an infectious organism such as a bacterium, virus, fungus or protozoan. In some embodiments, the bioagent is a bacterium identified in a biological product.

[341] In some embodiments, the molecular mass of an amplification product containing the codon being analyzed is measured by mass spectrometry. The mass spectrometry can be either electrospray (ESI) mass spectrometry or matrix-assisted laser desorption ionization (MALDI) mass spectrometry.

Time-of-flight (TOF) is an example of one mode of mass spectrometry compatible with the analyses of the present invention.

[342] The methods of the present invention can also be employed to determine the relative abundance of drug resistant strains of the organism being analyzed. Relative abundances can be calculated from amplitudes of mass spectral signals with relation to internal calibrants. In some embodiments, known quantities of internal amplification calibrants can be included in the amplification reactions and abundances of analyte amplification product estimated in relation to the known quantities of the calibrants.

[343] In some embodiments, upon identification of one or more drug-resistant strains of an infectious organism infecting an individual, one or more alternative treatments can be devised to treat the individual.

G. Determination of the Quantity of a Bioagent

[344] In some embodiments, the identity and quantity of an unknown bioagent can be determined using the process illustrated in Figure 2. Primers (500) and a known quantity of a calibration polynucleotide (505) are added to a sample containing nucleic acid of an unknown bioagent. The total nucleic acid in the sample is then subjected to an amplification reaction (510) to obtain amplification products. The molecular masses of amplification products are determined (515) from which are obtained molecular mass and abundance data. The molecular mass of the bioagent identifying amplicon (520) provides the means for its identification (525) and the molecular mass of the calibration amplicon obtained from the calibration polynucleotide (530) provides the means for its identification (535). The abundance data of the bioagent identifying amplicon is recorded (540) and the abundance data for the calibration data is recorded (545), both of which are used in a calculation (550) which determines the quantity of unknown bioagent in the sample.

[345] A sample comprising an unknown bioagent is contacted with a pair of primers that provide the means for amplification of nucleic acid from the bioagent, and a known quantity of a polynucleotide that comprises a calibration sequence. The nucleic acids of the bioagent and of the calibration sequence are amplified and the rate of amplification is reasonably assumed to be similar for the nucleic acid of the bioagent and of the calibration sequence. The amplification reaction then produces two amplification products: a bioagent identifying amplicon and a calibration amplicon. The bioagent identifying amplicon and the calibration amplicon should be distinguishable by molecular mass while being amplified at essentially the same rate. Effecting differential molecular masses can be accomplished by choosing as a calibration sequence, a representative bioagent identifying amplicon (from a specific species of bioagent) and performing, for example, a 2-8 nucleobase deletion or insertion within the variable region between

the two priming sites. The amplified sample containing the bioagent identifying amplicon and the calibration amplicon is then subjected to molecular mass analysis by mass spectrometry, for example. The resulting molecular mass analysis of the nucleic acid of the bioagent and of the calibration sequence provides molecular mass data and abundance data for the nucleic acid of the bioagent and of the calibration sequence. The molecular mass data obtained for the nucleic acid of the bioagent enables identification of the unknown bioagent and the abundance data enables calculation of the quantity of the bioagent, based on the knowledge of the quantity of calibration polynucleotide contacted with the sample.

[346] In some embodiments, construction of a standard curve where the amount of calibration polynucleotide spiked into the sample is varied provides additional resolution and improved confidence for the determination of the quantity of bioagent in the sample. The use of standard curves for analytical determination of molecular quantities is well known to one with ordinary skill and can be performed without undue experimentation.

[347] In some embodiments, multiplex amplification is performed where multiple bioagent identifying amplicons are amplified with multiple primer pairs which also amplify the corresponding standard calibration sequences. In this or other embodiments, the standard calibration sequences are optionally included within a single vector which functions as the calibration polynucleotide. Multiplex amplification methods are well known to those with ordinary skill and can be performed without undue experimentation.

[348] In some embodiments, the calibrant polynucleotide is used as an internal positive control to confirm that amplification conditions and subsequent analysis steps are successful in producing a measurable amplicon. Even in the absence of copies of the genome of a bioagent, the calibration polynucleotide should give rise to a calibration amplicon. Failure to produce a measurable calibration amplicon indicates a failure of amplification or subsequent analysis step such as amplicon purification or molecular mass determination. Reaching a conclusion that such failures have occurred is in itself, a useful event.

[349] In some embodiments, the calibration sequence is comprised of DNA. In some embodiments, the calibration sequence is comprised of RNA.

[350] In some embodiments, the calibration sequence is inserted into a vector that itself functions as the calibration polynucleotide. In some embodiments, more than one calibration sequence is inserted into the vector that functions as the calibration polynucleotide. Such a calibration polynucleotide is herein termed a "combination calibration polynucleotide." The process of inserting polynucleotides into vectors

is routine to those skilled in the art and can be accomplished without undue experimentation. Thus, it should be recognized that the calibration method should not be limited to the embodiments described herein. The calibration method can be applied for determination of the quantity of any bioagent identifying amplicon when an appropriate standard calibrant polynucleotide sequence is designed and used. The process of choosing an appropriate vector for insertion of a calibrant is also a routine operation that can be accomplished by one with ordinary skill without undue experimentation.

H. Identification of Bacteria

[351] In other embodiments of the present invention, the primer pairs produce bioagent identifying amplicons within stable and highly conserved regions of bacteria. The advantage to characterization of an amplicon defined by priming regions that fall within a highly conserved region is that there is a low probability that the region will evolve past the point of primer recognition, in which case, the primer hybridization of the amplification step would fail. Such a primer set is thus useful as a broad range survey-type primer. In another embodiment of the present invention, the intelligent primers produce bioagent identifying amplicons including a region which evolves more quickly than the stable region described above. The advantage of characterization bioagent identifying amplicon corresponding to an evolving genomic region is that it is useful for distinguishing emerging strain variants or the presence of virulence genes, drug resistance genes, or codon mutations that induce drug resistance.

[352] The present invention also has significant advantages as a platform for identification of diseases caused by emerging bacterial strains such as, for example, drug-resistant strains of *Staphylococcus aureus*. The present invention eliminates the need for prior knowledge of bioagent sequence to generate hybridization probes. This is possible because the methods are not confounded by naturally occurring evolutionary variations occurring in the sequence acting as the template for production of the bioagent identifying amplicon. Measurement of molecular mass and determination of base composition is accomplished in an unbiased manner without sequence prejudice.

[353] Another embodiment of the present invention also provides a means of tracking the spread of a bacterium, such as a particular drug-resistant strain when a plurality of samples obtained from different locations are analyzed by the methods described above in an epidemiological setting. In one embodiment, a plurality of samples from a plurality of different locations is analyzed with primer pairs which produce bioagent identifying amplicons, a subset of which contains a specific drug-resistant bacterial strain. The corresponding locations of the members of the drug-resistant strain subset indicate the spread of the specific drug-resistant strain to the corresponding locations.

I. Kits

[354] The present invention also provides kits for carrying out the methods described herein. In some embodiments, the kit may comprise a sufficient quantity of one or more primer pairs to perform an amplification reaction on a target polynucleotide from a bioagent to form a bioagent identifying amplicon. In some embodiments, the kit may comprise from one to fifty primer pairs, from one to twenty primer pairs, from one to ten primer pairs, or from two to five primer pairs. In some embodiments, the kit may comprise one or more primer pairs recited in Table 2.

[355] In some embodiments, the kit comprises one or more broad range survey primer(s), division wide primer(s), or drill-down primer(s), or any combination thereof. If a given problem involves identification of a specific bioagent, the solution to the problem may require the selection of a particular combination of primers to provide the solution to the problem. A kit may be designed so as to comprise particular primer pairs for identification of a particular bioagent. A drill-down kit may be used, for example, to distinguish different genotypes or strains, drug-resistant, or otherwise. In some embodiments, the primer pair components of any of these kits may be additionally combined to comprise additional combinations of broad range survey primers and division-wide primers so as to be able to identify a bacterium.

[356] In some embodiments, the kit contains standardized calibration polynucleotides for use as internal amplification calibrants. Internal calibrants are described in commonly owned U.S. Patent Application Serial No: 60/545,425 which is incorporated herein by reference in its entirety.

[357] In some embodiments, the kit comprises a sufficient quantity of reverse transcriptase (if RNA is to be analyzed for example), a DNA polymerase, suitable nucleoside triphosphates (including alternative dNTPs such as inosine or modified dNTPs such as the 5-propynyl pyrimidines or any dNTP containing molecular mass-modifying tags such as those described above), a DNA ligase, and/or reaction buffer, or any combination thereof, for the amplification processes described above. A kit may further include instructions pertinent for the particular embodiment of the kit, such instructions describing the primer pairs and amplification conditions for operation of the method. A kit may also comprise amplification reaction containers such as microcentrifuge tubes and the like. A kit may also comprise reagents or other materials for isolating bioagent nucleic acid or bioagent identifying amplicons from amplification, including, for example, detergents, solvents, or ion exchange resins which may be linked to magnetic beads. A kit may also comprise a table of measured or calculated molecular masses and/or base compositions of bioagents using the primer pairs of the kit.

[358] Some embodiments are kits that contain one or more survey bacterial primer pairs represented by primer pair compositions wherein each member of each pair of primers has 70% to 100% sequence

identity with the corresponding member from the group of primer pairs represented by any of the primer pairs of Table 5. The survey primer pairs may include broad range primer pairs which hybridize to ribosomal RNA, and may also include division-wide primer pairs which hybridize to housekeeping genes such as rplB, tufB, rpoB, rpoC, valS, and infB, for example.

[359] In some embodiments, a kit may contain one or more survey bacterial primer pairs and one or more triangulation genotyping analysis primer pairs such as the primer pairs of Tables 8, 12, 14, 19, 21, 23, or 24. In some embodiments, the kit may represent a less expansive genotyping analysis but include triangulation genotyping analysis primer pairs for more than one genus or species of bacteria. For example, a kit for surveying nosocomial infections at a health care facility may include, for example, one or more broad range survey primer pairs, one or more division wide primer pairs, one or more *Acinetobacter baumannii* triangulation genotyping analysis primer pairs and one or more *Staphylococcus aureus* triangulation genotyping analysis primer pairs. One with ordinary skill will be capable of analyzing *in silico* amplification data to determine which primer pairs will be able to provide optimal identification resolution for the bacterial bioagents of interest.

[360] In some embodiments, a kit may be assembled for identification of strains of bacteria involved in contamination of food. An example of such a kit embodiment is a kit comprising one or more bacterial survey primer pairs of Table 5 with one or more triangulation genotyping analysis primer pairs of Table 12 which provide strain resolving capabilities for identification of specific strains of *Campylobacter jejuni*.

[361] Some embodiments of the kits are 96-well or 384-well plates with a plurality of wells containing any or all of the following components: dNTPs, buffer salts, Mg²⁺, betaine, and primer pairs. In some embodiments, a polymerase is also included in the plurality of wells of the 96-well or 384-well plates.

[362] Some embodiments of the kit contain instructions for PCR and mass spectrometry analysis of amplification products obtained using the primer pairs of the kits.

[363] Some embodiments of the kit include a barcode which uniquely identifies the kit and the components contained therein according to production lots and may also include any other information relative to the components such as concentrations, storage temperatures, etc. The barcode may also include analysis information to be read by optical barcode readers and sent to a computer controlling amplification, purification and mass spectrometric measurements. In some embodiments, the barcode provides access to a subset of base compositions in a base composition database which is in digital

communication with base composition analysis software such that a base composition measured with primer pairs from a given kit can be compared with known base compositions of bioagent identifying amplicons defined by the primer pairs of that kit.

[364] In some embodiments, the kit contains a database of base compositions of bioagent identifying amplicons defined by the primer pairs of the kit. The database is stored on a convenient computer readable medium such as a compact disk or USB drive, for example.

[365] In some embodiments, the kit includes a computer program stored on a computer formatted medium (such as a compact disk or portable USB disk drive, for example) comprising instructions which direct a processor to analyze data obtained from the use of the primer pairs of the present invention. The instructions of the software transform data related to amplification products into a molecular mass or base composition which is a useful concrete and tangible result used in identification and/or classification of bioagents. In some embodiments, the kits of the present invention contain all of the reagents sufficient to carry out one or more of the methods described herein.

[366] While the present invention has been described with specificity in accordance with certain of its embodiments, the following examples serve only to illustrate the invention and are not intended to limit the same. In order that the invention disclosed herein may be more efficiently understood, examples are provided below. It should be understood that these examples are for illustrative purposes only and are not to be construed as limiting the invention in any manner.

EXAMPLES

Example 1: Design and Validation of Primers that Define Bioagent Identifying Amplicons for Identification of Bacteria

[367] For design of primers that define bacterial bioagent identifying amplicons, a series of bacterial genome segment sequences were obtained, aligned and scanned for regions where pairs of PCR primers would amplify products of about 45 to about 150 nucleotides in length and distinguish subgroups and/or individual strains from each other by their molecular masses or base compositions. A typical process shown in Figure 1 is employed for this type of analysis.

[368] A database of expected base compositions for each primer region was generated using an *in silico* PCR search algorithm, such as (ePCR). An existing RNA structure search algorithm (Macke et al., Nucl. Acids Res., 2001, 29, 4724-4735, which is incorporated herein by reference in its entirety) has been modified to include PCR parameters such as hybridization conditions, mismatches, and thermodynamic calculations (SantaLucia, Proc. Natl. Acad. Sci. U.S.A., 1998, 95, 1460-1465, which is incorporated

herein by reference in its entirety). This also provides information on primer specificity of the selected primer pairs.

[369] Table 2 represents a collection of primers (sorted by primer pair number) designed to identify bacteria using the methods described herein. The primer pair number is an in-house database index number. Primer sites were identified on segments of genes, such as, for example, the 16S rRNA gene. The forward or reverse primer name shown in Table 2 indicates the gene region of the bacterial genome to which the primer hybridizes relative to a reference sequence. In Table 2, for example, the forward primer name 16S_EC_1077_1106_F indicates that the forward primer (F) hybridizes to residues 1077-1106 of the reference sequence represented by a sequence extraction of coordinates 4033120..4034661 from GenBank gi number 16127994 (as indicated in Table 3). As an additional example: the forward primer name BONTA_X52066_450_473 indicates that the primer hybridizes to residues 450-437 of the gene encoding *Clostridium botulinum* neurotoxin type A (BoNT/A) represented by GenBank Accession No. X52066 (primer pair name codes appearing in Table 2 are defined in Table 3. One with ordinary skill knows how to obtain individual gene sequences or portions thereof from genomic sequences present in GenBank. In Table 2, Tp = 5-propynyluracil; Cp = 5-propynylcytosine; * = phosphorothioate linkage; I = inosine. T. GenBank Accession Numbers for reference sequences of bacteria are shown in Table 3 (below). In some cases, the reference sequences are extractions from bacterial genomic sequences or complements thereof.

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO.: 10593)

Table 2: Primer Pairs for Identification of Bacteria

Primer Pair Number	Forward Primer Name	Forward Sequence	Forward SEQ ID NO.	Reverse Primer Name	Reverse Sequence	Reverse SEQ ID NO.
1	16S BC 1077 1106 F	GTGATGTGTGGGTTAATGTTCCCGTAA	134	16S BC 1175 1195 R	GAACTCATCCCAACACTCTCTC	809
2	16S BC 1082 1106 F	ATGTGGGTGTATGTCGCCAAGAG	38	16S BC 1175 1197 R	TTCAGTGTATCCCACTCTCTCTC	1398
3	16S BC 1090 1111 F	TTAATGTCGCAACATATCCGAA	651	16S BC 1175 1196 R	TGAAGTCTACCCCACTCTCTC	1159
4	16S BC 1222 1241 F	GCTACACAGTGTGTATATG	114	16S BC 1303 1323 R	CGAGTTGCACTGACGATACCG	787
5	16S BC 1332 1353 F	ATGTCCGATATGTGTATATG	10	16S BC 1389 1407 R	GACGGCGGTGTGTATGAG	806
6	16S BC 30 54 F	TGAGACCTGTGGGCTGCTTAAAC	429	16S BC 105 126 R	TTACTGACACCTGCGGCGCT	897
7	16S BC 38 64 F	GTGCTATGCTCTATATACATGCAATG	136	16S BC 101 120 R	TTACTGACACCTGCGGCGCT	1365
8	16S BC 43 68 F	TAACTACATGCAATGCAAG	152	16S BC 104 120 R	TTACTGACACCTGCGGCGCT	1364
9	16S BC 683 700 F	GTTGACGCGTGAATGCG	137	16S BC 774 795 R	GTATCTAATCTGTTTGTCTCC	839
10	16S BC 713 732 F	AGACACCGATGCGAAGGC	21	16S BC 789 809 R	CTGTGACATGCGGGTATCTA	798
11	16S BC 785 806 F	GGATTAGACATCTGTGTATCTC	118	16S BC 880 897 R	GCGCTACTCTCCCAAGCG	830
12	16S BC 785 810 F	GGATTAGACATCTGTGTATCTC	119	16S BC 880 897 2 R	GCGCTACTCTCCCAAGCG	830
13	16S BC 789 810 F	TAACTACATCTGTGTATCTC	206	16S BC 880 894 R	CTACTCTCCCAAGCG	796
14	16S BC 960 981 F	TTCGATGCAAGGAGAGAACT	672	16S BC 1054 1073 R	ACGACTCTCCCAAGCG	735
15	16S BC 969 985 F	ACGAGAGAACTTAC	19	16S BC 1061 1078 R	ACGACAGACTGACAC	734
16	23S BC 1926 1843 P	CTGACACTGCTCCCGTGC	80	23S BC 1906 1924 R	ACGCTGTATGTTACGCG	805
17	23S BC 2645 2669 F	TGTGCTCTGTATGACAGAGACCG	408	23S BC 2744 2761 R	TGCTTATGACTCTGAC	1252
18	23S BC 2645 2669 2	CTGTCTCTGTATGACAGAGACCG	83	23S BC 2751 2767 R	GTTCAGTCTGAGAGCTTTCAC	846
19	23S BC 493 518 F	GGGAGTGAAGAGATCTGACACCG	125	23S BC 551 571 R	ACAAAAGTACAGCTGACCG	717
20	23S BC 493 518 2 F	GGGAGTGAAGAGATCTGACACCG	125	23S BC 551 571 2 R	ACAAAAGTACAGCTGACCG	716
21	23S BC 971 992 F	CGAGAGGAAACATCTGATACCG	66	23S BC 1059 1077 R	TGGCTGCTGTAGCGAC	1282
22	CAPC BA 104 131 F	GTATTTAGCACTGCTTTTATACG	139	CAPC BA 180 205 R	TGAACTGTGAAGACATGAGTACG	1150
23	CAPC BA 114 132 F	ACTGTGTTTATATGACCG	20	CAPC BA 185 205 R	TGAATCTGAAACACACGAC	1149
24	CAPC BA 274 301 F	GATTATGTTTATCTGATGATGCAAT	109	CAPC BA 349 376 R	GTAACCTGTCTTGAATATGATTCG	837
25	CAPC BA 276 296 F	TATATGTATCTGATGATGCAAT	663	CAPC BA 356 377 R	GTAACCTGTCTTGAATATGATTCG	834
26	CAPC BA 281 301 F	GTATATCTGATGATGATGCAAT	138	CAPC BA 361 378 R	TGGTACACCTGTCTTGAAT	1288
27	CAPC BA 315 334 F	CGATGATATGATGATGATG	59	CAPC BA 361 378 R	TGGTACACCTGTCTTGAAT	1288
28	CYA BA 1055 1072 F	GAAAGGTTTGGATGCG	92	CYA BA 1112 1130 R	TGTTGACATGCTGCTCTAC	1352

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

29	CYA BA 1349 1370 F	ACACQAGCTACATCATCAGC	12	CYA BA 1447 1426 R	CTTCTACATTTTACCATAC	800
30	CYA BA 1355 1379 F	CGAATCTACATCAGACAGAGAG	64	CYA BA 1448 1467 R	TGTTAAGGCTTACAGACC	1342
31	CYA BA 1359 1379 F	ACATATCAGACAGACAGAGAG	13	CYA BA 1447 1467 R	CGGCTTTACAGACC	794
32	CYA BA 914 937 F	CGAGTTTATCTCAGACAGAGAG	53	CYA BA 999 1026 R	ACACTTTTAAAGTTTGTACACAC	728
33	CYA BA 916 935 F	GGTTTATCTCAGACAGAGAG	131	CYA BA 1003 1025 R	CACTTTTATTAAGTTTGTAC	768
34	INF B EC 1365 1393 F	TTA	524	INF B EC 1439 1467 R	TGCTGCTTCCATGTTTATCTTCA	1248
35	LEF BA 1033 1052 F	TCAGAGAGAGAGAGC	254	LEF BA 1119 1135 R	AGATTAACATTTTAC	803
36	LEF BA 1036 1066 F	CAGAGAGAGAGAGCTTAAAG	44	LEF BA 1119 1149 R	AGATTAAGATCAGATACATATTT	745
37	LEF BA 756 781 F	AGCTTTGCTATTAATCAGC	26	LEF BA 843 872 R	TGTTCCAGAGATAGATTTATTTCTGTT	1135
38	LEF BA 759 778 F	CTTTTGCTATTAATCAGC	90	LEF BA 843 865 R	AGATTAAGATTTATTTCTGTTG	748
39	LEF BA 755 813 F	TTTACAGCTTATCAGC	700	LEF BA 883 900 R	TGTTCCAGATCCGTTG	1140
40	LEF BA 843 859 F	CACGATCTTCCAG	43	LEF BA 939 958 R	CAATTAAGATCCCTCAG	762
41	PAG BA 122 142 F	CGATATCAGTCCAGAG	49	PAG BA 190 209 R	CTGTATGATAGAGATAC	781
42	PAG BA 123 145 F	AGATCAGTTTCCAGAGGTTAC	22	PAG BA 377 210 R	CCGTATGATAGAGATGATACAG	774
43	PAG BA 269 287 F	AACTGCTATTTGCTGAG	11	PAG BA 325 344 R	TGATATCAGAGAGTAC	1186
44	PAG BA 655 675 F	GAGATATCAGTGTATGATC	93	PAG BA 755 772 R	CGTGTCTCATTTTAC	778
45	PAG BA 753 772 F	TCCTGAGAAATGAGACAG	341	PAG BA 849 868 R	TGATTAAGCTCCACAG	1089
46	PAG BA 763 781 F	TGAGACAGCTTCTATC	552	PAG BA 849 868 R	TGAGAGGCTATTTTGTAAACAC	1089
47	REOC EC 1018 1045 F	CT	39	REOC EC 1095 1124 R	AT	959
48	REOC EC 1018 1045 2 F	CT	39	REOC EC 1095 1124 2 R	AT	958
49	REOC EC 114 140 F	TAGAGGCGAGACATCAACTAC	158	REOC EC 213 232 R	GGGCTTGTACTTACGAC	831
50	REOC EC 2178 2196 F	TGATTTCTGCTCCGTTG	478	REOC EC 2225 2246 R	TGCTCATCAGACAGACATAC	1414
51	REOC EC 2178 2196 2 F	TGATTTCTGCTCCGTTG	477	REOC EC 2225 2246 2 R	TGCTCATCAGACAGACATAC	1413
52	REOC EC 2218 2241 F	CTTGAGATGCTGCTGCTGATG	81	REOC EC 2313 2337 R	CGACCTGCTGCTGATGATGATAC	790
53	REOC EC 2218 2241 2 F	CTTGCTATGCTGCTGCTGATG	86	REOC EC 2313 2337 2 R	CGACCTGCTGCTGATGATGATAC	789
54	REOC EC 808 833 F	CGTGGGTATTAACCTACACAG	75	REOC EC 865 889 R	GTCTTCTGCTGCTGATGATGATC	847
55	REOC EC 808 833 2 F	CGTGGGTATTAACCTACACAG	76	REOC EC 865 891 R	ACCTTTTCTGCTGCTGATGATGATC	741
56	REOC EC 993 1019 F	A	41	REOC EC 1036 1059 R	CAGAGGCTGATGATGATGATGATC	785
57	REOC EC 993 1019 2 F	A	40	REOC EC 1036 1059 2 R	CAGAGGCTGATGATGATGATGATC	784

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

58	SPPE BA 115 137 F	CAAGCTAAGCCGACATCAAGAC	45	SPPE BA 137 222 R	TCCAGCTGCTGTTCACTTGCAGATTC	1201
59	TUPE EC 239 259 F	TGAGCTGCGAGAGACAGAGCTG	204	TUPE EC 283 303 R	GCCCTCCATCTGAGCAGACAC	815
60	TUPE EC 239 259 2 F	TGAGCTGCGAGAGACAGAGCTG	678	TUPE EC 303 303 2 R	GCCCTCCATCTGAGCAGACAC	816
61	TUPE EC 976 1000 F	AACTACAGTCGCGAGCTTCACTTCC	4	TUPE EC 1045 1068 R	GTTCGCGCGGACGACATTTTC	845
62	TUPE EC 976 1000 2 F	AACTACAGTCGCGAGCTTCACTTCC	5	TUPE EC 1045 1068 2 R	GTTCGCGCGAGCATTACCATTTTC	844
63	TUPE EC 985 1012 F	CCAGCTTCTCTCTCTGCTACATCGA	56	TUPE EC 1033 1062 R	TCCAGGCACTCACTTCTGACTCTCT	1006
64	TUPE EC 985 1012 F	CCAGCTTCTCTCTCTGCTACATCGA	56	TUPE EC 1033 1062 R	TCCAGGCACTCACTTCTGACTCTCT	1006
65	REPE EC 650 679 F	GACCTACGTAGAGGTTCTGTAATG	98	REPE EC 739 762 R	TCCAGGCTGCTGTATCCCATCATG	999
66	REPE EC 688 710 F	CACTCACAAGTGCTGTGTAGAGG	54	REPE EC 735 757 R	GTCTGCTTTACCCATCATGAGT	842
67	REPE EC 688 710 F	CACTCACAAGTGCTGTGTAGAGG	54	REPE EC 735 757 R	GTCTGCTTTACCCATCATGAGT	842
68	RPOC EC 1036 1060 F	CGTGTGACATTCGCGGCGTTCAG	78	RPOC EC 1097 1126 R	ATTACAGACCATTTCTTTGCTAAAC	754
69	RPOB EC 3782 3790 F	TGCACACTCTTTGCGGTGTAAGCTC	248	RPOB EC 3836 3865 R	TTTCTTGAAGGTATGAGCTGCTCGTA	1435
70	REPE EC 688 710 F	CACTCACAAGTGCTGTGTAGAGG	54	REPE EC 743 771 R	TGTTTGTATATCCAGTGCTGTGTTACCC	1356
71	VAPE EC 1105 1124 F	CGTGGCGCTGTGTTGTGGA	77	VAPE EC 1195 1218 R	CGTGTGACATCTGAGATCTGCGCTT	795
72	RPOB EC 1845 1865 F	TATGTGTGCGAGCACTCGAC	233	RPOB EC 1909 1939 R	GCTGCGATTCGCTTTCTCAG	825
73	REPE EC 669 698 F	TGATATGACCTTATATACATCTC	623	REPE EC 735 761 R	CGATGTGCTGTTTACCCATCTGAGTA	767
74	REPE EC 671 700 F	TGATATGACCTTATATACATCTC	169	REPE EC 737 762 R	TCCAGGCTGCTGTGTTACCCATCTGAG	1000
75	SP101 SP2T11 1 29 F	AACTCAATTGTGAAAGAACACAGAA	2	SP101 SP2T11 92 116 R	CCCTACCCAAAGTCTCACAGAGGAG	779
76	SP101 SP2T11 118 14 F	CGTGTGAAATTAACCAAGATGTGT	115	SP101 SP2T11 213 238 R	TGTCGCGAGATCTCACAGCTGCTGCT	1340
77	SP101 SP2T11 216 24 F	ACGAGTGATGTAATCCGACATGGA	24	SP101 SP2T11 308 333 R	TGCACTGTTGACAACTCTCTGTTGCTG	1209
78	SP101 SP2T11 266 29 F	CTGTGATCTGTGCTGACAGCGCTG	89	SP101 SP2T11 355 380 R	CGTGTGTTGATGCTGATTCCTCTTC	824
79	SP101 SP2T11 312 34 F	GTCAATGTGACAGCTTTATGCG	132	SP101 SP2T11 423 441 R	ATCCCTGCTCTTCTCTGCG	753
80	SP101 SP2T11 358 38 F	GTCAATGTGACAGCTTTATGCG	126	SP101 SP2T11 448 473 R	CGAATCTTTCCACACAGATCTGAG	766
81	SP101 SP2T11 600 62 F	CGTCACTGCAAGCTTGAATCTTTG	62	SP101 SP2T11 686 714 R	CCATTTCTTCTGACATCTGAGAAAT	772
82	SP101 SP2T11 658 68 F	GAGATGTGATCACTACAGATAGAGA	127	SP101 SP2T11 756 784 R	GATGTGCAATGAACTGATTTCTTAA	813
83	SP101 SP2T11 776 80 F	TCCCACTCAAACTGAGGATGCG	364	SP101 SP2T11 871 896 R	GCCCAAGAGAGACTAGACGATAA	814
84	SP101 SP2T11 893 92 F	GCGAAGCAGCAAGATTCGCAATGC	123	SP101 SP2T11 988 1012 R	CATGACGCGAAGACTCTCACCAAC	763

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

85	SPI01_SPE111_1154_1	CAATACACACAGCGTGGCTGGG	47	SPI01_SPE111_1251_1277_R	GAACCCACCTGGCGCTTTTGTGTGTA AAACTATTTTTCAGTACTATCTGAACA C	804
86	SPI01_SPE111_1314_1	CCGCAAAATACCGCTATTATGCG	68	SPI01_SPE111_1403_1431_R	GGAATATTTTTCAGTACTATCTGAACA AG	711
87	SPI01_SPE111_1408_1	CGAGTATAGTAAATAAGATTAT GACA	67	SPI01_SPE111_1486_1515_R	GGAATATTTTTCAGTACTATCTGAACA AG	828
88	SPI01_SPE111_1688_1	CCATATTAATGGTTTTCAGAACTG GCT	60	SPI01_SPE111_1783_1808_R	ATATGATATATCTGAACTGGCGCG G	762
89	SPI01_SPE111_1711_1	CTGCTGTAATCTTGGCGAGTGG	82	SPI01_SPE111_1809_1835_R	GGTGGAGGACTCTTGAATGTAATCA A	821
90	SPI01_SPE111_1807_1	ATGATATGATATCAAGAGAGTCTCA CGC	33	SPI01_SPE111_1901_1927_R	TTGGACTGTATGAGTGTGATCTGG G	1412
91	SPI01_SPE111_1967_1	TTAGCGTCTATGCGCCAGATGGG	155	SPI01_SPE111_2062_2083_R	ATATGCCAGAAATCAATCATC	755
92	SPI01_SPE111_2260_2	CAGAGAGCGTTTATCTCTATCAGC	50	SPI01_SPE111_2375_2397_R	TCCTGGGTGACCTGGTGTATGA	1131
93	SPI01_SPE111_2375_2	TCCTAAACACACCGTCAACCGAG	390	SPI01_SPE111_2470_2497_R	AGCTCTGATGATGAGCTCTGCGCTGGC	747
94	SPI01_SPE111_2468_2	ATGCGCTGCGAGAGCTCA	35	SPI01_SPE111_2543_2570_R	CCATATGATCAACGTCACATCAAGC	770
95	SPI01_SPE111_2961_2	ACCATGACGAGAGCGATTTTGACA	15	SPI01_SPE111_3023_3045_R	GGATTTTACGCGATGAGACGCG	827
96	SPI01_SPE111_3075_3	GATGATCTTTTATGATGATGTCAGC AGC	108	SPI01_SPE111_3168_3196_R	ATGCGAGGACCTCTGGAGAGTTCT C	715
97	SPI01_SPE111_3386_3	AGCGTAAAGGTGAACCTT	25	SPI01_SPE111_3480_3506_R	CGACGAGTACTCTGCTCATCTTTG	769
98	SPI01_SPE111_3511_3	AGCTGTGAGATCAATGATGACAG GCA	116	SPI01_SPE111_3605_3629_R	GGTCTTACCTGCACTCTGCTTAC GC	832
111	RPOB_EC_3775_1803_F	CGATGCGTCTGCTCTAC	87	RPOB_EC_3829_3858_R	CGATATAGCTCAACATGAGCTTGAAT GC	797
112	VALS_EC_1433_1850_F	CGATGCGTCTGCTCTAC	65	VALS_EC_1920_1943_R	GGTTCACAGCTTTGTCAGAG	822
113	RPOB_EC_1936_1353_F	CGATGCGTCTGCTCTAC	97	RPOB_EC_1438_1455_R	TTGCTCTGCGCTGGCG	1386
114	TUFB_EC_225_251_F	CGATGCGTCTGCTCTAC	111	TUFB_EC_284_309_R	TATAGCACATCATCTGAGCGGAC	930
115	DNAX_EC_428_445_F	CGATGCGTCTGCTCTAC	72	DNAX_EC_503_522_R	CGGCTGCGCTGTTGATGA	792
116	VALS_EC_1920_1943_F	CGATGCGTCTGCTCTAC	85	VALS_EC_1948_1970_R	TGCGAGTCTGATGAGAGAGG	1075
117	TUFB_EC_757_774_F	AGAGAGCGTCTGCTGGCG	6	TUFB_EC_849_867_R	GGGCTCAATCTCTTCAAGC	819
118	238_EC_2646_2667_F	CGTCTCTATGATGAGAGGAC	84	238_EC_2745_2765_R	TTGCTGCTGATGATGCTTTTAC	1389
119	169_EC_969_985_1P_F	AGCGAGAGCTCTTAC	19	169_EC_1061_1078_2P_R	AGGACAGAGCTGTGGAGAC	733
120	169_EC_972_985_2P_F	CGAGAGAGCTCTTAC	63	169_EC_1064_1075_2P_R	AGCGAGAGCTGTGGAGAC	727
121	169_EC_972_985_2P_F	CGAGAGAGCTCTTAC	63	169_EC_1064_1075_2P_R	AGCGAGAGCTGTGGAGAC	727
122	DNAX_EC_32_50_2_F	CGATGATGAGGTGAGTGG	61	238_EC_40_59_R	AGGCTCTGATGCTGCTCTCA	740

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

123	235 EC -7 15 F	GTGTGGGTTGTAAGCACTTACG	140	235 EC 430 450 R	CTATCGTCTGCTGCAAGATAT	799
124	235 EC -7 15 F	GTGTGGGTTGTAAGCACTTACG	141	235 EC 891 910 R	TTCGTCTGGGTTGGTAGTTC	1403
125	235 EC 430 450 F	ATATCTCTGTAAGCACTTACG	30	235 EC 1424 1442 R	AACTATGACTCTTCGTCTC	712
126	235 EC 891 910 F	GAATCTACCAAGCACTTACG	100	235 EC 1908 1931 R	TACTCTAAGACGGTATAGTACG	893
127	235 EC 1424 1442 F	GGACCAAGACCACTTACG	117	235 EC 2475 2494 R	CGAAACACCGCGGTGATAT	765
128	235 EC 1908 1931 F	GTCTACTATACGGTCTTACGATTA	73	235 EC 2833 2853 R	GTTTACACACCGCGGTATAC	826
129	235 EC 2475 2494 F	ATATCGACCGGGTGTTCG	31	165 EC 23 41.2 R	GGTGTGACGAGGATATTC	820
131	165 EC -60 -59 F	AGTCTGAAGGTGAGCACTTACG	28	165 EC 508 525 R	GTCTGTGACGAGGATTA	823
132	165 EC 376 345 F	GACCGGTCTCGAGCTCTTAC	95	165 EC 1041 1059 R	CAATCGACCACTCTCTC	771
133	165 EC 705 724 F	GATCTGAGAGTACCGGTG	107	165 EC 1493 1512 R	ACGGTACCTTGTATACGACT	739
134	165 EC 1268 1287 F	GAGCGACACCGCACTCTATA	101	TRNA ALA--	CTCTCTGGTGTGCAAGC	780
135	165 EC 969 985 F	ACGCGAGACCACTTAC	19	165 EC 1061 1078.2 R	ACACACGAGCTGTGAGAC	719
137	165 EC 969 985 F	ACGCGAGACCACTTAC	19	165 EC 1061 1078.2 114 R	ACACACGAGCTGTGAGAC	721
138	165 EC 969 985 F	ACGCGAGACCACTTAC	19	165 EC 1061 1078.2 112 R	ACACACGAGCTGTGAGAC	718
139	165 EC 969 985 F	ACGCGAGACCACTTAC	19	165 EC 1061 1078.2 111 R	ACACACGAGCTGTGAGAC	722
140	165 EC 969 985 F	ACGCGAGACCACTTAC	19	165 EC 1061 1078.2 116 R	ACACACGAGCTGTGAGAC	720
141	165 EC 969 985 F	ACGCGAGACCACTTAC	19	165 EC 1061 1078.2 21 R	ACACACGAGCTGTGAGAC	723
142	165 EC 969 985 F	ACGCGAGACCACTTAC	19	165 EC 1061 1078.2 31 R	ACACACGAGCTGTGAGAC	724
143	165 EC 969 985 F	ACGCGAGACCACTTAC	19	165 EC 1061 1078.2 41 R	ACACACGAGCTGTGAGAC	725
147	235 EC 2452 2669 F	CTATGACGCTGAAATGCG	79	235 EC 2741 2760 R	ACTTAGAGTGTCTTACGGGT	743
156	165 EC 693 700 F	GTATGACGCTGAAATGCG	137	165 EC 880 894 R	CTACTATCCGAGGCG	796
159	165 EC 1100 1115 F	CGACGAGCTGAAATGCG	42	165 EC 1174 1188 R	TCCGACCTCTCTCC	1019
215	SSP8 BA 121 137 F	ACGCGACATGCAATGCG	3	SSP8 BA 197 216 R	TGTTGTGTTGTGAAATTC	1132
220	GR0L EC 941 959 F	TGGAAGATCTGGTGTGCG	544	GR0L EC 1039 1060 R	GAATCTGTGAGAGTGTGAGC	759
221	INF8 EC 1103 1124 F	GTCTGTGAAACGAGACCTGAGA	133	INF8 EC 1174 1191 R	CATGTATGTCTACACCGG	764
222	HF18 EC 1082 1102 F	TGCGACACTCTGTACAGACGA	569	HF18 EC 1144 1168 R	CTTTCGGTCTTCTCAACTCAACAT	802
223	INF8 EC 1969 1994 F	CTTCACTGTAAATCTCTGAACTTAA	74	INF8 EC 2038 2058 R	AACTTCGGCTTGGTGTGATGTT	713
224	GR0L EC 215 243 F	GTGTGAAGATCTGGTGTGAGAC	128	GR0L EC 328 350 R	TCTAGGCTGACGGTGTGAGCC	1377
225	VALUS EC 1105 1124 F	CTTTCGGCTGTGTGTATGCA	77	VALUS EC 1195 1214 R	ACCACTGTGATGTGCGGTT	732
226	165 EC 556 575 F	CGGAATATCTGGTGTGAGAC	70	165 EC 693 700 R	CGATTTATACCGCTGAGC	791
227	RPOC EC 1256 1277 F	ACCGAGTCTGTGAGACCTGTC	16	RPOC EC 1295 1315 R	CGATTTATACCGCTGAGC	843
228	165 EC 774 795 F	GGAGACACAGAGATGTATAC	122	165 EC 880 894 R	GTACTTCCCGAGGCG	796
229	RPOC EC 1584 1604 F	TGCGCCGAAAGAGATGTAGAG	567	RPOC EC 1623 1643 R	ACGCGGCACTGAGAGATGCG	737
230	165 EC 1082 1100 F	AGTTGGGTGTAGTATGAG	37	165 EC 1177 1196 R	TGAGTGTATCCCACTCTTC	1158
231	165 EC 1389 1407 F	CTTGTACACATCCCGCTTC	88	165 EC 1525 1541 R	AAAGAGGTGATCCAGCC	714

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

222	165 BC 1303 1323 F	CGHNTGAGTCTGCTACTCG	71	165 BC 1389 1407 R	GACGGGCGTGTGTACAG	808
223	235 BC 21 37 20 F	GGTGTGATCGCTCTTSC	129	235 BC 115 130 R	GGTTCCTCCGATTCGG	833
224	235 BC 137 207 F	GGGAACTGAAATCTCTAGTA	121	235 BC 242 256 R	TTCTCTCTCCGCTGTC	1385
225	235 BC 1602 1620 F	TACCCAGAACAGACGAG3	184	235 BC 1686 1703 R	CTCTCTCCGAGATTCG	762
226	235 BC 1683 1703 F	CGCTACTCTCGGAAAG3	58	235 BC 1828 1842 R	CACCGGCGAGCGCTC	760
227	235 BC 1827 1843 F	CGCTACTCTCGGAAAG3	99	235 BC 1919 1949 R	CGACAGAGAAITTCGCTAC	775
228	235 BC 2434 2456 F	GAGCTCTCCGCTGTC	9	235 BC 2490 2511 R	ACGACACATGAGTTCGAC	746
229	235 BC 2559 2616 F	GAGCTCTCCGCTGTC	96	235 BC 2653 2669 R	CGCTCTCTCTGCTACT	777
230	235 BC 2653 2669 F	GAGCTCTCCGCTGTC	227	235 BC 2737 2758 R	TTGAGTCTCTCAGCTCTATC	1369
241	235 BC -68 -44 F	TAGTACAGAGACGAGACAC	1	235 BC 5 21 R	GTGCGCTCTTCTATC	841
242	165 BC 8 27 F	AGATTTGATCTAGCTGCTAG	23	165 BC 342 358 R	ACTGCTCTCTCCGCTAG	742
243	165 BC 314 332 F	CTCTGGAATCTAGAG3	43	165 BC 556 575 R	CTTTACGCGAGTATTCG	801
244	165 BC 518 536 F	CGGAGCGGCGGATATC	57	165 BC 774 795 R	GTATCTACCTCTGTTTCTCC	839
245	165 BC 683 700 F	GTGTAGCGGTGAAATTCG	137	165 BC 967 985 R	GGTAAAGTCTCTCGCTG	835
246	165 BC 937 954 F	TAAGCGGTGAGCATTCG	7	165 BC 1220 1240 R	ATTGTAGCGGTGTGTACGC	714
247	165 BC 1195 1213 F	CAGTCTATCAAGCGCTTAT	46	165 BC 1525 1541 R	AAAGGAGTATCGAGCC	714
248	165 BC 8 27 F	AGATTTGATCTAGCTGCTAG	23	165 BC 1919 1936 R	TCGCTACTCTTGGACGCT	1080
249	235 BC 1831 1849 F	ACTCTCGATCTCTGAG	18	235 BC 1494 1513 R	CAGGCTACTCTGTATACAC	761
250	165 BC 1387 1407 F	CTCTGTACACAGCTCTCTCTC	112	165 BC 1486 1505 R	CTCTGTGTGAGCTCTACCC	783
251	165 BC 1390 1411 F	TTGTACACAGCGCGGTGTATC	693	165 BC 1486 1506 R	ACTCTGTGTGAGCTCTACCC	731
252	165 BC 1367 1387 F	TACGTGAAATCTGTTCCCG66	191	165 BC 909 929 R	CGCTGTGAGCTCTCTCTGCT	772
253	165 BC 804 822 F	ACACGCGGTAAAGATAT	14	165 BC 886 904 R	CGCTGTGAGCTCTCTCTGCT	817
254	165 BC 791 812 F	GATACCTGTGTGTCTACG3	206	165 BC 882 899 R	GCACCTGACTCTCTCTGCT	818
255	165 BC 789 810 F	TAGTACCTGTGTGTCTACG3	106	165 BC 1174 1195 R	GAGCTCATCCGCTCTCTCTGCT	810
256	165 BC 1092 1109 F	TAGTACCTGTGTGTCTACG3	228	235 BC 2658 2677 R	ATAGCCATCTGTGTGTCTATC	749
257	235 BC 2186 2607 F	TTGAGTGTGTGTGTGTGTGTG	203	235 BC 359 379 R	ATAGCCATCTGTGTGTCTATC	750
258	235 BC 31 49 F	GAGGAAATCTGTGTGTCTAC	103	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
259	235 BC 31 49 F	GAGGAAATCTGTGTGTCTAC	103	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
260	235 BC 31 49 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
261	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
262	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
263	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
264	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
265	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
266	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
267	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
268	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
269	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
270	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
271	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
272	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
273	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
274	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
275	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
276	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
277	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
278	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
279	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
280	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
281	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
282	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
283	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
284	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
285	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
286	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
287	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
288	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
289	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
290	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
291	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
292	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
293	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
294	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
295	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
296	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
297	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
298	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
299	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
300	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
301	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
302	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
303	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
304	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
305	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
306	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
307	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
308	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
309	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
310	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
311	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
312	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
313	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
314	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
315	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
316	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
317	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
318	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
319	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
320	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
321	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
322	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
323	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
324	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
325	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
326	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
327	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
328	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
329	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
330	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
331	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
332	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
333	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
334	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
335	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
336	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
337	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
338	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
339	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
340	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
341	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
342	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
343	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
344	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
345	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
346	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
347	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
348	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
349	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
350	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
351	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
352	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
353	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
354	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
355	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
356	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
357	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
358	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 363 384 R	ATAGCCATCTGTGTGTCTATC	838
359	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 358 379 R	ATAGCCATCTGTGTGTCTATC	750
360	235 BC 43 61 F	GAGGAAATCTGTGTGTCTAC	104	235 BC 345 367 R	ATAGCCATCTGTGTGTCTATC	751
361	235 BC 43 61 F	GAG				

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

260	RNASP EC 61_77_F	GAGGAAGTCAGCGCTC	105	RNASP EC 345_362_R	ATAGCCCGGTTCTGTCG	751
262	RNASP SA 31_49_F	GAGGAAGTCAGCTCTCC	103	RNASP SA 358_379_R	ATAGCCCGGTTCTGCTCTC	750
263	16S EC 1082_1100_F	ATGTGGGTAACTGCTCCG	37	16S EC 1535_1541_R	AGAGAGGTATCCGCGC	714
264	16S EC 556_575_F	CGGATTAATCGGTGGGTAAAG	70	16S EC 774_795_R	GTATCTAATCTGTTGTCTCC	839
265	16S EC 1082_1100_F	ATGTGGGTAACTGCTCCG	37	16S EC 1177_1196_100_R	TGAGCTATCCGCACTCTCC	1160
266	16S EC 1082_1100_F	ATGTGGGTAACTGCTCCG	37	16S EC 1177_1196_100_110_R	TGAGCTATCCGCACTCTCC	1161
268	VARF_EC 511_532_F_N	GGTGTAAATAGACCTGGG	130	THRA_ALA		744
269	MOD	16S EC 1082_1100_F	37	RHH EC 30_49_F MOD	AGACTCTCTGCTGCAAGC	1158
270	MOD	23S EC 2586_2607_F	203	16S EC 1177_1196_R MOD	TGAGCTATCCGCACTCTCC	749
272	16S EC 969_985_F	TAGACGTGCGAGCGAGTTCC	19	23S EC 2656_2677_R MOD	AGTCTCTCCGCTCTCTCCG	807
273	16S EC 683_700_F	ACGGAGAGACTCTTACC	137	16S EC 1189_1407_R	GAGGCGGCTGTGTACAG	788
274	16S EC 49_68_F	GTGTAGCGTGAATTCGC	152	16S EC 1303_1323_R	CGAGTTGCAAGTGCAGTCG	796
275	16S EC 49_68_F	TACACATGCACTGTAAG	152	16S EC 880_894_R	CGTATCTCCGAGCG	734
277	CYA BA 1349_1370_F	TACACATGCACTGTAAG	12	16S EC 1061_1078_R	ACACAGAGACTTACAC	800
278	F	ACACAGAGACTTACAC	650	CYA BA 1426_1447_R	CTTCTACTTTTACGCACTAC	1159
279	16S EC 405_432_F	TAAATCTCCGACAGCGACAA	464	16S EC 1175_1196_R	TGAGCTATCCGCACTCTCTC	793
280	GHOL EC 496_516_F	TAGATGTGACAGCTTAGGGTGA	34	16S EC 507_527_R	CGCTCTCTGCGCAAGATFAG	914
281	GHOL EC 511_536_F	ATGCGACAGTGTGCGACAGAG	8	GHOL EC 577_596_R	TACCCGCGCTCGAATTCAT	776
288	RPDB EC 3802_3821_F	AUGAGGGCGGTATCCCTTTGAGA	51	GHOL EC 577_593_R	CGCCGCTGGAATCTGAGCTCTC	786
289	RPDB EC 3799_3821_F	CGCGTCTCGGCAATGGA	124	RPDB EC 3862_3885_R	CGACTTGAAGTTAAGATTCCTG	840
290	RPDB EC 2146_2174_F	CGGAGCTGTTGCGGAAATGGA	52	RPDB EC 3862_3888_R	GTCCGACTTGAAGTGCATCTCTCG	736
291	APPS EC 405_422_F	CHGAGTCTGTAATCTGATCTACT	110	RPDB EC 2227_2245_R	ACGCTATCGACGACGCT	738
292	APPS EC 1374_1393_F	GCACACACTTACGCTCTCG	69	APPS EC 521_538_R	ACGCGACAGGTATFATGCG	811
293	TUPE EC 957_979_F	CGCGAGCTCTGTAATCTGATCT	55	RPDB EC 1437_1455_R	GACATCGAGCTGCTGCTCT	829
294	16S EC 7_33_F	CACAGCTGCTCTGCTGTAACGA	102	TUPE EC 1034_1056_R	GGCATCACTATTCTCTCTCTCTCG	1345
295	VALS EC 610_649_F	GAGATTTGACCTTGCTGTAACGA	17	16S EC 101_122_R	TGTTACTACACCGCTCTGCGACT	929
344	16S EC 971_980_F	ACCGAGAGAGAGACGACG	113	VALS EC 705_727_R	TATNACCACTGCTCGAGGTGA	736
345	16S EC 713_732_TWOD	GGGAGAACTTACCGGTC	202	16S EC 1043_1062_R	ACACACTGCACTCTCTCTC	1110
347	16S EC 785_806_TWOD	TAGAACACGATGCGGAGGC	560	16S EC 789_809_TWOD_R	TGCTTGACCTACACGGGTATCTA	1278
348	16S EC 960_981_TWOD	TGATTAGAGACCTGCTGATGTC	706	16S EC 880_897_TWOD_R	TGCGCTACTCTCTCTCTCTCTC	895
		TTTGATTCGACGCGAGAACT		16S EC 1054_1073_TWOD_R	TACGAGCTGACGACGCGCTG	

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

[illegible]

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

428	SP101_SPT11_1688_1 716_TMOD_F	TCTGCTATATATGCTTTACGAAACT GGCT	334	SP101_SPT11_1783_1808_TM OD_R	TATATGATATATCAATGAACTGCGCG	932
429	SP101_SPT11_1711_1 733_TMOD_F	TCTGCTAAATATGTCGCAAGGT TATGATATCAATCAAGAGAGGCTC	406	SP101_SPT11_1808_1835_TM OD_R	TCTGTAAGAGACTCTTGAATGTATTC A	1239
430	SP101_SPT11_1807_1 835_TMOD_F	ACGC	235	SP101_SPT11_1901_1927_TM OD_R	TTTGGACTGTATATGCTGATATCTGG	1439
431	SP101_SPT11_1967_1 991_TMOD_F	TTTACGCTTATCTAGCGCCGATGCG TTCGAGCTGTGAAATCGCGCAAG	649	SP101_SPT11_2062_2083_TM OD_R	TATGCGCCGAAATCAAAATCATC	940
432	SP101_SPT11_216_24 3_TMOD_F	ATT	210	SP101_SPT11_308_333_TMOD R	TTTCCACTTTTACACTCTCTGTCTGTC	1404
433	SP101_SPT11_2260_2 283_TMOD_F	TCTGAGAGCCGTTTATCTATCACG	272	SP101_SPT11_2375_2397_TM OD_R	TTCTGSGAGCTGCGGTGTATGA	1393
434	SP101_SPT11_2375_2 399_TMOD_F	TTCTTAACACACGCTGACCGAGAG	675	SP101_SPT11_2470_2497_TM OD_R	TACTCTGAAATGAGCTTCTGACGAGGC C	918
435	SP101_SPT11_2468_2 487_TMOD_F	TATGAGCCATGCGGAGCTCA	238	SP101_SPT11_2543_2570_TM OD_R	CTCTGTAAGGTCACTGACCATTCGAAG C	1007
436	SP101_SPT11_266_29 5_TMOD_F	TCTGTATCTTGGCTGACAGGCTG	417	SP101_SPT11_355_380_TMOD R	TCTGCTTTGATGCTGTAATCCCTTTC	1249
437	SP101_SPT11_2961_2 984_TMOD_F	TATCATGACAGAGGCACTTTTGACA	183	SP101_SPT11_3023_3045_TM OD_R	TGGAATTTACCGCACTATGACGCC	1264
438	SP101_SPT11_3075_3 103_TMOD_F	TGATGACTTTTATGCTATGCTGAG	473	SP101_SPT11_3168_3196_TM OD_R	TATGTGAGACCACTTGGAAAGATTC TC	875
439	SP101_SPT11_322_34 4_TMOD_F	TGTCMAAGTGCCTATTTATCTGCG	631	SP101_SPT11_423_441_TMOD R	TATCCCTGCTCTCTGCTGCG	934
440	SP101_SPT11_3386_3 403_TMOD_F	TAGCTATAGGTGATACCTT	215	SP101_SPT11_3480_3506_TM OD_R	TCCGACGATTTACTGCTGCTCATCTTG	1005
441	SP101_SPT11_3811_3 TACTCTGAGATCACTGATGAGCGAG	531	SP101_SPT11_3605_3629_TM OD_R	TOGCTCTACACCTGACTGCTGATAC	1294	
442	SP101_SPT11_358_38 7_TMOD_F	TGAGATTTCACTCACTGATGATCTTTT	588	SP101_SPT11_448_473_TMOD R	TCTGAACCTTTTCAACAGAGATCAGC	998
443	SP101_SPT11_600_62 9_TMOD_F	TCTGTATCTGAACTATGATCTTTT	348	SP101_SPT11_686_714_TMOD R	TCCCATTTTTCACGATCTGGAATA TC	1018
444	SP101_SPT11_658_68 4_TMOD_F	TCTGATATGATATCACTGATGAGAG	589	SP101_SPT11_756_784_TMOD R	TGATTTGGCATTAAGTATTTTCTAA AA	1189
445	SP101_SPT11_776_80 1_TMOD_F	TTTGCCTCACTCAAACTATGAGGATGG C	673	SP101_SPT11_871_896_TMOD R	TCCCAACAGAAAGCTATGAGCATTA	1217
446	SP101_SPT11_1_29_T MOD_F	TACTCTATTATTGAAAGAACCTAG AAGT	154	SP101_SPT11_92_116_TMOD R	TCTTACCACTGCTTCCAGAGCGAG	1044
447	SP101_SPT11_364_38 5_F	TGAGCCATCTGAGCGCATTTG	276	SP101_SPT11_448_471_R	TACTTTTTCACACAGAAATCAGC	894
448	SP101_SPT11_3085_3 104_F	TAGTATATGCTGCGCGGCC	216	SP101_SPT11_3170_3194_R	TCTGAGAGCACTTCTGGAAGATTTTC	1066
449	REF_EC_590_710_F	TCTACAGCTGTGTGTGAGAG	309	REF_EC_737_758_R	TGTGCTGTTCACCGGAGAG	1336
481	BOWTA_X52066_536_55 2_F	TATGGCTCTACTGCA	239	BOWTA_X52066_647_660_R	TGTTACTCTGCGAT	1346

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO.: 10593)

482	BOWT_X52066_538_55	BA*TP*SGC*TP*CP*TPA*CP*TP*C BA	143	BOWT_X52066_647_660P_R	TP*TPA*CP*TPA*CP*TP*SGCT	1146
483	BOWT_X52066_701_72	GAATGACANTATCAAAAT	94	BOWT_X52066_759_775_R	TPACTCTPACCOMCTC	1367
484	BOWT_X52066_701_72	GA*TPAG*CPAA*TP*TPAA*TP*C P*CPAAAT	91	BOWT_X52066_759_775P_R	TPA*CP*TP*CP*TPAA*CP*CP*C PA*CP*TPC	1359
485	BOWT_X52066_450_47	TCRTGTAAATATGACCTGACC	393	BOWT_X52066_517_533_R	TPACCAATTCGCTGAGTCAA	859
486	BOWT_X52066_450_47	TP*CP*TPAG*CPAA*TP*TPAA*TP*C P*CPAAAT	142	BOWT_X52066_517_539P_R	TPACCA*TP*TP*TP*CPGCTGAGAT	857
487	BOWT_X52066_591_62	TCRTGTAAATATGACCTGACC	463	BOWT_X52066_644_671_R	P*TP*CPAA	992
608	SSPE BA 156 168P_F	TCRTGTAAATATGACCTGACC	616	SSPE BA 243_255P_R	TCRTGTAAATATGACCTGACCTG	1241
609	SSPE BA 75 89P_F	TPCQAGAGTTP*TP*CPGACC	192	SSPE BA 163_177P_R	TCQAGCTGAGTTPCT	1336
610	SSPE BA 150 168P_F	TCRTGTAAATATGACCTGACC	533	SSPE BA 243_264P_R	TCQAGCTTP*TP*CPAA*TP*CPCT	1191
611	SSPE BA 72 89P_F	TCRTGTAAATATGACCTGACC	602	SSPE BA 163_182P_R	TCQAGCTTP*TP*CPAA*TP*CPCT	995
612	SSPE BA 114 137P_F	TCQAGCTTP*TP*CPGACC	255	SSPE BA 196_222P_R	TCQAGCTTP*TP*CPGACCTGACCT	1401
699	SSPE BA 123 153_F	TCQAGCTTP*TP*CPGACC	488	SSPE BA 202_231_R	TCQAGCTTP*TP*CPGACCTGACCT	1431
700	SSPE BA 156 168_F	TCQAGCTTP*TP*CPGACC	612	SSPE BA 243_255_R	TCQAGCTTP*TP*CPGACCTGACCT	1302
701	SSPE BA 75 89_F	TCQAGCTTP*TP*CPGACC	179	SSPE BA 163_177_R	TCQAGCTTP*TP*CPGACCTGACCT	1338
702	SSPE BA 150 168_F	TCQAGCTTP*TP*CPGACC	533	SSPE BA 243_264_R	TCQAGCTTP*TP*CPGACCTGACCT	1190
703	SSPE BA 72 89_F	TCQAGCTTP*TP*CPGACC	600	SSPE BA 163_182_R	TCQAGCTTP*TP*CPGACCTGACCT	995
704	SSPE BA 146 163_F	TCQAGCTTP*TP*CPGACC	484	SSPE BA 242_267_R	TCQAGCTTP*TP*CPGACCTGACCT	1421
705	SSPE BA 63 89_F	TCQAGCTTP*TP*CPGACC	518	SSPE BA 163_191_R	TCQAGCTTP*TP*CPGACCTGACCT	986
706	SSPE BA 114 137_F	TCQAGCTTP*TP*CPGACC	255	SSPE BA 196_222_R	TCQAGCTTP*TP*CPGACCTGACCT	1402
770	PLA_A053945_7377_7	TCQAGCTTP*TP*CPGACC	442	PLA_A053945_7434_7462_R	TCQAGCTTP*TP*CPGACCTGACCT	1313
771	PLA_A053945_7380_7	TCQAGCTTP*TP*CPGACC	327	PLA_A053945_7482_7502_R	TCQAGCTTP*TP*CPGACCTGACCT	1304
772	PLA_A053945_7481_7	TCQAGCTTP*TP*CPGACC	481	PLA_A053945_7539_7562_R	TCQAGCTTP*TP*CPGACCTGACCT	943
773	PLA_A053945_7486_7	TCQAGCTTP*TP*CPGACC	657	PLA_A053945_7577_7600_R	TCQAGCTTP*TP*CPGACCTGACCT	879
774	CAPI_A053947_33407	TCQAGCTTP*TP*CPGACC	292	CAPI_A053947_33494_33511	TCQAGCTTP*TP*CPGACCTGACCT	1235
775	CAPI_A053947_33515	TCQAGCTTP*TP*CPGACC	270	CAPI_A053947_33595_33621	TCQAGCTTP*TP*CPGACCTGACCT	1053
776	CAPI_A053947_33435	TCQAGCTTP*TP*CPGACC	542	CAPI_A053947_33495_33517	TCQAGCTTP*TP*CPGACCTGACCT	1183

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

777	CNFI_A003947_33607 33716_F	TGAGAGTGAATTAACCAACATCA CTGC	286	CNFI_A003947_33755_33782 R	TGAGGTCTCTACCGTTTACTTACGAG	962
778	INV_U22457_515_539_	TGCTCCCTGTATATGACTCTGCTTC	573	INV_U22457_571_598_R	TGTTAAGTGTGTGCGGTGTCTTTAT	1343
779	INV_U22457_689_724_	TGCTGAGGCTGAGCACTATATTAC	525	INV_U22457_753_776_R	TGCGCGAGGTGCAATCAATG	976
780	INV_U22457_834_850_	TTATTACTCTGACTTCCCAACTG	664	INV_U22457_942_966_R	TGACCAAGCTGAAGCTTACTG	1154
781	INV_U22457_1058_1088_	TGTTAACTGAGCTTTTATAGGCA	597	INV_U22457_1619_1643_R	TGCTGTTCAGATTAATCTTACAA	1408
782	LL_NC003143_2366996	TGAGCGCTGAGCACTACATCTC	637	LL_NC003143_2367073_23670	TCTATCCGATTAATACGCACTGA	1123
783	LL_NC003143_2367172	TGAGCGCACTGAGCACTCTAC	550	LL_NC003143_2367249_23672	TGCGCAAGCTCTACACCTTTGG	1272
784	RELB_EC_649_679_F	TGCGCTGCTGTTGGGCTCTGTAAT	620	RELB_EC_733_762_TM0D_R	TTCCAGTGTCTGTTTACCCCTGG	1380
785	RELB_EC_642_6799_F	TGCGCTGTGTGTGTGCTGCTGTT	646	RELB_EC_733_762_TM0D_R	TTCCAGTGTCTGTTTACCCCTGG	1380
786	MECA_Y14051_3311_3	TGAGCTGTGTAATGAC	653	MECA_Y14051_3367_3393_R	TGCGATGCGAGGTGAGAGGTGA	1333
787	MECA_Y14051_3774_38	TGAGCACTATGCTGAGCACTGACTG	144	MECA_Y14051_3828_3854_R	TCCCACTGACTTCCACATACATCT	1015
788	MECA_Y14051_3774_38	TGAGCACTATGCTGAGCACTGACTG	144	MECA_Y14051_3690_3719_R	TGATCTGTAATGTTTATCTTTACGC	1181
789	MECA_Y14051_3645_36	TGAGCACTATGCTGAGCACTGACTG	434	MECA_Y14051_4555_4581_R	TGCTATGAGCTCATATGAGGTGCT	1269
790	MECA_Y14051_4507_45	TGAGCACTATGCTGAGCACTGACTG	288	MECA_Y14051_4586_4610_R	TATCTCTGCTGACTGATCCATCA	939
791	MECA_Y14051_4510_45	TGAGCACTATGCTGAGCACTGACTG	626	MECA_Y14051_4765_4793_R	TACACCCCGAGATTAATCTTTTCC	858
792	MECA_Y14051_4669_46	TGAGCACTATGCTGAGCACTGACTG	262	MECA_Y14051_4520_4600_R	TGCTGTGCTGCTGCTGCTGCT	1357
793	MECA_Y14051_4520_45	TGAGCACTATGCTGAGCACTGACTG	389	MECA_Y14051_4600_4610_R	TGCTGTGCTGCTGCTGCTGCT	1358
794	MECA_Y14051_4520_45	TGAGCACTATGCTGAGCACTGACTG	389	TRPE_AY094355_1551_1580_R	TGCGCGAGCTTTTATTTGGGTTC	1231
795	TRPE_AY094355_1467_	TGAGCACTATGCTGAGCACTGACTG	36	TRPE_AY094355_1592_1592_R	TATTTGGTTTCTTCACTGAGATCT	944
796	TRPE_AY094355_1445_	TGAGCACTATGCTGAGCACTGACTG	557	TRPE_AY094355_1581_1580_R	TGCTCTTTTCCAGGCTCTACTCATC	1040
797	TRPE_AY094355_1471_F	TGAGCACTATGCTGAGCACTGACTG	247	TRPE_AY094355_1392_1418_R	TGAGGTGTGCTGCTGCTGCTGCT	885
798	TRPE_AY094355_1278_	TGAGCACTATGCTGAGCACTGACTG	357	TRPE_AY094355_769_791_R	TGCGAGGTGCTGCTGCTGCTGCT	1372
799	TRPE_AY094355_1064_	TGAGCACTATGCTGAGCACTGACTG	135	TRPE_AY094355_864_883_R	TGCGAGGTGCTGCTGCTGCTGCT	1218
800	TRPE_AY094355_666_	TGAGCACTATGCTGAGCACTGACTG	483			

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

1055	CJST_CJ_2869_2895_F	TKAGCTGTGTCTTACAGAGAGCTTC	432	TGACCTGTGTCTTACAGAGAGCTTC	1045	TCGCTGTGTCTTACAGAGAGCTTC	A
1056	CJST_CJ_1880_1910_F	KCCGATATATCTGCGATTTTCCCA	317	KCCGATATATCTGCGATTTTCCCA	1309	TCGCTGTGTCTTACAGAGAGCTTC	CCA
1057	CJST_CJ_2185_2212_F	TKAGNAGAGAGAGAGAGAGAGCTAT	208	TKAGNAGAGAGAGAGAGAGAGCTAT	1152	TATAGT	TATAGT
1058	CJST_CJ_1643_1670_F	TTATGCTTTGTGAGAGCTATGCTTAT	660	TTATGCTTTGTGAGAGCTATGCTTAT	1198	TGCAATGTGTGTATGTACAGAGAGAG	F
1059	CJST_CJ_2165_2194_F	TKGCGATCTTGTGTGTGTGTAGAG	511	TKGCGATCTTGTGTGTGTGTAGAG	1002	TCGACAGCTGTGTATGTATGTCTTGT	CCGCT
1060	CJST_CJ_599_632_F	TGAAAAATATCCAGAGAGAGAGAGCA	424	TGAAAAATATCCAGAGAGAGAGAGCA	1024	TCGACAGAGAGAGAGAGAGAGAGAT	TTTAC
1061	CJST_CJ_360_393_F	TCGCTGTATCTCGATAGATATATATTC	345	TCGCTGTATCTCGATAGATATATATTC	882	TACAGATCTTGTGCTTCATACAGAGCT	AATTTTC
1062	CJST_CJ_2678_2703_F	TCGCGAGAGAGAGAGAGAGAGAGATTC	321	TCGCGAGAGAGAGAGAGAGAGAGATTC	1319	TGTCGCTTTGTGCTTCATACAGAGCT	ATTTTC
1063	CJST_CJ_1268_1299_F	AGTTATATATATATATATATATATATAT	29	AGTTATATATATATATATATATATATAT	1096	TCGCTTTTACGCTATGATGATGAGAG	ATA
1064	CJST_CJ_1580_1713_F	TCGCTGTATATATATATATATATATAT	479	TCGCTGTATATATATATATATATATAT	938	TATGTGTATGTATGCTTACTACAGAGCT	TCGCTTCAAAGAGATTTTACATTTTC
1065	CJST_CJ_2857_2887_F	TCGCTGTATATATATATATATATATAT	565	TCGCTGTATATATATATATATATATAT	1253	TTAAG	TTAAG
1070	RNASEP_BKM_560_599_F	TCGCTGTATATATATATATATATATAT	512	TCGCTGTATATATATATATATATATAT	1034	TCGAGAGAGAGAGAGAGAGAGAGAG	TCGAGAGAGAGAGAGAGAGAGAGAG
1071	RNASEP_BKM_615_637_F	TCGCTGTATATATATATATATATATAT	333	TCGCTGTATATATATATATATATATAT	1222	TCGCTGTATATATATATATATATATAT	TCGCTGTATATATATATATATATATAT
1072	RNASEP_BDP_574_592_F	TCGCTGTATATATATATATATATATAT	561	TCGCTGTATATATATATATATATATAT	1115	TCGCTGTATATATATATATATATATAT	TCGCTGTATATATATATATATATATAT
1073	23S_BKM_1110_1129_F	TCGCTGTATATATATATATATATATAT	510	TCGCTGTATATATATATATATATATAT	1074	TCGCTGTATATATATATATATATATAT	TCGCTGTATATATATATATATATATAT
1074	23S_BKM_515_536_F	TCGCTGTATATATATATATATATATAT	496	TCGCTGTATATATATATATATATATAT	1088	TCGCTGTATATATATATATATATATAT	TCGCTGTATATATATATATATATATAT
1075	RNASEP_CLB_459_487_F	TCGCTGTATATATATATATATATATAT	162	TCGCTGTATATATATATATATATATAT	1247	TCGCTGTATATATATATATATATATAT	TCGCTGTATATATATATATATATATAT
1076	RNASEP_CLB_459_487_F	TCGCTGTATATATATATATATATATAT	162	TCGCTGTATATATATATATATATATAT	1426	TCGCTGTATATATATATATATATATAT	TCGCTGTATATATATATATATATATAT
1077	ICD_CXB_93_120_F	TCGCTGTATATATATATATATATATAT	343	TCGCTGTATATATATATATATATATAT	921	TATAGATTTTTCAGAGAGAGAGAGAG	TATAGATTTTTCAGAGAGAGAGAGAG
1078	ICD_CXB_93_120_F	TCGCTGTATATATATATATATATATAT	671	TCGCTGTATATATATATATATATATAT	921	TATAGATTTTTCAGAGAGAGAGAGAG	TATAGATTTTTCAGAGAGAGAGAGAG
1079	ICD_CXB_176_196_F	TCGCTGTATATATATATATATATATAT	369	TCGCTGTATATATATATATATATATAT	916	TATAGATTTTTCAGAGAGAGAGAGAG	TATAGATTTTTCAGAGAGAGAGAGAG
1080	ICB1111A_NC002971_60	TCGCTGTATATATATATATATATATAT	290	TCGCTGTATATATATATATATATATAT	848	TATAGATTTTTCAGAGAGAGAGAGAG	TATAGATTTTTCAGAGAGAGAGAGAG
1081	ICB1111A_NC002971_74	TCGCTGTATATATATATATATATATAT	594	TCGCTGTATATATATATATATATATAT	952	TATAGATTTTTCAGAGAGAGAGAGAG	TATAGATTTTTCAGAGAGAGAGAGAG
1082	RNASEP_RKP_419_448_F	TCGCTGTATATATATATATATATATAT	599	TCGCTGTATATATATATATATATATAT	957	TATAGATTTTTCAGAGAGAGAGAGAG	TATAGATTTTTCAGAGAGAGAGAGAG

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

1083	RNASP_RFP_422_443_F	TACAGGCGACGCGTAACTGGG	159	RNASP_RFP_542_565_R	TGAGCGATCTACCGCATACAA	957
1084	RNASP_RFP_466_491_F	TCCACDAGCGAGATCATATAGCC	310	RNASP_RFP_542_565_R	TGAGCGATCTACCGCATACAA	957
1085	RNASP_RFP_264_287_F	TCTAAATGCTGTGAGTGGTG	391	RNASP_RFP_295_321_R	TCTTAGAGTCGGACTTCTCGTGA	1119
1086	RNASP_RFP_426_448_F	TGCATACCGTAAAGTGCGACAA	497	RNASP_RFP_542_565_R	TGAGCGATCTACCGCATACAA	957
1087	OMP_RFP_860_890_F	TTACAGGAGGTTTGAAGTGAATCTA	654	OMP_RFP_972_996_R	TCTCGAGCTCTACTCTCTCATTA	1051
1088	OMP_RFP_1192_1221_F	TCTACTGATTTTGGTAATCTTCGAC	392	OMP_RFP_1288_1315_R	TAGCAGCAAAAGTATACGACTGCGAT	910
1089	OMP_RFP_3417_3440_F	TGCAGTGGTATCTCAACNCGGG	485	OMP_RFP_3520_3550_R	TGGTGTGATCTCTGTAGTGTGGCAT	1310
1090	GLTA_RFP_1043_1072_F	TGGACTTGAAGCTATCTCTCTTAA	576	GLTA_RFP_1138_1162_R	TGACATTTGCGAGGCTATACCCCT	1147
1091	GLTA_RFP_400_428_F	TCTCTCTCTCTCTCTCTCTCTCTCT	413	GLTA_RFP_499_529_R	TGAGTGGTATCTTAGAGCATCATCTAAT	1305
1092	GLTA_RFP_1023_1055_F	TGAGCTTGAAGCTATGCTCTTAA	330	GLTA_RFP_1129_1156_R	TTGCGAGCGTATACCCCTATATATA	1415
1093	GLTA_RFP_1043_1072_F	TGAGCTTGAAGCTATGCTCTTAA	553	GLTA_RFP_1138_1162_R	TGACATTTGCGAGGCTATACCCCT	1147
1094	GLTA_RFP_1043_1072_F	TGAGCTTGAAGCTATGCTCTTAA	543	GLTA_RFP_1138_1162_R	TGAGCATTTGCGAGGCTATACCCCT	1330
1095	GLTA_RFP_400_428_F	TCTCTCTCTCTCTCTCTCTCTCTCT	413	GLTA_RFP_505_534_R	TGGAGTGTAGTATCTTAGCGATCAT	1230
1096	CTTA_VBC_117_142_F	TCTTATCTCGAGGAGGAGAGTGT	410	CTTA_VBC_194_218_R	TGCTTACACATCTCTGTAGTTC	1226
1097	CTTA_VBC_351_377_F	TGATATAGGCGATGACGTCTCATC	630	CTTA_VBC_441_466_R	TGTCATCAAGCACCCCAATATCACT	1324
1098	RNASP_VBC_331_349_F	TGCGGCGATCTGACTGGT	325	RNASP_VBC_388_414_R	TGACTTCTCTCTCTCTCTCTCTCT	1163
1099	TOX_RFC_135_158_F	TCTTATAGGCGAGCGAGAGCGG	362	TOX_RFC_221_246_R	TGACACCTCTCTCTCTCTCTCTCT	1370
1100	ASD_FRT_1_29_F	TCTTATAGGCGAGCGAGAGCGG	690	ASD_FRT_86_116_R	TGAGATCTCGAAGAGAGGTGGCAAA	1156
1101	ASD_FRT_43_76_F	TGAGATCTGAGTCTCTCTCTCTCT	295	ASD_FRT_129_156_R	TGCTTATGTCGTATACACTGTGGC	1009
1102	GALE_FRT_168_199_F	TTATGCTAGACTTTTAGTAAAG	658	GALE_FRT_241_269_R	TGACTTACAGCTTAAAGCAGCAAT	973
1103	GALE_FRT_834_865_F	TTATGCTAGACTTTTAGTAAAG	245	GALE_FRT_901_925_R	TGACTTGGCAGCTATCGCAAACT	915
1104	GALE_FRT_308_339_F	TGAGATCTGACTTAACTTACTTAA	306	GALE_FRT_390_422_R	TCTCTGTAAGAGGTGGTATATCTTA	1136
1105	TPAH_SGF_258_277_F	TGAGAGCTCTCTCTCTCTCTCT	458	TPAH_SGF_301_327_R	TCTCTCTCTCTCTCTCTCTCTCTCT	1085
1106	TPAH_SGF_113_134_F	TCTTTCACGCTCTTCTCGTATC	350	TPAH_SGF_172_191_R	TTTTTCACGCTCTTCTCGTATC	1441

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

1107	IPAH SQF 462 486 F	TGACGACGCTTCGGGAGGACTT	271	IPAH SQF 522 540 R	TGTCNCTCCGACGACGCA	1322
1111	RNASEP BRM 467 486 F	TGACCCCTTCGGGAGGACGCA	147	RNASEP BRM 542 561 R	TGCTTCGGGAGGACTCCCG	1227
1112	RNASEP BRM 325 347 F	TACCCGAGGAGGAGTCCGACGA	185	RNASEP BRM 402 428 R	TGCTTACCCGAGGCTTTCACGCTAC	1125
1128	HUPB CJ 113 134 F	TAGTGTCTCCACGAGCTGGGCT	230	HUPB CJ 157 188 R	TGCTTATGTGTGAGGATCTGACGACG	1028
1129	HUPB CJ 76 102 F	AT	324	HUPB CJ 157 188 R	TGCGTAATGATGAACTTACGTCGACG	1028
1130	HUPB CJ 76 102 F	AT	324	HUPB CJ 114 135 R	TAGCCGACGCTGTTTGACGACT	913
1151	AB MUST-11- OIF007 62 91 F	TGAGATTTGCTGACCTTTTATCTGAA	454	AB MUST-11- OIF007 169 203 R	TGTACGATTTGAAACCAATGATCATGACA	1418
1152	OIF007 195 214 F	TATGTTTTCATATGTCAGAGTGAG	243	AB MUST-11- OIF007 291 324 R	TGACAGGTTTCTCTTCTTCTATCATATTTTC	969
1153	AB MUST-11- OIF007 230 289 F	TGCGACGTTTACAGGTGCGGCAAAA	541	AB MUST-11- OIF007 364 393 R	TGCGATGCAATATCATATTTGACGACG	1400
1154	AB MUST-11- OIF007 206 239 F	TGAGTGTGCTGATGATCATGATGAC	436	AB MUST-11- OIF007 318 344 R	TGCGCBAABACTCCCTTTTCACGAG	1036
1155	AB MUST-11- OIF007 522 552 F	TGAGTGTGATGAAAGAACTATTGCT	378	AB MUST-11- OIF007 597 610 R	TTCNCTTGAGGAAATGTGCTGG	1392
1156	AB MUST-11- OIF007 547 571 F	CAACC	250	AB MUST-11- OIF007 656 686 R	TAGTCTTACGAAATTTCTTCTACGATAC	902
1157	OIF007 601 627 F	TGACGAGAGGCTTTTGAGAGGAGG	256	AB MUST-11- OIF007 710 736 R	TATGCGGAGTGTGACATCATCTTCTC	881
1158	AB MUST-11- OIF007 1202 1235 F	TGTCGCCGCAATTTGCTNPAAGC	384	AB MUST-11- OIF007 1266 1296 R	TAG	878
1159	AB MUST-11- OIF007 1212 1235 P	TGTCGCCGCAATTTGCTNPAAGC	384	AB MUST-11- OIF007 1299 1315 R	TGCACTGTGCTGTGACGAG	1199
1160	AB MUST-11- OIF007 1214 1264 F	TGTCGCCGCAATTTGCTNPAAGC	694	AB MUST-11- OIF007 1335 1362 R	TGCACTGTGCTGTGACGAG	1215
1161	AB MUST-11- OIF007 1327 1356 F	TAGGTTTACGTCATGCTGCTGAT	225	AB MUST-11- OIF007 1422 1448 R	TGCGAGTTTCCGATTTTCACTGCTGTG	1212
1162	AB MUST-11- OIF007 1345 1359 F	ATGG	393	AB MUST-11- OIF007 1470 1494 R	TGCGATGTGATGTGATGATGCTGG	1083
1163	AB MUST-11- OIF007 1351 1375 F	TGTCGCCGCAATTTGCTNPAAGC	662	AB MUST-11- OIF007 1470 1494 R	TGCGATGTGATGTGATGATGCTGG	1083
1164	AB MUST-11- OIF007 1387 1412 F	TGTCGCCGCAATTTGCTNPAAGC	422	AB MUST-11- OIF007 1470 1494 R	TGCGATGTGATGTGATGATGCTGG	1083
1165	AB MUST-11- OIF007 1542 1569 F	TGTCGCCGCAATTTGCTNPAAGC	194	AB MUST-11- OIF007 1656 1680 R	TGAGTGTGATGTGATGATGCTGG	1173
1166	AB MUST-11- OIF007 1546 1593 F	TGTCGCCGCAATTTGCTNPAAGC	684	AB MUST-11- OIF007 1656 1680 R	TGAGTGTGATGTGATGATGCTGG	1173
1167	AB MUST-11- OIF007 1611 1638 F	TGTCGCCGCAATTTGCTNPAAGC	375	AB MUST-11- OIF007 1731 1757 R	TGAGTGTGATGTGATGATGCTGG	890

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

47_F	OMPU NC002505_931_9	TTACTACTTACGCGAAGCTCG	193	OMPU NC002505_1033_1053_R	TTAGAGTCTTAACCTGAGCC	1368
2018	53_F	TTACTACTTCTTACGCGAAGCTTC	197	OMPU NC002505_1033_1054_R	TGTTAGAGTCTGTAAGCTGAGCC	1307
2019	238_F	TCGCTGATCTTACGCGAAGCTTC	269	TCFA NC002505_148_170_R	TTCTGGAGTCTTAAGCTGCGCTG	1391
2020	238_F	TCGCTGATCTTACGCGAAGCTTC	574	TDH NC004605_357_386_R	TTCTGGAGTCTTACGCGAAGCTTC	1351
2021	TDH NC004605_265_28	TCGCTGATCTTACGCGAAGCTTC	412	VYFA NC004460_852_886_R	TGCTGAGTCTTACGCGAAGCTTC	887
2022	42_F	TCGCTGATCTTACGCGAAGCTTC	508	235_EC_2746_2770_R	TGCTGAGTCTTACGCGAAGCTTC	1297
2023	235_EC_2643_2667_F	TGCTGATCTTACGCGAAGCTTC	202	165_EC_789_811_R	TGCTGAGTCTTACGCGAAGCTTC	1240
2024	165_EC_713_732_TM0D	TGCTGATCTTACGCGAAGCTTC	560	165_EC_880_897_TM0D_R	TGCTGAGTCTTACGCGAAGCTTC	1278
2025	165_EC_784_806_F	TGCTGATCTTACGCGAAGCTTC	634	165_EC_1052_1074_R	TGCTGAGTCTTACGCGAAGCTTC	896
2026	165_EC_959_981_F	TGCTGATCTTACGCGAAGCTTC	489	TUFB_EC_1034_1058_2_R	TGCTGAGTCTTACGCGAAGCTTC	1204
2027	TUFB_EC_365_379_F	TGCTGATCTTACGCGAAGCTTC	284	RPOC_EC_2227_2249_R	TGCTGAGTCTTACGCGAAGCTTC	1244
2028	MOD_F_2146_2174_F	TGCTGATCTTACGCGAAGCTTC	617	RPOB_EC_1909_1929_TM0D_R	TGCTGAGTCTTACGCGAAGCTTC	1250
2029	1841_1866_F	TGCTGATCTTACGCGAAGCTTC	385	RPFB_EC_739_763_R	TGCTGAGTCTTACGCGAAGCTTC	1208
2030	RPFB_EC_550_579_TM0	TGCTGATCTTACGCGAAGCTTC	449	RPFB_EC_737_760_R	TGCTGAGTCTTACGCGAAGCTTC	1295
2031	D_F	TGCTGATCTTACGCGAAGCTTC	309	INFB_EC_1439_1469_R	TGCTGAGTCTTACGCGAAGCTTC	1335
2032	INFB_EC_590_710_F	TGCTGATCTTACGCGAAGCTTC	397	VALS_EC_1195_1219_R	TGCTGAGTCTTACGCGAAGCTTC	1292
2033	VALS_EC_1105_1124_F	TGCTGATCTTACGCGAAGCTTC	482	SSFE_BA_197_222_TM0D_R	TGCTGAGTCTTACGCGAAGCTTC	1402
2034	SSFE_BA_113_137_F	TGCTGATCTTACGCGAAGCTTC	405	RPOC_EC_2313_2338_R	TGCTGAGTCTTACGCGAAGCTTC	1273
2035	MOD_F_2218_2241_F	TGCTGATCTTACGCGAAGCTTC	698	MRCT-R_NC003923-41798-41609_86_113_R	TGCTGAGTCTTACGCGAAGCTTC	1420
2056	MRCT-R_NC003923-41798-41609_33_60_F	TGCTGATCTTACGCGAAGCTTC	263	AGR-III_NC003923-2108074-2109507_36_73_R	TGCTGAGTCTTACGCGAAGCTTC	730
2057	AGR-III_NC003923-2108074-	TGCTGATCTTACGCGAAGCTTC	457	AGR-III_NC003923-2108074-2109507_622_653_R	TGCTGAGTCTTACGCGAAGCTTC	906
2058	AGR-III_NC003923-2108074-	TGCTGATCTTACGCGAAGCTTC	701	AGR-III_NC003923-2108074-2109507_1070_1098_R	TGCTGAGTCTTACGCGAAGCTTC	1319
2059	AGR-III_NC003923-2108074-	TGCTGATCTTACGCGAAGCTTC	610	AGR-I_AJ617706_694_726_R	TGCTGAGTCTTACGCGAAGCTTC	1021
2060	AGR-I_AJ617706_622_651_R	TGCTGATCTTACGCGAAGCTTC				

DOCKET NO.: DIBIS-0083USI (COUNSEL DOCKET NO: 10593)

F	F	AGR- I_AJ617706_580_611_	TGGATATTTAAACAACTGTGTGACAA TGGCAG	579	AGR-I_AJ617706_626_655_R	TGGTACTTACCTTCATCCATATGAG TC	1302
2061	AGR- II_AJ617711_250_283	AGR-II_NC002745- 2079448- 2080879_620_651_F	TCTTCGCGAGTTTATTTGATGACCC TAAAGT	415	AGR-II_NC002745-2079448- 2080879_700_731_R	TGTTTATATGTTTCCATATGATCAGAC TTTC	1424
2062	AGR- II_AJ617711_250_283	AGR-II_NC002745- 2079448- 2080879_649_679_F	TGTATCCGCTGAATTAACGATTTAT ACGAC	624	AGR-II_NC002745-2079448- 2080879_715_745_R	TTCGACGACGTAGTTGTTATATGTTTC CAT	1077
2063	AGR- IV_AJ617711_931_961	TGTTATCTATTTTGTCTGATATGAC CTCCG		566	AGR- IV_AJ617711_1004_1035_R	TGCTCATATGACGATTTTACGACATAT GTGA	1233
2064	AGR- IV_AJ617711_250_283	TGGCATCTCTGCTTTTAAATTTAGTA AACATCA		562	AGR-IV_AJ617711_309_335_R	TCCATACCTATGCGATTAAGTTCAT TCC	1017
2065	BLAZ_NC002952(19138 27...1914672)_68_68_	TCCACTTATGCGAATGCGAATTA GCAA		312	BLAZ_NC002952(1913827...19 14672)_68_68_R	TGGCCACTTTTATGACGACCTTACAT C	1277
2066	BLAZ_NC002952(19138 27...1914672)_68_68_	TGCATCTATGCGAATGCGAATTA GCAA		494	BLAZ_NC002952(1913827...19 14672)_68_68_2_R	TAGCTTTTGGACGACGCTTTTAATTA AGAT	926
2067	BLAZ_NC002952(19138 27...1914672)_68_68_	TGGATCTCTGCGCTGCTGCTTC TGGATCTCTGCGCTGCTGCTTC		467	BLAZ_NC002952(1913827...19 14672)_68_68_3_R	TGGACGACGCTCTTAAATTAAGTATCT CC	1263
2068	BLAZ_NC002952(19138 27...1914672)_68_68_	TGATCTCTGCGCTGCTGCTTC TGGATCTCTGCGCTGCTGCTTC		232	BLAZ_NC002952(1913827...19 14672)_68_68_4_R	TCTTTCTTTGCTTAATTTCCATTTGCG GAT	1145
2069	BLAZ_NC002952(19138 27...1914672)_68_68_	TGATCTCTGCGCTGCTGCTTC TGGATCTCTGCGCTGCTGCTTC		487	BLAZ_NC002952(1913827...19 14672)_34_67_R	TATCTCTGCGCTGCTTATGATCTAA AGCATTA	1366
2070	BLAZ_NC002952(19138 27...1914672)_68_68_	TGATCTCTGCGCTGCTGCTTC TGGATCTCTGCGCTGCTGCTTC		351	BLAZ_NC002952(1913827...19 14672)_40_68_R	TGGGACCTCTTACGACCTTTTATGATC TAA	1289
2071	BSA-A_NC003923- 1304065-	TGCGGATGTGCGCTTACTTCAAT T		214	BSA-A_NC003923-1304065- 1303589_253_278_R	TGACGGGACCTTAGATTTACAAACCC T	1197
2072	BSA-A_NC003923- 1304065-	TGCGGATGTGCGCTTACTTCAAT T		32	BSA-A_NC003923-1304065- 1303589_253_278_R	TGATAGGAAGGATTAACCATAGTT T	1203
2073	BSA-A_NC003923- 1304065-	TGCGGATGTGCGCTTACTTCAAT T		679	BSA-A_NC003923-1304065- 1303589_253_278_R	TAAACAGCTTACTTTCGATCCATTA TAA	856
2074	BSA-A_NC003923- 1304065-	TGCGGATGTGCGCTTACTTCAAT T		519	BSA-A_NC003923-1304065- 1303589_253_278_R	TGTGTGCGCGCTGATTAATTTAATA TAA	1353
2075	BSA-A_NC003923- 1304065-	TGCGGATGTGCGCTTACTTCAAT T		209	BSA-B_NC003923-1917149- 1303589_253_278_R	TGTGTGCGCGCTGATTAATTTAATA TAA	1331

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

[illegible]

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2094	EPWB_Y13600-625- 1362_551 645 R	TRAGCTCCGAGGAGTCTTTC	161	EPWB_Y13600-625- 1362_551 645 R	TGATGTGCTATCGCGGTGAGTT	989
2095	EPWB_NC003923- 1529595- 1531285 688 713 F	TGAGCTGCATCAGCTGATGAGAG	456	PVLUK_NC003923-1529595- 1531285 775 804 R	TGGAAATCTCATGAAATTAAGAGTGAAG GA	1261
2096	PVLUK_NC003923- 1529595- 1531285 1039 1068 F	TGGAAATTAATCTCTCGGATTTT GACT	539	PVLUK_NC003923-1529595- 1531285 1095 1125 R	TGATAGTAAATATCTCGACATGATC CAA	993
2097	PVLUK_NC003923- 1529595- 1531285 908 936 F	TGAGTACATCATATTTCTGCGATA CGT	461	PVLUK_NC003923-1529595- 1531285 950 978 R	TTCATGAAAAAGCTCTGAGAGTACAA G	1124
2098	PVLUK_NC003923- 1529595- 1531285 610 633 F	TCGGAATCTGATTTGTCAGTTGTT	373	PVLUK_NC003923-1529595- 1531285 654 682 R	TCCACACTGTGATGAGAGAAAGTTTCA T	968
2099	SA442_NC003923- 2538576- 2538831 11 35 P	TGTCGTGACAGATATTTCTGACAA	635	SA442_NC003923-2538576- 2538831 98 124 R	TTTCGATCGACCTGATGAGATTTCA	1433
2100	SA442_NC003923- 2538576- 2538831 103 126 F	TGAAATCTCATCTGTCATCGGAA A	427	SA442_NC003923-2538576- 2538831 163 188 R	TGATATGACCACCTCTGCTACTACTA T	1098
2101	SA442_NC003923- 2538576- 2538831 166 188 P	TGTCATTAATCTGCTCTGCGAACA	395	SA442_NC003923-2538576- 2538831 231 257 R	TTTATGACACGCTCTGCTACTACTAA T	1428
2102	SEA_NC003923- 2052219- 2051456 115 135 P	TGATACCAATGCTGCTATACAA	226	SEA_NC003923-2052219- 2051456 173 200 R	TGATATGAGAGGAAACTTTTTCAGC	1179
2103	SEA_NC003923- 2052219- 2051456 572 598 P	TCGAGGACACGCTTTAGCCA	495	SEA_NC003923-2052219- 2051456 621 651 R	TGATCTGTAATCTCTTTATTTGATTT ATG	1070
2104	SEA_NC003923- 2052219- 2051456 382 411 F	TGATCTGATCTTTTTCATGCGAAG T	156	SEA_NC003923-2052219- 2051456 461 492 R	TGATATTAACCAAGCTCTCTGAGACT ATG	1315
2105	SEA_NC003923- 2052219- 2051456 377 406 F	TGATATGCTGCTACATCTACATAT ATATATC	629	SEA_NC003923-2052219- 2051456 455 492 R	TGACCTGTTCCAAAGCTCTGATTTTG T	861
2106	SEA_NC002758- 2135540- 2135140 208 237 F	TGATATGATGCTGCTGCTACATAC ATGAA	695	SEA_NC003923-2052219- 2051456 273 298 R	TGACCGTTTCGAAAGCTCTGATTTTG TTTACC	862
2107	SEA_NC002758- 2135540- 2135140 206 235 F	TTTTCATGATATATTTGATATTCG CACT	702	SEA_NC002758-2135540- 2135140 402 402 R	TGATCTGTTAGGACTCTGTTGACT T	988
2108	SEA_NC002758- 2135540- 2135140 206 235 F	TGATCTGCTGCTGATGAGAGGAT ATA	244	SEA_NC002758-2135540- 2135140 283 304 R	TGCACTCATCTGTTTAGAGCTCT T	1194
2109	SEA_NC002758- 2135540- 2135140 206 235 F	TGATCTGCTGCTGATGAGAGGAT ATA	151	SEA_NC002758-2135540- 2135140 283 304 R	TGTCAGGACATCTCTCTCTACAA T	1334

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2110	SEC_NCO02758- 2135540-402_402_2_F A	TGTATATATAGTGATGATACGAC	596	SEC_NCO02758-2135540- 2135140_402_402_2_R	TTACATCTTCMAATACCCAGCTAA	1361
2111	SEC_NCO03923- 852768_546_575_F ATGG	TTACATCTTCAGAGAAACACCTTTGATA	648	SEC_NCO03923-851678- 852768_620_647_R	TGAGTTCACCTTCMAAGAAATGTGT	1177
2112	SEC_NCO03923- 851678- 852768_537_566_F ACTT	TGGATATCAAAATGATGAAGAAAC	546	SEC_NCO03923-851678- 852768_619_647_R	TCAGTTCACCTTCMAAGAAATGTGT	985
2113	SEC_NCO03923- 851678- 852768_720_749_F CAGG	TGAGTTCACCTTCACCAATAGAA	466	SEC_NCO03923-851678- 852768_794_815_R	TGCGCTGTCGCGCATAT	1078
2114	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2115	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2116	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2117	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2118	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2119	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2120	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2121	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2122	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2123	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2124	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133
2125	SEC_NCO03923- 851678- 852768_757_810_F P_NCO_M28521_687_682	TGATATATATATATGCTTCACCA	504	SEC_NCO03923-851678- 852768_853_886_R	TTCTACACTTTTATAGATACACGTTTT	1133

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2126	SEQ_NC002758- 1955100- 1954171.633.651.F SEQ_NC002758- 1955100- 1954171.633.651.F	TGACACATGACGACATCTTGATTT ACA	548	SEQ_NC002758-1955100- 1954171.671.702.R	TGCTTTGTAATCTAGTTCTGTGATGTA ACCA	1260
2127	SEQ_NC002758- 1955100- 1954171.640.564.F	TGAGGTTTGTTGATGATGTTGTT	555	SEQ_NC002758-1955100- 1954171.607.635.R	TGCTATGTTGCAATGTTTACTGTACAG T	1329
2128	SEQ_NC002758- 1955100- 1954171.654.718.F	TACAGACAGACATGCTGCTACTA	173	SEQ_NC002758-1955100- 1954171.735.762.R	TGATCTAATGAGAACATCAATCAATGCG GA	1187
2129	SEQ_NC002758- 1955100- 1954171.649.472.F	TTCACATCTCTGATTTTACTCTGAGA	662	SEQ_NC002758-1955100- 1954171.735.762.R	TAGTGTGTGACTCAATGACATCTCA GA	827
2130	SEQ_NC002758- 1955100- 1954171.649.472.F	TAGACATCAAGGTGATGATGGACATG	201	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGTCTAATCTAGACATGCTGA T	1390
2131	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC ACTA	400	SEQ_NC002758-1955100- 1954171.607.635.R	TACCATCTACCCCAACATTAGACACAA	888
2132	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CACT	677	SEQ_NC002758-1955100- 1954171.607.635.R	TAGCAACATCACTCTTCTCTGT	909
2133	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	253	SEQ_NC002758-1955100- 1954171.607.635.R	TCAACAGACATCAATCAATCAATGACAA	966
2134	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	666	SEQ_NC002758-1955100- 1954171.607.635.R	TGTACAGACATCAATCAATCAATGACAA	1316
2135	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	471	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGCTCTCCCAACATGATGATGAG	1129
2136	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	394	SEQ_NC002758-1955100- 1954171.607.635.R	TGCTGAGGCTCTCCCAACATGATGATGAG	1293
2137	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	637	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGCTCTCCCAACATGATGATGAG	1119
2138	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	211	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGCTCTCCCAACATGATGATGAG	1049
2139	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	153	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGCTCTCCCAACATGATGATGAG	925
2140	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	301	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGCTCTCCCAACATGATGATGAG	984
2141	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	619	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGCTCTCCCAACATGATGATGAG	1312
2142	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	514	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGCTCTCCCAACATGATGATGAG	1221
2143	SEQ_NC002758- 1955100- 1954171.649.472.F	TCTGATCTCTGATGATGGATGCTACAC CA	304	SEQ_NC002758-1955100- 1954171.607.635.R	TCTGAGGCTCTCCCAACATGATGATGAG	907

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2144	TSST_NCO02758- 2138255 397 325 P 2138255 347 373 R	TCCTTTTCCAAAGGGGAAAAGTTGA CTT	423	TSST_NCO02758-2137564- 2138255 347 373 R	TAACTCTTCTCTAGTAACTTGGCTT	974
2145	ARCC_NCO03923- 2725050- 2724595 37 58 F 2724595 97 128 R	TGCGCGCAATCCCACTGGTGA	368	2724595 37 58 F 2724595 97 128 R	TCGATTAATTCACATTCATTTCTGTTT	1175
2146	ARCC_NCO03923- 2725050- 2724595 131 161 F 2724595 218 249 F	TGAATATGCTAAGAACTGTAGGCA ATCGT	437	ARCC_NCO03923-2725050- 2724595 218 249 F 2724595 214 245 R	TCCTCTCTTCTGTAATAAAGGCCAA TTGG	1137
2147	ARCC_NCO03923- 2725050- 2724595 154 176 F 2724595 271 393 P	TTGGTCTCTTTTATACGAAGAGAA GTTGAA	691	ARCC_NCO03923-2725050- 2724595 271 393 P 2724595 322 353 R	TGCTGTTCTAGTATAGATTCAGGTAGT GTGA	1306
2148	ARCC_NCO03923- 1674726- 1674377 30 62 F 1674726 154 176 F	TTGCGATAGAGAGATGCTCGT	686	ARCC_NCO03923-1674726- 1674377 30 62 F 1674726 154 176 F	TGAAATGCTGTAATTAATCTTTCGCAAT CT	1064
2149	ARCC_NCO03923- 1674726- 1674377 204 232 F 1674726 154 176 F	TGGGCTTTTAAATATCCCAATGAAG ATTTCGA	590	ARCC_NCO03923-1674726- 1674377 204 232 F 1674726 154 176 F	TAACTGCTGTAATTAATCTTTCGCAAT CT	891
2150	GLFP_NCO03923- 1296927- 1297391 270 301 F 1296927 81 100 R	TGATGGCGATGGATAGGGTAAATA CAG	474	ARCC_NCO03923-1674726- 1296927 81 100 R 1297391 323 359 R	TAACTGCTGTAATTAATCTTTCGCAAT CT	869
2151	GLFP_NCO03923- 1296927- 1297391 27 51 F 1296927 81 100 R	TGCAACCGCTATTAAGCAATATCTTGG CGAAT	491	GLFP_NCO03923-1296927- 1297391 27 51 F 1296927 81 100 R	TGCAACCGCTATTAAGCAATATCTTGG CGAAT	1193
2152	GLFP_NCO03923- 1296927- 1297391 239 260 F 1296927 81 100 R	TGATGGGCTATTAAGCAATATCTTGG CGAAT	558	GLFP_NCO03923-1296927- 1297391 239 260 F 1296927 81 100 R	TAACTGCTGTAATTAATCTTTCGCAAT CT	850
2153	OMK_NCO03923- 1191334 91 122 F 1191334 301 329 P	TAGCTGGCGCAATTAAGGTGT	218	OMK_NCO03923-1190906- 1191334 91 122 F 1191334 301 329 P	TAACTGCTGTAATTAATCTTTCGCAAT CT	972
2154	GLFP_NCO03923- 1190906- 1191334 240 267 F 1190906 81 100 R	TAGCTGGCGCAATTAAGGTGT	200	OMK_NCO03923-1190906- 1191334 240 267 F 1190906 81 100 R	TAACTGCTGTAATTAATCTTTCGCAAT CT	1180
2155	GLFP_NCO03923- 1190906- 1191334 301 329 P 1190906 81 100 R	TAGCTGGCGCAATTAAGGTGT	435	OMK_NCO03923-1190906- 1191334 301 329 P 1190906 81 100 R	TAACTGCTGTAATTAATCTTTCGCAAT CT	1082
2156	PTA_NCO03923- 623885- 623935 237 263 F 623885 314 345 R	TAGCTGGCGCAATTAAGGTGT	268	PTA_NCO03923-623885- 623935 237 263 F 623885 314 345 R	TAACTGCTGTAATTAATCTTTCGCAAT CT	1284
2157	PTA_NCO03923- 623885- 623935 237 263 F 623885 314 345 R	TAGCTGGCGCAATTAAGGTGT	418	PTA_NCO03923-623885- 623935 237 263 F 623885 314 345 R	TAACTGCTGTAATTAATCTTTCGCAAT CT	1301
2158	PTA_NCO03923- 623885- 623935 237 263 F 623885 314 345 R	TAGCTGGCGCAATTAAGGTGT	439	PTA_NCO03923-623885- 623935 237 263 F 623885 314 345 R	TAACTGCTGTAATTAATCTTTCGCAAT CT	1207

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2159	PTA_NCO03923- 629555 328 356 F	TCCBACCGGGGTATCAGACATC AGG	303	PTA_NCO03923-629555 393 422 R	TGTCTGGATTGATGACACATCAGCA AG	1349
2160	TFI_NCO03923- 830671- 831072 131 160 F	TGAGATTAGAAAGCTGTTCAGGT TTAT	486	TFI_NCO03923-830671- 831072 209 239 R	TGAGATTGATGATTTACCACTTCCGA TTG	1165
2161	TFI_NCO03923- 830671- 831072 1 34 F	TCCACGACAGATGAGAAATTA CAAAAG	318	TFI_NCO03923-830671- 831072 97 129 R	TGGTCAACATCTGACCTTACCACTT TACG	1300
2162	TFI_NCO03923- 830671- 831072 199 227 F	TCCAACTGGGCAATCGGAATGTTAA ATC	246	TFI_NCO03923-830671- 831072 253 286 R	TGGTCAACATTTGACCTGACAAATGC ACACAT	1275
2163	VOI_NCO03923- 379431 142 167 F	TGATTTCTGCTGATGAAAGTGGCTT T	440	VOI_NCO03923-378916- 379431 259 284 R	TCCCACTTAGCAGATGATCAATTC T	1076
2164	VOI_NCO03923- 379431 44 77 F	TACACATATATTAAAGAACGGGT TTGATTC	175	VOI_NCO03923-378916- 379431 120 145 R	TTCTGCTGGATTTGCTGTTGCTT T	1388
2165	VOI_NCO03923- 378916- 379431 135 160 F	TCCACGACGATTTGCTGCTGTAAG T	314	VOI_NCO03923-378916- 379431 193 221 R	TCCAAACCGACACGACATCTATTCA C	987
2166	VOI_NCO03923- 378916- 379431 235 300 F	TACTGCGGCTGAGGAAATATGCT T	219	VOI_NCO03923-378916- 379431 364 396 R	TCCACGCTTAAACCTCATATACCAAG CTATC	1013
2167	BLAZ_NCO03923- 1914672 546 575 F	TCCACATCTCGAAATGGAATTA GCA	312	BLAZ_NCO03923-1914672- 55 653 R	TGGCCACTTTATCAGACACCTTAGCT C	1277
2168	BLAZ_NCO03923- 1914672 546 575 2 F	TCCACATCTCGAAATGGAATTA GCA	494	BLAZ_NCO03923-1914672- 28 659 R	TAGTCTTTTGGACACGCTTTTATTA AACT	926
2169	BLAZ_NCO03923- 1914672 507 531 F	TGATCTCTGACGCTCTCTCTTTC T	467	BLAZ_NCO03923-1914672- 22 651 R	TGGACACGCTCTTTTAAATTAAGTACT CC	1263
2170	BLAZ_NCO03923- 1914672 508 531 F	TATCTCTCAAGCCTCTGCTTTC T	232	BLAZ_NCO03923-1914672- 53 593 R	TCTTTCTTTGCTTAAATTTGCAATGAC GAT	1145
2171	BLAZ_NCO03923- 1914672 546 575 F	TCCATGCTCTTAGTTTATGATGCT GTGATTC	487	BLAZ_NCO03923-1914672- 21 651 R	GAATCTCTTACCACTTTTATGATCAAA AAAT	1366
2172	BLAZ_NCO03923- 1914672 546 575 2 F	TCTCTGCTTTTATGTTTAAAGCAATG AATCAAA	351	BLAZ_NCO03923-1914672- 27 157 R	TGGGCACTTCTTACCACTTTTATGATC TAA	1289
2173	BLAZ_NCO02952- 1914672 546 575 F	TCCATCTTCCGAATGGAATTA GCA	312	BLAZ_NCO02952-1913827- 1914672 655 683 R	TGCCCACTTTATGCAACCACTTCAAGT C	1277
2174	BLAZ_NCO02952- 1913827- 1914672 546 575 2 F	TCCATCTTCCGAATGGAATTA GCA	494	BLAZ_NCO02952-1913827- 1914672 628 659 R	TAGCTTTTGGACACGCTCTTTATTA AACT	926
2175	BLAZ_NCO02952- 1913827- 1914672 507 531 F	TGATCTCTGACGCTCTGCTTTC T	467	BLAZ_NCO02952-1913827- 1914672 622 651 R	TGGACACGCTCTTTTAAATTAAGTACT CC	1263

DOCKET NO.: DIBIS-0083USI (COUNSEL DOCKET NO: 10593)

2176	BLAZ NC002952-1913827	232	BLAZ NC002952-1913827-1914672 553 583 R	TCTTTCTTCTTCTTAACTTTCTTCAATTCG GAT	1145
2177	BLAZ NC002952-1913827 24 56 F	487	BLAZ NC002952-1913827-1914672 121 154 R	TACTCTCTTACCACTTTTATATATCA AAGATA	1366
2178	BLAZ NC002952-1913827 26 58 F	351	BLAZ NC002952-1913827-1914672 127 157 R	TGSGACCTCTTACCACTTTTATATATC TRA	1289
2247	TUPE NC002758-615038 693 721 F	643	TUPE NC002758-615038-616222 793 820 R	TCTCACCGCTTTCAGCTTATCTAATAA	1321
2248	TUPE NC002758-615038 693 721 F	386	TUPE NC002758-615038-616222 793 820 R	TCTCACCGCTTTCAGCTTATCTAATAA	1321
2249	TUPE NC002758-615038 693 721 F	430	TUPE NC002758-615038-616222 793 820 R	TCTCACCGCTTTCAGCTTATCTAATAA	1321
2250	TUPE NC002758-615038 693 721 F	320	TUPE NC002758-615038-616222 601 630 R	TGCTTTGTTCAGATCACTCTTCTGAGATT GG	1311
2251	TUPE NC002758-615038 693 721 F	433	TUPE NC002758-615038-616222 1030 1060 R	TAGGCTAACCACTTTCAGCTTCTTCTGG TRA	922
2252	TUPE NC002758-615038 693 721 F	307	TUPE NC002758-615038-616222 424 459 R	TTCCTATTCAGCTTATCTTATATATCTT TCAATCTC	1382
2253	TUPE NC002758-615038 693 721 F	342	TUPE NC002758-615038-616222 424 459 R	TAGCTTACCACTTCCCACTTATCTA TCA	899
2254	TUPE NC002758-615038 693 721 F	349	TUPE NC002758-615038-616222 424 459 R	TGCTTTGTTCAGCTTCTTCTGAGCTA A	1354
2255	TUPE NC002758-615038 693 721 F	273	TUPE NC002758-615038-616222 424 459 R	TAGCTTACCACTTCCCACTTATCTA TCA	928
2256	TUPE NC002758-615038 693 721 F	174	TUPE NC002758-615038-616222 424 459 R	TAGCTTACCACTTCCCACTTATCTA TCA	853
2270	TUPE NC002758-615038 693 721 F	566	TUPE NC002758-615038-616222 424 459 R	TAGCTTACCACTTCCCACTTATCTA TCA	979
2271	TUPE NC002758-615038 693 721 F	294	TUPE NC002758-615038-616222 424 459 R	TAGCTTACCACTTCCCACTTATCTA TCA	1107
2272	TUPE NC002758-615038 693 721 F	294	TUPE NC002758-615038-616222 424 459 R	TAGCTTACCACTTCCCACTTATCTA TCA	1102

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

[illegible]

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO.: 10593)

2327	3772_1555_1581_P INV 122457-74- 3772_1558_1585_P	A TGGTAACTAGAGGCTTTAGAGGCGATA TG	598	3772_1619_1647_R INV 122457-74- 3772_1622_1652_R	A TACATAGAGGTTGCGTTCGACGATATCTTT TAC	987
2328	ASD NC006570- 439714- 439608_3_37_F	TGAGCGTTTATCTTAAAGTGGTT TTATATGGTT	459	ASD NC006570-439714- 439608_51_84_R	TGATTCGATACATACAGACATATAAAGCT GAG	1188
2329	ASD NC006570- 439714- 439608_18_45_F	TAAAGTGGTTTATATGGTTGCGCG GA	149	ASD NC006570-439714- 439608_66_95_R	TCGAAATCTTTGATCTCATCTACAG AC	948
2330	ASD NC006570- 439714- 439608_17_45_F	TTAAAGTGGTTTATATGGTTGCGCG GGA	647	ASD NC006570-439714- 439608_67_95_R	TCGCAATCTTTGATCTCATCTACAG TCC	1016
2331	ASD NC006570- 439714- 439608_9_40_F	TTTATAGCTTAAAGTGGTTTATATG GTGCG	709	ASD NC006570-439714- 439608_107_114_R	TCGCTGCGAGATGTCGAAAGACGTTG ATG	1128
2332	GALE AF513299_171_2 00_F	TGAGCTAGACCTTTAGAGTAAGCTA ACCT	280	GALE AF513299_241_271_R	TACGATGTAAATTCGCAAGACTTTG GCTTAG	900
2333	GALE AF513299_168_1 99_F	TTATCACTAGACCTTTAGGTAAG CTAAGC	658	GALE AF513299_245_271_R	TCGCACTACAGCTTTAAGCCGCGA ACAG	1121
2334	GALE AF513299_168_1 99_F	TGCACTAGACCTTTAGGTAAG CTAAGC	658	GALE AF513299_233_264_R	TACGCTTTAAGAGCGCGAATGATTT ACAG	883
2335	GALE AF513299_169_1 99_F	TGCACTAGACCTTTAGGTAAG CTAAGC	319	GALE AF513299_252_279_R	TCGCACTCTCAGCTACAGCTTTAAG GCTTAG	1374
2336	PIA AF053945_7377_7 403_F	TTGAGAGAGATCCCGCTACGTTAT TGGGTA	680	PIA AF053945_7434_7468_R	TACGATGTAAATTCGCAAGACTTTG GCTTAG	900
2337	PIA AF053945_7377_7 403_F	TGCACTCCGCTCAAGTTATATAGTT AC	443	PIA AF053945_7428_7455_R	TCGCAAGAGCTTTGCTTTAGGTGGA TGT	1035
2338	PIA AF053945_7377_7 404_F	TCGTTATCGCATGTCATATTTGG AACT	444	PIA AF053945_7430_7460_R	TGATTCGCGAGAGATTTGCGATTAGG AT	884
2339	CAF AF053947_33412_1 33458_F	TCGTTATCGCATGTCATATTTGG AACT	329	CAF AF053947_33412_1 33458_F	TGAGGTGATGCGGCTGCTTTGACCA TGT	866
2340	CAF AF053947_33412_1 33458_F	TGCTTTATCGCATGTCATATTTGG AACT	499	CAF AF053947_33483_33507_R	TGCTTCACAGAGATTTGCGTSCA TGT	1308
2341	CAF AF053947_33407_1 33429_F	TCAGTTCGTTTATCGCATGCA TGT	291	CAF AF053947_33483_33504_R	TCGCAAGAGATTTGCGTTCGA TGT	1373
2342	CAF AF053947_33407_1 33429_F	TCAGTTCGTTTATCGCATGCA TGT	293	CAF AF053947_33494_33517_R	TGCTTCGCGATGTTTACAGAGAG AT	1184
2344	GALE NC_002505_1_28 F_1	TCAGTTCGTTTATCGCATGCA TGT	260	GALE NC_002505_29_58_R_1	TCCTTATGCAAGTTGCTATCAAGGA AT	1060
2472	OMPA NC000117_68_89 F	TGCTGTAAGGATCTCTCTGA TGT	507	OMPA NC000117_145_167_R	TCGCAAGAGATTTGCGAGGATC C	967
2473	OMPA NC000117_798_8 F	TGCTTACGATGATGCGACGACG TGT	475	OMPA NC000117_865_893_R	TCGCAATCTCTTACAGACCTTTAATC C	947
2474	OMPA NC000117_645_6 F	TGCTTACGATGATGCGACGACG TGT	521	OMPA NC000117_757_777_R	TGCTTCAGAGATTTGCTTCGCG TGT	1328

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2475	OWPA_N0000117_947_9 73 F	TACGTGATGAGACACCTCTTACTCA	157	OWPA_N0000117_1011_040_R AG	TGACGAGCAGCACTCTGAGAGTCTG AG	1153
2476	OWPA_N0000117_774_7 2476	TACTGAGACAGAGCTGAGACC	196	OWPA_N0000117_871_894_R	TGAGAGAGTTGCTGAGAGCACTG	1371
2477	OWPA_N0000117_457_4 83 F	TCTCTCTGCTTGCTGTATTCGAGT T	676	OWPA_N0000117_511_534_R	TAAAGAGAGCTTGCTGTAGTCTCATTCG CA	851
2478	OWPA_N0000117_687_7 10 F	TAGCCGAGCAGCAATTTGATGATCA	212	OWPA_N0000117_787_815_R	TGACCATCTCACTGGTATTTAAGTGTAGCA	1406
2479	OWPA_N0000117_540_5 56 F	TGCGGTAGTAGACATTTATTCAGACA C	571	OWPA_N0000117_649_672_R	TCTCTGAGACGAGGTTTGATATG	1395
2480	OWPA_N0000117_338_3 50 F	TGCACTATCCGAGATGTTGATCA	492	OWPA_N0000117_417_444_R	TCTTTAAATATACGCTTAGTACTCTT	1058
2481	OWPA_N0000117_19_40 F	TATGACCAACTCATCATGAGCAG	234	OWPA_N0000117_71_91_R	TCCGCTGCGGAAATATACCTG	1025
2482	OWPA_N0000117_354_3 82 F	TGCTGAGGTGAGGATCTCTTATCTCA TTG	516	OWPA_N0000117_445_471_R	TGATCTCATCTGCTTACGAACTCAGCTTC	1270
2483	OWPA_N0000117_1297_1 1319 F	TGGAAGAGGTGTGAGCTACTCA	537	OWPA_N0000117_1396_1419_R	TACGTTTGCTCTCTCTGAGAAC	903
2484	OWPA_N0000117_1465 1453 F	TCTGCTCAGCAAGAGAGAGATAC AGG	407	OWPA_N0000117_1541_1569_R	TGCTTTCAAGTTGACAGAAACTCTACA G	1052
2485	OWPA_N0000117_44_66 F	TGAGCATCTTCGCGTGAATGT	450	OWPA_N0000117_120_148_R	TGTGAGCTAGCTATTAAGTTTGTAGA G	1323
2486	OWPA_N0000117_166_1 90 F	TGACAGCAGAGAGAGTTAGACTGTC C	441	OWPA_N0000117_240_261_R	TGACATGCTTCCTCTTCTCAAG	1386
2487	OWPA_N0000117_514_5 36 F	TGAGCATCTTCGCTTGAGATGAGC	287	OWPA_N0000117_640_660_R	TGCTGTGAGGAATCTAGGAGCC	1251
2488	OWPA_N0000117_801_8 27 F	TGTGATATATCTAGATGATGATGAGC A	636	OWPA_N0000117_871_893_R	TGCTGAGACTCTCTGAGAGATC	1419
2489	OWPA_N0000117_219_2 42 F	TGTATGAGTATATATCACTCTCA	632	OWPA_N0000117_319_345_R	TGATCTGATGAGCAAGTGTACTCTG	1010
2490	OWPA_N0000117_964_9 83 F	TGAGCATCTTCGAGTCTCACTCA	176	OWPA_N0000117_1024_1041_R	TGCGAGCTGCTGCTGAGC	1073
2491	OWPA_N0000117_1505 1520 F	TGCGCCGCGAGAGCT	366	OWPA_N0000117_1546_1562_R	TGCTGCTCTTGGGATTA	1416
2492	OWPA_N0000117_216_2 39 F	TGAGCATCTTCGAGTCTCACTCTG	279	OWPA_N0000117_124_143_R	TGCGATGAGCACTGAGCTGAG	1279
2493	OWPA_N0000117_219_2 42 F	TGAGTCTGAGTATGATGATGAGC	452	OWPA_N0000117_313_333_R	TGCGAGTGTGCTGTGCTGTC	1032
2494	OWPA_N0000117_115_1 41 F	TGAGTCTGAGTATGATGATGAGC C	625	OWPA_N0000117_308_330_R	TGAGTCTGCTGCGGCTCAACA	873
2495	OWPA_N0000117_517_5 39 F	TGAGCATCTTCGAGTCTCACTCTG	453	OWPA_N0000117_220_242_R	TGCGGATGATCTTACGAGTAC	1296
2496	OWPA_N0000117_273_2 380	TGAGTCTGAGTATGATGATGAGC	287	OWPA_N0000117_643_663_R	TGCTGTGAGGGAATCTAGGAGCC	1251
2497	OWPA_N0000117_273_2 380	TGAGTCTGAGTATGATGATGAGC	380	OWPA_N0000117_338_360_R	TGCGGATGATCTTACGAGTAC	1296

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

[illegible]

DOCKET NO.: DBIS-0083US1 (COUNSEL DOCKET NO: 10593)

112166-	112647_80_113_F	GATTGAG				112647_346_171_R		
2566	URCA_NC002163- 112166-	TAACTGATGATGAGTCGGGTCCTT T	170			URCA_NC002163-112166- 112647_394_329_R	TGGATATACATCTGGATGGAATATAGCA GAGGCTC	1205
2567	FGM_NC002163- 327773-	TCTTGAATCTCTTAATATGGAGATA AATATCT	414			FGM_NC002163-327773- 328270_365_396_R	TCTATCCGACGCTTTTTCGATGATCGCTA AAAA	1012
2568	URCA_NC002163- 1569415-	TTAGAGAGCTGTCTTTTAGACGAGAC ATCA	661			TKT_NC002163-1569415- 1569873_350_383_R	TGAAAAGGCAATTTTATCACTCTCTGTAA AGGCTA	946
2570	GLFA_NC002163- 1604930-	TGCTCTTTTGATCTTCTCCCTATA ATCC	381			GLFA_NC002163-1604930- 1604529_109_142_R	TGTTCACTGTTTAAATGATCAGCTAAAA AGCACT	1347
2571	TKT_NC002163- 1569803_31_62_F	TGATCTTAAATTTCCGCGACATTC ATTC	472			TKT_NC002163-1569415- 1569903_139_162_R	TGCCATAGCAAGCTCTACAGCAAT T	1214
2572	GLFA_NC002163- 1569415-	TGAGCTTATATGCTTTTGTGAGATG GGGATTT	164			TKT_NC002163-1569415- 1569903_313_345_R	TACATCTCTCTGATAGAAATTTATG CTCT	886
2573	URCA_NC002163- 1569403_207_239_F	TAGCTCTTACAGCAATGTAAATATG CGTTTAA	213			TKT_NC002163-1569415- 1569903_449_481_R	TAGACAGGCTTTGTGATGTTTTTAGC TTGAT	865
2574	TKT_NC002163- 1569415-	TTCAAATCTACAGGCACTCTGAAAA TTTCAC	665			TKT_NC002163-1569415- 1569903_139_163_R	TTGCCATAGCAAGGCTACAGCAAT T	1405
2575	GLFA_NC002163- 1604930-	TGCTCTTTTATATCTTCTCCCTGATA ATGCTC	382			GLFA_NC002163-1604930- 1604529_139_168_R	TGCCATTTCCATATGATCTTCTTACAA TT	1216
2576	URCA_NC002163- 367572-	TGACATTTTCTTCAGGATATCAGAG TGG	281			GLFA_NC002163-367572- 368079_476_508_R	ATGCTCTCTTACTGCTTNGCATAAAT TTCCA	756
2577	GLFA_NC002163- 368079_148_174_F	TGTGTGAGATGCTCTATCTCTGTATA T	611			GLFA_NC002163-367572- 368079_242_270_R	TGCTGCTCTCTACAGCAAGTCCGCA T	1246
2578	TKT_NC002163- 367572-	TGTATGATCTTACACGACCAAACT TAGC	622			GLFA_NC002163-367572- 368079_381_416_R	TTCCATCTTGTATACCTGGAAAAATGC TGAAT	1381
2579	GLFA_NC002163- 367572-	TGTGTGACATTTAACCATGCTGCA A	614			GLFA_NC002163-367572- 368079_52_81_R	TGACGCTCTGACCTATAAAAAGCTCT CA	961
2580	FGM_NC002163- 327746-	TGAGCATCGGGCTTTGAAGAAATTT TTATAT	455			FGM_NC002163-327746- 328270_356_379_R	TTTGCTCTCGCAAGTTTCCAC TTTCTCTCTGATGATGAGTAC	1438
2581	FGM_NC002163- 327746-	TGAABAAGGATGATGATCAATGAG ATAG	425			FGM_NC002163-327746- 328270_241_267_R	TGCCCCATCTCTGATGATGAGTAC TTTCTCTCTGCAAGTTTCCAC	1219

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2582	328270_13_182_F	FGM_NC002163-327746- 327746- 328270_19_50_F	TGCGCTTAATGCGCTTAATCAATGA AAATATG	568	FGM_NC002163-327746- 328270_19_102_R	TGCGCAAGAGCTTACTCAGC	1200
2583	UNCA_NC002163- 112647_114_141_F	UNCA_NC002163- 112647_114_141_F	TBAGCTGCTGCTGCTTATCTGAGAA TG	160	UNCA_NC002163-112166- 112647_196_225_R	TGCGCTTCTAAAGCTCTTGAGTGAGA TA	1220
2584	UNCA_NC002163- 112166- 112647_3_25_F	UNCA_NC002163- 112166- 112647_3_25_F	TGCTTGGATCAGCAGCACTTCAAT A	532	UNCA_NC002163-112166- 112647_88_123_R	TGCTGCTCTCAATCAATCATCAATAACA ATTAAAGC	1206
2585	ASPA_NC002163- 96692- 97166_304_335_F	ASPA_NC002163- 96692- 97166_304_335_F	TTATTTTGCCTAAATGCAACAGGT AG	652	ASPA_NC002163-96692- 97166_403_432_R	TGCAAAATGAGGCTTACATCTGCTCA AT	1192
2586	GLNA_NC002163- 96692- 97166_228_258_F	GLNA_NC002163- 96692- 97166_228_258_F	TGCGTTGCGACCAACTTCTAAG TATGT	370	ASPA_NC002163-96692- 97166_316_346_R	TGATGATGAGCACTGCTGCTGCTTTT TGG	991
2587	GLNA_NC002163- 96692- 97166_244_275_F	GLNA_NC002163- 96692- 97166_244_275_F	TGCGATCATGCTAAGAGTTTCCGGA TAGCTA	547	GLNA_NC002163-658085- 957609_340_371_R	TGAGTTTGACCATTTCCAGAGCAAT CTAC	1176
2588	TKT_NC002163- 1569415- 1569903_107_130_F	TKT_NC002163- 1569415- 1569903_107_130_F	TGCGTACAGGCCCTTAGGCAAG TGATTA	371	TKT_NC002163-1569415- 1569903_212_236_R	TGCGCATTTCCGCAAGCAATATA TCTT	1020
2589	GLNA_NC002163- 367572- 368095_214_246_F	GLNA_NC002163- 367572- 368095_214_246_F	TGTTTCTTTAGCAGCACTTCCAAACT TGATTA	642	TKT_NC002163-1569415- 1569903_361_393_R	TGCTTGTGCTTCAAAAGCATTTTACA TTTC	1057
2590	GLNA_NC002163- 367572- 368095_415_444_F	GLNA_NC002163- 367572- 368095_415_444_F	TGCTACTCTTTTGTGATATGAC AAATGTC	505	GLNA_NC002163-367572- 368095_317_340_R	TGCTTGTGGCCGCGAAGTTTT TGCT	1047
2591	FGM_NC002163_21_54_F	FGM_NC002163_21_54_F	TGCTTTGATGCTATGATATGCTCA ATATGGA	353	GLNA_NC002163-367572- 368095_485_516_R	TGCTGAGCATGCTGCTTACTTGTGTTG CMTA	1141
2592	GLNA_NC002163- 658085- 657609_149_17_F	GLNA_NC002163- 658085- 657609_149_17_F	TGCTTGTGATGCTATGATATGCTCA ATATGGA	332	FGM_NC002163_116_142_R	TGCAAGATCCGCTCATGCTCAAGAG AAT	949
2593	GLNA_NC002163- 658085- 657609_149_17_F	GLNA_NC002163- 658085- 657609_149_17_F	TGCTTGTGATGCTATGATATGCTCA ATATGGA	207	FGM_NC002163_247_277_R	TGCGCTTTAAGCACTCATCTCATAT AAT	1023
2594	ASPA_NC002163- 96692- 97166_367_402_F	ASPA_NC002163- 96692- 97166_367_402_F	TGCTTGTGATGCTATGATATGCTCA ATATGGA	633	GLNA_NC002163-658085- 657609_148_179_R	TGCAAAACAAAGATTAATTTTGTGCTC CAAA	945
2595	ASPA_NC002163- 96692- 97166_367_402_F	ASPA_NC002163- 96692- 97166_367_402_F	TGCTTGTGATGCTATGATATGCTCA ATATGGA	347	ASPA_NC002163-96685- 97196_467_497_R	TGAGCTTATCTACACTGCTGCTCA AAG	960
2596	ASPA_NC002163- 96685- 97196_1_13_F	ASPA_NC002163- 96685- 97196_1_13_F	TGCTTGTGATGCTATGATATGCTCA ATATGGA	502	ASPA_NC002163-96685- 97196_95_127_R	TGAGCTTGTGATTAATGAGTGAAGA TMAAT	880
2597	ASPA_NC002163- 96685- 97196_188_210_R	ASPA_NC002163- 96685- 97196_188_210_R	TGCTTGTGATGCTATGATATGCTCA ATATGGA	540	ASPA_NC002163-96685- 97196_188_210_R	TGAGCTTGTGATTAATGAGTGAAGA TMAAT	872

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO.: 10593)

2598	97196_85_117_P PBM_NC002163- 327746- 328270_165 195 P	TGCGAGCTAGATATAGTACTGATTAATC GATAC	563	PBM_NC002163-327746- 328270_230 261 R	TGCGAGCTAAATTTGGATGAGCGGTAG GAAA	975
2599	PBM_NC002163- 327746- 328270_252 286 P	TGCGCTCGGTGTTTACGGAATATTC TTATATATG	593	PBM_NC002163-327746- 328270_353 381 R	TTTGTCTCGATCTCGATGAAGCTAA A	1443
2600	PBM_NC002163- 327746- 328270_1 30 P	TGCGAATGAAAGAGGTTCTTTTATCC ATGA	577	PBM_NC002163-327746- 328270_95 123 R	TGATAAAGCACTAAGCGATGAAACAG C	1178
2601	PBM_NC002163- 328270_220 250 P	TBAACAGCGCTTTCCTATGCGCTATC CAAT	146	PBM_NC002163-327746- 328270_314 345 R	TCAAGTCTTTTACTTCTATGATTTAA GCTC	963
2602	UNCA_NC002163- 112166- 112647_123 152 P	TGAGCTATGCGGAATGCTTTGGA TTTT	628	UNCA_NC002163-112166- 112647_199 229 R	TGCTTGTCTTTTGAGCGAGCTTGGAAG AAG	1258
2603	UNCA_NC002163- 112647_333 365 P	TGCGAGTGGACAAATTTCTTGAGAAA GCAATTT	313	UNCA_NC002163-112166- 112647_430 461 R	TCCGAAGCTGTTTGTAGCTTTAATTT GAGC	1031
2734	GYBA_AY291534_237_2 64_P	TGACGCTCATGTAATTCAGCTGTTT AT	265	GYBA_AY291534_260 288 R	TTGCGCTACTGACTACATCGT TGCAATGACTACATCGTTTCTTAACA GC	1407
2735	GYBA_AY291534_224_2 52_P	TATTCGGTAATGATGACCTCATGGT GAT	167	GYBA_AY291534_256 285 R	TTGCGCTACTGACTACATCGT TTGCGCTACTGACTACATCGT	1407
2736	GYBA_AY291534_170_1 98_P	TAGGAATACGCTGATTAAGCGGTAT AAA	221	GYBA_AY291534_268 288 R	TTGCGCTACTGACTACATCGT TTGCGCTACTGACTACATCGT	935
2737	GYBA_AY291534_234_2 52_P	TATTCGGTAATGATGACCTCATGGT GAT	167	GYBA_AY291534_319 346 R	TTTTCGAGATCCAAAGTTACATCGCC GAT	1142
2738	GYBA_NC002953-7005- 8668_221 249 P	TGAGGTATGACGGAATTAATCAATA GAT	163	GYBA_NC002953-7005- 8668_265 287 R	TCTTGAGCATACCTACATCGATGC GAT	912
2739	GYBA_NC002953-7005- 8668_221 249 P	TGAGGTATGACGGAATTAATCAATA GAT	171	GYBA_NC002953-7005- 8668_253 283 R	TATCCATTAACCAAAAGTTACTTGGCC AGA	1142
2740	GYBA_NC002953-7005- 8668_221 249 P	TGAGGTATGACGGAATTAATCAATA GAT	171	GYBA_NC002953-7005- 8668_265 287 R	TATCCATTAACCAAAAGTTACTTGGCC AGA	1142
2741	GYBA_NC002953-7005- 8668_224 261 P	TGACCTCATGTTGTAATCACTATTT AT	264	GYBA_NC002953-7005- 8668_265 287 R	TTGCGCATACGACTACATCGC GAT	1299
2842	CAPC_AF188935- 56074- 56074_271 304 P	TGGAATATATTTATCTGTTATATCC ATTGAGA	578	CAPC_AF188935-56074- 55628_346 378 R	TGCGTACCTGCTCTTGATATGATTT GCA	1314
2843	CAPC_AF188935- 56074- 56074_271 302 P	TGATATATATATCTGTTATATCCP ATPTATGAG	476	CAPC_AF188935-56074- 55628_349 377 R	TTTTCGATCTGCTCTTGATATGATTT GAT	1344
2844	CAPC_AF188935- 56074- 56074_268 303 P	TGATATATATATCTGTTATATCCP GCAATGAG	331	CAPC_AF188935-56074- 55628_349 384 R	TGATATATGATTAACCTCTGCTGATTT GATATTC	860
2845	CAPC_AF188935- 56074- 56074_268 303 P	TGATATATATATCTGTTATATCCP GCAATGAG	331	CAPC_AF188935-56074- 55628_349 384 R	TGATATATGATTAACCTCTGCTGATTT GATATTC	860

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

56074- 35626_366_303_F	GCCTTTTGG		55628_337_375_R	ATTATCTCG	
2846	PABC_X956119_33_58_P TCCGAAATACGCGTACGATGG TACTTGGTAAATACCCACATGGT GA	302	PABC_X956119_121_153_R	TAAAGATAGCGTACTTAATGCTGGA GCAT	852
2847	PABC_X956119_65_92_P	199	PABC_X958119_157_178_R	TACCCAGTTTCCCTGACCTTC TACGACATGATGATGATGCTTACAC ATGC	889
2848	PABC_X956119_69_93_F PABC_NC0033997- 3362578- 3365001_181_205_P	596	PABC_X958119_97_128_R	TACCCAGTTTCCCTGACCTTC TACGACATGATGATGATGCTTACAC ATGC	1169
2849	PABC_NC0033997- 3362578- 3365001_217_240_P	668	PABC_NC0033997-3362578- 3365001_256_283_R	TCCAGTTTACCTTAAAGTACACATGC TCTC	1001
2850	PABC_NC0033997- 3362578- 3365001_217_240_P	621	PABC_NC0033997-3362578- 3365001_304_335_R	TCCAGTTTACCTTAAAGTACACATGC TCTC	1099
2851	GYPA_AY642140_26_54 F	621	PABC_NC0033997-3362578- 3365001_244_275_R	TGACTTAAACCTACCATCTCTCTATTA CAGA	1162
2852	GYPA_AY642140_26_54 F	150	GYPA_AY642140_71_100_R	TGCTTAATGCTTTACAGTACACAT GC	1242
2853	GYPA_AY642140_26_54 F	166	GYPA_AY642140_121_146_R	TGATGACACCATGATTTACCTTACCC CAGC	1069
2854	CYA_AY065404_1348_1 P	166	GYPA_AY642140_58_89_R	TGACGATACGACATGATTTTCTATAA CAGC	1168
2860	LEF_BA_AY065404_751 P	305	CYA_AY065404_1448_1472_R	TGATGCTTACGATGCTTACAGCC TACGATGCTTACGATGCTTACAGCC	983
2861	LEF_BA_AY065404_762 P	354	LEF_BA_AY065404_843_881_R	TGCTTATGATGCGAGGATGATTA TTTCTTTGTCG	1144
2862	MUTS_AY698802_106_1 P	498	LEF_BA_AY065404_843_881_R	TGCTTATGATGCGAGGATGATTA TTTCTTTGTCG	1144
2917	MUTS_AY698802_172_1 P	326	MUTS_AY698802_172_193_R	TGCGCTCTGGCGACATAGATTA T	1237
2918	MUTS_AY698802_228_2 P	187	MUTS_AY698802_228_252_R	TGATCTTGACTTTTTCGCGCTA T	965
2919	MUTS_AY698802_315_3 P	186	MUTS_AY698802_314_342_R	TGCTTATGATGCTTTCGCGATAGAG T	1097
2920	MUTS_AY698802_394_4 P	419	MUTS_AY698802_413_433_R	TCCAGACGACACATCTGTA T	1210
2921	AD_MLST-11- OIF007_991_1018_P	585	MUTS_AY698802_497_519_R	TCCGCTTACCTTGGCTGACCTGTA T	1040
2922	GAFA_NC002505_694_7 P	583	GAFA_NC002505_769_798_2_R	TAGTATACACATGACACCGCTGCT AT	923
2927	GAFA_NC002505_694_7 P	259	GAFA_NC002505_769_798_2_R	TGCTTATGCTTGGTATGATGACGGA AT	1060
2928	GAFA_NC002505_694_7 P	361	GAFA_NC002505_769_798_2_R	TGCTTATGCTTGGTATGATGACGGA AT	1061

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2929	2929	GRF_NCO02505_694_7 21_2_F	TCGRGAAAGACCAACCAAGTATGA	361	GRF_NCO02505_769_798_3 R	TCCTTATGACACTTATGATACACCGA	1059
2932	2932	INF8_EC_1364_1394_P	TGCTCTCTGCTGTCACAAATACAGAT	688	INF8_EC_1439_1468_R	TGCTCTCTGCTGTCACAAATACAGAT	1410
2933	2933	INF8_EC_1364_1394_2 F	TGCTCTCTGCTGTCACAAATACAGAT	689	INF8_EC_1439_1468_R	TGCTCTCTGCTGTCACAAATACAGAT	1410
2934	2934	INF8_EC_80_110_F	TGCTCTCTGCTGTCACAAATACAGAT	685	INF8_EC_1439_1468_R	TGCTCTCTGCTGTCACAAATACAGAT	1410
2949	2949	ACS_NCO02516- 970624-	TCGCGCCCTGCTCTGATGA	376	ACS_NCO02516-970624- 971013_364_383_R	TGCGCACGCGCAGAGACG	1265
2950	2950	ARO_NCO02516-26883- 27380_8_26_F	TCACCTGCTGCTGTCACAAATACAGAT	267	ARO_NCO02516-26883- 27380_8_26_F	TGCTCTCTGCTGTCACAAATACAGAT	1341
2951	2951	ARO_NCO02516-26883- 27380_8_26_F	TCACCTGCTGCTGTCACAAATACAGAT	705	ARO_NCO02516-26883- 27380_8_26_F	TGCTCTCTGCTGTCACAAATACAGAT	1056
2952	2952	GRF_NCO02516- 422654-	TGCGCTCTGCTGCTGCTGCTGCTG	551	GRF_NCO02516-422654- 4226174_127_146_R	TGCGCTCTGCTGCTGCTGCTGCTG	1091
2953	2953	GRF_NCO02516- 422654-	TGCGCTCTGCTGCTGCTGCTGCTG	448	GRF_NCO02516-422654- 4226174_214_233_R	TGCTCTCTGCTGCTGCTGCTGCTG	1256
2954	2954	GRF_NCO02516- 422654-	TTTGAAAGTATGCTGCTGCTGCTG	710	GRF_NCO02516-422654- 4226174_255_287_R	TGCTCTCTGCTGCTGCTGCTGCTG	1259
2955	2955	GRF_NCO02516- 422654-	TTTGAAAGTATGCTGCTGCTGCTG	670	GRF_NCO02516-422654- 4226174_288_309_R	TGCTCTCTGCTGCTGCTGCTGCTG	1239
2956	2956	GRF_NCO02516- 422654-	TTTGAAAGTATGCTGCTGCTGCTG	374	GRF_NCO02516-422654- 4226174_355_371_R	TGCTCTCTGCTGCTGCTGCTGCTG	1111
2957	2957	MUT_NCO02516- 5551158-	TCGCGCTCTGCTGCTGCTGCTGCTG	545	MUT_NCO02516-5551158- 5550717_99_116_R	TGCTCTCTGCTGCTGCTGCTGCTG	978
2958	2958	MUT_NCO02516- 5551158-	TCGCGCTCTGCTGCTGCTGCTGCTG	358	MUT_NCO02516-5551158- 5550717_256_277_R	TGCTCTCTGCTGCTGCTGCTGCTG	971
2959	2959	GRF_NCO02516- 2984589- 2984589_8_26_F	TCGCGCTCTGCTGCTGCTGCTGCTG	249	GRF_NCO02516-2984589- 2984589_8_26_F	TGCTCTCTGCTGCTGCTGCTGCTG	1095
2960	2960	GRF_NCO02516- 2984589- 2984589_218_233_F	TCGCGCTCTGCTGCTGCTGCTGCTG	195	GRF_NCO02516-2984589- 2984589_218_233_F	TGCTCTCTGCTGCTGCTGCTGCTG	1376
2961	2961	PPS_NCO02516- 1915014- 1915383_44_63_F	TCGCGCTCTGCTGCTGCTGCTGCTG	311	PPS_NCO02516-1915014- 1915383_44_63_F	TGCTCTCTGCTGCTGCTGCTGCTG	1014

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

2962	PPS_NC002516-1915014-1915383 240 258 F	TCGCATCTGTCACACCG	365	PPS_NC002516-1915014-1915383 341 360 R	TCCTCCGATCTCTGAGGAT	1052
2963	TRP_NC002516-671831-672273 24 42 F	TCCTCGTATCGSGTCGCGGA	527	TRP_NC002516-671831-672273 131 150 R	TCGAACTCTCTTGSGCTGCGA	1071
2964	TRP_NC002516-671831-672273 261 292 F	TCGCATCTGTCACACCGTC	490	TRP_NC002516-671831-672273 282 303 R	TCGATCTCTGCGCGCGGATCT	1182
2972	OTF007_1007_1034 F	TCGCATCTGTCACACCGTC	592	OTF007_1126 1153 R	TCGATCTCTGCGCGGATCT	924
2993	OMPU_NC002505-674528-675880 428 455 F	TCGCATCTGTCACACCGTC	667	OMPU_NC002505 544 567 R	TCGTCGTCGAAAACGTCGCTTC	1094
2994	GAPA_NC002505-506780-507937 691 721 F	TCGCATCTGTCACACCGTC	335	GAPA_NC002505-506780-507937 769 802 R	TTTTCTCTTTGTCGACTTCGATCTAC	1442
2995	GAPA_NC002505-506780-507937 691 721 2 F	TCGCATCTGTCACACCGTC	339	GAPA_NC002505-506780-507937 769 803 R	TCGTCGCTTTGTCGACTTCGATCTAC	1008
2996	GAPA_NC002505-506780-507937 682 721 F	TCGCATCTGTCACACCGTC	396	GAPA_NC002505-506780-507937 785 817 R	TCGTCGCTTTGTCGACTTCGATCTAC	1005
2997	GAPA_NC002505-506780-507937 691 721 3 F	TCGCATCTGTCACACCGTC	337	GAPA_NC002505-506780-507937 785 817 R	TCGTCGCTTTGTCGACTTCGATCTAC	1085
2998	GAPA_NC002505-506780-507937 691 721 4 F	TCGCATCTGTCACACCGTC	336	GAPA_NC002505-506780-507937 784 817 R	TCGTCGCTTTGTCGACTTCGATCTAC	1087
2999	GAPA_NC002505-506780-507937 691 721 5 F	TCGCATCTGTCACACCGTC	340	GAPA_NC002505-506780-507937 784 817 2 R	TCGTCGCTTTGTCGACTTCGATCTAC	1086
3000	GAPA_NC002505-506780-507937 691 721 6 F	TCGCATCTGTCACACCGTC	338	GAPA_NC002505-506780-507937 769 805 R	TCGTCGCTTTGTCGACTTCGATCTAC	1430
3001	CTMX_NC002505-1565967-1567341 46 71 F	TCGCATCTGTCACACCGTC	275	CTMX_NC002505-1565967-1567341 139 163 R	TCGTCGCTTTGTCGACTTCGATCTAC	1025
3002	CTMX_NC002505-1565967-1567341 46 70 F	TCGCATCTGTCACACCGTC	274	CTMX_NC002505-1565967-1567341 132 162 R	TCGTCGCTTTGTCGACTTCGATCTAC	1018
3003	CTMX_NC002505-1565967-1567341 46 70 F	TCGCATCTGTCACACCGTC	274	CTMX_NC002505-1565967-1567341 118 350 R	TCGTCGCTTTGTCGACTTCGATCTAC	1225
3004	TUPR_NC002758-615038-	TCGCATCTGTCACACCGTC	180	TUPR_NC002758-615038-	TCGTCGCTTTGTCGACTTCGATCTAC	982

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

616222_694_P	TUPE_NCO02758- 616222_694_P	TGCGTGTGAGAGTGTGCAAT	503	TUPE_NCO02758-615038- 616222_783_813_R	TGCTTCAGCGTAGCTAATTAATTAAGG AAC	1255
3005	TUPE_NCO02758- 615038- 616222_700_726_P	TGCTGTCAATCAATGTTGGAGAGA A	638	TUPE_NCO02758-615038- 616222_778_807_R	TGCTAGCTCTAATTAATTAATTAAGCAACT TC	1238
3006	TUPE_NCO02758- 615038- 616222_702_726_P	TGCTCAATCAATGTTGGAGAGA	607	TUPE_NCO02758-615038- 616222_778_807_R	TGCTAGCTCTAATTAATTAATTAAGCAACT TC	1238
3007	TUPE_NCO02758- 615038- 616222_696_726_P	TGAGCTGTGTCAATCAATGTTGGT AGAGA	431	TUPE_NCO02758-615038- 616222_785_818_R	TCACACCTTCAGCGTAGCTCTAATTAAT TAGAGA	970
3008	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3009	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3010	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3011	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3012	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3013	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3014	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3015	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3016	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3017	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3018	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3019	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3020	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3021	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134
3022	TUPE_NCO02758- 615038- 616222_700_716_P	TGCTGTGTGAGCTGTGTCAATCAAG T	386	TUPE_NCO02758-615038- 616222_778_812_R	TGCTACACCTTACTAATTAATTAATTAAGGA AACTTC	1134

DOCKET NO.: DIBIS-0083US1 (COUNSEL DOCKET NO: 10593)

[illegible]

[370] Primer pair name codes and reference sequences are shown in Table 3. The primer name code typically represents the gene to which the given primer pair is targeted. The primer pair name may include specific coordinates with respect to a reference sequence defined by an extraction of a section of sequence or defined by a GenBank gi number, or the corresponding complementary sequence of the extraction, or the entire GenBank gi number as indicated by the label "no extraction." Where "no extraction" is indicated for a reference sequence, the coordinates of a primer pair named to the reference sequence are with respect to the GenBank gi listing. Gene abbreviations are shown in bold type in the "Gene Name" column.

[371] To determine the exact primer hybridization coordinates of a given pair of primers on a given bioagent nucleic acid sequence and to determine the sequences, molecular masses and base compositions of an amplification product to be obtained upon amplification of nucleic acid of a known bioagent with known sequence information in the region of interest with a given pair of primers, one with ordinary skill in bioinformatics is capable of obtaining alignments of the primers of the present invention with the GenBank gi number of the relevant nucleic acid sequence of the known bioagent. For example, the reference sequence GenBank gi numbers (Table 3) provide the identities of the sequences which can be obtained from GenBank. Alignments can be done using a bioinformatics tool such as BLASTn provided to the public by NCBI (Bethesda, MD). Alternatively, a relevant GenBank sequence may be downloaded and imported into custom programmed or commercially available bioinformatics programs wherein the alignment can be carried out to determine the primer hybridization coordinates and the sequences, molecular masses and base compositions of the amplification product. For example, to obtain the hybridization coordinates of primer pair number 2095 (SEQ ID NOs: 456:1261), First the forward primer (SEQ ID NO: 456) is subjected to a BLASTn search on the publicly available NCBI BLAST website. "RefSeq_Genomic" is chosen as the BLAST database since the gi numbers refer to genomic sequences. The BLAST query is then performed. Among the top results returned is a match to GenBank gi number 21281729 (Accession Number NC_003923). The result shown below, indicates that the forward primer hybridizes to positions 1530282..1530307 of the genomic sequence of *Staphylococcus aureus* subsp. *aureus* MW2 (represented by gi number 21281729).

Staphylococcus aureus subsp. *aureus* MW2, complete genome
Length=2820462

Features in this part of subject sequence:

Panton-Valentine leukocidin chain F precursor

Score = 52.0 bits (26), Expect = 2e-05
Identities = 26/26 (100%), Gaps = 0/26 (0%)

Strand=Plus/Plus

```

Query   1          TGAGCTGCATCAACTGTATTGGATAG   26
          |||||
Sbjct   1530282   TGAGCTGCATCAACTGTATTGGATAG   1530307

```

[372] The hybridization coordinates of the reverse primer (SEQ ID NO: 1261) can be determined in a similar manner and thus, the bioagent identifying amplicon can be defined in terms of genomic coordinates. The query/subject arrangement of the result would be presented in Strand = Plus/Minus format because the reverse strand hybridizes to the reverse complement of the genomic sequence. If the preceding sequence analyses are well known to one with ordinary skill in bioinformatics and thus, Table 3 contains sufficient information to determine the primer hybridization coordinates of any of the primers of Table 2 to the applicable reference sequences described therein.

Table 3: Primer Name Codes and Reference Sequence

Primer name code	Gene Name	Organism	Reference GenBank gi number
16S_EC	16S rRNA (16S ribosomal RNA gene)	<i>Escherichia coli</i>	16127994
23S_EC	23S rRNA (23S ribosomal RNA gene)	<i>Escherichia coli</i>	16127994
CAPC_BA	capC (capsule biosynthesis gene)	<i>Bacillus anthracis</i>	6470151
CYA_BA	cya (cyclic AMP gene)	<i>Bacillus anthracis</i>	4894216
DNAK_EC	dnaK (chaperone dnaK gene)	<i>Escherichia coli</i>	16127994
GROL_EC	groL (chaperonin groL)	<i>Escherichia coli</i>	16127994
HFLB_EC	hflB (cell division protein peptidase ftsH)	<i>Escherichia coli</i>	16127994
INF3_EC	inf3 (protein chain initiation factor inf3 gene)	<i>Escherichia coli</i>	16127994
LEF_BA	lef (lethal factor)	<i>Bacillus anthracis</i>	21392688
PAG_BA	pag (protective antigen)	<i>Bacillus anthracis</i>	21392688
RPLB_EC	rplB (50S ribosomal protein L2)	<i>Escherichia coli</i>	16127994
RPOB_EC	rpoB (DNA-directed RNA polymerase beta chain)	<i>Escherichia coli</i>	6127994
RPOC_EC	rpoC (DNA-directed RNA polymerase beta' chain)	<i>Escherichia coli</i>	16127994
SP101FT_SPET_11	Artificial Sequence Concatenation comprising: gki (glucose kinase) gtr (glutamine transporter protein) murI (glutamate racemase) mutS (DNA mismatch repair protein) xpt (xanthine phosphoribosyl transferase) yqiL (acetyl-CoA-acetyl transferase) tkt (transketolase)	Artificial Sequence* - partial gene sequences of <i>Streptococcus pyogenes</i>	15674250
SSPE_BA	sspE (small acid-soluble spore protein)	<i>Bacillus anthracis</i>	30253828
TUFD_EC	tufB (Elongation factor Tu)	<i>Escherichia coli</i>	16127994

VALS EC	valS (Valyl-tRNA synthetase)	<i>Escherichia coli</i>	16127994
ASPS EC	aspS (Aspartyl-tRNA synthetase)	<i>Escherichia coli</i>	16127994
CAF1 AF053947	cafI (capsular protein cafI)	<i>Yersinia pestis</i>	2996286
INV U22457	inv (invasin)	<i>Yersinia pestis</i>	1256565
LL NC003143	<i>Y. pestis</i> specific chromosomal genes - difference region	<i>Yersinia pestis</i>	16120353
BONTA X52066	BoNT/A (neurotoxin type A)	<i>Clostridium botulinum</i>	40381
MECA Y14051	mecA methicillin resistance gene	<i>Staphylococcus aureus</i>	2791983
TRPE AY094355	trpE (anthranilate synthase (large component))	<i>Acinetobacter baumannii</i>	20853695
RECA AF251469	recA (recombinase A)	<i>Acinetobacter baumannii</i>	9965210
GYRA AF100557	gyrA (DNA gyrase subunit A)	<i>Acinetobacter baumannii</i>	4240540
GYRB AB008700	gyrB (DNA gyrase subunit B)	<i>Acinetobacter baumannii</i>	4514436
WAA3 Z96925	waaA (3-deoxy-D-manno-octulosonic-acid transferase)	<i>Acinetobacter baumannii</i>	2765828
CJST_CJ	Artificial Sequence Concatenation comprising: tkk (transketolase) glyA (serine hydroxymethyltransferase) gltA (citrate synthase) aspA (aspartate ammonia lyase) glnA (glutamine synthase) pgm (phosphoglycerate mutase) uncA (ATP synthetase alpha chain)	Artificial Sequence* - partial gene sequences of <i>Campylobacter jejuni</i>	15791399
RNASEP BDP	RNase P (ribonuclease P)	<i>Bordetella pertussis</i>	33591275
RNASEP BKM	RNase P (ribonuclease P)	<i>Burkholderia mallei</i>	53723370
RNASEP BS	RNase P (ribonuclease P)	<i>Bacillus subtilis</i>	16077068
RNASEP CLB	RNase P (ribonuclease P)	<i>Clostridium perfringens</i>	18308982
RNASEP EC	RNase P (ribonuclease P)	<i>Escherichia coli</i>	16127994
RNASEP RKP	RNase P (ribonuclease P)	<i>Rickettsia prowazekii</i>	15603881
RNASEP SA	RNase P (ribonuclease P)	<i>Staphylococcus aureus</i>	15922990
RNASEP VBC	RNase P (ribonuclease P)	<i>Vibrio cholerae</i>	15640032
ICD CXB	icd (isocitrate dehydrogenase)	<i>Coxiella burnetii</i>	29732244
IS1111A	multi-locus IS1111A insertion element	<i>Acinetobacter baumannii</i>	29732244
OMPA AY485227	ompA (outer membrane protein A)	<i>Rickettsia prowazekii</i>	40287451
OMP B RKP	ompB (outer membrane protein B)	<i>Rickettsia prowazekii</i>	15603881
GLTA RKP	gltA (citrate synthase)	<i>Vibrio cholerae</i>	15603881
TOXR VBC	toxR (transcription regulator toxR)	<i>Francisella tularensis</i>	15640032
ASD FRT	asd (Aspartate semialdehyde dehydrogenase)	<i>Francisella tularensis</i>	56707187
GALE FRT	galE (UDP-glucose 4-epimerase)	<i>Shigella flexneri</i>	56707187
IPAH SGF	ipah (invasion plasmid antigen)	<i>Campylobacter jejuni</i>	30061571
HUPB CJ	hupB (DNA-binding protein Hu-beta)	<i>Coxiella burnetii</i>	15791399
AB MLST	Artificial Sequence Concatenation comprising: trpE (anthranilate synthase component)	Artificial Sequence* - partial gene	Sequenced in-house (SEQ ID NO: 1444)

	I)) adk (adenylate kinase) mutY (adenine glycosylase) fumC (fumarate hydratase) efp (elongation factor p) ppa (pyrophosphate phospho- hydratase)	sequences of <i>Acinetobacter baumannii</i>	
MUPR X75439	mupB (mupirocin resistance gene)	<i>Staphylococcus aureus</i>	438226
PARC X95819	parC (topoisomerase IV)	<i>Acinetobacter baumannii</i>	1212748
SED M28521	sed (enterotoxin D)	<i>Staphylococcus aureus</i>	1492109
PLA AF053945	pla (plasminogen activator)	<i>Yersinia pestis</i>	2996216
SEJ AF053140	sej (enterotoxin J)	<i>Staphylococcus aureus</i>	3372540
GYRA NC000912	gyrA (DNA gyrase subunit A)	<i>Mycoplasma pneumoniae</i>	13507739
ACS NC002516	acsA (Acetyl CoA Synthase)	<i>Pseudomonas aeruginosa</i>	15595198
ARO NC002516	aroE (shikimate 5-dehydrogenase)	<i>Pseudomonas aeruginosa</i>	15595198
GUA NC002516	guaA (GMP synthase)	<i>Pseudomonas aeruginosa</i>	15595198
MUT NC002516	mutL (DNA mismatch repair protein)	<i>Pseudomonas aeruginosa</i>	15595198
NUO NC002516	nuoD (NADH dehydrogenase I chain C, D)	<i>Pseudomonas aeruginosa</i>	15595198
PPS NC002516	ppsA (Phosphoenolpyruvate synthase)	<i>Pseudomonas aeruginosa</i>	15595198
TRP NC002516	trpE (Anthranilate synthetase component I)	<i>Pseudomonas aeruginosa</i>	15595198
OMP2 NC000117	ompB (outer membrane protein B)	<i>Chlamydia trachomatis</i>	15604717
OMPA NC000117	ompA (outer membrane protein B)	<i>Chlamydia trachomatis</i>	15604717
GYRA NC000117	gyrA (DNA gyrase subunit A)	<i>Chlamydia trachomatis</i>	15604717
CTXA NC002505	ctxA (Cholera toxin A subunit)	<i>Vibrio cholerae</i>	15640032
CTXB NC002505	ctxB (Cholera toxin B subunit)	<i>Vibrio cholerae</i>	15640032
FUR NC002505	fur (ferric uptake regulator protein)	<i>Vibrio cholerae</i>	15640032
GAPA NC 002505	gapA (glyceraldehyde-3-phosphate dehydrogenase)	<i>Vibrio cholerae</i>	15640032
GYRB NC002505	gyrB (DNA gyrase subunit B)	<i>Vibrio cholerae</i>	15640032
OMPU NC002505	ompU (outer membrane protein)	<i>Vibrio cholerae</i>	15640032
TCPA NC002505	tcpA (toxin-coregulated pilus)	<i>Vibrio cholerae</i>	15640032
ASPA NC002163	aspA (aspartate ammonia lyase)	<i>Campylobacter jejuni</i>	15791399
GLNA NC002163	glnA (glutamine synthetase)	<i>Campylobacter jejuni</i>	15791399
GLTA NC002163	gltA (glutamate synthase)	<i>Campylobacter jejuni</i>	15791399
GLYA NC002163	glyA (serine hydroxymethyltransferase)	<i>Campylobacter jejuni</i>	15791399
PGM NC002163	pgm (phosphoglyceromutase)	<i>Campylobacter jejuni</i>	15791399
TKT NC002163	tkt (transketolase)	<i>Campylobacter jejuni</i>	15791399
UNCA NC002163	uncA (ATP synthetase alpha chain)	<i>Campylobacter jejuni</i>	15791399
AGR-III NC003923	agr-III (accessory gene regulator-III)	<i>Staphylococcus aureus</i>	21281729
AROC NC003923	aroC (carbamate kinase)	<i>Staphylococcus aureus</i>	21281729
AROE NC003923	aroE (shikimate 5-dehydrogenase)	<i>Staphylococcus</i>	21281729

		aureus	
BSA-A NC003923	bsa-a (glutathione peroxidase)	Staphylococcus aureus	21281729
BSA-B NC003923	bsa-b (epidermin biosynthesis protein EpiB)	Staphylococcus aureus	21281729
GLPF NC003923	glpF (glycerol transporter)	Staphylococcus aureus	21281729
GMK NC003923	gmK (guanylate kinase)	Staphylococcus aureus	21281729
MECI-R NC003923	meclR (truncated methicillin resistance protein)	Staphylococcus aureus	21281729
PTA NC003923	pta (phosphate acetyltransferase)	Staphylococcus aureus	21281729
PVLUK NC003923	pvluk (Panton-Valentine leukocidin chain F precursor)	Staphylococcus aureus	21281729
SA442 NC003923	sa442 gene	Staphylococcus aureus	21281729
SEA NC003923	sea (staphylococcal enterotoxin A precursor)	Staphylococcus aureus	21281729
SEC NC003923	sec4 (enterotoxin type C precursor)	Staphylococcus aureus	21281729
TPI NC003923	tpi (triosephosphate isomerase)	Staphylococcus aureus	21281729
YQI NC003923	yqi (acetyl-CoA C-acetyltransferase homologue)	Staphylococcus aureus	21281729
GALE AF513299	galE (galactose epimerase)	Francisella tularensis	23506418
VVHA NC004460	vvha (cytotoxin, cytotoxin precursor)	Vibrio vulnificus	27366463
TDH NC004605	tdh (thermostable direct hemolysin A)	Vibrio parahaemolyticus	28899855
AGR-II NC002745	agr-II (accessory gene regulator-II)	Staphylococcus aureus	29165615
PARC NC003997	parC (topoisomerase IV)	Bacillus anthracis	30260195
GYRA AY291534	gyrA (DNA gyrase subunit A)	Bacillus anthracis	31323274
AGR-I AJ617706	agr-I (accessory gene regulator-I)	Staphylococcus aureus	46019543
AGR-IV AJ617711	agr-IV (accessory gene regulator-III)	Staphylococcus aureus	46019563
BLA3 NC002952	bla3 (beta lactamase III)	Staphylococcus aureus	49482253
ERMA NC002952	erma (rRNA methyltransferase A)	Staphylococcus aureus	49482253
ERMB Y13600	ermB (rRNA methyltransferase B)	Staphylococcus aureus	49482253
SEA-SEE NC002952	sea (staphylococcal enterotoxin A precursor)	Staphylococcus aureus	49482253
SEA-SEE NC002952	sea (staphylococcal enterotoxin A precursor)	Staphylococcus aureus	49482253
SEE NC002952	sea (staphylococcal enterotoxin A precursor)	Staphylococcus aureus	49482253
SEH NC002953	seh (staphylococcal enterotoxin H)	Staphylococcus aureus	49484912
ERMC NC005908	ermC (rRNA methyltransferase C)	Staphylococcus aureus	49489772
MUTS AY698802	mutS (DNA mismatch repair protein)	Shigella boydii	52698233
NUC NC002758	nuc (staphylococcal nuclease)	Staphylococcus aureus	57634611
SEB NC002758	seb (enterotoxin type B precursor)	Staphylococcus aureus	57634611
SEG NC002758	seg (staphylococcal enterotoxin G)	Staphylococcus aureus	57634611
SEI NC002758	sei (staphylococcal enterotoxin I)	Staphylococcus aureus	57634611
TSST NC002758	tsst (toxic shock syndrome toxin-1)	Staphylococcus aureus	57634611
TUFB NC002758	tufB (Elongation factor Tu)	Staphylococcus aureus	57634611

[373] Note: artificial reference sequences represent concatenations of partial gene extractions from the indicated reference gi number. Partial sequences were used to create the concatenated sequence because complete gene sequences were not necessary for primer design.

Example 2: Sample Preparation and PCR

[374] Genomic DNA was prepared from samples using the DNeasy Tissue Kit (Qiagen, Valencia, CA) according to the manufacturer's protocols.

[375] All PCR reactions were assembled in 50 μ L reaction volumes in a 96-well microtiter plate format using a Packard MP11 liquid handling robotic platform and M.J. Dyad thermocyclers (MJ research, Waltham, MA) or Eppendorf Mastercycler thermocyclers (Eppendorf, Westbury, NY). The PCR reaction mixture consisted of 4 units of Amplitaq Gold, 1x buffer II (Applied Biosystems, Foster City, CA), 1.5 mM MgCl₂, 0.4 M betaine, 800 μ M dNTP mixture and 250 nM of each primer. The following typical PCR conditions were used: 95°C for 10 min followed by 8 cycles of 95°C for 30 seconds, 48°C for 30 seconds, and 72°C 30 seconds with the 48°C annealing temperature increasing 0.9°C with each of the eight cycles. The PCR was then continued for 37 additional cycles of 95°C for 15 seconds, 56°C for 20 seconds, and 72°C 20 seconds.

Example 3: Purification of PCR Products for Mass Spectrometry with Ion Exchange Resin-Magnetic Beads

[376] For solution capture of nucleic acids with ion exchange resin linked to magnetic beads, 25 μ L of a 2.5 mg/mL suspension of BioClone amine terminated superparamagnetic beads were added to 25 μ L of a PCR (or RT-PCR) reaction containing approximately 10 pM of a typical PCR amplification product. The above suspension was mixed for approximately 5 minutes by vortexing or pipetting, after which the liquid was removed after using a magnetic separator. The beads containing bound PCR amplification product were then washed three times with 50mM ammonium bicarbonate/50% MeOH or 100mM ammonium bicarbonate/50% MeOH, followed by three more washes with 50% MeOH. The bound PCR amplicon was eluted with a solution of 25mM piperidine, 25mM imidazole, 35% MeOH which included peptide calibration standards.

Example 4: Mass Spectrometry and Base Composition Analysis

[377] The ESI-FTICR mass spectrometer is based on a Bruker Daltonics (Billerica, MA) Apex II 70e electrospray ionization Fourier transform ion cyclotron resonance mass spectrometer that employs an actively shielded 7 Tesla superconducting magnet. The active shielding constrains the majority of the fringing magnetic field from the superconducting magnet to a relatively small volume. Thus,

components that might be adversely affected by stray magnetic fields, such as CRT monitors, robotic components, and other electronics, can operate in close proximity to the FTICR spectrometer. All aspects of pulse sequence control and data acquisition were performed on a 600 MHz Pentium II data station running Bruker's Xmass software under Windows NT 4.0 operating system. Sample aliquots, typically 15 μ l, were extracted directly from 96-well microtiter plates using a CTC HTS PAL autosampler (LEAP Technologies, Carboro, NC) triggered by the FTICR data station. Samples were injected directly into a 10 μ l sample loop integrated with a fluidics handling system that supplies the 100 μ l/hr flow rate to the ESI source. Ions were formed via electrospray ionization in a modified Analytica (Branford, CT) source employing an off axis, grounded electrospray probe positioned approximately 1.5 cm from the metalized terminus of a glass desolvation capillary. The atmospheric pressure end of the glass capillary was biased at 6000 V relative to the ESI needle during data acquisition. A counter-current flow of dry N_2 was employed to assist in the desolvation process. Ions were accumulated in an external ion reservoir comprised of an rf-only hexapole, a skimmer cone, and an auxiliary gate electrode, prior to injection into the trapped ion cell where they were mass analyzed. Ionization duty cycles greater than 99% were achieved by simultaneously accumulating ions in the external ion reservoir during ion detection. Each detection event consisted of 1M data points digitized over 2.3 s. To improve the signal-to-noise ratio (S/N), 32 scans were co-added for a total data acquisition time of 74 s.

[378] The ESI-TOF mass spectrometer is based on a Bruker Daltonics MicroTOF™. Ions from the ESI source undergo orthogonal ion extraction and are focused in a reflectron prior to detection. The TOF and FTICR are equipped with the same automated sample handling and fluidics described above. Ions are formed in the standard MicroTOF™ ESI source that is equipped with the same off-axis sprayer and glass capillary as the FTICR ESI source. Consequently, source conditions were the same as those described above. External ion accumulation was also employed to improve ionization duty cycle during data acquisition. Each detection event on the TOF was comprised of 75,000 data points digitized over 75 μ s.

[379] The sample delivery scheme allows sample aliquots to be rapidly injected into the electrospray source at high flow rate and subsequently be electrosprayed at a much lower flow rate for improved ESI sensitivity. Prior to injecting a sample, a bolus of buffer was injected at a high flow rate to rinse the transfer line and spray needle to avoid sample contamination/carryover. Following the rinse step, the autosampler injected the next sample and the flow rate was switched to low flow. Following a brief equilibration delay, data acquisition commenced. As spectra were co-added, the autosampler continued rinsing the syringe and picking up buffer to rinse the injector and sample transfer line. In general, two syringe rinses and one injector rinse were required to minimize sample carryover. During a

routine screening protocol a new sample mixture was injected every 106 seconds. More recently a fast wash station for the syringe needle has been implemented which, when combined with shorter acquisition times, facilitates the acquisition of mass spectra at a rate of just under one spectrum/minute.

[380] Raw mass spectra were post-calibrated with an internal mass standard and deconvoluted to monoisotopic molecular masses. Unambiguous base compositions were derived from the exact mass measurements of the complementary single-stranded oligonucleotides. Quantitative results are obtained by comparing the peak heights with an internal PCR calibration standard present in every PCR well at 500 molecules per well. Calibration methods are commonly owned and disclosed in U.S. Provisional Patent Application Serial No. 60/545,425 which is incorporated herein by reference in entirety.

Example 5: *De Novo* Determination of Base Composition of Amplification Products using Molecular Mass Modified Deoxynucleotide Triphosphates

[381] Because the molecular masses of the four natural nucleobases have a relatively narrow molecular mass range (A = 313.058, G = 329.052, C = 289.046, T = 304.046 – See Table 4), a persistent source of ambiguity in assignment of base composition can occur as follows: two nucleic acid strands having different base composition may have a difference of about 1 Da when the base composition difference between the two strands is $G \leftrightarrow A$ (-15.994) combined with $C \leftrightarrow T$ (+15.000). For example, one 99-mer nucleic acid strand having a base composition of $A_{27}G_{30}C_{21}T_{21}$ has a theoretical molecular mass of 30779.058 while another 99-mer nucleic acid strand having a base composition of $A_{26}G_{31}C_{22}T_{20}$ has a theoretical molecular mass of 30780.052. A 1 Da difference in molecular mass may be within the experimental error of a molecular mass measurement and thus, the relatively narrow molecular mass range of the four natural nucleobases imposes an uncertainty factor.

[382] The present invention provides for a means for removing this theoretical 1 Da uncertainty factor through amplification of a nucleic acid with one mass-tagged nucleobase and three natural nucleobases. The term “nucleobase” as used herein is synonymous with other terms in use in the art including “nucleotide,” “deoxynucleotide,” “nucleotide residue,” “deoxynucleotide residue,” “nucleotide triphosphate (NTP),” or deoxynucleotide triphosphate (dNTP).

[383] Addition of significant mass to one of the 4 nucleobases (dNTPs) in an amplification reaction, or in the primers themselves, will result in a significant difference in mass of the resulting amplification product (significantly greater than 1 Da) arising from ambiguities arising from the $G \leftrightarrow A$ combined with $C \leftrightarrow T$ event (Table 4). Thus, the same the $G \leftrightarrow A$ (-15.994) event combined with 5-Iodo-C \leftrightarrow T (-110.900) event would result in a molecular mass difference of 126.894. If the molecular mass of the

base composition $A_{27}G_{30}5\text{-Iodo-C}_{21}T_{21}$ (33422.958) is compared with $A_{26}G_{31}5\text{-Iodo-C}_{22}T_{20}$, (33549.852) the theoretical molecular mass difference is +126.894. The experimental error of a molecular mass measurement is not significant with regard to this molecular mass difference. Furthermore, the only base composition consistent with a measured molecular mass of the 99-mer nucleic acid is $A_{27}G_{30}5\text{-Iodo-C}_{21}T_{21}$. In contrast, the analogous amplification without the mass tag has 18 possible base compositions.

Table 4: Molecular Masses of Natural Nucleobases and the Mass-Modified Nucleobase 5-Iodo-C and Molecular Mass Differences Resulting from Transitions

Nucleobase	Molecular Mass	Transition	Δ Molecular Mass
A	313.058	A \rightarrow T	-9.012
A	313.058	A \rightarrow C	-24.012
A	313.058	A \rightarrow 5-Iodo-C	101.888
A	313.058	A \rightarrow G	15.994
T	304.046	T \rightarrow A	9.012
T	304.046	T \rightarrow C	-15.000
T	304.046	T \rightarrow 5-Iodo-C	110.900
T	304.046	T \rightarrow G	25.006
C	289.046	C \rightarrow A	24.012
C	289.046	C \rightarrow T	15.000
C	289.046	C \rightarrow G	40.006
5-Iodo-C	414.946	5-Iodo-C \rightarrow A	-101.888
5-Iodo-C	414.946	5-Iodo-C \rightarrow T	-110.900
5-Iodo-C	414.946	5-Iodo-C \rightarrow G	-85.894
G	329.052	G \rightarrow A	-15.994
G	329.052	G \rightarrow T	-25.006
G	329.052	G \rightarrow C	-40.006
G	329.052	G \rightarrow 5-Iodo-C	85.894

[384] Mass spectra of bioagent-identifying amplicons were analyzed independently using a maximum-likelihood processor, such as is widely used in radar signal processing. This processor, referred to as GenX, first makes maximum likelihood estimates of the input to the mass spectrometer for each primer by running matched filters for each base composition aggregate on the input data. This includes the GenX response to a calibrant for each primer.

[385] The algorithm emphasizes performance predictions culminating in probability-of-detection versus probability-of-false-alarm plots for conditions involving complex backgrounds of naturally occurring organisms and environmental contaminants. Matched filters consist of *a priori* expectations of signal values given the set of primers used for each of the bioagents. A genomic sequence database is

used to define the mass base count matched filters. The database contains the sequences of known bacterial bioagents and includes threat organisms as well as benign background organisms. The latter is used to estimate and subtract the spectral signature produced by the background organisms. A maximum likelihood detection of known background organisms is implemented using matched filters and a running-sum estimate of the noise covariance. Background signal strengths are estimated and used along with the matched filters to form signatures which are then subtracted. The maximum likelihood process is applied to this "cleaned up" data in a similar manner employing matched filters for the organisms and a running-sum estimate of the noise-covariance for the cleaned up data.

[386] The amplitudes of all base compositions of bioagent-identifying amplicons for each primer are calibrated and a final maximum likelihood amplitude estimate per organism is made based upon the multiple single primer estimates. Models of all system noise are factored into this two-stage maximum likelihood calculation. The processor reports the number of molecules of each base composition contained in the spectra. The quantity of amplification product corresponding to the appropriate primer set is reported as well as the quantities of primers remaining upon completion of the amplification reaction.

[387] Base count blurring can be carried out as follows. "Electronic PCR" can be conducted on nucleotide sequences of the desired bioagents to obtain the different expected base counts that could be obtained for each primer pair. See for example, ncbi.nlm.nih.gov/sutils/e-pcr/; Schuler, *Genome Res.* 7:541-50, 1997. In one illustrative embodiment, one or more spreadsheets, such as Microsoft Excel workbooks contain a plurality of worksheets. First in this example, there is a worksheet with a name similar to the workbook name; this worksheet contains the raw electronic PCR data. Second, there is a worksheet named "filtered bioagents base count" that contains bioagent name and base count; there is a separate record for each strain after removing sequences that are not identified with a genus and species and removing all sequences for bioagents with less than 10 strains. Third, there is a worksheet, "Sheet1" that contains the frequency of substitutions, insertions, or deletions for this primer pair. This data is generated by first creating a pivot table from the data in the "filtered bioagents base count" worksheet and then executing an Excel VBA macro. The macro creates a table of differences in base counts for bioagents of the same species, but different strains. One of ordinary skill in the art may understand additional pathways for obtaining similar table differences without undo experimentation.

[388] Application of an exemplary script, involves the user defining a threshold that specifies the fraction of the strains that are represented by the reference set of base counts for each bioagent. The reference set of base counts for each bioagent may contain as many different base counts as are needed

to meet or exceed the threshold. The set of reference base counts is defined by taking the most abundant strain's base type composition and adding it to the reference set and then the next most abundant strain's base type composition is added until the threshold is met or exceeded. The current set of data was obtained using a threshold of 55%, which was obtained empirically.

[389] For each base count not included in the reference base count set for that bioagent, the script then proceeds to determine the manner in which the current base count differs from each of the base counts in the reference set. This difference may be represented as a combination of substitutions, $S_i=X_i$, and insertions, $I_i=Y_i$, or deletions, $D_i=Z_i$. If there is more than one reference base count, then the reported difference is chosen using rules that aim to minimize the number of changes and, in instances with the same number of changes, minimize the number of insertions or deletions. Therefore, the primary rule is to identify the difference with the minimum sum (X_i+Y_i) or (X_i+Z_i), e.g., one insertion rather than two substitutions. If there are two or more differences with the minimum sum, then the one that will be reported is the one that contains the most substitutions.

[390] Differences between a base count and a reference composition are categorized as one, two, or more substitutions, one, two, or more insertions, one, two, or more deletions, and combinations of substitutions and insertions or deletions. The different classes of nucleobase changes and their probabilities of occurrence have been delineated in U.S. Patent Application Publication No. 2004209260 (U.S. Application Serial No. 10/418,514) which is incorporated herein by reference in entirety.

Example 6: Use of Broad Range Survey and Division Wide Primer Pairs for Identification of Bacteria in an Epidemic Surveillance Investigation

[391] This investigation employed a set of 16 primer pairs which is herein designated the "surveillance primer set" and comprises broad range survey primer pairs, division wide primer pairs and a single *Bacillus* clade primer pair. The surveillance primer set is shown in Table 5 and consists of primer pairs originally listed in Table 2. This surveillance set comprises primers with T modifications (note TMOD designation in primer names) which constitutes a functional improvement with regard to prevention of non-templated adenylation (*vide supra*) relative to originally selected primers which are displayed below in the same row. Primer pair 449 (non-T modified) has been modified twice. Its predecessors are primer pairs 70 and 357, displayed below in the same row. Primer pair 360 has also been modified twice and its predecessors are primer pairs 17 and 118.

Table 5: Bacterial Primer Pairs of the Surveillance Primer Set

Primer Pair	Forward Primer Name	Forward Primer	Reverse Primer Name	Reverse Primer	Target Gene
-------------	---------------------	----------------	---------------------	----------------	-------------

No.		(SEQ ID NO:)		(SEQ ID NO:)	
346	16S_EC_713_732_TM0D_F	202	16S_EC_789_809_TM0D_R	1110	16S rRNA
10	16S_EC_713_732_F	21	16S_EC_789_809	798	16S rRNA
347	16S_EC_785_806_TM0D_F	560	16S_EC_880_897_TM0D_R	1278	16S rRNA
11	16S_EC_785_806_F	118	16S_EC_880_897_R	830	16S rRNA
348	16S_EC_960_981_TM0D_F	706	16S_EC_1054_1073_TM0D_R	895	16S rRNA
14	16S_EC_960_981_F	672	16S_EC_1054_1073_R	735	16S rRNA
349	23S_EC_1826_1843_TM0D_F	401	23S_EC_1906_1924_TM0D_R	1156	23S rRNA
16	23S_EC_1826_1843_F	80	23S_EC_1906_1924_R	805	23S rRNA
352	INFB_EC_1365_1393_TM0D_F	687	INFB_EC_1439_1467_TM0D_R	1411	infB
34	INFB_EC_1365_1393_F	524	INFB_EC_1439_1467_R	1248	infB
354	RPOC_EC_2218_2241_TM0D_F	405	RPOC_EC_2313_2337_TM0D_R	1072	rpoC
52	RPOC_EC_2218_2241_F	81	RPOC_EC_2313_2337_R	790	rpoC
355	SSPE_BA_115_137_TM0D_F	255	SSPE_BA_197_222_TM0D_R	1402	sspE
58	SSPE_BA_115_137_F	45	SSPE_BA_197_222_R	1201	sspE
356	RPLB_EC_650_679_TM0D_F	232	RPLB_EC_739_762_TM0D_R	592	rplB
66	RPLB_EC_650_679_F	98	RPLB_EC_739_762_R	999	rplB
358	VALS_EC_1105_1124_TM0D_F	385	VALS_EC_1195_1218_TM0D_R	1093	vals
71	VALS_EC_1105_1124_F	77	VALS_EC_1195_1218_R	795	vals
359	RPOB_EC_1845_1866_TM0D_F	659	RPOB_EC_1909_1929_TM0D_R	1250	rpoB
72	RPOB_EC_1845_1866_F	233	RPOB_EC_1909_1929_R	825	rpoB
360	23S_EC_2646_2667_TM0D_F	409	23S_EC_2745_2765_TM0D_R	1434	23S rRNA
118	23S_EC_2646_2667_F	84	23S_EC_2745_2765_R	1389	23S rRNA
17	23S_EC_2645_2669_F	408	23S_EC_2744_2761_R	1252	23S rRNA
361	16S_EC_1090_1111_2_TM0D_F	697	16S_EC_1175_1196_TM0D_R	1398	16S rRNA
3	16S_EC_1090_1111_2_F	651	16S_EC_1175_1196_R	1359	16S rRNA
362	RPOB_EC_3799_3821_TM0D_F	581	RPOB_EC_3862_3888_TM0D_R	1325	rpoB
289	RPOB_EC_3799_3821_F	124	RPOB_EC_3862_3888_R	840	rpoB
363	RPOC_EC_2146_2174_TM0D_F	284	RPOC_EC_2227_2245_TM0D_R	898	rpoC
290	RPOC_EC_2146_2174_F	52	RPOC_EC_2227_2245_R	736	rpoC
367	TUFB_EC_957_979_TM0D_F	308	TUFB_EC_1034_1058_TM0D_R	1276	tufB
293	TUFB_EC_957_979_F	55	TUFB_EC_1034_1058_R	829	tufB
449	RPLB_EC_690_710_F	309	RPLB_EC_737_758_R	1336	rplB
357	RPLB_EC_688_710_TM0D_F	296	RPLB_EC_736_757_TM0D_R	1337	rplB
67	RPLB_EC_688_710_F	54	RPLB_EC_736_757_R	842	rplB

[392] The 16 primer pairs of the surveillance set are used to produce bioagent identifying amplicons whose base compositions are sufficiently different amongst all known bacteria at the species level to identify, at a reasonable confidence level, any given bacterium at the species level. As shown in Tables 6A-E, common respiratory bacterial pathogens can be distinguished by the base compositions of bioagent identifying amplicons obtained using the 16 primer pairs of the surveillance set. In some cases, triangulation identification improves the confidence level for species assignment. For example, nucleic acid from *Streptococcus pyogenes* can be amplified by nine of the sixteen surveillance primer pairs and *Streptococcus pneumoniae* can be amplified by ten of the sixteen surveillance primer pairs. The base

compositions of the bioagent identifying amplicons are identical for only one of the analogous bioagent identifying amplicons and differ in all of the remaining analogous bioagent identifying amplicons by up to four bases per bioagent identifying amplicon. The resolving power of the surveillance set was confirmed by determination of base compositions for 120 isolates of respiratory pathogens representing 70 different bacterial species and the results indicated that natural variations (usually only one or two base substitutions per bioagent identifying amplicon) amongst multiple isolates of the same species did not prevent correct identification of major pathogenic organisms at the species level.

[393] *Bacillus anthracis* is a well known biological warfare agent which has emerged in domestic terrorism in recent years. Since it was envisioned to produce bioagent identifying amplicons for identification of *Bacillus anthracis*, additional drill-down analysis primers were designed to target genes present on virulence plasmids of *Bacillus anthracis* so that additional confidence could be reached in positive identification of this pathogenic organism. Three drill-down analysis primers were designed and are listed in Tables 2 and 6. In Table 6, the drill-down set comprises primers with T modifications (note TMOD designation in primer names) which constitutes a functional improvement with regard to prevention of non-templated adenylation (*vide supra*) relative to originally selected primers which are displayed below in the same row.

Table 6: Drill-Down Primer Pairs for Confirmation of Identification of *Bacillus anthracis*

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO.)	Reverse Primer Name	Reverse Primer (SEQ ID NO.)	Target Gene
350	CAPC_BA_274_303_TMOD_F	476	CAPC_BA_349_376_TMOD_R	1314	capC
24	CAPC_BA_274_303_F	109	CAPC_BA_349_376_R	837	capC
351	CYA_BA_1353_1379_TMOD_F	355	CYA_BA_1448_1467_TMOD_R	1423	cyA
30	CYA_BA_1353_1379_F	64	CYA_BA_1448_1467_R	1342	cyA
353	LEF_BA_756_781_TMOD_F	220	LEF_BA_843_872_TMOD_R	1394	lef
37	LEF_BA_756_781_F	26	LEF_BA_843_872_R	1135	lef

[394] Phylogenetic coverage of bacterial space of the sixteen surveillance primers of Table 5 and the three *Bacillus anthracis* drill-down primers of Table 6 is shown in Figure 3 which lists common pathogenic bacteria. Figure 3 is not meant to be comprehensive in illustrating all species identified by the primers. Only pathogenic bacteria are listed as representative examples of the bacterial species that can be identified by the primers and methods of the present invention. Nucleic acid of groups of bacteria enclosed within the polygons of Figure 3 can be amplified to obtain bioagent identifying amplicons using the primer pair numbers listed in the upper right hand corner of each polygon. Primer coverage for polygons within polygons is additive. As an illustrative example, bioagent identifying amplicons can be obtained for *Chlamydia trachomatis* by amplification with, for example, primer pairs 346-349, 360 and

361, but not with any of the remaining primers of the surveillance primer set. On the other hand, bioagent identifying amplicons can be obtained from nucleic acid originating from *Bacillus anthracis* (located within 5 successive polygons) using, for example, any of the following primer pairs: 346-349, 360, 361 (base polygon), 356, 449 (second polygon), 352 (third polygon), 355 (fourth polygon), 350, 351 and 353 (fifth polygon). Multiple coverage of a given organism with multiple primers provides for increased confidence level in identification of the organism as a result of enabling broad triangulation identification.

[395] In Tables 7A-E, base compositions of respiratory pathogens for primer target regions are shown. Two entries in a cell, represent variation in ribosomal DNA operons. The most predominant base composition is shown first and the minor (frequently a single operon) is indicated by an asterisk (*). Entries with NO DATA mean that the primer would not be expected to prime this species due to mismatches between the primer and target region, as determined by theoretical PCR.

Table 7A – Base Compositions of Common Respiratory Pathogens for Bioagent Identifying Amplicons Corresponding to Primer Pair Nos: 346, 347 and 348

Organism	Strain	Primer 346 [A G C T]	Primer 347 [A G C T]	Primer 348 [A G C T]
<i>Klebsiella pneumoniae</i>	MGH78578	[29 32 25 13] [29 31 25 13]*	[23 38 28 26] [23 37 28 26]*	[26 32 28 30] [26 31 28 30]*
<i>Yersinia pestis</i>	CO-92 Biovar Orientalis	[29 32 25 13]	[22 39 28 26]	[29 30 28 29] [30 30 27 29]*
<i>Yersinia pestis</i>	KIMS P12 (Biovar Mediaevalis)	[29 32 25 13]	[22 39 28 26]	[29 30 28 29] [29 30 28 29]
<i>Yersinia pestis</i>	91001	[29 32 25 13]	[22 39 28 26]	[30 30 27 29]*
<i>Haemophilus influenzae</i>	KW20	[28 31 23 17]	[24 37 25 27]	[29 30 28 29]
<i>Pseudomonas aeruginosa</i>	PAO1	[30 31 23 15]	[26 36 29 24] [27 36 29 23]*	[26 32 29 29] [25 32 29 29]
<i>Pseudomonas fluorescens</i>	PF0-1	[30 31 23 15]	[26 35 29 25]	[28 31 28 29]
<i>Pseudomonas putida</i>	KF2440	[30 31 23 15]	[28 33 27 27]	[27 32 29 28]
<i>Legionella pneumophila</i>	Philadelphia-1	[30 30 24 15]	[33 33 23 27]	[29 28 28 31]
<i>Francisella tularensis</i>	schu 4	[32 29 22 16]	[28 38 26 26]	[25 32 28 31]
<i>Bordetella pertussis</i>	Tohama I	[30 29 24 16]	[23 37 30 24]	[30 32 30 26] [27 36 31 24] [20 42 35 19]*
<i>Burkholderia cepacia</i>	J2315	[29 29 27 14]	[27 32 26 29]	[27 36 31 24] [20 42 35 19]*
<i>Burkholderia pseudomallei</i>	K96243	[29 29 27 14]	[27 32 26 29]	[27 36 31 24]
<i>Neisseria gonorrhoeae</i>	FA 1090, ATCC 700825	[29 28 24 18]	[27 34 26 28]	[24 36 29 27]
<i>Neisseria meningitidis</i>	MC58 (serogroup B)	[29 28 26 16]	[27 34 27 27]	[25 35 30 26]
<i>Neisseria meningitidis</i>	serogroup C, FAM18	[29 28 26 16]	[27 34 27 27]	[25 35 30 26]
<i>Neisseria meningitidis</i>	Z2491 (serogroup A)	[29 28 26 16]	[27 34 27 27]	[25 35 30 26]
<i>Chlamydia pneumoniae</i>	TW-183	[31 27 22 19]	NO DATA	[32 27 27 29]

<i>Chlamydia pneumoniae</i>	AR39	[31 27 22 19]	NO DATA	[32 27 27 29]
<i>Chlamydia pneumoniae</i>	CWL029	[31 27 22 19]	NO DATA	[32 27 27 29]
<i>Chlamydia pneumoniae</i>	J138	[31 27 22 19]	NO DATA	[32 27 27 29]
<i>Corynebacterium diphtheriae</i>	NCTC13129	[29 34 21 15]	[22 38 31 25]	[22 33 25 34]
<i>Mycobacterium avium</i>	k10	[27 36 21 15]	[22 37 30 28]	[21 36 27 30]
<i>Mycobacterium avium</i>	104	[27 36 21 15]	[22 37 30 28]	[21 36 27 30]
<i>Mycobacterium tuberculosis</i>	CSU#93	[27 36 21 15]	[22 37 30 28]	[21 36 27 30]
<i>Mycobacterium tuberculosis</i>	CDC 1551	[27 36 21 15]	[22 37 30 28]	[21 36 27 30]
<i>Mycobacterium tuberculosis</i>	H37Rv (lab strain)	[27 36 21 15]	[22 37 30 28]	[21 36 27 30]
<i>Mycoplasma pneumoniae</i>	M129	[31 29 19 20]	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	MRSA252	[27 30 21 21]	[25 35 30 26]	[30 29 30 29] [29 31 30 29] *
<i>Staphylococcus aureus</i>	MSSA476	[27 30 21 21]	[25 35 30 26]	[30 29 30 29] [30 29 29 30] *
<i>Staphylococcus aureus</i>	COL	[27 30 21 21]	[25 35 30 26]	[30 29 30 29] [30 29 29 30] *
<i>Staphylococcus aureus</i>	Mu50	[27 30 21 21]	[25 35 30 26]	[30 29 30 29] [30 29 29 30] *
<i>Staphylococcus aureus</i>	MW2	[27 30 21 21]	[25 35 30 26]	[30 29 30 29] [30 29 29 30] *
<i>Staphylococcus aureus</i>	N315	[27 30 21 21]	[25 35 30 26]	[30 29 30 29] [30 29 29 30] *
<i>Staphylococcus aureus</i>	NCTC 8325	[27 30 21 21]	[25 35 30 26] [25 35 31 26] *	[30 29 30 29] [30 29 29 30]
<i>Streptococcus agalactiae</i>	NEM316	[26 32 23 18]	[24 36 31 25] [24 36 30 26] *	[25 32 29 30]
<i>Streptococcus equi</i>	NC 002955	[26 32 23 18]	[23 37 31 25]	[29 30 25 32]
<i>Streptococcus pyogenes</i>	MGAS8232	[26 32 23 18]	[24 37 30 25]	[25 31 29 31]
<i>Streptococcus pyogenes</i>	MGAS315	[26 32 23 18]	[24 37 30 25]	[25 31 29 31]
<i>Streptococcus pyogenes</i>	SSI-1	[26 32 23 18]	[24 37 30 25]	[25 31 29 31]
<i>Streptococcus pyogenes</i>	MGAS10394	[26 32 23 18]	[24 37 30 25]	[25 31 29 31]
<i>Streptococcus pyogenes</i>	Manfredo (M5)	[26 32 23 18]	[24 37 30 25]	[25 31 29 31]
<i>Streptococcus pyogenes</i>	SF370 (M1)	[26 32 23 18]	[24 37 30 25]	[25 31 29 31]
<i>Streptococcus pneumoniae</i>	670	[26 32 23 18]	[25 35 28 28]	[25 32 29 30]
<i>Streptococcus pneumoniae</i>	R6	[26 32 23 18]	[25 35 28 28]	[25 32 29 30]
<i>Streptococcus pneumoniae</i>	TIGR4	[26 32 23 18]	[25 35 28 28]	[25 32 30 29]
<i>Streptococcus gordonii</i>	NCTC7868	[25 33 23 18]	[24 36 31 25]	[25 31 29 31]
<i>Streptococcus mitis</i>	NCTC 12261	[26 32 23 18]	[25 35 30 26]	[25 32 29 30] [24 31 35 29] *
<i>Streptococcus mutans</i>	UA159	[24 32 24 19]	[25 37 30 24]	[28 31 26 31]

Table 7B – Base Compositions of Common Respiratory Pathogens for Bioagent Identifying

Amplicons Corresponding to Primer Pair Nos: 349, 360, and 356

Organism	Strain	Primer 349	Primer 360	Primer 356
----------	--------	------------	------------	------------

		[A G C T]	[A G C T]	[A G C T]
<i>Klebsiella pneumoniae</i>	MGH78578	[25 31 25 22]	[33 37 25 27]	NO DATA
<i>Yersinia pestis</i>	CO-92 Biovar Orientalis	[25 31 27 20]		
		[25 32 26 20]*	[34 35 25 28]	NO DATA
<i>Yersinia pestis</i>	KIMS P12 (Biovar Mediaevalis)	[25 31 27 20]		
		[25 32 26 20]*	[34 35 25 28]	NO DATA
<i>Yersinia pestis</i>	91001	[25 31 27 20]	[34 35 25 28]	NO DATA
<i>Haemophilus influenzae</i>	KW20	[28 28 25 20]	[32 38 25 27]	NO DATA
<i>Pseudomonas aeruginosa</i>	PAO1	[24 31 26 20]	[31 36 27 27]	
			[31 36 27 28]*	NO DATA
<i>Pseudomonas fluorescens</i>	PF0-1	NO DATA	[30 37 27 28]	
			[30 37 27 28]	NO DATA
<i>Pseudomonas putida</i>	KT2440	[24 31 26 20]	[30 37 27 28]	NO DATA
<i>Legionella pneumophila</i>	Philadelphia-1	[23 30 25 23]	[30 39 29 24]	NO DATA
<i>Francisella tularensis</i>	schu 4	[26 31 25 19]	[32 36 27 27]	NO DATA
<i>Bordetella pertussis</i>	Tohama I	[21 29 24 18]	[33 36 26 27]	NO DATA
<i>Burkholderia cepacia</i>	J2315	[23 27 22 20]	[31 37 28 26]	NO DATA
<i>Burkholderia pseudomallei</i>	K96243	[23 27 22 20]	[31 37 28 26]	NO DATA
<i>Neisseria gonorrhoeae</i>	FA 1090, ATCC 700825	[24 27 24 17]	[34 37 25 26]	NO DATA
<i>Neisseria meningitidis</i>	MC58 (serogroup B)	[25 27 22 18]	[34 37 25 26]	NO DATA
<i>Neisseria meningitidis</i>	serogroup C, FAM18	[25 26 23 18]	[34 37 25 26]	NO DATA
<i>Neisseria meningitidis</i>	Z2491 (serogroup A)	[25 26 23 18]	[34 37 25 26]	NO DATA
<i>Chlamydia pneumoniae</i>	TW-183	[30 28 27 18]	NO DATA	NO DATA
<i>Chlamydia pneumoniae</i>	AR39	[30 28 27 18]	NO DATA	NO DATA
<i>Chlamydia pneumoniae</i>	CWL029	[30 28 27 18]	NO DATA	NO DATA
<i>Chlamydia pneumoniae</i>	J138	[30 28 27 18]	NO DATA	NO DATA
<i>Corynebacterium diphtheriae</i>	NCTC13129	NO DATA	[29 40 28 25]	NO DATA
<i>Mycobacterium avium</i>	x10	NO DATA	[33 35 32 22]	NO DATA
<i>Mycobacterium avium</i>	104	NO DATA	[33 35 32 22]	NO DATA
<i>Mycobacterium tuberculosis</i>	CSU#93	NO DATA	[30 36 34 22]	NO DATA
<i>Mycobacterium tuberculosis</i>	CDC 1551	NO DATA	[30 36 34 22]	NO DATA
<i>Mycobacterium tuberculosis</i>	H37Rv (lab strain)	NO DATA	[30 36 34 22]	NO DATA
<i>Mycoplasma pneumoniae</i>	MI29	[28 30 24 19]	[34 31 29 28]	NO DATA
<i>Staphylococcus aureus</i>	MRSA252	[26 30 25 20]	[31 38 24 29]	[33 30 31 27]
<i>Staphylococcus aureus</i>	MRSA476	[26 30 25 20]	[31 38 24 29]	[33 30 31 27]
<i>Staphylococcus aureus</i>	COL	[26 30 25 20]	[31 38 24 29]	[33 30 31 27]
<i>Staphylococcus aureus</i>	Mu50	[26 30 25 20]	[31 38 24 29]	[33 30 31 27]
<i>Staphylococcus aureus</i>	MW2	[26 30 25 20]	[31 38 24 29]	[33 30 31 27]
<i>Staphylococcus aureus</i>	N315	[26 30 25 20]	[31 38 24 29]	[33 30 31 27]
<i>Staphylococcus</i>	NCTC 8325	[26 30 25 20]	[31 38 24 29]	[33 30 31 27]

<i>aureus</i>				
<i>Streptococcus agalactiae</i>	NEM316	[28 31 22 20]	[33 37 24 28]	[37 30 28 26]
<i>Streptococcus equi</i>	NC 002955	[28 31 23 19]	[33 38 24 27]	[37 31 28 25]
<i>Streptococcus pyogenes</i>	MGAS8232	[28 31 23 19]	[33 37 24 28]	[38 31 29 23]
<i>Streptococcus pyogenes</i>	MGAS315	[28 31 23 19]	[33 37 24 28]	[38 31 29 23]
<i>Streptococcus pyogenes</i>	SSI-1	[28 31 23 19]	[33 37 24 28]	[38 31 29 23]
<i>Streptococcus pyogenes</i>	MGAS10394	[28 31 23 19]	[33 37 24 28]	[38 31 29 23]
<i>Streptococcus pyogenes</i>	Manfredo (M5)	[28 31 23 19]	[33 37 24 28]	[38 31 29 23]
<i>Streptococcus pyogenes</i>	SF370 (M1)	[28 31 22 20] *	[33 37 24 28]	[38 31 29 23]
<i>Streptococcus pneumoniae</i>	670	[28 31 22 20]	[34 36 24 28]	[37 30 29 25]
<i>Streptococcus pneumoniae</i>	R6	[28 31 22 20]	[34 36 24 28]	[37 30 29 25]
<i>Streptococcus pneumoniae</i>	TIGR4	[28 31 22 20]	[34 36 24 28]	[37 30 29 25]
<i>Streptococcus gordonii</i>	NCTC7868	[28 32 23 20]	[34 36 24 28]	[36 31 29 25]
<i>Streptococcus mitis</i>	NCTC 12261	[28 31 22 20] [29 30 22 20] *	[34 36 24 28]	[37 30 29 25]
<i>Streptococcus mutans</i>	UA159	[26 32 23 22]	[34 37 24 27]	NO DATA

Table 7C – Base Compositions of Common Respiratory Pathogens for Bioagent Identifying Amplicons Corresponding to Primer Pair Nos: 449, 354, and 352

Organism	Strain	Primer 449 [A G C T]	Primer 354 [A G C T]	Primer 352 [A G C T]
<i>Klebsiella pneumoniae</i>	MGH78578	NO DATA	[27 33 36 26]	NO DATA
<i>Yersinia pestis</i>	CO-92 Biovar Orientalis	NO DATA	[29 31 33 29]	[32 28 20 25]
<i>Yersinia pestis</i>	KIM5 F12 (Biovar Mediaevalis)	NO DATA	[29 31 33 29]	[32 28 20 25]
<i>Yersinia pestis</i>	91001	NO DATA	[29 31 33 29]	NO DATA
<i>Haemophilus influenzae</i>	KW20	NO DATA	[30 29 31 32]	NO DATA
<i>Pseudomonas aeruginosa</i>	PAO1	NO DATA	[26 33 39 24]	NO DATA
<i>Pseudomonas fluorescens</i>	Pf0-1	NO DATA	[26 33 34 29]	NO DATA
<i>Pseudomonas putida</i>	KT2440	NO DATA	[25 34 36 27]	NO DATA
<i>Legionella pneumophila</i>	Philadelphia-1	NO DATA	NO DATA	NO DATA
<i>Francisella tularensis</i>	schu 4	NO DATA	[33 32 25 32]	NO DATA
<i>Bordetella pertussis</i>	Tohama I	NO DATA	[26 33 39 24]	NO DATA
<i>Burkholderia cepacia</i>	J2315	NO DATA	[25 37 33 27]	NO DATA
<i>Burkholderia pseudomallei</i>	K96243	NO DATA	[25 37 34 26]	NO DATA
<i>Neisseria gonorrhoeae</i>	FA 1090, ATCC 700825	[17 23 22 10]	[29 31 32 30]	NO DATA
<i>Neisseria meningitidis</i>	MC58 (serogroup B)	NO DATA	[29 30 32 31]	NO DATA
<i>Neisseria meningitidis</i>	serogroup C, FAM18	NO DATA	[29 30 32 31]	NO DATA

<i>Neisseria meningitidis</i>	EZ491 (serogroup A)	NO DATA	[29 30 32 31]	NO DATA
<i>Chlamydomonas pneumoniae</i>	TW-183	NO DATA	NO DATA	NO DATA
<i>Chlamydomonas pneumoniae</i>	AR39	NO DATA	NO DATA	NO DATA
<i>Chlamydomonas pneumoniae</i>	CWL029	NO DATA	NO DATA	NO DATA
<i>Chlamydomonas pneumoniae</i>	J138	NO DATA	NO DATA	NO DATA
<i>Corynebacterium diphtheriae</i>	NCTC13129	NO DATA	NO DATA	NO DATA
<i>Mycobacterium avium</i>	k10	NO DATA	NO DATA	NO DATA
<i>Mycobacterium avium</i>	104	NO DATA	NO DATA	NO DATA
<i>Mycobacterium tuberculosis</i>	CSU#93	NO DATA	NO DATA	NO DATA
<i>Mycobacterium tuberculosis</i>	CDC 1551	NO DATA	NO DATA	NO DATA
<i>Mycobacterium tuberculosis</i>	H37Rv (lab strain)	NO DATA	NO DATA	NO DATA
<i>Mycoplasma pneumoniae</i>	M129	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	MRSA252	[17 20 21 17]	[30 27 30 35]	[36 24 19 26]
<i>Staphylococcus aureus</i>	MSSA476	[17 20 21 17]	[30 27 30 35]	[36 24 19 26]
<i>Staphylococcus aureus</i>	COL	[17 20 21 17]	[30 27 30 35]	[35 24 19 27]
<i>Staphylococcus aureus</i>	Mu50	[17 20 21 17]	[30 27 30 35]	[36 24 19 26]
<i>Staphylococcus aureus</i>	MW2	[17 20 21 17]	[30 27 30 35]	[36 24 19 26]
<i>Staphylococcus aureus</i>	N315	[17 20 21 17]	[30 27 30 35]	[36 24 19 26]
<i>Staphylococcus aureus</i>	NCTC 8325	[17 20 21 17]	[30 27 30 35]	[35 24 19 27]
<i>Streptococcus agalactiae</i>	NEM316	[22 20 19 14]	[26 31 27 38]	[29 26 22 28]
<i>Streptococcus equi</i>	NC 002955	[22 21 19 13]	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	MGAS8232	[23 21 19 12]	[24 32 30 36]	NO DATA
<i>Streptococcus pyogenes</i>	MGAS315	[23 21 19 12]	[24 32 30 36]	NO DATA
<i>Streptococcus pyogenes</i>	SSI-1	[23 21 19 12]	[24 32 30 36]	NO DATA
<i>Streptococcus pyogenes</i>	MGAS10394	[23 21 19 12]	[24 32 30 36]	NO DATA
<i>Streptococcus pyogenes</i>	Manfredo (M5)	[23 21 19 12]	[24 32 30 36]	NO DATA
<i>Streptococcus pyogenes</i>	SF370 (M1)	[23 21 19 12]	[24 32 30 36]	NO DATA
<i>Streptococcus pneumoniae</i>	670	[22 20 19 14]	[25 33 29 35]	[30 29 21 25]
<i>Streptococcus pneumoniae</i>	R5	[22 20 19 14]	[25 33 29 35]	[30 29 21 25]
<i>Streptococcus pneumoniae</i>	TIGRA	[22 20 19 14]	[25 33 29 35]	[30 29 21 25]
<i>Streptococcus gordonii</i>	NCTC7868	[21 21 19 14]	NO DATA	[29 26 22 28]
<i>Streptococcus mitis</i>	NCTC 12261	[22 20 19 14]	[26 30 32 34]	NO DATA
<i>Streptococcus mutans</i>	UA159	NO DATA	NO DATA	NO DATA

Table 7D – Base Compositions of Common Respiratory Pathogens for Bioagent Identifying Amplicons Corresponding to Primer Pair Nos: 355, 358, and 359

Organism	Strain	Primer 355 [A G C T]	Primer 358 [A G C T]	Primer 359 [A G C T]
<i>Klebsiella pneumoniae</i>	MGH78578	NO DATA	[24 39 33 20]	[25 21 24 17]
<i>Yersinia pestis</i>	CO-92 Biovar Orientalis	NO DATA	[26 34 35 21]	[23 23 19 22]
<i>Yersinia pestis</i>	KHS P12 (Biovar Mediævalis)	NO DATA	[26 34 35 21]	[23 23 19 22]
<i>Yersinia pestis</i>	91001	NO DATA	[26 34 35 21]	[23 23 19 22]
<i>Haemophilus influenzae</i>	KN20	NO DATA	NO DATA	NO DATA
<i>Pseudomonas aeruginosa</i>	PAO1	NO DATA	NO DATA	NO DATA
<i>Pseudomonas fluorescens</i>	Pf0-1	NO DATA	NO DATA	NO DATA
<i>Pseudomonas putida</i>	KT2440	NO DATA	[21 37 37 21]	NO DATA
<i>Legionella pneumophila</i>	Philadelphia-1	NO DATA	NO DATA	NO DATA
<i>Francisella tularensis</i>	schu 4	NO DATA	NO DATA	NO DATA
<i>Bordetella pertussis</i>	Tohama I	NO DATA	NO DATA	NO DATA
<i>Burkholderia cepacia</i>	J2315	NO DATA	NO DATA	NO DATA
<i>Burkholderia pseudomallei</i>	K96243	NO DATA	NO DATA	NO DATA
<i>Neisseria gonorrhoeae</i>	FA 1090, ATCC 700825	NO DATA	NO DATA	NO DATA
<i>Neisseria meningitidis</i>	MCS8 (serogroup B)	NO DATA	NO DATA	NO DATA
<i>Neisseria meningitidis</i>	serogroup C, FAM18	NO DATA	NO DATA	NO DATA
<i>Neisseria meningitidis</i>	Z2491 (serogroup A)	NO DATA	NO DATA	NO DATA
<i>Chlamydia pneumoniae</i>	TW-183	NO DATA	NO DATA	NO DATA
<i>Chlamydia pneumoniae</i>	AR39	NO DATA	NO DATA	NO DATA
<i>Chlamydia pneumoniae</i>	CWL029	NO DATA	NO DATA	NO DATA
<i>Chlamydia pneumoniae</i>	J138	NO DATA	NO DATA	NO DATA
<i>Corynebacterium diphtheriae</i>	NCTC13129	NO DATA	NO DATA	NO DATA
<i>Mycobacterium avium</i>	k10	NO DATA	NO DATA	NO DATA
<i>Mycobacterium avium</i>	104	NO DATA	NO DATA	NO DATA
<i>Mycobacterium tuberculosis</i>	CSU#93	NO DATA	NO DATA	NO DATA
<i>Mycobacterium tuberculosis</i>	CDC 1551	NO DATA	NO DATA	NO DATA
<i>Mycobacterium tuberculosis</i>	H37Rv (lab strain)	NO DATA	NO DATA	NO DATA
<i>Mycoplasma pneumoniae</i>	M129	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	MRSA252	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	MSSA476	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	COL	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	Mu50	NO DATA	NO DATA	NO DATA

<i>Staphylococcus aureus</i>	MF2	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	N315	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	NCTC 8325	NO DATA	NO DATA	NO DATA
<i>Streptococcus agalactiae</i>	NEM316	NO DATA	NO DATA	NO DATA
<i>Streptococcus equi</i>	NC 002955	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	MGAS8232	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	MGAS315	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	SSI-1	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	MGAS10394	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	Manfredo (MS)	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	SP370 (M1)	NO DATA	NO DATA	NO DATA
<i>Streptococcus pneumoniae</i>	670	NO DATA	NO DATA	NO DATA
<i>Streptococcus pneumoniae</i>	R6	NO DATA	NO DATA	NO DATA
<i>Streptococcus pneumoniae</i>	TIGR4	NO DATA	NO DATA	NO DATA
<i>Streptococcus gordonii</i>	NCTC7868	NO DATA	NO DATA	NO DATA
<i>Streptococcus mitis</i>	NCTC 12261	NO DATA	NO DATA	NO DATA
<i>Streptococcus mutans</i>	UA159	NO DATA	NO DATA	NO DATA

Table 7E – Base Compositions of Common Respiratory Pathogens for Bioagent Identifying Amplicons Corresponding to Primer Pair Nos: 362, 363, and 367

Organism	Strain	Primer 362 [A G C T]	Primer 363 [A G C T]	Primer 367 [A G C T]
<i>Klebsiella pneumoniae</i>	MGH78578	[21 33 22 16]	[16 34 26 26]	NO DATA
<i>Yersinia pestis</i>	CO-92 Biovar orientalis	[20 34 18 20]	NO DATA	NO DATA
<i>Yersinia pestis</i>	KIM5 P12 (Biovar Mediaevalis)	[20 34 18 20]	NO DATA	NO DATA
<i>Yersinia pestis</i>	91001	[20 34 18 20]	NO DATA	NO DATA
<i>Haemophilus influenzae</i>	KW20	NO DATA	NO DATA	NO DATA
<i>Pseudomonas aeruginosa</i>	PAO1	[19 35 21 17]	[16 36 28 22]	NO DATA
<i>Pseudomonas fluorescens</i>	PF0-1	NO DATA	[18 35 26 23]	NO DATA
<i>Pseudomonas putida</i>	KT2440	NO DATA	[16 35 28 23]	NO DATA
<i>Legionella pneumophila</i>	Philadelphia-1	NO DATA	NO DATA	NO DATA
<i>Francisella tularensis</i>	schu 4	NO DATA	NO DATA	NO DATA
<i>Bordetella pertussis</i>	Tohama I	[20 31 24 17]	[15 34 32 21]	[26 25 34 19]
<i>Burkholderia cepacia</i>	J2315	[20 33 21 18]	[15 36 26 25]	[25 27 32 20]
<i>Burkholderia pseudomallei</i>	K96243	[19 34 19 20]	[15 37 28 22]	[25 27 32 20]
<i>Neisseria gonorrhoeae</i>	FA 1090, ATCC 700825	NO DATA	NO DATA	NO DATA

<i>Neisseria meningitidis</i>	MCS8 (serogroup B)	NO DATA	NO DATA	NO DATA
<i>Neisseria meningitidis</i>	serogroup C, FAM18	NO DATA	NO DATA	NO DATA
<i>Neisseria meningitidis</i>	E2491 (serogroup A)	NO DATA	NO DATA	NO DATA
<i>Chlamydomonas pneumoniae</i>	TW-183	NO DATA	NO DATA	NO DATA
<i>Chlamydomonas pneumoniae</i>	AR39	NO DATA	NO DATA	NO DATA
<i>Chlamydomonas pneumoniae</i>	CWL029	NO DATA	NO DATA	NO DATA
<i>Chlamydomonas pneumoniae</i>	J138	NO DATA	NO DATA	NO DATA
<i>Corynebacterium diphtheriae</i>	NCTC13129	NO DATA	NO DATA	NO DATA
<i>Mycobacterium avium</i>	K10	[19 34 23 16]	NO DATA	[24 26 35 19]
<i>Mycobacterium avium</i>	104	[19 34 23 16]	NO DATA	[24 26 35 19]
<i>Mycobacterium tuberculosis</i>	CSU#93	[19 31 25 17]	NO DATA	[25 25 34 20]
<i>Mycobacterium tuberculosis</i>	CDC 1551	[19 31 24 18]	NO DATA	[25 25 34 20]
<i>Mycobacterium tuberculosis</i>	H37Rv (lab strain)	[19 31 24 18]	NO DATA	[25 25 34 20]
<i>Mycoplasma pneumoniae</i>	M129	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	MRSA252	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	MSSA476	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	COL	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	Mu50	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	MW2	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	N315	NO DATA	NO DATA	NO DATA
<i>Staphylococcus aureus</i>	NCTC 8325	NO DATA	NO DATA	NO DATA
<i>Streptococcus agalactiae</i>	NEM316	NO DATA	NO DATA	NO DATA
<i>Streptococcus equi</i>	NC 002955	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	MGAS8232	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	MGAS315	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	SSI-1	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	MGAS10394	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	Manfredo (M5)	NO DATA	NO DATA	NO DATA
<i>Streptococcus pyogenes</i>	SF370 (M1)	NO DATA	NO DATA	NO DATA
<i>Streptococcus pneumoniae</i>	670	NO DATA	NO DATA	NO DATA
<i>Streptococcus pneumoniae</i>	R6	[20 30 19 23]	NO DATA	NO DATA
<i>Streptococcus pneumoniae</i>	TIGR4	[20 30 19 23]	NO DATA	NO DATA
<i>Streptococcus gordonii</i>	NCTC7868	NO DATA	NO DATA	NO DATA
<i>Streptococcus mitis</i>	NCTC 12261	NO DATA	NO DATA	NO DATA
<i>Streptococcus</i>	UA159	NO DATA	NO DATA	NO DATA

<i>mutans</i>				
---------------	--	--	--	--

[396] Four sets of throat samples from military recruits at different military facilities taken at different time points were analyzed using the primers of the present invention. The first set was collected at a military training center from November 1 to December 20, 2002 during one of the most severe outbreaks of pneumonia associated with group A *Streptococcus* in the United States since 1968. During this outbreak, fifty-one throat swabs were taken from both healthy and hospitalized recruits and plated on blood agar for selection of putative group A *Streptococcus* colonies. A second set of 15 original patient specimens was taken during the height of this group A *Streptococcus* -associated respiratory disease outbreak. The third set were historical samples, including twenty-seven isolates of group A *Streptococcus*, from disease outbreaks at this and other military training facilities during previous years. The fourth set of samples was collected from five geographically separated military facilities in the continental U.S. in the winter immediately following the severe November/December 2002 outbreak.

[397] Pure colonies isolated from group A *Streptococcus*-selective media from all four collection periods were analyzed with the surveillance primer set. All samples showed base compositions that precisely matched the four completely sequenced strains of *Streptococcus pyogenes*. Shown in Figure 4 is a 3D diagram of base composition (axes A, G and C) of bioagent identifying amplicons obtained with primer pair number 14 (a precursor of primer pair number 348 which targets 16S rRNA). The diagram indicates that the experimentally determined base compositions of the clinical samples closely match the base compositions expected for *Streptococcus pyogenes* and are distinct from the expected base compositions of other organisms.

[398] In addition to the identification of *Streptococcus pyogenes*, other potentially pathogenic organisms were identified concurrently. Mass spectral analysis of a sample whose nucleic acid was amplified by primer pair number 349 (SEQ ID NOs: 401:1156) exhibited signals of bioagent identifying amplicons with molecular masses that were found to correspond to analogous base compositions of bioagent identifying amplicons of *Streptococcus pyogenes* (A27 G32 C24 T18), *Neisseria meningitidis* (A25 G27 C22 T18), and *Haemophilus influenzae* (A28 G28 C25 T20) (see Figure 5 and Table 7B). These organisms were present in a ratio of 4:5:20 as determined by comparison of peak heights with peak height of an internal PCR calibration standard as described in commonly owned U.S. Patent Application Serial No: 60/545,425 which is incorporated herein by reference in its entirety.

[399] Since certain division-wide primers that target housekeeping genes are designed to provide coverage of specific divisions of bacteria to increase the confidence level for identification of bacterial species, they are not expected to yield bioagent identifying amplicons for organisms outside of the specific divisions. For example, primer pair number 356 (SEQ ID NOs: 449:1380) primarily amplifies the nucleic acid of members of the classes *Bacilli* and *Clostridia* and is not expected to amplify proteobacteria such as *Neisseria meningitidis* and *Haemophilus influenzae*. As expected, analysis of the mass spectrum of amplification products obtained with primer pair number 356 does not indicate the presence of *Neisseria meningitidis* and *Haemophilus influenzae* but does indicate the presence of *Streptococcus pyogenes* (Figures 3 and 6, Table 7B). Thus, these primers or types of primers can confirm the absence of particular bioagents from a sample.

[400] The 15 throat swabs from military recruits were found to contain a relatively small set of microbes in high abundance. The most common were *Haemophilus influenza*, *Neisseria meningitidis*, and *Streptococcus pyogenes*. *Staphylococcus epidermidis*, *Moraxella catarrhalis*, *Corynebacterium pseudodiphtheriticum*, and *Staphylococcus aureus* were present in fewer samples. An equal number of samples from healthy volunteers from three different geographic locations, were identically analyzed. Results indicated that the healthy volunteers have bacterial flora dominated by multiple, commensal non-beta-hemolytic *Streptococcal* species, including the viridans group *streptococci* (*S. parasanguis*, *S. vestibularis*, *S. mitis*, *S. oralis* and *S. pneumoniae*; data not shown), and none of the organisms found in the military recruits were found in the healthy controls at concentrations detectable by mass spectrometry. Thus, the military recruits in the midst of a respiratory disease outbreak had a dramatically different microbial population than that experienced by the general population in the absence of epidemic disease.

Example 7: Triangulation Genotyping Analysis for Determination of emm-Type of *Streptococcus pyogenes* in Epidemic Surveillance

[401] As a continuation of the epidemic surveillance investigation of Example 6, determination of sub-species characteristics (genotyping) of *Streptococcus pyogenes*, was carried out based on a strategy that generates strain-specific signatures according to the rationale of Multi-Locus Sequence Typing (MLST). In classic MLST analysis, internal fragments of several housekeeping genes are amplified and sequenced (Enright et al. *Infection and Immunity*, 2001, 69, 2416-2427). In classic MLST analysis, internal fragments of several housekeeping genes are amplified and sequenced. In the present investigation, bioagent identifying amplicons from housekeeping genes were produced using drill-down primers and analyzed by mass spectrometry. Since mass spectral analysis results in molecular mass,

from which base composition can be determined, the challenge was to determine whether resolution of *emm* classification of strains of *Streptococcus pyogenes* could be determined.

[402] For the purpose of development of a triangulation genotyping assay, an alignment was constructed of concatenated alleles of seven MLST housekeeping genes (glucose kinase (*gki*), glutamine transporter protein (*gtr*), glutamate racemase (*murI*), DNA mismatch repair protein (*mutS*), xanthine phosphoribosyl transferase (*xpt*), and acetyl-CoA acetyl transferase (*yqiL*)) from each of the 212 previously *emm*-typed strains of *Streptococcus pyogenes*. From this alignment, the number and location of primer pairs that would maximize strain identification via base composition was determined. As a result, 6 primer pairs were chosen as standard drill-down primers for determination of *emm*-type of *Streptococcus pyogenes*. These six primer pairs are displayed in Table 8. This drill-down set comprises primers with T modifications (note TMOD designation in primer names) which constitutes a functional improvement with regard to prevention of non-templated adenylation (*vide supra*) relative to originally selected primers which are displayed below in the same row.

Table 8: Triangulation Genotyping Analysis Primer Pairs for Group A *Streptococcus* Drill-Down

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO:)	Reverse Primer Name	Reverse Primer (SEQ ID NO:)	Target Gene
442	SP101_SPET11_358_387_TMOD_F	588	SP101_SPET11_448_473_TMOD_R	998	<i>gki</i>
80	SP101_SPET11_358_387_F	126	SP101_SPET11_448_473_TMOD_R	766	<i>gki</i>
443	SP101_SPET11_600_629_TMOD_F	348	SP101_SPET11_686_714_TMOD_R	1018	<i>gtr</i>
81	SP101_SPET11_600_629_F	62	SP101_SPET11_686_714_R	772	<i>gtr</i>
426	SP101_SPET11_1314_1336_TMOD_F	363	SP101_SPET11_1403_1431_TMOD_R	849	<i>murI</i>
86	SP101_SPET11_1314_1336_F	68	SP101_SPET11_1403_1431_R	711	<i>murI</i>
430	SP101_SPET11_1807_1835_TMOD_F	235	SP101_SPET11_1901_1927_TMOD_R	1439	<i>mutS</i>
90	SP101_SPET11_1807_1835_F	33	SP101_SPET11_1901_1927_R	1412	<i>mutS</i>
438	SP101_SPET11_3075_3103_TMOD_F	473	SP101_SPET11_3168_3196_TMOD_R	875	<i>xpt</i>
96	SP101_SPET11_3075_3103_F	108	SP101_SPET11_3168_3196_R	715	<i>xpt</i>
441	SP101_SPET11_3511_3535_TMOD_F	531	SP101_SPET11_3605_3629_TMOD_R	1294	<i>yqiL</i>
98	SP101_SPET11_3511_3535_F	116	SP101_SPET11_3605_3629_R	832	<i>yqiL</i>

[403] The primers of Table 8 were used to produce bioagent identifying amplicons from nucleic acid present in the clinical samples. The bioagent identifying amplicons which were subsequently analyzed by mass spectrometry and base compositions corresponding to the molecular masses were calculated.

[404] Of the 51 samples taken during the peak of the November/December 2002 epidemic (Table 9A-C rows 1-3), all except three samples were found to represent *emm3*, a Group A *Streptococcus* genotype previously associated with high respiratory virulence. The three outliers were from samples obtained from healthy individuals and probably represent non-epidemic strains. Archived samples (Tables 9A-C, rows 5-13) from historical collections showed a greater heterogeneity of base compositions and *emm* types as would be expected from different epidemics occurring at different places and dates. The results of the mass spectrometry analysis and *emm* gene sequencing were found to be concordant for the epidemic and historical samples.

Table 9A: Base Composition Analysis of Bioagent Identifying Amplicons of Group A *Streptococcus* samples from Six Military Installations Obtained with Primer Pair Nos. 426 and 430

# of Instances	emm-type by Mass Spectrometry	emm-Gene Sequencing	Location (sample)	Year	murI (Primer Pair No. 426)	mutS (Primer Pair No. 430)
48	3	3	MCRD San Diego (Cultured)	2002	A39 G25 C20 T34	A38 G27 C23 T33
2	6	6			A40 G24 C20 T34	A38 G27 C23 T33
1	28	28			A39 G25 C20 T34	A38 G27 C23 T33
15	3	ND			A39 G25 C20 T34	A38 G27 C23 T33
6	3	3			A39 G25 C20 T34	A38 G27 C23 T33
3	5,58	5	NHRC San Diego-Archive (Cultured)	2003	A40 G24 C20 T34	A38 G27 C23 T33
6	6	6			A40 G24 C20 T34	A38 G27 C23 T33
1	11	11			A39 G25 C20 T34	A38 G27 C23 T33
3	12	12			A40 G24 C20 T34	A38 G26 C24 T33
1	22	22			A39 G25 C20 T34	A38 G27 C23 T33
3	25,75	75	Ft. Leonard Wood (Cultured)	2003	A39 G25 C20 T34	A38 G27 C23 T33
4	44/61,82,9	44/61			A40 G24 C20 T34	A38 G26 C24 T33
2	53,91	91			A39 G25 C20 T34	A38 G27 C23 T33
1	2	2			A39 G25 C20 T34	A38 G27 C24 T32
2	3	3			A39 G25 C20 T34	A38 G27 C23 T32
1	4	4	Ft. Sill (Cultured)	2003	A39 G25 C20 T34	A38 G27 C23 T33
1	6	6			A40 G24 C20 T34	A38 G27 C23 T33
11	25 or 75	75			A39 G25 C20 T34	A38 G27 C23 T33
1	25,75, 33, 34,4,52,84	75			A39 G25 C20 T34	A38 G27 C23 T33
1	44/61 or 82 or 9	44/61			A40 G24 C20 T34	A38 G26 C24 T33
2	5 or 58	5	Ft. Sill (Cultured)	2003	A40 G24 C20 T34	A38 G27 C23 T33
3	1	1			A40 G24 C20 T34	A38 G27 C23 T33
2	3	3			A39 G25 C20 T34	A38 G27 C23 T33
1	4	4			A39 G25 C20 T34	A38 G27 C23 T33
1	28	28			A39 G25 C20 T34	A38 G27 C23 T33
1	3	3	Ft.	2003	A39 G25 C20 T34	A38 G27 C23 T33

1	4	4	Benning (Cultured)		A39 G25 C20 T34	A38 G27 C23 T33
3	6	6			A40 G24 C20 T34	A38 G27 C23 T33
1	11	11			A39 G25 C20 T34	A38 G27 C23 T33
1	13	94**			A40 G24 C20 T34	A38 G27 C23 T33
1	44/61 or 82 or 9	82	Lackland AFB (Throat Swabs)	2003	A40 G24 C20 T34	A38 G26 C24 T33
1	5 or 58	58			A40 G24 C20 T34	A38 G27 C23 T33
1	78 or 89	89			A39 G25 C20 T34	A38 G27 C23 T33
2	5 or 58				A40 G24 C20 T34	A38 G27 C23 T33
1	2				A39 G25 C20 T34	A38 G27 C24 T32
1	81 or 90	ND			A40 G24 C20 T34	A38 G27 C23 T33
1	78				A38 G26 C20 T34	A38 G27 C23 T33
3***	No detection				No detection	No detection
7	3	ND			A39 G25 C20 T34	A38 G27 C23 T33
1	3	ND			No detection	A38 G27 C23 T33
1	3	ND	MCRD San Diego (Throat Swabs)	2002	No detection	No detection
1	3	ND			No detection	No detection
2	3	ND			No detection	A38 G27 C23 T33
3	No detection	ND			No detection	No detection

**Table 9B: Base Composition Analysis of Bioagent Identifying Amplicons of Group A
Streptococcus samples from Six Military Installations Obtained with Primer Pair Nos. 438 and 441**

# of Instances	emm-type by Mass Spectrometry	emm-Gene Sequencing	Location (sample)	Year	xpt (Primer Pair No. 438)	yqilL (Primer Pair No. 441)
48	3	3	MCRD San Diego (Cultured)	2002	A30 G36 C20 T36	A40 G29 C19 T31
2	6	6			A30 G36 C20 T36	A40 G29 C19 T31
1	28	28			A30 G36 C20 T36	A41 G28 C18 T32
15	3	ND			A30 G36 C20 T36	A40 G29 C19 T31
6	3	3	NHRC San Diego- Archive (Cultured)	2003	A30 G36 C20 T36	A40 G29 C19 T31
3	5,58	5			A30 G36 C20 T36	A40 G29 C19 T31
6	6	6			A30 G36 C20 T36	A40 G29 C19 T31
1	11	11			A30 G36 C20 T36	A40 G29 C19 T31
3	12	12			A30 G36 C19 T37	A40 G29 C19 T31
1	22	22			A30 G36 C20 T36	A40 G29 C19 T31
3	25,75	75			A30 G36 C20 T36	A40 G29 C19 T31
4	44/61,82,9	44/61			A30 G36 C20 T36	A41 G28 C19 T31
2	53,91	91			A30 G36 C19 T37	A40 G29 C19 T31
1	2	2			A30 G36 C20 T36	A40 G29 C19 T31
2	3	3	Ft. Leonard Wood (Cultured)	2003	A30 G36 C20 T36	A40 G29 C19 T31
1	4	4			A30 G36 C19 T37	A41 G28 C19 T31
1	6	6			A30 G36 C20 T36	A40 G29 C19 T31
11	25 or 75	75			A30 G36 C20 T36	A40 G29 C19 T31
1	25,75, 33, 34,4,52,84	75	Ft. Sill (Cultured)	2003	A30 G36 C19 T37	A40 G29 C19 T31
1	44/61 or 82 or 9	44/61			A30 G36 C20 T36	A41 G28 C19 T31
2	5 or 58	5			A30 G36 C20 T36	A40 G29 C19 T31
3	1	1			A30 G36 C19 T37	A40 G29 C19 T31
2	3	3	(Cultured)	2003	A30 G36 C20 T36	A40 G29 C19 T31
1	4	4			A30 G36 C19 T37	A41 G28 C19 T31

1	28	28	Ft. Benning (Cultured)	2003	A30 G36 C20 T36	A41 G28 C18 T32
1	3	3			A30 G36 C20 T36	A40 G29 C19 T31
1	4	4			A30 G36 C19 T37	A41 G28 C19 T31
3	6	6			A30 G36 C20 T36	A40 G29 C19 T31
1	11	11			A30 G36 C20 T36	A40 G29 C19 T31
1	13	94**			A30 G36 C20 T36	A41 G28 C19 T31
1	44/61 or 82 or 9	82			A30 G36 C20 T36	A41 G28 C19 T31
1	5 or 58	58			A30 G36 C20 T36	A40 G29 C19 T31
1	78 or 89	89			A30 G36 C20 T36	A41 G28 C19 T31
2	5 or 58				A30 G36 C20 T36	A40 G29 C19 T31
1	2	ND	Lackland AFB (Throat Swabs)	2003	A30 G36 C20 T36	A40 G29 C19 T31
1	81 or 90				A30 G36 C20 T36	A40 G29 C19 T31
1	78				A30 G36 C20 T36	A41 G28 C19 T31
3***	No detection				No detection	No detection
7	3	ND	MCRD San Diego (Throat Swabs)	2002	A30 G36 C20 T36	A40 G29 C19 T31
1	3	ND			A30 G36 C20 T36	A40 G29 C19 T31
1	3	ND			A30 G36 C20 T36	No detection
1	3	ND			No detection	A40 G29 C19 T31
2	3	ND			A30 G36 C20 T36	A40 G29 C19 T31
3	No detection	ND			No detection	No detection

Table 9C: Base Composition Analysis of Bioagent Identifying Amplicons of Group A *Streptococcus* samples from Six Military Installations Obtained with Primer Pair Nos. 438 and 441

# of Instances	emm-type by Mass Spectrometry	emm-Genes Sequencing	Location (sample)	Year	gk1 (Primer Pair No. 442)	gtr ((Primer Pair No. 443)
48	3	3	MCRD San Diego (Cultured)	2002	A32 G35 C17 T32	A39 G28 C16 T32
2	6	6			A31 G35 C17 T33	A39 G28 C15 T33
1	28	28			A30 G36 C17 T33	A39 G28 C16 T32
15	3	ND			A32 G35 C17 T32	A39 G28 C16 T32
6	3	3			A32 G35 C17 T32	A39 G28 C16 T32
3	5,58	5	NHRC San Diego-Archive (Cultured)	2003	A30 G36 C20 T30	A39 G28 C15 T33
6	6	6			A31 G35 C17 T33	A39 G28 C15 T33
1	11	11			A30 G36 C20 T30	A39 G28 C16 T32
3	12	12			A31 G35 C17 T33	A39 G28 C15 T33
1	22	22			A31 G35 C17 T33	A38 G29 C15 T33
3	25,75	75			A30 G36 C17 T33	A39 G28 C15 T32
4	44/61,82,9	44/61			A30 G36 C18 T32	A39 G28 C15 T33
2	53,91	91			A32 G35 C17 T32	A39 G28 C16 T32
1	2	2			A30 G36 C17 T33	A39 G28 C15 T33
2	3	3			A32 G35 C17 T32	A39 G28 C16 T32
1	4	4	Ft. Leonard Wood (Cultured)	2003	A31 G35 C17 T33	A39 G28 C15 T33
1	6	6			A31 G35 C17 T33	A39 G28 C15 T33
11	25 or 75	75			A30 G36 C17 T33	A39 G28 C15 T33
1	25,75, 33, 34,4,52,84	75			A30 G36 C17 T33	A39 G28 C15 T33
1	44/61 or 82 or 9	44/61			A30 G36 C18 T32	A39 G28 C15 T33
2	5 or 58	5	Ft. Sill	2003	A30 G36 C20 T30	A39 G28 C15 T33
3	1	1			A30 G36 C18 T32	A39 G28 C15 T33
2	3	3			A32 G35 C17 T32	A39 G28 C16 T32

1	4	4	(Cultured)		A31 G35 C17 T33	A39 G28 C15 T33
1	28	28			A30 G36 C17 T33	A39 G28 C16 T32
1	3	3			A32 G35 C17 T32	A39 G28 C16 T32
1	4	4			A31 G35 C17 T33	A39 G28 C15 T33
3	6	6			A31 G35 C17 T33	A39 G28 C15 T33
1	11	11			A30 G36 C20 T30	A39 G28 C16 T32
1	13	94**			A30 G36 C19 T31	A39 G28 C15 T33
1	44/61 or 82 or 9	82			A30 G36 C18 T32	A39 G28 C15 T33
1	5 or 58	58			A30 G36 C20 T30	A39 G28 C15 T33
1	78 or 89	89			A30 G36 C18 T32	A39 G28 C15 T33
2	5 or 58				A30 G36 C20 T30	A39 G28 C15 T33
1	2				A30 G36 C17 T33	A39 G28 C15 T33
1	81 or 90	ND			A30 G36 C17 T33	A39 G28 C15 T33
1	78				A30 G36 C18 T32	A39 G28 C15 T33
3***	No detection				No detection	No detection
7	3	ND			A32 G35 C17 T32	A39 G28 C16 T32
1	3	ND			No detection	No detection
1	3	ND			A32 G35 C17 T32	A39 G28 C16 T32
1	3	ND			A32 G35 C17 T32	No detection
2	3	ND			A32 G35 C17 T32	No detection
3	No detection	ND			No detection	No detection

Example 8: Design of Calibrant Polynucleotides based on Bioagent Identifying Amplicons for Identification of Species of Bacteria (Bacterial Bioagent Identifying Amplicons)

[405] This example describes the design of 19 calibrant polynucleotides based on bacterial bioagent identifying amplicons corresponding to the primers of the broad surveillance set (Table 5) and the *Bacillus anthracis* drill-down set (Table 6).

[406] Calibration sequences were designed to simulate bacterial bioagent identifying amplicons produced by the T modified primer pairs shown in Tables 5 and 6 (primer names have the designation "TMOD"). The calibration sequences were chosen as a representative member of the section of bacterial genome from specific bacterial species which would be amplified by a given primer pair. The model bacterial species upon which the calibration sequences are based are also shown in Table 10. For example, the calibration sequence chosen to correspond to an amplicon produced by primer pair no. 361 is SEQ ID NO: 1445. In Table 10, the forward (_F) or reverse (_R) primer name indicates the coordinates of an extraction representing a gene of a standard reference bacterial genome to which the primer hybridizes e.g.: the forward primer name 16S_EC_713_732_TMOD_F indicates that the forward primer hybridizes to residues 713-732 of the gene encoding 16S ribosomal RNA in an *E. coli* reference sequence (in this case, the reference sequence is an extraction consisting of residues 4033120-4034661 of the genomic sequence of *E. coli* K12 (GenBank gi number 16127994). Additional gene coordinate reference information is shown in Table 11. The designation "TMOD" in the primer names indicates that the 5' end of the primer has been modified with a non-matched template T residue which

prevents the PCR polymerase from adding non-templated adenosine residues to the 5' end of the amplification product, an occurrence which may result in miscalculation of base composition from molecular mass data (*vide supra*).

[0143] The 19 calibration sequences described in Tables 10 and 11 were combined into a single calibration polynucleotide sequence (SEQ ID NO: 1464 - which is herein designated a "combination calibration polynucleotide") which was then cloned into a pCR®-Blunt vector (Invitrogen, Carlsbad, CA). This combination calibration polynucleotide can be used in conjunction with the primers of Tables 5 or 6 as an internal standard to produce calibration amplicons for use in determination of the quantity of any bacterial bioagent. Thus, for example, when the combination calibration polynucleotide vector is present in an amplification reaction mixture, a calibration amplicon based on primer pair 346 (16S rRNA) will be produced in an amplification reaction with primer pair 346 and a calibration amplicon based on primer pair 363 (rpoC) will be produced with primer pair 363. Coordinates of each of the 19 calibration sequences within the calibration polynucleotide (SEQ ID NO: 1464) are indicated in Table 11.

Table 10: Bacterial Primer Pairs for Production of Bacterial Bioagent Identifying Amplicons and Corresponding Representative Calibration Sequences

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO)	Reverse Primer Name	Reverse Primer (SEQ ID NO)	Calibration Sequence Model Species	Calibration Sequence (SEQ ID NO)
361	16S_EC_1090_1111_2_TMOD_F	697	16S_EC_1175_1196_TMOD_R	1398	<i>Bacillus anthracis</i>	1445
346	16S_EC_713_732_TMOD_F	202	16S_EC_789_809_TMOD_R	1110	<i>Bacillus anthracis</i>	1446
347	16S_EC_785_806_TMOD_F	560	16S_EC_880_897_TMOD_R	1278	<i>Bacillus anthracis</i>	1447
348	16S_EC_960_981_TMOD_F	706	16S_EC_1054_1073_TMOD_R	895	<i>Bacillus anthracis</i>	1448
349	23S_EC_1826_1843_TMOD_F	401	23S_EC_1906_1924_TMOD_R	1156	<i>Bacillus anthracis</i>	1449
360	23S_EC_2646_2667_TMOD_F	409	23S_EC_2745_2765_TMOD_R	1434	<i>Bacillus anthracis</i>	1450
350	CAPC_BA_274_303_TMOD_F	476	CAPC_BA_349_376_TMOD_R	1314	<i>Bacillus anthracis</i>	1451
351	CYA_BA_1353_1379_TMOD_F	355	CYA_BA_1448_1467_TMOD_R	1423	<i>Bacillus anthracis</i>	1452
352	INF8_EC_1365_1393_TMOD_F	687	INF8_EC_1439_1467_TMOD_R	1411	<i>Bacillus anthracis</i>	1453
353	LEF_BA_756_781_TMOD_F	220	LEF_BA_843_872_TMOD_R	1394	<i>Bacillus anthracis</i>	1454
356	RPLB_EC_650_679_TMOD_F	449	RPLB_EC_739_762_TMOD_R	1380	<i>Clostridium botulinum</i>	1455
449	RPLB_EC_680_710_F	309	RPLB_EC_737_758_R	1336	<i>Clostridium botulinum</i>	1456
359	RPOB_EC_1845_1866_TMOD_F	659	RPOB_EC_1909_1929_TMOD_R	1250	<i>Yersinia Pestis</i>	1457
362	RPOB_EC_3759_3821_TMOD_F	581	RPOB_EC_3862_3888_TMOD_R	1325	<i>Burkholderia mallei</i>	1458
363	RPOC_EC_2146_2174_TMOD_F	284	RPOC_EC_2227_2245_TMOD_R	898	<i>Burkholderia mallei</i>	1459

354	RPOC_EC_2218_2241_TM OD_F	405	RPOC_EC_2313_2337_TM OD_R	1072	<i>Bacillus anthracis</i>	1460
355	SSPE_BA_115_137_TM _F	255	SSPE_BA_197_222_TM _R	1402	<i>Bacillus anthracis</i>	1461
367	TUFB_EC_957_979_TM _F	308	TUFB_EC_1034_1058_TM _R	1276	<i>Sarkoidia mallei</i>	1462
358	VALS_EC_1105_1124_TM OD_F	385	VALS_EC_1195_1218_TM _R	1093	<i>Yersinia Pestis</i>	1463

Table 11: Primer Pair Gene Coordinate References and Calibration Polynucleotide Sequence
Coordinates within the Combination Calibration Polynucleotide

Bacterial Gene and Species	Gene Extraction Coordinates of Genomic or Plasmid Sequence	Reference GenBank GI No. of Genomic (G) or Plasmid (P) Sequence	Primer Pair No.	Coordinates of Calibration Sequence in Combination Calibration polynucleotide (SEQ ID NO: 1464)
16S <i>E. coli</i>	4033120..4034661	16127994 (G)	346	16..109
16S <i>E. coli</i>	4033120..4034661	16127994 (G)	347	83..190
16S <i>E. coli</i>	4033120..4034661	16127994 (G)	348	246..353
23S <i>E. coli</i>	4166220..4169123	16127994 (G)	361	368..469
23S <i>E. coli</i>	4166220..4169123	16127994 (G)	349	743..837
23S <i>E. coli</i>	4166220..4169123	16127994 (G)	360	865..981
rcpB <i>E. coli</i>	4178823..4182851 (complement strand)	16127994 (G)	359	1591..1672
rcpB <i>E. coli</i>	4178823..4182851 (complement strand)	16127994 (G)	362	2081..2167
rcpC <i>E. coli</i>	4182928..4187151	16127994 (G)	354	1810..1926
rcpC <i>E. coli</i>	4182928..4187151	16127994 (G)	363	2183..2279
infB <i>E. coli</i>	3311655..3310983 (complement strand)	16127994 (G)	352	1692..1751
tuftB <i>E. coli</i>	4173523..4174707	16127994 (G)	367	2400..2498
rp1B <i>E. coli</i>	3449001..3448180	16127994 (G)	356	1945..2060
rp1B <i>E. coli</i>	3449001..3448180	16127994 (G)	449	1986..2055
valS <i>E. coli</i>	4481405..4478550 (complement strand)	16127994 (G)	358	1462..1572
capC <i>B. anthracis</i>	56074..55628 (complement strand)	6470151 (P)	350	2517..2616
cya <i>B. anthracis</i>	156626..154288 (complement strand)	4894216 (P)	351	1338..1449
lef <i>B. anthracis</i>	127442..129921	4894216 (P)	353	1321..1234
sapB <i>B. anthracis</i>	226496..226783	30253828 (P)	355	1007..1104

Example 9: Use of a Calibration Polynucleotide for Determining the Quantity of *Bacillus Anthracis* in a Sample Containing a Mixture of Microbes

[407] The process described in this example is shown in Figure 2. The capC gene is a gene involved in capsule synthesis which resides on the pX02 plasmid of *Bacillus anthracis*. Primer pair number 350 (see Tables 10 and 11) was designed to identify *Bacillus anthracis* via production of a bacterial bioagent identifying amplicon. Known quantities of the combination calibration polynucleotide vector described in Example 8 were added to amplification mixtures containing bacterial bioagent nucleic acid from a mixture of microbes which included the Ames strain of *Bacillus anthracis*. Upon amplification of the bacterial bioagent nucleic acid and the combination calibration polynucleotide vector with primer pair no. 350, bacterial bioagent identifying amplicons and calibration amplicons were obtained and characterized by mass spectrometry. A mass spectrum measured for the amplification reaction is shown in Figure 7. The molecular masses of the bioagent identifying amplicons provided the means for identification of the bioagent from which they were obtained (Ames strain of *Bacillus*

anthracis) and the molecular masses of the calibration amplicons provided the means for their identification as well. The relationship between the abundance (peak height) of the calibration amplicon signals and the bacterial bioagent identifying amplicon signals provides the means of calculation of the copies of the pX02 plasmid of the Ames strain of *Bacillus anthracis*. Methods of calculating quantities of molecules based on internal calibration procedures are well known to those of ordinary skill in the art.

[408] Averaging the results of 10 repetitions of the experiment described above, enabled a calculation that indicated that the quantity of Ames strain of *Bacillus anthracis* present in the sample corresponds to approximately 10 copies of pX02 plasmid.

Example 10: Triangulation Genotyping Analysis of *Campylobacter* Species

[409] A series of triangulation genotyping analysis primers were designed as described in Example 1 with the objective of identification of different strains of *Campylobacter jejuni*. The primers are listed in Table 12 with the designation "CJST_CJ." Housekeeping genes to which the primers hybridize and produce bioagent identifying amplicons include: tkt (transketolase), glyA (serine hydroxymethyltransferase), gltA (citrate synthase), aspA (aspartate ammonia lyase), glnA (glutamine synthase), pgm (phosphoglycerate mutase), and uncA (ATP synthetase alpha chain).

Table 12: *Campylobacter* Genotyping Primer Pairs

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO:)	Reverse Primer Name	Reverse Primer (SEQ ID NO:)	Target Gene
1053	CJST_CJ_1080_1110_F	681	CJST_CJ_1166_1198_R	1022	gltA
1047	CJST_CJ_584_616_F	315	CJST_CJ_663_692_R	1379	glnA
1048	CJST_CJ_360_394_F	346	CJST_CJ_442_476_R	955	aspA
1049	CJST_CJ_2636_2668_F	504	CJST_CJ_2753_2777_R	1409	tkt
1054	CJST_CJ_2060_2090_F	323	CJST_CJ_2148_2174_R	1068	pgm
1064	CJST_CJ_1680_1713_F	479	CJST_CJ_1795_1822_R	938	glyA

[410] The primers were used to amplify nucleic acid from 50 food product samples provided by the USDA, 25 of which contained *Campylobacter jejuni* and 25 of which contained *Campylobacter coli*. Primers used in this study were developed primarily for the discrimination of *Campylobacter jejuni* clonal complexes and for distinguishing *Campylobacter jejuni* from *Campylobacter coli*. Finer discrimination between *Campylobacter coli* types is also possible by using specific primers targeted to loci where closely-related *Campylobacter coli* isolates demonstrate polymorphisms between strains. The conclusions of the comparison of base composition analysis with sequence analysis are shown in Tables 13A-C.

Table 13A – Results of Base Composition Analysis of 50 *Campylobacter* Samples with Drill-down
MLST Primer Pair Nos: 1048 and 1047

Group	Species	Isolate origin	MLST type or Clonal Complex by Base Composition analysis	MLST Type or Clonal Complex by Sequence analysis	Strain	Base Composition of Bioagent Identifying Amplicon Obtained with Primer Pair No: 1048 (ampH)	Base Composition of Bioagent Identifying Amplicon Obtained with Primer Pair No: 1047 (glnA)
J-1	<i>C. jejuni</i>	Goose	ST 690 /692/707/991	ST 991	RM3673	A30 G25 C16 T46	A47 G21 C16 T25
J-2	<i>C. jejuni</i>	Human	Complex 296/48/353	ST 356, complex 353	RM4192	A30 G25 C16 T46	A48 G21 C17 T23
J-3	<i>C. jejuni</i>	Human	Complex 354/179	ST 436	RM4194	A30 G25 C15 T47	A48 G21 C18 T22
J-4	<i>C. jejuni</i>	Human	Complex 257	ST 257, complex 257	RM4197	A30 G25 C16 T46	A48 G21 C18 T22
J-5	<i>C. jejuni</i>	Human	Complex 52	ST 52, complex 52	RM4277	A30 G25 C16 T46	A48 G21 C17 T23
J-6	<i>C. jejuni</i>	Human	Complex 443	ST 51, complex 443	RM4275 RM4279	A30 G25 C15 T47 A30 G25 C15 T47	A48 G21 C17 T23 A48 G21 C17 T23
J-7	<i>C. jejuni</i>	Human	Complex 42	ST 604, complex 42	RM1064	A30 G25 C15 T47	A48 G21 C18 T22
J-8	<i>C. jejuni</i>	Human	Complex 42/49/362	ST 362, complex 362	RM3193	A30 G25 C15 T47	A48 G21 C18 T22
J-9	<i>C. jejuni</i>	Human	Complex 45/203	ST 147, Complex 45	RM3203	A30 G25 C15 T47	A47 G21 C18 T23
C-1	<i>C. jejuni</i>	Human	Consistent with 74 closely related sequence types (none belong to a clonal complex)	ST 828	RM4103	A31 G27 C20 T39	A48 G21 C16 T24
	ST 832			RM1169	A31 G27 C20 T39	A48 G21 C16 T24	
	ST 1056			RM1057	A31 G27 C20 T39	A48 G21 C16 T24	
	ST 889			RM1166	A31 G27 C20 T39	A48 G21 C16 T24	
	ST 829			RM1182	A31 G27 C20 T39	A48 G21 C16 T24	
	ST 1050			RM1518	A31 G27 C20 T39	A48 G21 C16 T24	
	ST 1051			RM1521	A31 G27 C20 T39	A48 G21 C16 T24	
	ST 1053			RM1523	A31 G27 C20 T39	A48 G21 C16 T24	
	ST 1055			RM1527	A31 G27 C20 T39	A48 G21 C16 T24	
	ST 1017			RM1529	A31 G27 C20 T39	A48 G21 C16 T24	
	Poultry	ST 860		RM1840	A31 G27 C20 T39	A48 G21 C16 T24	
		ST 1063		RM2219	A31 G27 C20 T39	A48 G21 C16 T24	
		ST 1066		RM2241	A31 G27 C20 T39	A48 G21 C16 T24	
		ST 1067		RM2243	A31 G27 C20 T39	A48 G21 C16 T24	
		ST 1068		RM2439	A31 G27 C20 T39	A48 G21 C16 T24	
		ST 1016		RM3230	A31 G27 C20 T39	A48 G21 C16 T24	
		ST 1069		RM3231	A31 G27 C20 T39	A48 G21 C16 T24	
		ST 1061		RM1904	A31 G27 C20 T39	A48 G21 C16 T24	
		Unknown		ST 825	RM1534	A31 G27 C20 T39	A48 G21 C16 T24
				ST 901	RM1505	A31 G27 C20 T39	A48 G21 C16 T24
C-2	<i>C. coli</i>	Human	ST 895	ST 895	RM1532	A31 G27 C19 T40	A48 G21 C16 T24
C-3	<i>C. coli</i>	Poultry	Consistent	ST 1064	RM2223	A31 G27 C20 T39	A48 G21 C16 T24

			with 63 closely related sequence types (none belong to a clonal complex)	ST 1082	RM1178	A31 G27 C20 T39	A48 G21 C16 T24
				ST 1054	RM1525	A31 G27 C20 T39	A48 G21 C16 T24
				ST 1049	RM1517	A31 G27 C20 T39	A48 G21 C16 T24
		Marmoset		ST 891	RM1531	A31 G27 C20 T39	A48 G21 C16 T24

Table 13B – Results of Base Composition Analysis of 50 *Campylobacter* Samples with Drill-down MLST Primer Pair Nos: 1053 and 1064

Group	Species	Isolate origin	MLST type or Clonal Complex by Base Composition analysis	MLST Type or Clonal Complex by Sequence analysis	Strain	Base Composition of Bioagent Identifying Amplicon Obtained with Primer Pair No: 1053 (gltA)	Base Composition of Bioagent Identifying Amplicon Obtained with Primer Pair No: 1064 (glyA)
J-1	<i>C. jejuni</i>	Goose	ST 690 /692/707/991	ST 991	RM3673	A24 G25 C23 T47	A40 G29 C29 T45
J-2	<i>C. jejuni</i>	Human	Complex 206/48/353	ST 356, complex 353	RM4192	A24 G25 C23 T47	A40 G29 C29 T45
J-3	<i>C. jejuni</i>	Human	Complex 354/172	ST 436	RM4194	A24 G25 C23 T47	A40 G29 C29 T45
J-4	<i>C. jejuni</i>	Human	Complex 257	ST 257, complex 257	RM4197	A24 G25 C23 T47	A40 G29 C29 T45
J-5	<i>C. jejuni</i>	Human	Complex 52	ST 52, complex 52	RM4277	A24 G25 C23 T47	A39 G30 C26 T48
J-6	<i>C. jejuni</i>	Human	Complex 443	ST 51, complex 443	RM4275 RM4279	A24 G25 C23 T47	A39 G30 C28 T46 A39 G30 C28 T46
J-7	<i>C. jejuni</i>	Human	Complex 42	ST 604, complex 42	RM1864	A24 G25 C23 T47	A39 G30 C26 T48
J-8	<i>C. jejuni</i>	Human	Complex 42/49/362	ST 362, complex 362	RM3193	A24 G25 C23 T47	A38 G31 C28 T46
J-9	<i>C. jejuni</i>	Human	Complex 45/283	ST 147, Complex 45	RM3203	A24 G25 C23 T47	A38 G31 C28 T46
C-1	<i>C. jejuni</i>	Human	Consistent with 74 closely related sequence types (none belong to a clonal complex)	ST 828	RM4183	A23 G24 C26 T46	A39 G30 C27 T47
	ST 832			RM1169	A23 G24 C26 T46	A39 G30 C27 T47	
	ST 1055			RM1857	A23 G24 C26 T46	A39 G30 C27 T47	
	ST 889			RM1166	A23 G24 C26 T46	A39 G30 C27 T47	
	ST 829			RM1182	A23 G24 C26 T46	A39 G30 C27 T47	
	ST 1050			RM1518	A23 G24 C26 T46	A39 G30 C27 T47	
	ST 1051			RM1521	A23 G24 C26 T46	A39 G30 C27 T47	
	ST 1053			RM1523	A23 G24 C26 T46	A39 G30 C27 T47	
	ST 1055			RM1527	A23 G24 C26 T46	A39 G30 C27 T47	
	ST 1017			RM1529	A23 G24 C26 T46	A39 G30 C27 T47	
	Poultry			ST 860	RM1840	A23 G24 C26 T46	A39 G30 C27 T47
				ST 1063	RM2219	A23 G24 C26 T46	A39 G30 C27 T47
				ST 1066	RM2241	A23 G24 C26 T46	A39 G30 C27 T47
				ST 1067	RM2243	A23 G24 C26 T46	A39 G30 C27 T47
				ST 1068	RM2439	A23 G24 C26 T46	A39 G30 C27 T47
				ST 1016	RM3230	A23 G24 C26 T46	A39 G30 C27 T47
	Swine						

		Unknown		ST 1069	RM3231	A23 G24 C26 T46	NO DATA
				ST 1061	RM1904	A23 G24 C26 T46	A39 G30 C27 T47
				ST 825	RM1534	A23 G24 C26 T46	A39 G30 C27 T47
				ST 901	RM1505	A23 G24 C26 T46	A39 G30 C27 T47
C-2	C. coli	Human	ST 895	ST 895	RM1532	A23 G24 C26 T46	A39 G30 C27 T47
C-3	C. coli	Poultry	Consistent with 63 closely related sequence types (none belong to a clonal complex)	ST 1064	RM2223	A23 G24 C26 T46	A39 G30 C27 T47
				ST 1082	RM1178	A23 G24 C26 T46	A39 G30 C27 T47
				ST 1054	RM1525	A23 G24 C25 T47	A39 G30 C27 T47
				ST 1049	RM1517	A23 G24 C26 T46	A39 G30 C27 T47
		Marmoset		ST 891	RM1531	A23 G24 C26 T46	A39 G30 C27 T47

Table 13C – Results of Base Composition Analysis of 50 *Campylobacter* Samples with Drill-down MLST Primer Pair Nos: 1054 and 1049

Group	Species	Isolate origin	MLST type or Clonal Complex by Base Composition analysis	MLST Type or Clonal Complex by Sequence analysis	Strain	Base Composition of Bioagent Identifying Amplicon Obtained with Primer Pair No: 1054 (pgm)	Base Composition of Bioagent Identifying Amplicon Obtained with Primer Pair No: 1049 (tkbt)
J-1	C. jejuni	Goose	ST 690 /692/707/991	ST 991	RM3673	A26 G33 C19 T38	A41 G28 C35 T38
J-2	C. jejuni	Human	Complex 206/48/353	ST 356, complex 353	RM4192	A26 G33 C19 T37	A41 G28 C36 T37
J-3	C. jejuni	Human	Complex 354/179	ST 436	RM4194	A27 G32 C19 T37	A42 G28 C36 T36
J-4	C. jejuni	Human	Complex 257	ST 257, complex 257	RM4197	A27 G32 C19 T37	A41 G29 C35 T37
J-5	C. jejuni	Human	Complex 52	ST 52, complex 52	RM4277	A26 G33 C18 T38	A41 G28 C36 T37
J-6	C. jejuni	Human	Complex 443	ST 51, complex 443	RM4275	A27 G31 C19 T38	A41 G28 C36 T37
					RM4279	A27 G31 C19 T38	A41 G28 C36 T37
J-7	C. jejuni	Human	Complex 42	ST 604, complex 42	RM1864	A27 G32 C19 T37	A42 G28 C35 T37
J-8	C. jejuni	Human	Complex 42/49/362	ST 362, complex 362	RM3193	A26 G33 C19 T37	A42 G28 C35 T37
J-9	C. jejuni	Human	Complex 45/283	ST 147, Complex 45	RM3203	A28 G31 C19 T37	A43 G28 C36 T35
	C. jejuni	Human <					

				ST 860	RM1840	A27 G30 C19 T39	A46 G28 C32 T36
				ST 1063	RM2219	A27 G30 C19 T39	A46 G28 C32 T36
				ST 1066	RM2241	A27 G30 C19 T39	A46 G28 C32 T36
				ST 1067	RM2243	A27 G30 C19 T39	A46 G28 C32 T36
				ST 1068	RM2439	A27 G30 C19 T39	A46 G28 C32 T36
				ST 1016	RM3230	A27 G30 C19 T39	A46 G28 C32 T36
				ST 1069	RM3231	A27 G30 C19 T39	A46 G28 C32 T36
				ST 1061	RM1904	A27 G30 C19 T39	A46 G28 C32 T36
				ST 825	RM1534	A27 G30 C19 T39	A46 G28 C32 T36
				ST 901	RM1505	A27 G30 C19 T39	A46 G28 C32 T36
				ST 895	RM1532	A27 G30 C19 T39	A46 G29 C32 T36
				ST 1064	RM2223	A27 G30 C19 T39	A46 G29 C32 T36
ST 1082	RM1178	A27 G30 C19 T39	A46 G29 C32 T36				
ST 1054	RM1525	A27 G30 C19 T39	A46 G29 C32 T36				
ST 1049	RM1517	A27 G30 C19 T39	A46 G29 C32 T36				
ST 891	RM1531	A27 G30 C19 T39	A46 G29 C32 T36				
C-2	C. coli	Human	ST 895	RM1532	A27 G30 C19 T39	A46 G29 C32 T36	
C-3	C. coli	Poultry	Consistent with 63 closely related sequence types (none belong to a clonal complex)	ST 1064	RM2223	A27 G30 C19 T39	A46 G29 C32 T36
				ST 1082	RM1178	A27 G30 C19 T39	A46 G29 C32 T36
				ST 1054	RM1525	A27 G30 C19 T39	A46 G29 C32 T36
				ST 1049	RM1517	A27 G30 C19 T39	A46 G29 C32 T36
		Marmoset		ST 891	RM1531	A27 G30 C19 T39	A46 G29 C32 T36

[411] The base composition analysis method was successful in identification of 12 different strain groups. *Campylobacter jejuni* and *Campylobacter coli* are generally differentiated by all loci. Ten clearly differentiated *Campylobacter jejuni* isolates and 2 major *Campylobacter coli* groups were identified even though the primers were designed for strain typing of *Campylobacter jejuni*. One isolate (RM4183) which was designated as *Campylobacter jejuni* was found to group with *Campylobacter coli* and also appears to actually be *Campylobacter coli* by full MLST sequencing.

Example 11: Identification of *Acinetobacter baumannii* Using Broad Range Survey and Division-Wide Primers in Epidemiological Surveillance

[412] To test the capability of the broad range survey and division-wide primer sets of Table 5 in identification of *Acinetobacter* species, 183 clinical samples were obtained from individuals participating in, or in contact with individuals participating in Operation Iraqi Freedom (including US service personnel, US civilian patients at the Walter Reed Army Institute of Research (WRAIR), medical staff, Iraqi civilians and enemy prisoners. In addition, 34 environmental samples were obtained from hospitals in Iraq, Kuwait, Germany, the United States and the USNS Comfort, a hospital ship.

[413] Upon amplification of nucleic acid obtained from the clinical samples, primer pairs 346-349, 360, 361, 354, 362 and 363 (Table 5) all produced bacterial bioagent amplicons which identified *Acinetobacter baumannii* in 215 of 217 samples. The organism *Klebsiella pneumoniae* was identified in the remaining two samples. In addition, 14 different strain types (containing single nucleotide polymorphisms relative to a reference strain of *Acinetobacter baumannii*) were identified and assigned

arbitrary numbers from 1 to 14. Strain type 1 was found in 134 of the sample isolates and strains 3 and 7 were found in 46 and 9 of the isolates respectively.

[414] The epidemiology of strain type 7 of *Acinetobacter baumannii* was investigated. Strain 7 was found in 4 patients and 5 environmental samples (from field hospitals in Iraq and Kuwait). The index patient infected with strain 7 was a pre-war patient who had a traumatic amputation in March of 2003 and was treated at a Kuwaiti hospital. The patient was subsequently transferred to a hospital in Germany and then to WRAIR. Two other patients from Kuwait infected with strain 7 were found to be non-infectious and were not further monitored. The fourth patient was diagnosed with a strain 7 infection in September of 2003 at WRAIR. Since the fourth patient was not related involved in Operation Iraqi Freedom, it was inferred that the fourth patient was the subject of a nosocomial infection acquired at WRAIR as a result of the spread of strain 7 from the index patient.

[415] The epidemiology of strain type 3 of *Acinetobacter baumannii* was also investigated. Strain type 3 was found in 46 samples, all of which were from patients (US service members, Iraqi civilians and enemy prisoners) who were treated on the USNS Comfort hospital ship and subsequently returned to Iraq or Kuwait. The occurrence of strain type 3 in a single locale may provide evidence that at least some of the infections at that locale were a result of nosocomial infections.

[416] This example thus illustrates an embodiment of the present invention wherein the methods of analysis of bacterial bioagent identifying amplicons provide the means for epidemiological surveillance.

Example 12: Selection and Use of Triangulation Genotyping Analysis Primer Pairs for *Acinetobacter baumannii*

[417] To combine the power of high-throughput mass spectrometric analysis of bioagent identifying amplicons with the sub-species characteristic resolving power provided by triangulation genotyping analysis, an additional 21 primer pairs were selected based on analysis of housekeeping genes of the genus *Acinetobacter*. Genes to which the drill-down triangulation genotyping analysis primers hybridize for production of bacterial bioagent identifying amplicons include anthranilate synthase component I (trpE), adenylate kinase (adk), adenine glycosylase (mutY), fumarate hydratase (fumC), and pyrophosphate phospho-hydratase (ppa). These 21 primer pairs are indicated with reference to sequence listings in Table 14. Primer pair numbers 1151-1154 hybridize to and amplify segments of trpE. Primer pair numbers 1155-1157 hybridize to and amplify segments of adk. Primer pair numbers 1158-1164 hybridize to and amplify segments of mutY. Primer pair numbers 1165-1170 hybridize to and amplify segments of fumC. Primer pair number 1171 hybridizes to and amplifies a segment of ppa.

Primer pair numbers: 2846-2848 hybridize to and amplify segments of the *parC* gene of DNA topoisomerase which include a codon known to confer quinolone drug resistance upon sub-types of *Acinetobacter baumannii*. Primer pair numbers 2852-2854 hybridize to and amplify segments of the *gyrA* gene of DNA gyrase which include a codon known to confer quinolone drug resistance upon sub-types of *Acinetobacter baumannii*. Primer pair numbers 2922 and 2972 are speciating primers which are useful for identifying different species members of the genus *Acinetobacter*. The primer names given in Table 14A (with the exception of primer pair numbers 2846-2848, 2852-2854) indicate the coordinates to which the primers hybridize to a reference sequence which comprises a concatenation of the genes *TrpE*, *efp* (elongation factor p), *adk*, *mutT*, *fumC*, and *ppa*. For example, the forward primer of primer pair 1151 is named AB_MLST-11-OIF007_62_91_F because it hybridizes to the *Acinetobacter* primer reference sequence of strain type 11 in sample 007 of Operation Iraqi Freedom (OIF) at positions 62 to 91. DNA was sequenced from strain type 11 and from this sequence data and an artificial concatenated sequence of partial gene extractions was assembled for use in design of the triangulation genotyping analysis primers. The stretches of arbitrary residues "N"s in the concatenated sequence were added for the convenience of separation of the partial gene extractions (40N for AB_MLST (SEQ ID NO: 1444)).

[418] The hybridization coordinates of primer pair numbers 2846-2848 are with respect to GenBank Accession number X95819. The hybridization coordinates of primer pair numbers 2852-2854 are with respect to GenBank Accession number AY642140. Sequence residue "T" appearing in the forward and reverse primers of primer pair number 2972 represents inosine.

Table 14A: Triangulation Genotyping Analysis Primer Pairs for Identification of Sub-species characteristics (Strain Type) of Members of the Bacterial Genus *Acinetobacter*

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO:)	Reverse Primer Name	Reverse Primer (SEQ ID NO:)
1151	AB_MLST-11-OIF007_62_91_F	454	AB_MLST-11-OIF007_169_203_R	1418
1152	AB_MLST-11-OIF007_185_214_F	243	AB_MLST-11-OIF007_291_324_R	959
1153	AB_MLST-11-OIF007_260_289_F	541	AB_MLST-11-OIF007_364_393_R	1400
1154	AB_MLST-11-OIF007_206_239_F	436	AB_MLST-11-OIF007_318_344_R	1016
1155	AB_MLST-11-OIF007_522_552_F	378	AB_MLST-11-OIF007_587_610_R	1392
1156	AB_MLST-11-OIF007_547_571_F	250	AB_MLST-11-OIF007_656_686_R	902
1157	AB_MLST-11-OIF007_601_627_F	256	AB_MLST-11-OIF007_710_736_R	881
1158	AB_MLST-11-OIF007_1202_1225_F	384	AB_MLST-11-OIF007_1266_1296_R	878
1159	AB_MLST-11-OIF007_1202_1225_F	384	AB_MLST-11-OIF007_1299_1316_R	1199
1160	AB_MLST-11-OIF007_1234_1264_F	694	AB_MLST-11-OIF007_1335_1362_R	1215

1161	AB MLST-11-OIF007 1327 1356 F	225	AB MLST-11-OIF007 1422 1448 R	1212
1162	AB MLST-11-OIF007 1345 1369 F	383	AB MLST-11-OIF007 1470 1494 R	1083
1163	AB MLST-11-OIF007 1351 1375 F	662	AB MLST-11-OIF007 1470 1494 R	1083
1164	AB MLST-11-OIF007 1387 1412 F	422	AB MLST-11-OIF007 1470 1494 R	1083
1165	AB MLST-11-OIF007 1542 1569 F	194	AB MLST-11-OIF007 1656 1680 R	1173
1166	AB MLST-11-OIF007 1566 1593 F	684	AB MLST-11-OIF007 1656 1680 R	1173
1167	AB MLST-11-OIF007 1611 1638 F	375	AB MLST-11-OIF007 1731 1757 R	890
1168	AB MLST-11-OIF007 1726 1752 F	182	AB MLST-11-OIF007 1790 1821 R	1195
1169	AB MLST-11-OIF007 1792 1826 F	656	AB MLST-11-OIF007 1876 1909 R	1151
1170	AB MLST-11-OIF007 1792 1826 F	656	AB MLST-11-OIF007 1895 1927 R	1224
1171	AB MLST-11-OIF007 1970 2002 F	618	AB MLST-11-OIF007 2097 2118 R	1157
2846	PARC X95819 33 58 F	302	PARC X95819 121 153 R	852
2847	PARC X95819 33 58 F	199	PARC X95819 157 178 R	889
2848	PARC X95819 33 58 F	596	PARC X95819 97 128 R	1169
2852	GYRA AY642140 -1 24 F	150	GYRA AY642140 71 100 R	1242
2853	GYRA AY642140 25 54 F	166	GYRA AY642140 121 146 R	1069
2854	GYRA AY642140 26 54 F	166	GYRA AY642140 58 89 R	1168
2922	AB MLST-11-OIF007 991 1018 F	583	AB MLST-11-OIF007 1110 1137 R	923
2972	AB MLST-11-OIF007 1007 1034 F	592	AB MLST-11-OIF007 1126 1153 R	924

Table 14B: Triangulation Genotyping Analysis Primer Pairs for Identification of Sub-species characteristics (Strain Type) of Members of the Bacterial Genus *Acinetobacter*

Primer Pair No.	Forward Primer (SEQ ID NO:)	SEQUENCE	Reverse Primer (SEQ ID NO:)	SEQUENCE
1151	454	TGAGATTCTGAACATTTAATGCTGATTGA	1418	TTGTACATTTGAACAATATGCATGACATGTGAAT
1152	243	TATGTGTTCAAACTGTACAAGTGGAAGTGCG	969	TCAACAGGTTCTACATCATCAATATTCATTTGC
1153	541	TGCAGCTTATCAGTGTCCCGAATAATTCG	1400	TTGCAATGACATATTCATTCACCATGCG
1154	436	TGAAGTGGCTGATGATATCATGCACTTGATGTA	1036	TCCGCCAAAACCTCCCTTTTCCACAG
1155	378	TGCGTTTAGTAAAGAACGATATTGCTCAACC	1392	TTCTGCTTGAGGATAATGTGCTGCG
1156	250	TCAACTGATGTCGCTGAATGGTGT	903	TACGTTCTACATTTCTCATCAAGTACATC
1157	256	TCAACGAGAGCTTTGGAAAGAGAGG	881	TACAACTGATAAACAGGACAGAAAGC
1158	384	TGCTGCCCAATTTGCATAAAGC	878	TAAAGCCGGGTAGTGCAATCATTTCTTAG
1159	384	TGCTGCCCAATTTGCATAAAGC	1199	TGCACCTGCGTTGAGCG
1160	694	TTTATGACAGCAAGGCAATTTCTGAAAC	1215	TGCCATCATCAATCAAGCCATATGACG
1161	225	TAGGTTTAACTCAGTATGGCGTATTATGG	1212	TGCCAGTTTCCATATTCACGTTCTGG
1162	383	TGCTGATATGGAATGCAACGTGAA	1083	TGCTTGAGTGTAGTCATGATTGCG

1163	662	TTATGAGTGGCAAGTGAAGACCGT	1083	TGCGTTGAGTGTATCATGATTCGG
1164	422	TCCTTGCCATGAGATGACTTAAGC	1083	TGCGTTGAGTGTATCATGATTCGG
1165	194	TACTAGCGGTAACTTAACAAGATTGC	1173	TGAGTCGCGTTCACTTTAAGTGGCA
1166	684	TTGCGATGATATTGCTTGGTAGCAG	1173	TGAGTCGCGTTCACTTTAAGTGGCA
1167	375	TGGCGAAATCGTATTCCTGAAGAATGA	890	TACCGAAGCAGCAGCAGCATTAAAG
1168	182	TACCACTATTAAATGTCGCTGCTTC	1195	TGCACTGAAATGAGATTGAGTAAGTATAAGC
1169	656	TTATACTTACTGCAATCTATTCAGTTGCTGCTG	1151	TGATTTATGCAAGATGATCAATTTCTCAGCA
1170	656	TTATACTTACTGCAATCTATTCAGTTGCTGCTG	1224	TGCGTTAATTAACATAGAGAAATATGCAAGA
1171	618	TGGTTATGTACCAATACCTTGTCTGAAGATGG	1157	TGAGCGCATGATACCAAGCTC
2846	302	TCCAAAAAATCAGCGTACAGTGG	852	TAAAGCATAGCGGTAACTAAATGGCTGAGCCAT
2847	199	TACTTGGTAAATACCAACCAATGCTGA	883	TACCCGATTTCCCTTGACCTTC
2848	596	TGGTAAATACCAACCAATGCTGAC	1169	TGAGCCATGATACCATGGCTTCATACATGC
2852	190	TAAATCTGCGGTGTGCTGCTGAC	1242	TGCTAAAGTCTTGAGCCATACGAAATAG
2853	166	TAACTCGTAAATATCACCCGATGCTGAC	1069	TGCTATGACGAGATTAACCTGACC
2854	166	TAACTCGTAAATATCACCCGATGCTGAC	1168	TGAGCCATACGATCAATGCTTATATAACAGC
2922	583	TGGCGATGCTGCGAAATGTTAAAGA	923	TAGTATCACCAAGTACACCCGATCAGT
2972	592	TGGCGATGCTGCGAAATGTTAAAGA	924	TAGTATCACCAAGTACACCCGATCAGT

[419] Analysis of bioagent identifying amplicons obtained using the primers of Table 14B for over 200 samples from Operation Iraqi Freedom resulted in the identification of 50 distinct strain type clusters. The largest cluster, designated strain type 11 (ST11) includes 42 sample isolates, all of which were obtained from US service personnel and Iraqi civilians treated at the 28th Combat Support Hospital in Baghdad. Several of these individuals were also treated on the hospital ship USNS Comfort. These observations are indicative of significant epidemiological correlation/linkage.

[420] All of the sample isolates were tested against a broad panel of antibiotics to characterize their antibiotic resistance profiles. As an example of a representative result from antibiotic susceptibility testing, ST11 was found to consist of four different clusters of isolates, each with a varying degree of sensitivity/resistance to the various antibiotics tested which included penicillins, extended spectrum penicillins, cephalosporins, carbapenem, protein synthesis inhibitors, nucleic acid synthesis inhibitors, anti-metabolites, and anti-cell membrane antibiotics. Thus, the genotyping power of bacterial bioagent identifying amplicons, particularly drill-down bacterial bioagent identifying amplicons, has the potential to increase the understanding of the transmission of infections in combat casualties, to identify the source of infection in the environment, to track hospital transmission of nosocomial infections, and to

rapidly characterize drug-resistance profiles which enable development of effective infection control measures on a time-scale previously not achievable.

Example 13: Triangulation Genotyping Analysis and Codon Analysis of *Acinetobacter baumannii* Samples from Two Health Care Facilities

[421] In this investigation, 88 clinical samples were obtained from Walter Reed Hospital and 95 clinical samples were obtained from Northwestern Medical Center. All samples from both healthcare facilities were suspected of containing sub-types of *Acinetobacter baumannii*, at least some of which were expected to be resistant to quinolone drugs. Each of the 183 samples was analyzed by the method of the present invention. DNA was extracted from each of the samples and amplified with eight triangulation genotyping analysis primer pairs represented by primer pair numbers: 1151, 1156, 1158, 1160, 1165, 1167, 1170, and 1171. The DNA was also amplified with speciating primer pair number 2922 and codon analysis primer pair numbers 2846-2848 which interrogate a codon present in the *parC* gene, and primer pair numbers 2852-2854 which bracket a codon present in the *gyrA* gene. The *parC* and *gyrA* codon mutations are both responsible for causing drug resistance in *Acinetobacter baumannii*. During evolution of drug resistant strains, the *gyrA* mutation usually occurs before the *parC* mutation. Amplification products were measured by ESI-TOF mass spectrometry as indicated in Example 4. The base compositions of the amplification products were calculated from the average molecular masses of the amplification products and are shown in Tables 15-18. The entries in each of the tables are grouped according to strain type number, which is an arbitrary number assigned to *Acinetobacter baumannii* strains in the order of observance beginning from the triangulation genotyping analysis OIF genotyping study described in Example 12. For example, strain type 11 which appears in samples from the Walter Reed Hospital is the same strain as the strain type 11 mentioned in Example 12. Ibis# refers to the order in which each sample was analyzed. Isolate refers to the original sample isolate numbering system used at the location from which the samples were obtained (either Walter Reed Hospital or Northwestern Medical Center). ST = strain type. ND = not detected. Base compositions highlighted with **bold** type indicate that the base composition is a unique base composition for the amplification product obtained with the pair of primers indicated.

Table 15A: Base Compositions of Amplification Products of 88 *A. baumannii* Samples Obtained from Walter Reed Hospital and Amplified with Codon Analysis Primer Pairs Targeting the *gyrA*

Gene

Species	Ibis#	Isolate	ST	PP No: 2852 <i>gyrA</i>	PP No: 2853 <i>gyrA</i>	PP No: 2854 <i>gyrA</i>
<i>A. baumannii</i>	20	1082	1	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. baumannii</i>	13	854	10	A25G23C21T32	A29G28C21T43	A17G13C13T21

A. baumannii	22	1162	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	27	1230	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	31	1367	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	37	1459	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	55	1700	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	64	1777	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	73	1861	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	74	1877	10	ND	A29G28C21T43	A17G13C13T21
A. baumannii	86	1972	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	3	684	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	6	720	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	7	726	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	19	1079	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	21	1123	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	23	1188	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	33	1417	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	34	1431	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	38	1496	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	40	1523	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	42	1640	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	50	1666	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	51	1668	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	52	1695	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	65	1781	11	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	44	1649	12	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	49A	1658.1	12	A25G23C22T31	A29G28C21T43	A17G13C13T21
A. baumannii	49B	1658.2	12	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	56	1707	12	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	80	1893	12	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	5	693	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	8	749	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	10	839	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	14	865	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	16	888	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	29	1326	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	35	1440	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	41	1524	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	46	1652	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	47	1653	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	48	1657	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	57	1709	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	61	1727	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	63	1762	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	67	1806	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	75	1881	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	77	1886	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	1	649	46	A25G23C21T32	A29G28C21T43	A17G13C13T21

<i>A. baumannii</i>	2	653	46	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	39	1497	16	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	24	1198	15	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	28	1243	15	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	43	1648	15	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	62	1746	15	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	4	689	15	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	68	1822	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	69	1823A	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	70	1823B	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	71	1826	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	72	1860	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	81	1924	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	82	1929	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	85	1966	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	11	841	3	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. baumannii</i>	32	1415	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	45	1651	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	54	1697	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	58	1712	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	60	1725	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	66	1802	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	76	1883	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	78	1891	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	79	1892	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	83	1947	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	84	1964	24	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	53	1696	24	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. baumannii</i>	36	1458	49	A25G23C21T32	A29G28C21T43	A17G13C13T21
<i>A. baumannii</i>	59	1716	9	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. baumannii</i>	9	805	30	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. baumannii</i>	18	967	39	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. baumannii</i>	30	1322	48	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. baumannii</i>	26	1218	50	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. sp. 13TU</i>	15	875	A1	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. sp. 13TU</i>	17	895	A1	A25G23C22T31	A29G28C22T42	A17G13C14T20
<i>A. sp. 3</i>	12	853	B7	A25G22C22T32	A30G29C22T40	A17G13C14T20
<i>A. johnsonii</i>	25	1202	NEW1	A25G22C22T32	A30G29C22T40	A17G13C14T20
<i>A. sp. 2082</i>	87	2082	NEW2	A25G22C22T32	A31G28C22T40	A17G13C14T20

Table 15B: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Walter Reed Hospital and Amplified with Codon Analysis Primer Pairs Targeting the parC Gene

Species	Ibis#	Isolate	ST	PP No: 2846 parC	PP No: 2847 parC	PP No: 2848 parC
<i>A. baumannii</i>	20	1082	1	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	13	854	10	A33G26C28T34	A29G28C25T32	A16G14C14T16

A. baumannii	22	1162	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	27	1230	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	31	1367	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	37	1459	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	55	1700	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	64	1777	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	73	1861	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	74	1877	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	86	1972	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	3	684	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	6	720	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	7	726	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	19	1079	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	21	1123	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	23	1188	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	33	1417	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	34	1431	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	38	1496	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	40	1523	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	42	1640	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	50	1666	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	51	1668	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	52	1695	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	65	1781	11	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	44	1649	12	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	49A	1658.1	12	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	49B	1658.2	12	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	56	1707	12	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	80	1893	12	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	5	693	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	8	749	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	10	839	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	14	865	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	16	888	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	29	1326	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	35	1440	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	41	1524	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	46	1652	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	47	1653	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	48	1657	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	57	1709	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	61	1727	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	63	1762	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	67	1806	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	75	1881	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	77	1886	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	1	649	46	A33G26C28T34	A29G28C25T32	A16G14C14T16

<i>A. baumannii</i>	2	653	46	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	39	1497	16	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	24	1198	15	A33G26C28T34	A29G29C23T33	A16G14C14T16
<i>A. baumannii</i>	28	1243	15	A33G26C28T34	A29G29C23T33	A16G14C14T16
<i>A. baumannii</i>	43	1648	15	A33G26C28T34	A29G29C23T33	A16G14C14T16
<i>A. baumannii</i>	62	1746	15	A33G26C28T34	A29G29C23T33	A16G14C14T16
<i>A. baumannii</i>	4	689	15	A34G25C29T33	A30G27C26T31	A16G14C15T15
<i>A. baumannii</i>	68	1822	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	69	1823A	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	70	1823B	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	71	1826	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	72	1860	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	81	1924	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	82	1929	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	85	1966	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	11	841	3	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	32	1415	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	45	1651	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	54	1697	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	58	1712	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	60	1725	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	66	1802	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	76	1883	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	78	1891	24	A34G25C29T33	A30G27C26T31	A16G14C15T15
<i>A. baumannii</i>	79	1892	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	83	1947	24	A34G25C29T33	A30G27C26T31	A16G14C15T15
<i>A. baumannii</i>	84	1964	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	53	1696	24	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	36	1458	49	A34G26C29T32	A30G28C24T32	A16G14C15T15
<i>A. baumannii</i>	59	1716	9	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	9	805	30	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	18	967	39	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	30	1322	48	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	26	1218	50	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. sp. 137U</i>	15	875	A1	A32G26C28T35	A28G28C24T34	A16G14C15T15
<i>A. sp. 137U</i>	17	895	A1	A32G26C28T35	A28G28C24T34	A16G14C15T15
<i>A. sp. 3</i>	12	853	B7	A29G26C27T39	A26G32C21T35	A16G14C15T15
<i>A. johnsonii</i>	25	1202	NEW1	A32G28C26T35	A29G29C22T34	A16G14C15T15
<i>A. sp. 2082</i>	87	2082	NEW2	A33G27C26T35	A31G28C20T35	A16G14C15T15

Table 16A: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Northwestern Medical Center and Amplified with Codon Analysis Primer Pairs Targeting the *gyrA* Gene

Species	Ibis#	Isolate	ST	PP No: 2852 <i>gyrA</i>	PP No: 2853 <i>gyrA</i>	PP No: 2854 <i>gyrA</i>
---------	-------	---------	----	----------------------------	----------------------------	----------------------------

A. baumannii	54	536	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	87	665	3	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	8	80	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	9	91	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	10	92	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	11	131	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	12	137	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	21	218	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	26	242	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	94	678	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	1	9	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	2	13	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	3	19	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	4	24	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	5	36	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	6	39	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	13	139	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	15	165	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	16	170	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	17	186	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	20	202	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	22	221	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	24	234	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	25	239	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	33	370	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	34	389	10	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	19	201	14	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	27	257	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	29	301	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	31	354	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	36	422	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	37	424	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	38	434	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	39	473	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	40	482	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	44	512	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	45	516	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	47	522	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	48	526	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	50	528	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	52	531	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	53	533	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	56	542	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	59	550	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	62	556	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	64	557	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	70	588	51	A25G23C21T32	A29G28C21T43	A17G13C13T21

A. baumannii	73	603	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	74	605	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	75	606	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	77	611	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	79	622	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	83	643	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	85	653	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	89	669	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	93	674	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	23	228	51	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	32	369	52	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	35	393	52	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	30	339	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	41	485	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	42	493	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	43	502	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	46	520	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	49	527	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	51	529	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	65	562	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	68	579	53	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	57	546	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	58	548	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	60	552	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	61	555	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	63	557	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	66	570	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	67	578	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	69	584	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	71	593	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	72	602	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	76	609	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	78	621	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	80	625	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	81	628	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	82	632	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	84	649	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	86	655	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	88	668	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	90	671	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	91	672	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	92	673	54	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	18	196	55	A25G23C22T31	A29G28C21T43	A17G13C13T21
A. baumannii	55	537	27	A25G23C21T32	A29G28C21T43	A17G13C13T21
A. baumannii	28	263	27	A25G23C22T31	A29G28C22T42	A17G13C14T20
A. sp. 3	14	164	B7	A25G22C22T32	A30G29C22T40	A17G13C14T20
mixture	7	71	-	ND	ND	A17G13C15T19

Table 16B: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Northwestern Medical Center and Amplified with Codon Analysis Primer Pairs Targeting the *parC* Gene

Species	Ibis#	Isolate	ST	PP No: 2846 parC	PP No: 2847 parC	PP No: 2848 parC
<i>A. baumannii</i>	54	536	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	87	665	3	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	8	80	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	9	91	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	10	92	10	A33G26C28T34	A29G28C25T32	ND
<i>A. baumannii</i>	11	131	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	12	137	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	21	218	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	26	242	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	94	678	10	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	1	9	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	2	13	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	3	19	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	4	24	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	5	36	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	6	39	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	13	139	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	15	165	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	16	170	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	17	186	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	20	202	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	22	221	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	24	234	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	25	239	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	33	370	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	34	389	10	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	19	201	14	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	27	257	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	29	301	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	31	354	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	36	422	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	37	424	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	38	434	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	39	473	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	40	482	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	44	512	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	45	516	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	47	522	51	A33G26C28T34	A29G28C25T32	A16G14C14T16

A. baumannii	48	526	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	50	528	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	52	531	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	53	533	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	56	542	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	59	550	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	62	556	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	64	557	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	70	588	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	73	603	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	74	605	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	75	606	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	77	611	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	79	622	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	83	643	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	85	653	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	89	669	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	93	674	51	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	23	228	51	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	32	369	52	A34G25C28T34	A30G27C25T32	A16G14C14T16
A. baumannii	35	393	52	A34G25C28T34	A30G27C25T32	A16G14C14T16
A. baumannii	30	339	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	41	485	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	42	493	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	43	502	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	46	520	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	49	527	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	51	529	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	65	562	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	68	579	53	A34G25C29T33	A30G27C26T31	A16G14C15T15
A. baumannii	57	546	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	58	548	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	60	552	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	61	555	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	63	557	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	66	570	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	67	578	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	69	584	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	71	593	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	72	602	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	76	609	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	78	621	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	80	625	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	81	628	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	82	632	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	84	649	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
A. baumannii	86	655	54	A33G26C28T34	A29G28C25T32	A16G14C14T16

<i>A. baumannii</i>	88	668	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	90	671	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	91	672	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	92	673	54	A33G26C28T34	A29G28C25T32	A16G14C14T16
<i>A. baumannii</i>	18	196	55	A33G27C28T33	A29G28C25T31	A15G14C15T16
<i>A. baumannii</i>	55	537	27	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. baumannii</i>	28	263	27	A33G26C29T33	A29G28C26T31	A16G14C15T15
<i>A. sp. 3</i>	14	164	B7	A35G25C29T32	A30G28C17T39	A16G14C15T15
mixture	7	71	-	ND	ND	A17G14C15T14

Table 17A: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Walter Reed Hospital and Amplified with Speciating Primer Pair No. 2922 and Triangulation Genotyping Analysis Primer Pair Nos. 1151 and 1156

Species	Ibis#	Isolate	ST	PP No: 2922 efp	PP No: 1151 trpE	PP No: 1156 Adk
<i>A. baumannii</i>	20	1082	1	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	13	854	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	22	1162	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	27	1230	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	31	1367	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	37	1459	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	55	1700	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	64	1777	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	73	1861	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	74	1877	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	86	1972	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	3	684	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	6	720	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	7	726	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	19	1079	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	21	1123	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	23	1188	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	33	1417	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	34	1431	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	38	1496	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	40	1523	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	42	1640	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	50	1666	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	51	1668	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	52	1695	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	65	1781	11	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	44	1649	12	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	49A	1658.1	12	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	49B	1658.2	12	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	56	1707	12	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	80	1893	12	A45G34C25T43	A44G35C21T42	A44G32C26T38

A. baumannii	5	693	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	8	749	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	10	839	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	14	865	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	16	888	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	29	1326	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	35	1440	14	A44G35C25T43	ND	A44G32C27T37
A. baumannii	41	1524	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	46	1652	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	47	1653	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	48	1657	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	57	1709	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	61	1727	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	63	1762	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	67	1806	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	75	1881	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	77	1886	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	1	649	46	A44G35C25T43	A44G35C22T41	A44G32C26T38
A. baumannii	2	653	46	A44G35C25T43	A44G35C22T41	A44G32C26T38
A. baumannii	39	1497	16	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	24	1198	15	A44G35C25T43	A44G35C22T41	A44G32C26T38
A. baumannii	28	1243	15	A44G35C25T43	A44G35C22T41	A44G32C26T38
A. baumannii	43	1648	15	A44G35C25T43	A44G35C22T41	A44G32C26T38
A. baumannii	62	1746	15	A44G35C25T43	A44G35C22T41	A44G32C26T38
A. baumannii	4	689	15	A44G35C25T43	A44G35C22T41	A44G32C26T38
A. baumannii	68	1822	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	69	1823A	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	70	1823B	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	71	1826	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	72	1860	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	81	1924	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	82	1929	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	85	1966	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	11	841	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
A. baumannii	32	1415	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	45	1651	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	54	1697	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	58	1712	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	60	1725	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	66	1802	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	76	1883	24	ND	A43G36C20T43	A44G32C27T37
A. baumannii	78	1891	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	79	1892	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	83	1947	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	84	1964	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	53	1696	24	A44G35C25T43	A43G36C20T43	A44G32C27T37
A. baumannii	36	1458	49	A44G35C25T43	A44G35C22T41	A44G32C27T37

<i>A. baumannii</i>	59	1716	9	A44G35C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	9	805	30	A44G35C25T43	A44G35C19T44	A44G32C27T37
<i>A. baumannii</i>	18	967	39	A45G34C25T43	A44G35C22T41	A44G32C26T38
<i>A. baumannii</i>	30	1322	48	A44G35C25T43	A43G36C20T43	A44G32C27T37
<i>A. baumannii</i>	26	1218	50	A44G35C25T43	A44G35C21T42	A44G32C26T38
<i>A. sp. 13TU</i>	15	875	A1	A47G33C24T43	A46G32C20T44	A44G33C27T36
<i>A. sp. 13TU</i>	17	895	A1	A47G33C24T43	A46G32C20T44	A44G33C27T36
<i>A. sp. 3</i>	12	853	B7	A46G35C24T42	A42G34C20T46	A43G33C24T40
<i>A. johnsonii</i>	25	1202	NEW1	A46G35C23T43	A42G35C21T44	A43G33C23T41
<i>A. sp. 2082</i>	87	2082	NEW2	A46G36C22T43	A42G32C20T48	A42G34C23T41

Table 17B: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Walter Reed Hospital and Amplified with Triangulation Genotyping Analysis Primer Pair Nos.

1158 and 1160 and 1165

Species	Ibis#	Isolate	ST	PP No: 1158 mutY	PP No: 1160 mutY	PP No: 1165 fumC
<i>A. baumannii</i>	20	1082	1	A27G21C25T22	A32G35C29T33	A40G33C30T36
<i>A. baumannii</i>	13	854	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	22	1162	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	27	1230	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	31	1367	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	37	1459	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	55	1700	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	64	1777	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	73	1861	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	74	1877	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	86	1972	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	3	684	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	6	720	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	7	726	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	19	1079	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	21	1123	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	23	1188	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	33	1417	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	34	1431	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	38	1496	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	40	1523	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	42	1640	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	50	1666	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	51	1668	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	52	1695	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	65	1781	11	A27G21C25T22	A32G34C28T35	A40G33C30T36
<i>A. baumannii</i>	44	1649	12	A27G21C26T21	A32G34C29T34	A40G33C30T36
<i>A. baumannii</i>	49A	1658.1	12	A27G21C26T21	A32G34C29T34	A40G33C30T36
<i>A. baumannii</i>	49B	1658.2	12	A27G21C26T21	A32G34C29T34	A40G33C30T36

A. baumannii	56	1707	12	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	80	1893	12	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	5	693	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	8	749	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	10	839	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	14	865	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	16	888	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	29	1326	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	35	1440	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	41	1524	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	46	1652	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	47	1653	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	48	1657	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	57	1709	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	61	1727	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	63	1762	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	67	1806	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	75	1881	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	77	1886	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	1	649	46	A29G19C26T21	A31G35C29T34	A40G33C29T37
A. baumannii	2	653	46	A29G19C26T21	A31G35C29T34	A40G33C29T37
A. baumannii	39	1497	16	A29G19C26T21	A31G35C29T34	A40G34C29T36
A. baumannii	24	1198	15	A29G19C26T21	A31G35C29T34	A40G33C29T37
A. baumannii	28	1243	15	A29G19C26T21	A31G35C29T34	A40G33C29T37
A. baumannii	43	1648	15	A29G19C26T21	A31G35C29T34	A40G33C29T37
A. baumannii	62	1746	15	A29G19C26T21	A31G35C29T34	A40G33C29T37
A. baumannii	4	689	15	A29G19C26T21	A31G35C29T34	A40G33C29T37
A. baumannii	68	1822	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	69	1823A	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	70	1823B	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	71	1826	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	72	1860	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	81	1924	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	82	1929	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	85	1966	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	11	841	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
A. baumannii	32	1415	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	45	1651	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	54	1697	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	58	1712	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	60	1725	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	66	1802	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	76	1883	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	78	1891	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	79	1892	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	83	1947	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	84	1964	24	A27G21C26T21	A32G35C28T34	A40G33C30T36

<i>A. baumannii</i>	53	1696	24	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	36	1458	49	A27G20C27T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	59	1716	9	A27G21C25T22	A32G35C28T34	A39G33C30T37
<i>A. baumannii</i>	9	805	30	A27G21C25T22	A32G35C28T34	A39G33C30T37
<i>A. baumannii</i>	18	967	39	A27G21C26T21	A32G35C28T34	A39G33C30T37
<i>A. baumannii</i>	30	1322	48	A28G21C24T22	A32G35C29T33	A40G33C30T36
<i>A. baumannii</i>	26	1218	50	A27G21C25T22	A31G36C28T34	A40G33C29T37
<i>A. sp.</i> 13TU	15	875	A1	A27G21C25T22	A30G36C26T37	A41G34C28T36
<i>A. sp.</i> 13TU	17	895	A1	A27G21C25T22	A30G36C26T37	A41G34C28T36
<i>A. sp.</i> 3	12	853	B7	A26G23C23T23	A30G36C27T36	A39G37C26T37
<i>A. johnsonii</i>	25	1202	NEW1	A25G23C24T23	A30G35C30T34	A38G37C26T38
<i>A. sp.</i> 2082	87	2082	NEW2	A26G22C24T23	A31G35C28T35	A42G34C27T36

Table 17C: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Walter Reed Hospital and Amplified with Triangulation Genotyping Analysis Primer Pair Nos. 1167 and 1170 and 1171

Species	Ibis#	Isolate	ST	PP No: 1167 fumC	PP No: 1170 fumC	PP No: 1171 ppa
<i>A. baumannii</i>	20	1082	1	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	13	854	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	22	1162	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	27	1230	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	31	1367	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	37	1459	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	55	1700	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	64	1777	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	73	1861	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	74	1877	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	86	1972	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	3	684	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	6	720	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	7	726	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	19	1079	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	21	1123	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	23	1188	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	33	1417	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	34	1431	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	38	1496	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	40	1523	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	42	1640	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	50	1666	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	51	1668	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	52	1695	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	65	1781	11	A41G34C34T38	A38G27C21T50	A35G37C33T44
<i>A. baumannii</i>	44	1649	12	A41G34C34T38	A38G27C21T50	A35G37C33T44

A. baumannii	49A	1658.1	12	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	49B	1658.2	12	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	56	1707	12	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	80	1893	12	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	5	693	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	8	749	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	10	839	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	14	865	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	16	888	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	29	1326	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	35	1440	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	41	1524	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	46	1652	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	47	1653	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	48	1657	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	57	1709	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	61	1727	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	63	1762	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	67	1806	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	75	1881	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	77	1886	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	1	649	46	A41G35C32T39	A37G28C20T51	A35G37C32T45
A. baumannii	2	653	46	A41G35C32T39	A37G28C20T51	A35G37C32T45
A. baumannii	39	1497	16	A41G35C32T39	A37G28C20T51	A35G37C30T47
A. baumannii	24	1198	15	A41G35C32T39	A37G28C20T51	A35G37C30T47
A. baumannii	28	1243	15	A41G35C32T39	A37G28C20T51	A35G37C30T47
A. baumannii	43	1648	15	A41G35C32T39	A37G28C20T51	A35G37C30T47
A. baumannii	62	1746	15	A41G35C32T39	A37G28C20T51	A35G37C30T47
A. baumannii	4	689	15	A41G35C32T39	A37G28C20T51	A35G37C30T47
A. baumannii	68	1822	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	69	1823A	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	70	1823B	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	71	1826	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	72	1860	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	81	1924	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	82	1929	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	85	1966	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	11	841	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	32	1415	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
A. baumannii	45	1651	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
A. baumannii	54	1697	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
A. baumannii	58	1712	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
A. baumannii	60	1725	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
A. baumannii	66	1802	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
A. baumannii	76	1883	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
A. baumannii	78	1891	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
A. baumannii	79	1892	24	A40G35C34T38	A39G26C22T49	A35G37C33T44

<i>A. baumannii</i>	83	1947	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
<i>A. baumannii</i>	84	1964	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
<i>A. baumannii</i>	53	1696	24	A40G35C34T38	A39G26C22T49	A35G37C33T44
<i>A. baumannii</i>	36	1458	49	A40G35C34T38	A39G26C22T49	A35G37C30T47
<i>A. baumannii</i>	59	1716	9	A40G35C32T40	A38G27C20T51	A36G35C31T47
<i>A. baumannii</i>	9	805	30	A40G35C32T40	A38G27C21T50	A35G36C29T49
<i>A. baumannii</i>	18	967	39	A40G35C33T39	A38G27C20T51	A35G37C30T47
<i>A. baumannii</i>	30	1322	48	A40G35C35T37	A38G27C21T50	A35G37C30T47
<i>A. baumannii</i>	26	1218	50	A40G35C34T38	A38G27C21T50	A35G37C33T44
<i>A. sp. 13TU</i>	15	875	A1	A41G39C31T36	A37G26C24T49	A34G38C31T46
<i>A. sp. 13TU</i>	17	895	A1	A41G39C31T36	A37G26C24T49	A34G38C31T46
<i>A. sp. 3</i>	12	853	B7	A43G37C30T37	A36G27C24T49	A34G37C31T47
<i>A. johnsonii</i>	25	1202	NEW1	A42G38C31T36	A40G27C19T50	A35G37C32T45
<i>A. sp. 2082</i>	87	2082	NEW2	A43G37C32T35	A37G26C21T52	A35G38C31T45

Table 18A: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Northwestern Medical Center and Amplified with Speciating Primer Pair No. 2922 and Triangulation Genotyping Analysis Primer Pair Nos. 1151 and 1156

Species	Ibis#	Isolate	ST	PP No: 2922 efp	PP No: 1151 trpE	PP No: 1156 adk
<i>A. baumannii</i>	54	536	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
<i>A. baumannii</i>	87	665	3	A44G35C24T44	A44G35C22T41	A44G32C26T38
<i>A. baumannii</i>	8	80	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	9	91	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	10	92	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	11	131	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	12	137	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	21	218	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	26	242	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	94	678	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	1	9	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	2	13	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	3	19	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	4	24	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	5	36	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	6	39	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	13	139	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	15	165	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	16	170	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	17	186	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	20	202	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	22	221	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	24	234	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	25	239	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
<i>A. baumannii</i>	33	370	10	A45G34C25T43	A44G35C21T42	A44G32C26T38

A. baumannii	34	389	10	A45G34C25T43	A44G35C21T42	A44G32C26T38
A. baumannii	19	201	14	A44G35C25T43	A44G35C22T41	A44G32C27T37
A. baumannii	27	257	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	29	301	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	31	354	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	36	422	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	37	424	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	38	434	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	39	473	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	40	482	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	44	512	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	45	516	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	47	522	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	48	526	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	50	528	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	52	531	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	53	533	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	56	542	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	59	550	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	62	556	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	64	557	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	70	588	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	73	603	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	74	605	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	75	606	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	77	611	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	79	622	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	83	643	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	85	653	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	89	669	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	93	674	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	23	228	51	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	32	369	52	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	35	393	52	A44G35C25T43	A43G36C20T43	A44G32C26T38
A. baumannii	30	339	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	41	485	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	42	493	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	43	502	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	46	520	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	49	527	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	51	529	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	65	562	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	68	579	53	A44G35C25T43	A44G35C19T44	A44G32C27T37
A. baumannii	57	546	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
A. baumannii	58	548	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
A. baumannii	60	552	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
A. baumannii	61	555	54	A44G35C25T43	A44G35C20T43	A44G32C26T38

<i>A. baumannii</i>	63	557	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	66	570	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	67	578	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	69	584	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	71	593	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	72	602	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	76	609	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	78	621	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	80	625	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	81	628	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	82	632	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	84	649	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	86	655	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	88	668	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	90	671	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	91	672	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	92	673	54	A44G35C25T43	A44G35C20T43	A44G32C26T38
<i>A. baumannii</i>	18	196	55	A44G35C25T43	A44G35C20T43	A44G32C27T37
<i>A. baumannii</i>	55	537	27	A44G35C25T43	A44G35C19T44	A44G32C27T37
<i>A. baumannii</i>	28	263	27	A44G35C25T43	A44G35C19T44	A44G32C27T37
<i>A. sp. 3</i>	14	164	B7	A46G35C24T42	A42G34C20T46	A43G33C24T40
mixture	7	71	?	mixture	ND	ND

Table 18B: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Northwestern Medical Center and Amplified with Triangulation Genotyping Analysis Primer Pair Nos. 1158, 1160 and 1165

Species	Ibis#	Isolate	ST	PP No: 1158 mutY	PP No: 1160 mutY	PP No: 1165 fumC
<i>A. baumannii</i>	54	536	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	87	665	3	A27G20C27T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	8	80	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	9	91	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	10	92	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	11	131	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	12	137	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	21	218	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	26	242	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	94	678	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	1	9	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	2	13	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	3	19	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	4	24	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	5	36	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	6	39	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	13	139	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
<i>A. baumannii</i>	15	165	10	A27G21C26T21	A32G35C28T34	A40G33C30T36

A. baumannii	16	170	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	17	186	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	20	202	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	22	221	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	24	234	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	25	239	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	33	370	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	34	389	10	A27G21C26T21	A32G35C28T34	A40G33C30T36
A. baumannii	19	201	14	A27G21C25T22	A31G36C28T34	A40G33C29T37
A. baumannii	27	257	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	29	301	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	31	354	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	36	422	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	37	424	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	38	434	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	39	473	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	40	482	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	44	512	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	45	516	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	47	522	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	48	526	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	50	528	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	52	531	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	53	533	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	56	542	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	59	550	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	62	556	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	64	557	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	70	588	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	73	603	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	74	605	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	75	606	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	77	611	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	79	622	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	83	643	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	85	653	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	89	669	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	93	674	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	23	228	51	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	32	369	52	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	35	393	52	A27G21C25T22	A32G35C28T34	A40G33C29T37
A. baumannii	30	339	53	A28G20C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	41	485	53	A28G20C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	42	493	53	A28G20C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	43	502	53	A28G20C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	46	520	53	A28G20C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	49	527	53	A28G20C26T21	A32G34C29T34	A40G33C30T36

A. baumannii	51	529	53	A28G20C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	65	562	53	A28G20C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	68	579	53	A28G20C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	57	546	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	58	548	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	60	552	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	61	555	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	63	557	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	66	570	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	67	578	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	69	584	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	71	593	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	72	602	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	76	609	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	78	621	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	80	625	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	81	628	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	82	632	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	84	649	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	86	655	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	88	668	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	90	671	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	91	672	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	92	673	54	A27G21C26T21	A32G34C29T34	A40G33C30T36
A. baumannii	18	196	55	A27G21C25T22	A31G36C27T35	A40G33C29T37
A. baumannii	55	537	27	A27G21C25T22	A32G35C28T34	A40G33C30T36
A. baumannii	28	263	27	A27G21C25T22	A32G35C28T34	A40G33C30T36
A. sp. 3	14	164	B7	A26G23C23T23	A30G36C27T36	A39G37C26T37
mixture	7	71	?	ND	ND	ND

Table 18C: Base Compositions Determined from *A. baumannii* DNA Samples Obtained from Northwestern Medical Center and Amplified with Triangulation Genotyping Analysis Primer Pair Nos. 1167, 1170 and 1171

Species	Ibis#	Isolate	ST	PP No: 1167 fumC	PP No: 1170 fumC	PP No: 1171 ppa
A. baumannii	54	536	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	87	665	3	A41G34C35T37	A38G27C20T51	A35G37C31T46
A. baumannii	8	80	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	9	91	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	10	92	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	11	131	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	12	137	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	21	218	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	26	242	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	94	678	10	A41G34C34T38	A38G27C21T50	A35G37C33T44

A. baumannii	1	9	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	2	13	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	3	19	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	4	24	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	5	36	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	6	39	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	13	139	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	15	165	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	16	170	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	17	186	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	20	202	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	22	221	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	24	234	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	25	239	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	33	370	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	34	389	10	A41G34C34T38	A38G27C21T50	A35G37C33T44
A. baumannii	19	201	14	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	27	257	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	29	301	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	31	354	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	36	422	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	37	424	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	38	434	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	39	473	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	40	482	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	44	512	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	45	516	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	47	522	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	48	526	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	50	528	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	52	531	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	53	533	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	56	542	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	59	550	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	62	556	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	64	557	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	70	588	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	73	603	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	74	605	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	75	606	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	77	611	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	79	622	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	83	643	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	85	653	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	89	669	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	93	674	51	A40G35C34T38	A38G27C21T50	A35G37C30T47
A. baumannii	23	228	51	A40G35C34T38	A38G27C21T50	A35G37C30T47

<i>A. baumannii</i>	32	369	52	A40G35C34T38	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	35	393	52	A40G35C34T38	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	30	339	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	41	485	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	42	493	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	43	502	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	46	520	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	49	527	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	51	529	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	65	562	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	68	579	53	A40G35C35T37	A38G27C21T50	A35G37C31T46
<i>A. baumannii</i>	57	546	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	58	548	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	60	552	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	61	555	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	63	557	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	66	570	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	67	578	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	69	584	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	71	593	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	72	602	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	76	609	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	78	621	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	80	625	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	81	628	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	82	632	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	84	649	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	86	655	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	88	668	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	90	671	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	91	672	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	92	673	54	A40G35C34T38	A39G26C22T49	A35G37C31T46
<i>A. baumannii</i>	18	196	55	A42G34C33T38	A38G27C20T51	A35G37C31T46
<i>A. baumannii</i>	55	537	27	A40G35C33T39	A38G27C20T51	A35G37C33T44
<i>A. baumannii</i>	28	263	27	A40G35C33T39	A38G27C20T51	A35G37C33T44
<i>A. sp. 3</i>	14	164	B7	A43G37C30T37	A36G27C24T49	A34G37C31T47
mixture	7	71	-	ND	ND	ND

[422] Base composition analysis of the samples obtained from Walter Reed hospital indicated that a majority of the strain types identified were the same strain types already characterized by the OIF study of Example 12. This is not surprising since at least some patients from which clinical samples were obtained in OIF were transferred to the Walter Reed Hospital (WRAIR). Examples of these common strain types include: ST10, ST11, ST12, ST14, ST15, ST16 and ST46. A strong correlation was noted between these strain types and the presence of mutations in the *gyrA* and *parC* which confer quinolone drug resistance.

[423] In contrast, the results of base composition analysis of samples obtained from Northwestern Medical Center indicate the presence of 4 major strain types: ST10, ST51, ST53 and ST54. All of these strain types have the *gyrA* quinolone resistance mutation and most also have the *parC* quinolone resistance mutation, with the exception of ST35. This observation is consistent with the current understanding that the *gyrA* mutation generally appears before the *parC* mutation and suggests that the acquisition of these drug resistance mutations is rather recent and that resistant isolates are taking over the wild-type isolates. Another interesting observation was that a single isolate of ST3 (isolate 841) displays a triangulation genotyping analysis pattern similar to other isolates of ST3, but the codon analysis amplification product base compositions indicate that this isolate has not yet undergone the quinolone resistance mutations in *gyrA* and *parC*.

[424] The six isolates that represent species other than *Acinetobacter baumannii* in the samples obtained from the Walter Reed Hospital were each found to not carry the drug resistance mutations.

[425] The results described above involved analysis of 183 samples using the methods and compositions of the present invention. Results were provided to collaborators at the Walter Reed hospital and Northwestern Medical center within a week of obtaining samples. This example highlights the rapid throughput characteristics of the analysis platform and the resolving power of triangulation genotyping analysis and codon analysis for identification of and determination of drug resistance in bacteria.

Example 14: Identification of Drug Resistance Genes and Virulence Factors in *Staphylococcus aureus*

[426] An eight primer pair panel was designed for identification of drug resistance genes and virulence factors of *Staphylococcus aureus* and is shown in Table 19. The primer sequences are found in Table 2 and are cross-referenced by the primer pair numbers, primer pair names or SEQ ID NOs listed in Table 19.

Table 19: Primer Pairs for Identification of Drug Resistance Genes and Virulence Factors in *Staphylococcus aureus*

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO:)	Reverse Primer Name	Reverse Primer (SEQ ID NO:)	Target Gene
879	MECA_Y14051_4507_4530_F	288	MECA_Y14051_4555_4581_R	1269	<i>mecA</i>
2056	MECI-R_NC003923-41798-41609_33_60_F	698	MECI-R_NC003923-41798-41609_86_113_R	1420	<i>MecI-R</i>
2081	ERMA_NC002952-55890-	217	ERMA_NC002952-55890-	1167	<i>ermA</i>

	56621_366_395_F		56621_438_465_R		
2086	ERM_C0005908-2004- 2738_85_116_F	399	ERM_C0005908-2004- 2738_173_206_R	1041	ermC
2095	PVL_C0003923-1529595- 1531285_688_713_F	456	PVL_C0003923-1529595- 1531285_775_804_R	1261	Pv-luk
2249	TUF_C0002758-615038- 616222_696_725_F	430	TUF_C0002758-615038- 616222_793_820_R	1321	tufB
2256	NUC_C0002758-894288- 894974_316_345_F	174	NUC_C0002758-894288- 894974_396_421_R	853	Nuc
2313	MUP_C00075439_2486_2516_F	172	MUP_C00075439_2548_2574_R	1360	mupR

[427] Primer pair numbers 2256 and 2249 are confirmation primers designed with the aim of high level identification of *Staphylococcus aureus*. The nuc gene is a *Staphylococcus aureus*-specific marker gene. The tufB gene is a universal housekeeping gene but the bioagent identifying amplicon defined by primer pair number 2249 provides a unique base composition (A43 G28 C19 T35) which distinguishes *Staphylococcus aureus* from other members of the genus *Staphylococcus*.

[428] High level methicillin resistance in a given strain of *Staphylococcus aureus* is indicated by bioagent identifying amplicons defined by primer pair numbers 879 and 2056. Analyses have indicated that primer pair number 879 is not expected to prime *S. sciuri* homolog or *Enterococcus faecalis/facium* ampicillin-resistant PBP5 homologs.

[429] Macrolide and erythromycin resistance in a given strain of *Staphylococcus aureus* is indicated by bioagent identifying amplicons defined by primer pair numbers 2081 and 2086.

[430] Resistance to mupirocin in a given strain of *Staphylococcus aureus* is indicated by bioagent identifying amplicons defined by primer pair number 2313.

[431] Virulence in a given strain of *Staphylococcus aureus* is indicated by bioagent identifying amplicons defined by primer pair number 2095. This primer pair can simultaneously identify the pvl (lukS-PV) gene and the lukD gene which encodes a homologous enterotoxin. A bioagent identifying amplicon of the lukD gene has a six nucleobase length difference relative to the lukS-PV gene.

[432] A total of 32 blinded samples of different strains of *Staphylococcus aureus* were provided by the Center for Disease Control (CDC). Each sample was analyzed by PCR amplification with the eight primer pair panel, followed by purification and measurement of molecular masses of the amplification products by mass spectrometry. Base compositions for the amplification products were calculated. The base compositions provided the information summarized above for each primer pair. The results are shown in Tables 20A and B. One result noted upon un-blinding of the samples is that each of the PVL+ identifications agreed with PVL+ identified in the same samples by standard PCR assays. These results

indicate that the panel of eight primer pairs is useful for identification of drug resistance and virulence sub-species characteristics for *Staphylococcus aureus*. It is expected that a kit comprising one or more of the members of this panel will be a useful embodiment of the present invention.

Table 20A: Drug Resistance and Virulence Identified in Blinded Samples of Various Strains of *Staphylococcus aureus* with Primer Pair Nos. 2081, 2086, 2095 and 2256

Sample Index No.	Primer Pair No. 2081 (exmA)	Primer Pair No. 2086 (exmC)	Primer Pair No. 2095 (pv-luk)	Primer Pair No. 2256 (nuc)
CDC0010	-	-	PVL-/lukD+	+
CDC0015	-	-	PVL+/lukD+	+
CDC0019	-	+	PVL-/lukD+	+
CDC0026	+	-	PVL-/lukD+	+
CDC0030	+	-	PVL-/lukD+	+
CDC004	-	-	PVL+/lukD+	+
CDC0014	-	+	PVL+/lukD+	+
CDC008	-	-	PVL-/lukD+	+
CDC001	+	-	PVL-/lukD+	+
CDC0022	+	-	PVL-/lukD+	+
CDC006	+	-	PVL-/lukD+	+
CDC007	-	-	PVL-/lukD+	+
CDCVRS1	+	-	PVL-/lukD+	+
CDCVRS2	+	+	PVL-/lukD+	+
CDC0011	+	-	PVL-/lukD+	+
CDC0012	-	-	<u>PVL+/lukD-</u>	+
CDC0021	+	-	PVL-/lukD+	+
CDC0023	+	-	PVL-/lukD+	+
CDC0025	+	-	PVL-/lukD+	+
CDC005	-	-	PVL-/lukD+	+
CDC0018	+	-	<u>PVL+/lukD-</u>	+
CDC002	-	-	PVL-/lukD+	+
CDC0028	+	-	PVL-/lukD+	+
CDC003	-	-	PVL-/lukD+	+
CDC0013	-	-	PVL+/lukD+	+
CDC0016	-	-	PVL-/lukD+	+
CDC0027	+	-	PVL-/lukD+	+
CDC0029	-	-	PVL+/lukD+	+

CDC0020	-	+	FVL-/lukD+	+
CDC0024	-	-	FVL-/lukD+	+
CDC0031	-	-	FVL-/lukD+	+

Table 20B: Drug Resistance and Virulence Identified in Blinded Samples of Various Strains of *Staphylococcus aureus* with Primer Pair Nos. 2249, 879, 2056, and 2313

Sample Index No.	Primer Pair No. 2249 (tufB)	Primer Pair No. 879 (mecA)	Primer Pair No. 2056 (mecI-R)	Primer Pair No. 2313 (mupR)
CDC0010	<i>Staphylococcus aureus</i>	+	+	-
CDC0015	<i>Staphylococcus aureus</i>	-	-	-
CDC0019	<i>Staphylococcus aureus</i>	+	+	-
CDC0026	<i>Staphylococcus aureus</i>	+	+	-
CDC0030	<i>Staphylococcus aureus</i>	+	+	-
CDC004	<i>Staphylococcus aureus</i>	+	+	-
CDC0014	<i>Staphylococcus aureus</i>	+	+	-
CDC008	<i>Staphylococcus aureus</i>	+	+	-
CDC001	<i>Staphylococcus aureus</i>	+	+	-
CDC0022	<i>Staphylococcus aureus</i>	+	+	-
CDC006	<i>Staphylococcus aureus</i>	+	+	+
CDC007	<i>Staphylococcus aureus</i>	+	+	-
CDCVRS1	<i>Staphylococcus aureus</i>	+	+	-
CDCVRS2	<i>Staphylococcus aureus</i>	+	+	-
CDC0011	<i>Staphylococcus aureus</i>	-	-	-
CDC0012	<i>Staphylococcus aureus</i>	+	+	-
CDC0021	<i>Staphylococcus aureus</i>	+	+	-
CDC0023	<i>Staphylococcus aureus</i>	+	+	-
CDC0025	<i>Staphylococcus aureus</i>	+	+	-
CDC005	<i>Staphylococcus aureus</i>	+	+	-
CDC0018	<i>Staphylococcus aureus</i>	+	+	-
CDC002	<i>Staphylococcus aureus</i>	+	+	-
CDC0028	<i>Staphylococcus aureus</i>	+	+	-
CDC003	<i>Staphylococcus aureus</i>	+	+	-
CDC0013	<i>Staphylococcus aureus</i>	+	+	-
CDC0016	<i>Staphylococcus aureus</i>	+	+	-
CDC0027	<i>Staphylococcus aureus</i>	+	+	-
CDC0029	<i>Staphylococcus aureus</i>	+	+	-

CDC0020	<i>Staphylococcus aureus</i>	-	-	-
CDC0024	<i>Staphylococcus aureus</i>	+	+	-
CDC0031	<i>Staphylococcus schleiferi</i>	-	-	-

Example 15: Selection and Use of Triangulation Genotyping Analysis Primer Pairs for *Staphylococcus aureus*

[433] To combine the power of high-throughput mass spectrometric analysis of bioagent identifying amplicons with the sub-species characteristic resolving power provided by triangulation genotyping analysis, a panel of eight triangulation genotyping analysis primer pairs was selected. The primer pairs are designed to produce bioagent identifying amplicons within six different housekeeping genes which are listed in Table 21. The primer sequences are found in Table 2 and are cross-referenced by the primer pair numbers, primer pair names or SEQ ID NOs listed in Table 21.

Table 21: Primer Pairs for Triangulation Genotyping Analysis of *Staphylococcus aureus*

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO:)	Reverse Primer Name	Reverse Primer (SEQ ID NO:)	Target Gene
2146	ARCC_NC003923-2725050-2724595 131 161 F	437	ARCC_NC003923-2725050-2724595 214 245 R	1137	arcC
2149	ARO_E_NC003923-1674726-1674277 30 62 F	530	ARO_E_NC003923-1674726-1674277 155 181 R	891	aroE
2150	ARO_E_NC003923-1674726-1674277 204 232 F	474	ARO_E_NC003923-1674726-1674277 308 335 R	869	aroE
2156	GMK_NC003923-1190906-1191334 301 329 F	268	GMK_NC003923-1190906-1191334 403 432 R	1284	gmK
2157	PTA_NC003923-628885-629355 237 263 F	418	PTA_NC003923-628885-629355 314 345 R	1301	pta
2161	TP1_NC003923-830671-831072 1 34 F	318	TP1_NC003923-830671-831072 87 129 R	1300	tpi
2163	YQI_NC003923-378916-379431 142 167 F	440	YQI_NC003923-378916-379431 259 284 R	1076	yqi
2166	YQI_NC003923-378916-379431 275 300 F	219	YQI_NC003923-378916-379431 364 396 R	1013	yqi

[434] The same samples analyzed for drug resistance and virulence in Example 14 were subjected to triangulation genotyping analysis. The primer pairs of Table 21 were used to produce amplification products by PCR, which were subsequently purified and measured by mass spectrometry. Base compositions were calculated from the molecular masses and are shown in Tables 22A and 22B.

Table 22A: Triangulation Genotyping Analysis of Blinded Samples of Various Strains of *Staphylococcus aureus* with Primer Pair Nos. 2146, 2149, 2150 and 2156

Sample Index No.	Strain	Primer Pair No. 2146 (arcC)	Primer Pair No. 2149 (aroE)	Primer Pair No. 2150 (aroE)	Primer Pair No. 2156 (gmK)
CDC0010	COL	A44 G24 C18 T25	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32

CDC0015	COL	A44 G24 C18 T29	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC0019	COL	A44 G24 C18 T29	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC0026	COL	A44 G24 C18 T29	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC0030	COL	A44 G24 C18 T29	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC004	COL	A44 G24 C18 T29	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC0014	COL	A44 G24 C18 T29	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC008	????	A44 G24 C18 T29	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC001	Mu50	A45 G23 C20 T27	A58 G24 C18 T52	A40 G36 C13 T43	A51 G29 C21 T31
CDC0022	Mu50	A45 G23 C20 T27	A58 G24 C18 T52	A40 G36 C13 T43	A51 G29 C21 T31
CDC006	Mu50	A45 G23 C20 T27	A58 G24 C18 T52	A40 G36 C13 T43	A51 G29 C21 T31
CDC0011	MRS252	A45 G24 C18 T28	A58 G24 C19 T51	A41 G36 C12 T43	A51 G29 C21 T31
CDC0012	MRS252	A45 G24 C18 T28	A58 G24 C19 T51	A41 G36 C12 T43	A51 G29 C21 T31
CDC0021	MRS252	A45 G24 C18 T28	A58 G24 C19 T51	A41 G36 C12 T43	A51 G29 C21 T31
CDC0023	ST:110	A45 G24 C18 T28	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC0025	ST:110	A45 G24 C18 T28	A59 G24 C18 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC005	ST:338	A44 G24 C18 T29	A59 G23 C19 T51	A40 G36 C14 T42	A51 G29 C21 T31
CDC0018	ST:338	A44 G24 C18 T29	A59 G23 C19 T51	A40 G36 C14 T42	A51 G29 C21 T31
CDC002	ST:108	A46 G23 C20 T26	A58 G24 C19 T51	A42 G36 C12 T42	A51 G29 C20 T32
CDC0028	ST:108	A46 G23 C20 T26	A58 G24 C19 T51	A42 G36 C12 T42	A51 G29 C20 T32
CDC003	ST:107	A45 G23 C20 T27	A58 G24 C18 T52	A40 G36 C13 T43	A51 G29 C21 T31
CDC0013	ST:112	ND	A59 G24 C18 T51	A40 G36 C13 T43	A51 G29 C21 T31
CDC0016	ST:120	A45 G23 C18 T29	A58 G24 C19 T51	A40 G37 C13 T42	A51 G29 C21 T31
CDC0027	ST:105	A45 G23 C20 T27	A58 G24 C18 T52	A40 G36 C13 T43	A51 G29 C21 T31
CDC0029	MRS2476	A45 G23 C20 T27	A58 G24 C19 T51	A40 G36 C13 T43	A50 G30 C20 T32
CDC0020	ST:15	A44 G23 C21 T27	A59 G23 C18 T52	A40 G36 C13 T43	A50 G30 C20 T32
CDC0024	ST:137	A45 G23 C20 T27	A57 G25 C19 T51	A40 G36 C13 T43	A51 G29 C22 T30
CDC0031	***	No product	No product	No product	No product

Table 22B: Triangulation Genotyping Analysis of Blinded Samples of Various Strains of *Staphylococcus aureus* with Primer Pair Nos. 2146, 2149, 2150 and 2156

Sample Index No.	Strain	Primer Pair No. 2157 (pta)	Primer Pair No. 2161 (tpi)	Primer Pair No. 2163 (yqi)	Primer Pair No. 2166 (yqi)
CDC0010	COL	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC0015	COL	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC0019	COL	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC0026	COL	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37

CDC0030	COL	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC004	COL	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC0014	COL	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC008	unknown	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC001	Mu50	A33 G25 C22 T29	A50 G28 C22 T29	A42 G36 C22 T43	A36 G31 C19 T36
CDC0022	Mu50	A33 G25 C22 T29	A50 G28 C22 T29	A42 G36 C22 T43	A36 G31 C19 T36
CDC006	Mu50	A33 G25 C22 T29	A50 G28 C22 T29	A42 G36 C22 T43	A36 G31 C19 T36
CDC0011	MRSA252	A32 G25 C23 T29	A50 G28 C22 T29	A42 G36 C22 T43	A37 G30 C18 T37
CDC0012	MRSA252	A32 G25 C23 T29	A50 G28 C22 T29	A42 G36 C22 T43	A37 G30 C18 T37
CDC0021	MRSA252	A32 G25 C23 T29	A50 G28 C22 T29	A42 G36 C22 T43	A37 G30 C18 T37
CDC0023	ST:110	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC0025	ST:110	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC005	ST:338	A32 G25 C24 T28	A51 G27 C21 T30	A42 G36 C22 T43	A37 G30 C18 T37
CDC0018	ST:338	A32 G25 C24 T28	A51 G27 C21 T30	A42 G36 C22 T43	A37 G30 C18 T37
CDC002	ST:108	A33 G25 C23 T28	A50 G28 C22 T29	A42 G36 C22 T43	A37 G30 C18 T37
CDC0028	ST:108	A33 G25 C23 T28	A50 G28 C22 T29	A42 G36 C22 T43	A37 G30 C18 T37
CDC003	ST:107	A32 G25 C23 T29	A51 G28 C22 T28	A41 G37 C22 T43	A37 G30 C18 T37
CDC0013	ST:12	A32 G25 C23 T29	A51 G28 C22 T28	A42 G36 C22 T43	A37 G30 C18 T37
CDC0016	ST:120	A32 G25 C24 T28	A50 G28 C21 T30	A42 G36 C22 T43	A37 G30 C18 T37
CDC0027	ST:105	A33 G25 C22 T29	A50 G28 C22 T29	A43 G36 C21 T43	A36 G31 C19 T36
CDC0029	MSSA476	A33 G25 C22 T29	A50 G28 C22 T29	A42 G36 C22 T43	A36 G31 C19 T36
CDC0020	ST:115	A33 G25 C22 T29	A50 G28 C21 T30	A42 G36 C22 T43	A36 G31 C18 T37
CDC0024	ST:137	A33 G25 C22 T29	A51 G28 C22 T28	A42 G36 C22 T43	A37 G30 C18 T37
CDC0031	***	A34 G25 C25 T25	A51 G27 C24 T27	No product	No product

[435] Note: *** The sample CDC0031 was identified as *Staphylococcus scleiferi* as indicated in Example 14. Thus, the triangulation genotyping primers designed for *Staphylococcus aureus* would generally not be expected to prime and produce amplification products of this organism. Tables 22A and 22B indicate that amplification products are obtained for this organism only with primer pair numbers 2157 and 2161.

[436] A total of thirteen different genotypes of *Staphylococcus aureus* were identified according to the unique combinations of base compositions across the eight different bioagent identifying amplicons obtained with the eight primer pairs. These results indicate that this eight primer pair panel is useful for analysis of unknown or newly emerging strains of *Staphylococcus aureus*. It is expected that a kit

comprising one or more of the members of this panel will be a useful embodiment of the present invention.

Example 16: Selection and Use of Triangulation Genotyping Analysis Primer Pairs for Members of the Bacterial Genus *Vibrio*

[437] To combine the power of high-throughput mass spectrometric analysis of bioagent identifying amplicons with the sub-species characteristic resolving power provided by triangulation genotyping analysis, a panel of eight triangulation genotyping analysis primer pairs was selected. The primer pairs are designed to produce bioagent identifying amplicons within seven different housekeeping genes which are listed in Table 23. The primer sequences are found in Table 2 and are cross-referenced by the primer pair numbers, primer pair names or SEQ ID NOs listed in Table 23.

Table 23: Primer Pairs for Triangulation Genotyping Analysis of Members of the Bacterial Genus *Vibrio*

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO:)	Reverse Primer Name	Reverse Primer (SEQ ID NO:)	Target Gene
1098	RVASEP VBC 331 349 F	325	RVASEP VBC 388 414 R	1163	RNAse P
2000	CTXB_NC002505 46 70 F	278	CTXB_NC002505 132 162 R	1039	ctxB
2001	FUR_NC002505 87 113 F	465	FUR_NC002505 205 229 R	1037	fur
2011	GYRB_NC002505_1161_1190 F	148	GYRB_NC002505 1255 1284 R	1172	gyrB
2012	OMP NC002505 85 110 F	190	OMP NC002505 154 180 R	1254	ompU
2014	OMP NC002505 431 455 F	266	OMP NC002505 544 567 R	1094	ompU
2323	CTXA_NC002505-1568114-1567341 122 149 F	508	CTXA_NC002505-1568114-1567341 186 214 R	1297	ctxA
2927	GAPA_NC002505 694 721 F	259	GAPA_NC002505 29 58 R	1060	gapA

[438] A group of 50 bacterial isolates containing multiple strains of both environmental and clinical isolates of *Vibrio cholerae*, 9 other *Vibrio* species, and 3 species of Photobacteria were tested using this panel of primer pairs. Base compositions of amplification products obtained with these 8 primer pairs were used to distinguish amongst various species tested, including sub-species differentiation within *Vibrio cholerae* isolates. For instance, the non-O1/non-O139 isolates were clearly resolved from the O1 and the O139 isolates, as were several of the environmental isolates of *Vibrio cholerae* from the clinical isolates.

[439] It is expected that a kit comprising one or more of the members of this panel will be a useful embodiment of the present invention.

Example 17: Selection and Use of Triangulation Genotyping Analysis Primer Pairs for Members of the Bacterial Genus *Pseudomonas*

[440] To combine the power of high-throughput mass spectrometric analysis of bioagent identifying amplicons with the sub-species characteristic resolving power provided by triangulation genotyping analysis, a panel of twelve triangulation genotyping analysis primer pairs was selected. The primer pairs are designed to produce bioagent identifying amplicons within seven different housekeeping genes which are listed in Table 24. The primer sequences are found in Table 2 and are cross-referenced by the primer pair numbers, primer pair names or SEQ ID NOs listed in Table 24.

Table 24: Primer Pairs for Triangulation Genotyping Analysis of Members of the Bacterial Genus *Pseudomonas*

Primer Pair No.	Forward Primer Name	Forward Primer (SEQ ID NO:1)	Reverse Primer Name	Reverse Primer (SEQ ID NO:1)	Target Gene
2949	ACS_NC002516-970624-971013 299 316 F	376	ACS_NC002516-970624-971013 364 383 R	1265	acsA
2950	ARO_NC002516-26883-27380 4 26 F	267	ARO_NC002516-26883-27380 111 128 R	1341	aroE
2951	ARO_NC002516-26883-27380 356 377 F	705	ARO_NC002516-26883-27380 459 484 R	1056	aroE
2954	GUA_NC002516-4226546-4226174 155 178 F	710	GUA_NC002516-4226546-4226174 265 287 R	1259	guaA
2956	GUA_NC002516-4226546-4226174 242 263 F	374	GUA_NC002516-4226546-4226174 355 371 R	1111	guaA
2957	MUT_NC002516-5551150-5550717 5 26 F	545	MUT_NC002516-5551150-5550717 99 116 R	978	mutL
2959	NUO_NC002516-2984589-2984954 8 26 F	249	NUO_NC002516-2984589-2984954 97 117 R	1095	nucD
2960	NUO_NC002516-2984589-2984954 218 239 F	195	NUO_NC002516-2984589-2984954 301 326 R	1376	nucD
2961	PPS_NC002516-1915014-1915383 44 63 F	311	PPS_NC002516-1915014-1915383 140 165 R	1014	pps
2962	PPS_NC002516-1915014-1915383 240 258 F	365	PPS_NC002516-1915014-1915383 341 360 R	1052	pps
2963	TRP_NC002516-671831-672273 24 42 F	527	TRP_NC002516-671831-672273 131 150 R	1071	trpE
2964	TRP_NC002516-671831-672273 261 282 F	490	TRP_NC002516-671831-672273 362 383 R	1182	trpE

[441] It is expected that a kit comprising one or more of the members of this panel will be a useful embodiment of the present invention.

[442] The present invention includes any combination of the various species and subgeneric groupings falling within the generic disclosure. This invention therefore includes the generic description of the invention with a proviso or negative limitation removing any subject matter from the genus, regardless of whether or not the excised material is specifically recited herein.

[443] While in accordance with the patent statutes, description of the various embodiments and examples have been provided, the scope of the invention is not to be limited thereto or thereby.

Modifications and alterations of the present invention will be apparent to those skilled in the art without departing from the scope and spirit of the present invention.

[444] Therefore, it will be appreciated that the scope of this invention is to be defined by the appended claims, rather than by the specific examples which have been presented by way of example.

[445] Each reference (including, but not limited to, journal articles, U.S. and non-U.S. patents, patent application publications, international patent application publications, gene bank gi or accession numbers, internet web sites, and the like) cited in the present application is incorporated herein by reference in its entirety.

CLAIMS

What is claimed is:

1. An oligonucleotide primer 14 to 35 nucleobases in length comprising at least 70% sequence identity with SEQ ID NO: 456.
2. An oligonucleotide primer 14 to 35 nucleobases in length comprising at least 70% sequence identity with SEQ ID NO: 1261.
3. A composition comprising the primer of claim 1.
4. The composition of claim 3 further comprising an oligonucleotide primer 14 to 35 nucleobases in length comprising at least 70% sequence identity with SEQ ID NO: 1261.
5. The composition of claim 4 wherein either or both of said primers comprises at least one modified nucleobase.
6. The composition of claim 5 wherein said modified nucleobase is 5-propynyluracil or 5-propynylcytosine.
7. The composition of claim 4 wherein either or both of said primers comprises at least one universal nucleobase.
8. The composition of claim 7 wherein said universal nucleobase is inosine.
9. The composition of claim 4 wherein either or both of said primers further comprises a non-templated T residue on the 5'-end.
10. The composition of claim 4 wherein either or both of said primers comprises at least one non-template tag.
11. The composition of claim 4 wherein either or both of said primers comprises at least one molecular mass modifying tag.
12. A kit comprising the composition of claim 4.

13. The kit of claim 12 further comprising one or more primer pairs wherein each member of said one or more primer pairs is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from the group of primer pairs represented by SEQ ID NOs: 288:1269, 698:1420, 217:1167, 399:1041, 430:1321, 174:853, and 172:1360.
14. The kit of claim 12 further comprising one or more calibration polynucleotides.
15. The kit of claim 12 further comprising at least one anion exchange functional group linked to a magnetic bead.
16. A method for identification of a strain of *Staphylococcus aureus* in a sample comprising:
amplifying nucleic acid from said strain of *Staphylococcus aureus* using the composition of claim 4 to obtain an amplification product;
determining the molecular mass of said amplification product;
optionally, determining the base composition of said amplification product from said molecular mass; and
comparing said molecular mass or said base composition with a plurality of molecular masses or base compositions of known amplification products of strains of *Staphylococcus aureus* defined by the composition of claim 4, wherein a match between said molecular mass or base composition and a member of said plurality of molecular masses or base compositions identifies said strain of *Staphylococcus aureus*.
17. The method of claim 16 further comprising repeating said amplifying, determining, optionally determining, and comparing steps using at least one additional primer pair, wherein each member of said at least one additional primer pair is of a length of 14 to 35 nucleobases and has 70% to 100% sequence identity with the corresponding member from the group of primer pairs represented by SEQ ID NOs: 288:1269, 698:1420, 217:1167, 399:1041, 430:1321, 174:853, and 172:1360.
18. The method of claim 16 wherein said strain of *Staphylococcus aureus* is a virulent strain.
19. The method of claim 18 wherein said strain of *Staphylococcus aureus* is a virulent strain.
20. A method for determination of the quantity of a strain of *Staphylococcus aureus* in a sample comprising:

contacting said sample with the composition of claim 4 and a known quantity of a calibration polynucleotide comprising a calibration sequence;

concurrently amplifying nucleic acid from said a strain of *Staphylococcus aureus* and nucleic acid from said calibration polynucleotide in said sample with the composition of claim 4 to obtain a first amplification product comprising a bacterial bioagent identifying amplicon and a second amplification product comprising a calibration amplicon;

determining the molecular mass and abundance for said bacterial bioagent identifying amplicon and said calibration amplicon; and

distinguishing said bacterial bioagent identifying amplicon from said calibration amplicon based on molecular mass, wherein comparison of bacterial bioagent identifying amplicon abundance and calibration amplicon abundance indicates the quantity of said strain of *Staphylococcus aureus* in said sample.

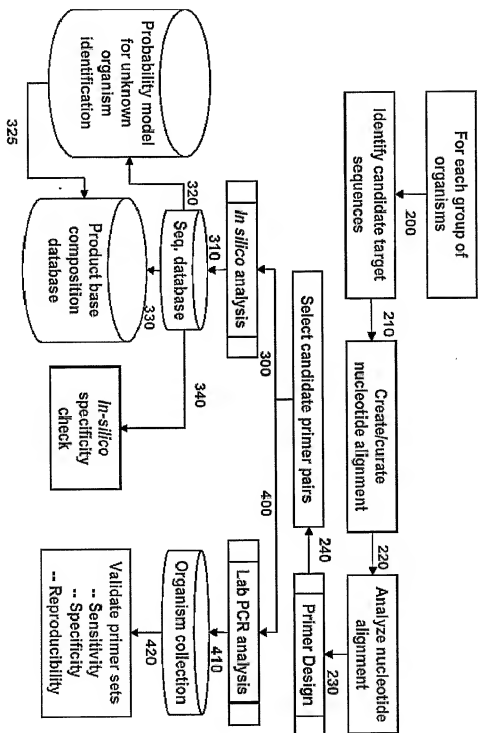


Figure 1

2/7

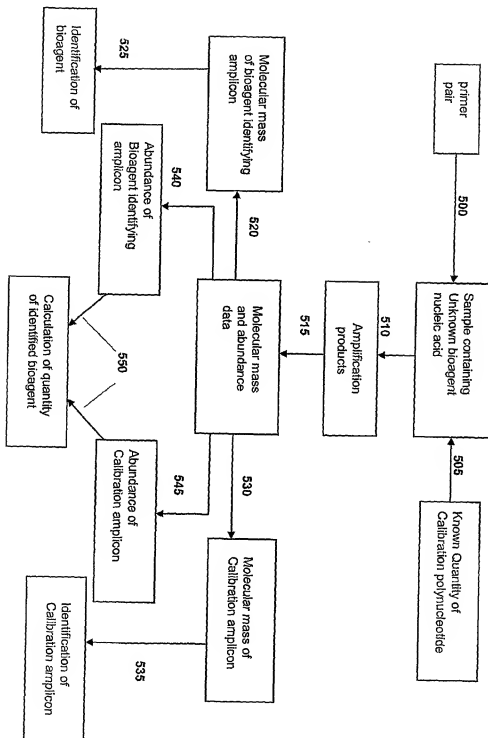


Figure 2

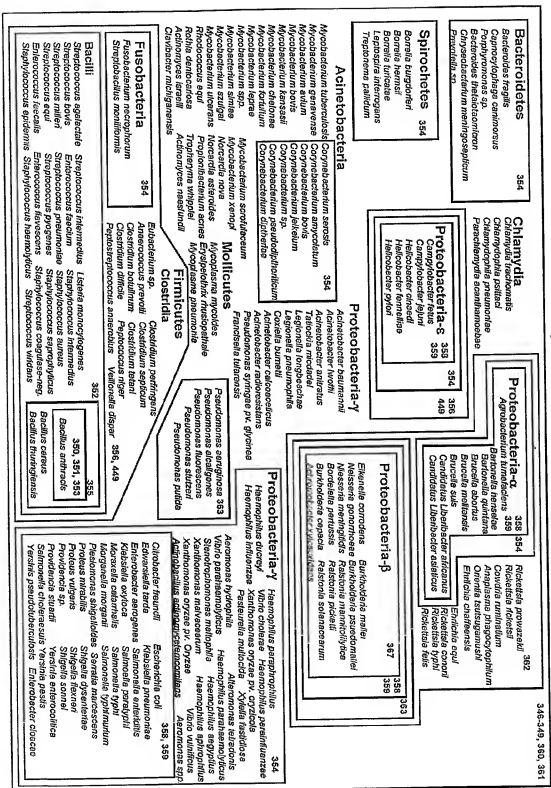


Figure 3

4/7

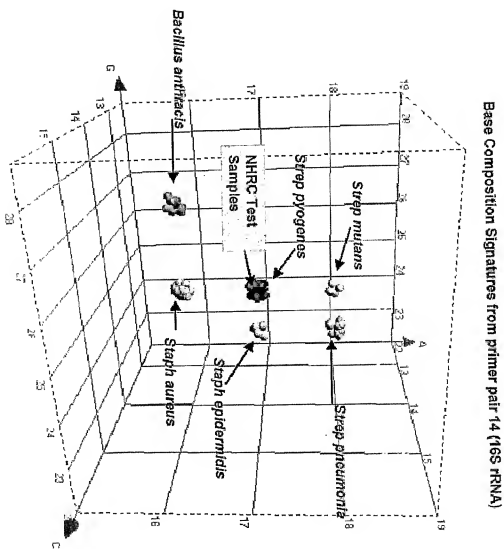


Figure 4

5/7

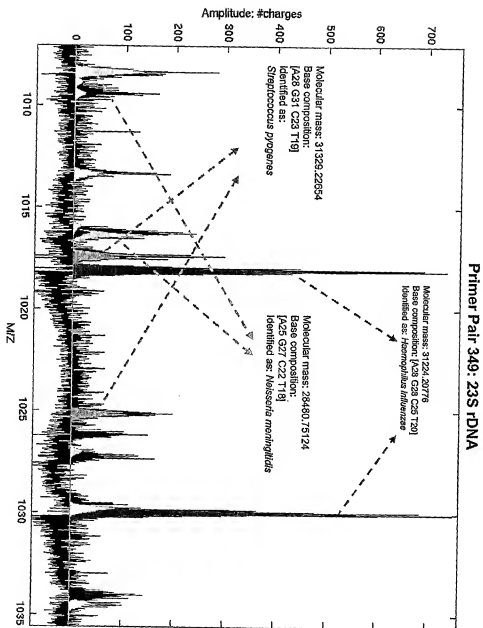
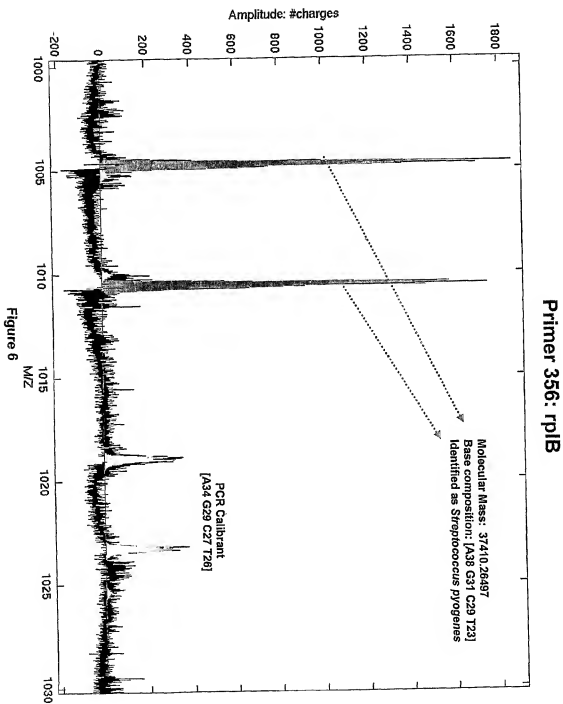


Figure 5

6/7



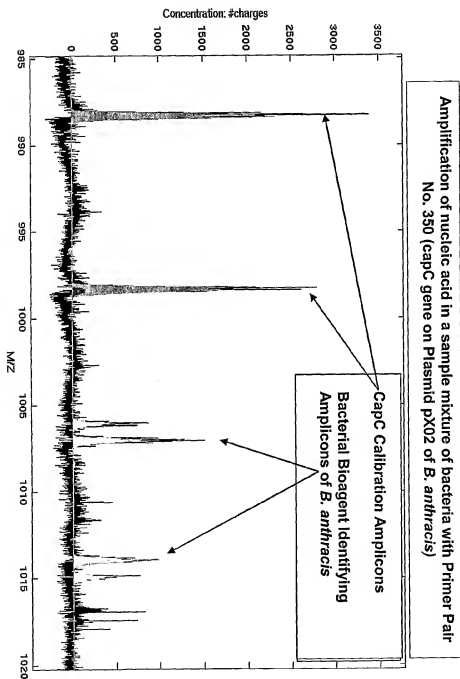


Figure 7

SEQUENCE LISTING

<110> Sampath, Rangarajan
Hall, Thomas A.
Eshoo, Mark W.

<120> COMPOSITIONS FOR IDENTIFICATION OF
BACTERIA

<130> DIBIS-0083WO

<150> 60/674,118
<151> 2005-04-21

<150> 60/705,631
<151> 2005-08-03

<150> 60/732,539
<151> 2005-11-01

<150> 60/773,124
<151> 2006-02-13

<160> 1464

<170> FastSEQ for Windows Version 4.0

<210> 1
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1
aaactagata acagtagaca tcac 24

<210> 2
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 2
aaccttaatt ggaaagaaac ccaagaagt 29

<210> 3
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 3
aacgcacaat cagaagc 17

<210> 4
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 4
aactaccgct cgcagttcta cttcc

25

<210> 5
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 5
aactaccgct ctcagttcta cttcc

25

<210> 6
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 6
aagacgacct gcacgggc

18

<210> 7
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 7
aagcggtgga gcatgtgg

18

<210> 8
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 8
aagggaaggcg tgatcaccgt tgaaga

26

<210> 9
<211> 23
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 9

aagggtactcc ggggataaca ggc

23

<210> 10

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 10

aagtcggaat cgctagtaat cg

22

<210> 11

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 11

aatctgctat ttggtcagg

19

<210> 12

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 12

acaacgaagt acaatacaag ac

22

<210> 13

<211> 21

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 13

acaatacaag acaaaagaag g

21

<210> 14

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 14

accacgcgt aaacgatga

19

<210> 15

<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 15
accatgacag aaggcatttt gaca 24

<210> 16
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 16
acccagtgtct gctgaacogt gc 22

<210> 17
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 17
accgagcaag gagaccagc 19

<210> 18
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 18
acctgcccag tgotggaag 19

<210> 19
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 19
acgcgaagaa ccttacc 17

<210> 20
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 20
actcgttttt aatcagcccg 20

<210> 21
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 21
agaacaccga tggcgaaggc 20

<210> 22
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 22
agaatcaagt tcccaggggt tac 23

<210> 23
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 23
agagtttgat catggctcag 20

<210> 24
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 24
agcagggtgt gaaatcggcc acatgatt 28

<210> 25
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 25
agcgtaaagg tgaacctt 18

<210> 26
<211> 26

<212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 26
 agctttttgca tattatatcg agccac 26

 <210> 27
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 27
 aggacagagt gactactttg accgaggt 28

 <210> 28
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 28
 agtctcaaga gtgaacacgt aa 22

 <210> 29
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 29
 agttataaac accgctttcc tatggcttat cc 32

 <210> 30
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 30
 atactectga ctgacogata g 21

 <210> 31
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 31
atatcgacgg cggtgtttgg 20

<210> 32
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 32
atcaatttgg tggccaagaa cctgg 25

<210> 33
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 33
atgattacaa ttcaagaagg tcgtcaogc 29

<210> 34
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 34
atggacaagg ttggcaagga agg 23

<210> 35
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 35
atggccatgg cagaagctca 20

<210> 36
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 36
atgtcgattg caatcggtac ttgtg 25

<210> 37
<211> 19
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 37

atggtggggtt aagtccgc

19

<210> 38

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 38

atggtggggtt aagtccgcga acgag

25

<210> 39

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 39

caaaacttat taggtaagcg tgttgact

28

<210> 40

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 40

caaaggtaag caaggacgtt tccgtca

27

<210> 41

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 41

caaaggtaag caaggtcggt tccgtca

27

<210> 42

<211> 17

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 42
caacgagcgc aaccctt 17
<210> 43
<211> 17
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 43
caacggatgc tggcaag 17
<210> 44
<211> 31
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 44
caagaagaaa aagagcttct aaaaagaata c 31
<210> 45
<211> 23
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 45
caagcaaacy cacaatcaga agc 23
<210> 46
<211> 19
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 46
caagtcatca tggccotta 19
<210> 47
<211> 26
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 47
caataacgca acagcggtgg cttggg 26
<210> 48
<211> 19

<212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 48
 cactggaact gagacacgg 19

<210> 49
 <211> 20
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 49
 cagaatcaag ttcccagggg 20

<210> 50
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 50
 cagagaccgt tttatcctat cagc 24

<210> 51
 <211> 20
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 51
 cagcggttcg gcgaatgga 20

<210> 52
 <211> 29
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 52
 caggagtctg tcaactcgat ctacatgat 29

<210> 53
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 53
cagggttttagt accagaacat gcag 24

<210> 54
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 54
catccacacg gtgggtggtga agg 23

<210> 55
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 55
ccacacgccg ttcttcaaca act 23

<210> 56
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 56
ccacagttct acttccgtac tactgacg 28

<210> 57
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 57
ccagcagccg cggtaatac 19

<210> 58
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 58
ccgtaacttc gggagaagg 19

<210> 59
<211> 20
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 59

ccgtggtatt ggagttattg

20

<210> 60

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 60

cctatattaa tcgtttacag aaactggct

29

<210> 61

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 61

cctgataagg gtgaggtcg

19

<210> 62

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 62

ccttacttcg aactatgaat cttttggaag

30

<210> 63

<211> 14

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 63

cgaagaacct tacc

14

<210> 64

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 64
cgaagtacaa tacaagacaa aagaagg 27

<210> 65
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 65
cgacgcgctg cgcttcac 18

<210> 66
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 66
cgagaggggaa acaaccoaga cc 22

<210> 67
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 67
cgagtatagc taaaaaata gtttatgaca 30

<210> 68
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 68
cgcaaaaaaa tccagctatt agc 23

<210> 69
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 69
cgccgacttc gacggtgacc 20

<210> 70
<211> 20

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 70
cggaattact gggcgtaaag 20

<210> 71
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 71
cggattggag tctgcaactc g 21

<210> 72
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 72
cggcgtactt caacgacagc ca 22

<210> 73
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 73
cgtaactata acggtcctaa ggta 24

<210> 74
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 74
cgtcagggtg aattccgtga agttaa 26

<210> 75
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 75
 cgtcgggtga ttaaccgtaa caaccg 26

 <210> 76
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 76
 cgtcgtgtaa ttaaccgtaa caaccg 26

 <210> 77
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 77
 cgtggcggcg tggttatcga 20

 <210> 78
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 78
 cgtgttgact attcggggcg ttcag 25

 <210> 79
 <211> 18
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 79
 ctagtacgag aggaccgg 18

 <210> 80
 <211> 18
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 80
 ctgacacctg cccggtgc 18

 <210> 81

<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 81
ctggcaggta tgcgtggtct gatg

<210> 82
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 82
ctggctaaaa ctttggcaac ggt

<210> 83
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 83
ctgtcccctag tacgagagga cggg

<210> 84
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 84
ctgttcttag tacgagagga cc

<210> 85
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 85
cttctgcaac aagctgtgga acgc

<210> 86
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

24

23

24

22

24

<400> 86
 cttgctggta tgcgtgggtct gatg 24

 <210> 87
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 87
 cttggaggta agtctcattt tggtagggca 29

 <210> 88
 <211> 19
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 88
 cttgtacaca ccgcccgtc 19

 <210> 89
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 89
 cttgtacttg tggctcacac ggctgtttgg 30

 <210> 90
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 90
 cttttgcata ttatatcgag c 21

 <210> 91
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 91
 gaatagcaat taatccaaat 20

 <210> 92
 <211> 18

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 92
gaaagagttc ggattggg

18

<210> 93
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 93
gaaggatata cggttgatgt c

21

<210> 94
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 94
gaatagcaat taatccaaat

20

<210> 95
<211> 20
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 95
gacacggtcc agactctac

20

<210> 96
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 96
gacagttcgg tcctatc

18

<210> 97
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 97
gaccacctcg gcaaccgt 18

<210> 98
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 98
gacctacagt aagaggttct gtaatgaacc 30

<210> 99
<211> 16
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 99
gacgcctgcc cgggtgc 16

<210> 100
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 100
gacttaccaa ccgatgcaa 20

<210> 101
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 101
gagagcaagc ggacctcata 20

<210> 102
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 102
gagagtttga tcctggctca gaacgaa 27

<210> 103

<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 103
gaggaaagtc catgctcac 19

<210> 104
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 104
gaggaaagtc catgctcgc 19

<210> 105
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 105
gaggaaagtc cgggctc 17

<210> 106
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 106
gataccctgg tagtccacac cg 22

<210> 107
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 107
gatctggagg aataccggtg 20

<210> 108
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 108 gatgactttt tagctaattg tcaggcagc	29
<210> 109 <211> 30 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 109 gattattgtt atcctgttat gccatttgag	30
<210> 110 <211> 18 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 110 gcacaacctg cggtcgcg	18
<210> 111 <211> 27 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 111 gcactatgca cacgtagatt gtctctgg	27
<210> 112 <211> 21 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 112 gccttgtaca cacctcccg c	21
<210> 113 <211> 20 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 113 gcgaagaacc ttaccaggtc	20
<210> 114 <211> 20	

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 114
gctacacacg tgctacaatg 20

<210> 115
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 115
gctgggtgaaa ataaccacaga tgtogtcttc 30

<210> 116
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 116
gcttcaggaa tcaatgatgg agcag 25

<210> 117
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 117
ggacggagaa ggctatggt 19

<210> 118
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 118
ggattagaga ccctggtagt cc 22

<210> 119
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 119	
ggattagata ccctggtagt ccacgc	26
<210> 120	
<211> 25	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 120	
ggctcagcca tttagttacc gctat	25
<210> 121	
<211> 21	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 121	
gggaactgaa acatctaagt a	21
<210> 122	
<211> 22	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 122	
gggagcaaac aggattagat ac	22
<210> 123	
<211> 29	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 123	
gggcaacagc agcggattgc gattgcgcg	29
<210> 124	
<211> 23	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 124	
gggcagcgtt tcggcgaaat gga	23
<210> 125	
<211> 26	
<212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 125

ggggagtgaa agagatccctg aaaccg

26

<210> 126

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 126

ggggattcag ccatcaaagc agctattgac

30

<210> 127

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 127

ggggattgat atcaccgata agaagaa

27

<210> 128

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 128

ggtgaaagaa gttgcctcta aagc

24

<210> 129

<211> 15

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 129

ggtggatgcc ttggc

15

<210> 130

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 130
 ggtgtttaat agcctggcag 20

 <210> 131
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 131
 ggtttagtac cagaacatgc 20

 <210> 132
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 132
 gtcaaatgg cacgtttact ggc 23

 <210> 133
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 133
 gtcgtgaaaa cgagctggaa ga 22

 <210> 134
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 134
 gtgagatgtt ggggttaagtc cagtaacgag 30

 <210> 135
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 135
 gtgcatgcgg atacagagca gag 23

 <210> 136
 <211> 26
 <212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 136

gtggcatgcc taatacatgc aagtcg

26

<210> 137

<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 137

gtgtagcggg gaaatgcg

18

<210> 138

<211> 21

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 138

ggtatcctgt tatgccattt g

21

<210> 139

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 139

gttatttagc actogttttt aatcagcc

28

<210> 140

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 140

gttgtagagg taagcgacta ag

22

<210> 141

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 141

gtttgtgaggt taagcgacta ag 22

<210> 142
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 142 24
 tctagtaata ataggaccct cagc

<210> 143
 <211> 15
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 143 15
 tatgggtccta ctcaa

<210> 144
 <211> 29
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 144 29
 taaaacaaac tacggtaaca ttgatcgca

<210> 145
 <211> 33
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 145 33
 taaaactttt gcogtaatga tgggtgaaga tat

<210> 146
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 146 31
 taaacacggc ttctctatgg cttatccaaa t

<210> 147
 <211> 28
 <212> DNA
 <213> Artificial Sequence

<220>
<223> Primer

<400> 147
taaaccatcg cgggagcaag accgaata 28

<210> 148
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 148
taaagcccggt gaaatgactc gtgctaaagg 30

<210> 149
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 149
taaagttggt tttattggtt ggcgcgga 28

<210> 150
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 150
taaatctgcc cgtgtcgttg gtgac 25

<210> 151
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 151
taacaactcg ccttatgaaa cgggatata 29

<210> 152
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 152
taacacatgc aagtogaacg 20

<210> 153
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 153
 taaaccattca agaactagat cttcaggca 29

 <210> 154
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 154
 taacottaat tggaagaaaa cccaagaagt 30

 <210> 155
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 155
 taacgggttat catggcccag atggg 25

 <210> 156
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 156
 taactctgat gtttttgatg ggaaggt 27

 <210> 157
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 157
 taactgcatg gaacccttct ttactag 27

 <210> 158
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
<223> Primer

<400> 158
taagaagccg gaaaccatca actaccg 27

<210> 159
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 159
taagagcgca ccggtagtt gg 22

<210> 160
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 160
taagcatgct gtggcttacc gtgaaatg 28

<210> 161
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 161
taagctgcc gccgaatgct ttc 23

<210> 162
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 162
taaggatagt gcaacagaga tataccgcc 29

<210> 163
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 163
taaggtatga caccggataa atcatataaa 30

<210> 164
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 164
taaggtttat tgtctttgtg gagatgggga ttt 33

<210> 165
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 165
taatcaagca ttggaagatg aaatgcatac c 31

<210> 166
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 166
taatcggtaa atatcaccog catggtgac 29

<210> 167
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 167
taatcggtaa gtatcaccct catggtgat 29

<210> 168
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 168
taatcggtga atacgggttt gcta 24

<210> 169
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
 <223> Primer

 <400> 169
 taatgaaccc taatgaccat ccacacggtg 30

 <210> 170
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 170
 taatgatgaa ttaggtgogg gttcttt 27

 <210> 171
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 171
 taatgggtaa atatcacctt catggtgac 29

 <210> 172
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 172
 taattgggct ctttctogct taaacacctt a 31

 <210> 173
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 173
 tacaaagcaa gacactggct cacta 25

 <210> 174
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 174
 tacaaaggtc aaccaatgac attcagacta 30

<210> 175
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 175
tacaacatat tattaagag acgggttga atcc 34

<210> 176
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 176
tacaagcact ccagctgca 20

<210> 177
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 177
tacaatgctt gtttatgctg gtaaagcag 29

<210> 178
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 178
tacacaaca tggcggtaaa gatgg 25

<210> 179
<211> 15
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 179
tacagagttt gcgac 15

<210> 180
<211> 21
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 180

tacaggccgt gttgaacgtg g

21

<210> 181

<211> 21

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 181

tacatgctag ccgcgtctta c

21

<210> 182

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 182

taccactatt aatgtcgctg gtgcttc

27

<210> 183

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 183

taccatgaca gaaggcattt tgaca

25

<210> 184

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 184

taccccaaac cgacacagg

19

<210> 185

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 185

tacccacagg aaagtgccac aga

23

<210> 186

<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 186
taccggcgca aaaagtcgag attgg 25

<210> 187
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 187
tacctatatg cgccagacg c 21

<210> 188
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 188
tacgatttca cttccgcagc cagatt 26

<210> 189
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 189
tacgcgtctt gaagcgtttc gttatga 27

<210> 190
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 190
tacgctgacg gaatcaacca aagcgg 26

<210> 191
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 191
tacggtgaat acgttcccg g 21

<210> 192
<211> 15
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 192
tacagagttt gcgac 15

<210> 193
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 193
tactacttca agccgaactt ccg 23

<210> 194
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 194
tactagcggg aagcttaaac aagattgc 28

<210> 195
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 195
tactctcggg ggagaagctc gc 22

<210> 196
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 196
tactggaaca aagtctgcga cc 22

<210> 197

<211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 197
 tacttactac ttcaagccga acttccg 27

 <210> 198
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 198
 tacttacttg agaatccaca agctgcaa 28

 <210> 199
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 199
 tacttggtga ataccaccca catggtga 28

 <210> 200
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 200
 tactttttta aaactaggga tgcgtttgaa gc 32

 <210> 201
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 201
 tagaaatcaa ggtgatagtg gcaatga 27

 <210> 202
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 202
tagaacaccg atggcgaagg c 21

<210> 203
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 203
tagaacgtcg cgagacagtt cg 22

<210> 204
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 204
tagactgccc aggacacgct g 21

<210> 205
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 205
tagataattg ggctctttct cgcttaaac 29

<210> 206
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 206
tagataccct ggtagtccac gc 22

<210> 207
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 207
tagatgaaaa aggcgaagtg gctaatgg 28

<210> 208
<211> 28

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 208
tagatgaaaa gggcgaagtg gctaattg 28

<210> 209
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 209
tagcaacaaaa tatattctgaa gcagcgtact 30

<210> 210
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 210
tagcagggtgg tgaaatcggc cacatgatt 29

<210> 211
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 211
tagcatcaga actgttgttc cgctag 26

<210> 212
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 212
tagcccagca caatttgtga ttca 24

<210> 213
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 213
 tagcctttaa cgaaaatgta aaatgcgtt ttga 34
 <210> 214
 <211> 27
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 214
 tagcgaatgt ggctttactt cacaatt 27
 <210> 215
 <211> 19
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 215
 tagcgtaaag gtgaacctt 19
 <210> 216
 <211> 20
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 216
 tagctaattg tcaggcagcc 20
 <210> 217
 <211> 30
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 217
 tagctatctt atcgttgaga agggatttgc 30
 <210> 218
 <211> 22
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 218
 tagctggcgc gaaattaggt gt 22
 <210> 219
 <211> 26
 <212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 219
tagctggcgg tatggagaat atgtct 26

<210> 220
<211> 27
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 220
tagcttttgc atattataac gagccac 27

<210> 221
<211> 29
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 221
taggaattac ggctgataaa gcgtataaa 29

<210> 222
<211> 35
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 222
taggcgaaga tatacaaaga gtattagaag ctaga 35

<210> 223
<211> 25
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 223
taggcgtgaa agcaagctac cgttt 25

<210> 224
<211> 25
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 224

tagtgcctgg ttacgcagat caaga	25
<210> 225	
<211> 30	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 225	
taggtttacg tcagtatggc gtgattatgg	30
<210> 226	
<211> 23	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 226	
tagtaccgaa gctggtcata cga	23
<210> 227	
<211> 17	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 227	
tagtaccgaga ggaccgg	17
<210> 228	
<211> 18	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 228	
tagtcccga acgagcgc	18
<210> 229	
<211> 27	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 229	
tagtgataga actgtaggca caatcgt	27
<210> 230	
<211> 22	
<212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 230

tagttgctca aacagctggg ct 22

<210> 231

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 231

tataagtggg taaacctga atatcgtgt 29

<210> 232

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 232

tataactcaa cgctgctgc ttcc 24

<210> 233

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 233

tatcgctcag gcgaactcca ac 22

<210> 234

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 234

tatgacaaa ctcacgac gag 23

<210> 235

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 235

tatgattaca attcaagaag gtcgtcacgc 30

<210> 236
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 236
 tatgcagtggt aacgatgggt tccaaga 27

<210> 237
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 237
 tatgctgacc gaccagtgggt acgt 24

<210> 238
 <211> 21
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 238
 tatggccatg gcagaagctc a 21

<210> 239
 <211> 15
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 239
 tatggctcta ctcaa 15

<210> 240
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 240
 tatgtccaag aagcatagca aaaaaagcaa t 31

<210> 241
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer
 <400> 241
 tattcaaggt ggctccttga tgcattgt 27
 <210> 242
 <211> 22
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 242
 tattggacaa cggctgcgc gg 22
 <210> 243
 <211> 30
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 243
 tattgtttca aatgtacaag gtgaagtgg 30
 <210> 244
 <211> 30
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 244
 tatttcacat gtaattttga tattgcact 30
 <210> 245
 <211> 32
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 245
 tcaaaaagcc ctaggtaag agattcata tc 32
 <210> 246
 <211> 29
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 246
 tcaaaactggg caatcggaac tggtaaatc 29

<210> 247
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 247
tcaaatgtac aaggtgaagt gcgtga

26

<210> 248
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 248
tcaacaacct cttggaggta aagctcagt

29

<210> 249
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 249
tcaacctcgg ccgaacca

19

<210> 250
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 250
tcaacctgac tgcgtgaatg gttgt

25

<210> 251
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 251
tcaacgaagg taaaaacat ctcaacy

27

<210> 252
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
 <223> Primer

 <400> 252
 tcaacggtaa cttctatggt actctctg 27

 <210> 253
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 253
 tcaactcgaa ttttcaacag gtacca 26

 <210> 254
 <211> 17
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 254
 tcaagaagaa aaagagc 17

 <210> 255
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 255
 tcaagcaaac gcacaatcag aagc 24

 <210> 256
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 256
 tcaagcagaa gctttggaag aagaagg 27

 <210> 257
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 257
 tcaagccgta cgtattatta ggtgctg 27

<210> 258
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 258
tcaataccgc aacagcgggtg gcttggg

27

<210> 259
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 259
tcaatgaacg accaacaagt gattgatg

28

<210> 260
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 260
tcaatgaacg atcaacaagt gattgatg

28

<210> 261
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 261
tcacatatcg tgagcaatga actg

24

<210> 262
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 262
tcaccagggt caactcaaaa aatattaaca

30

<210> 263
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 263
tcaccagttt gccacgtatc ttcaa 25

<210> 264
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 264
tcacctcat ggtgactcat ctatttat 28

<210> 265
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 265
tcacctcat ggtgattcag ctgtttat 28

<210> 266
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 266
tcaccgatat catggcttac cacgg 25

<210> 267
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 267
tcaccgtgcc gttaaggaa gag 23

<210> 268
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 268
tcacctcaa gtttagatca cttgagaga 29

<210> 269
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 269
tcacgataag aaaaccggtc aagagg 26

<210> 270
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 270
tcactcttac atataaggaa ggcgctc 27

<210> 271
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 271
tcagaccatg ctgcagaga aactt 25

<210> 272
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 272
tcagagaccg ttttatccta tcagc 25

<210> 273
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 273
tcagcaaatg catcacaaac agataa 26

<210> 274
<211> 25
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 274
tcagcatatg cacatggaac acctc

25

<210> 275
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 275
tcagcatatg cacatggaac acctca

26

<210> 276
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 276
tcagccatca aagcagctat tg

22

<210> 277
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 277
tcagcgcgta cagtgggtga t

21

<210> 278
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 278
tcagcgtatg cacatggaac tcctc

25

<210> 279
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 279
tcagctacat cgactatgcg atg

23

<210> 280

<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 280
tcagctagac cttttaggta aagctaagct 30

<210> 281
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 281
tcagctatatt ttccaggat ccaaggtgg 29

<210> 282
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 282
tcagctgtcg cagttcatgg acc 23

<210> 283
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 283
tcaggaag ggcattttac ccttg 25

<210> 284
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 284
tcaggagtcg ttcaactcga tctacatgat 30

<210> 285
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 285
tcaggagtcg ttcaactcga tctacatgat g 31
<210> 286
<211> 30
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 286
tcaggatgga aataaccacc aattcactac 30
<210> 287
<211> 23
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 287
tcaggcattg cgggtgggat ggc 23
<210> 288
<211> 24
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 288
tcaggtactg ctatccaccc tcaa 24
<210> 289
<211> 22
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 289
tcaggtggct tacacggcgt ag 22
<210> 290
<211> 26
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 290
tcagtatgta tccacogtag ccagtc 26
<210> 291

<211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 291
 tcagttccgt tatcgccatt gca 23

<210> 292
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 292
 tcagttccgt tatcgccatt gcat 24

<210> 293
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 293
 tcagttccgt tatcgccatt gcatt 25

<210> 294
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 294
 tcagttcggc ggtcagcgct tcgg 24

<210> 295
 <211> 34
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 295
 tcagtttttaa tgtctcgat gatcgaatca aaag 34

<210> 296
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 296
 tcattccacac ggtgggtgag aagg 24

 <210> 297
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 297
 tcattcctaag ccaagtgtag actctgta 26

 <210> 298
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 298
 tcattgataat atctttgaaa tcggctcagg a 31

 <210> 299
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 299
 tcattgttgag cttaaaccta tagaagtaaa agc 33

 <210> 300
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 300
 tcattatcat gcgccaatga gtgcaga 27

 <210> 301
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 301
 tcattcaaga actagatctt caggcaag 28

 <210> 302
 <211> 26

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 302
tccaaaaaaa tcagcgcgta cagtgg 26

<210> 303
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 303
tccaaaccag gtgtatcaag aacatcagg 29

<210> 304
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 304
tccaaataag tggcggttaca aatactgaa 29

<210> 305
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 305
tccaacgaag tacaatacaa gacaaaagaa gg 32

<210> 306
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 306
tccaaggtac actaaactta cttgagctaa tg 32

<210> 307
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 307
 tccaatgccca caaactcgtg aaca 24

 <210> 308
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 308
 tccacacgcc gttcttcaac aact 24

 <210> 309
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 309
 tccacacggt ggtggtgaag g 21

 <210> 310
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 310
 tccaccaaga gcaagatcaa ataggc 26

 <210> 311
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 311
 tccacggtca tggagcgcta 20

 <210> 312
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 312
 tccacttato gcaaatggaa aattaagcaa 30

 <210> 313
 <211> 33
 <212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 313 33
tccagatgga caaatcttct tagaaactga ttt

<210> 314
<211> 26
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 314 26
tccagcaoga attgctgcta tgaaag

<210> 315
<211> 33
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 315 33
tccaggacaa atgtatgaaa aatgtccaag aag

<210> 316
<211> 25
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 316 25
tccattgttc gtatgggtca agact

<210> 317
<211> 31
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 317 31
tcccaattaa ttctgccatt ttccaggta t

<210> 318
<211> 34
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 318

tccccagaaa	cagatgaaga	aattaacaaa	aaag	34
<210>	319			
<211>	30			
<212>	DNA			
<213>	Artificial Sequence			
<220>				
<223>	Primer			
<400>	319			
tcccagctag	accttttagg	taaagctaag		30
<210>	320			
<211>	26			
<212>	DNA			
<213>	Artificial Sequence			
<220>				
<223>	Primer			
<400>	320			
tcccaggtga	cgatgtacct	gtaatc		26
<210>	321			
<211>	26			
<212>	DNA			
<213>	Artificial Sequence			
<220>				
<223>	Primer			
<400>	321			
tccccaggac	accctgaaat	ttcaac		26
<210>	322			
<211>	34			
<212>	DNA			
<213>	Artificial Sequence			
<220>				
<223>	Primer			
<400>	322			
tccccacgc	tttaattggt	tatgatgatt	tgag	34
<210>	323			
<211>	31			
<212>	DNA			
<213>	Artificial Sequence			
<220>				
<223>	Primer			
<400>	323			
tcccggactt	aatatcaatg	aaaattgtgg	a	31
<210>	324			
<211>	28			
<212>	DNA			

<213> Artificial Sequence

<220>

<223> Primer

<400> 324

tcccgaggact ttatgacta aagcagat

28

<210> 325

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 325

tcccgaggagt tgactgggt

19

<210> 326

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 326

tccgctgaat ctgtcgccgc

20

<210> 327

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 327

tcgggtcac gttattatgg tac

23

<210> 328

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 328

tccgtacgta ttattagtg ctggtca

27

<210> 329

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 329

tccgttatcg ccattgcatt atttgggaact 30
 <210> 330
 <211> 33
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 330
 tccgttcttta caaatagcaa tagaacttga agc 33
 <210> 331
 <211> 36
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 331
 tccgttgatt attgttatcc tgttatgccca tttgag 36
 <210> 332
 <211> 34
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 332
 tcctaattgga cttaatatca atgaaaattg tgga 34
 <210> 333
 <211> 22
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 333
 tcctagagga atgggtgccca cg 22
 <210> 334
 <211> 30
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 334
 tcctatatatta atcgtttaca gaaactgggt 30
 <210> 335
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <220>
 <221> modified_base
 <222> (15)...(15)
 <223> I

 <400> 335
 tcctcaatga acgatcaaca agtgattgat g 31

 <210> 336
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 336
 tcctcaatga atgatcaaca agtgattgat g 31

 <210> 337
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <220>
 <221> modified_base
 <222> (15)...(15)
 <223> I

 <220>
 <221> modified_base
 <222> (24)...(24)
 <223> I

 <400> 337
 tcctogatga acgatcaaca agtgattgat g 31

 <210> 338
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <220>
 <221> modified_base
 <222> (15)...(15)
 <223> I

 <220>
 <221> modified_base
 <222> (24)...(24)
 <223> I

<400> 338
tcctcgatga atgatcaaca agtgattgat g

31

<210> 339
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

220>
<221> modified_base
<222> (6)...(6)
<223> I

<220>
<221> modified_base
<222> (15)...(15)
<223> I

<400> 339
tcctcgatga acgatcaaca agtgattgat g

31

<210> 340
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (6)...(6)
<223> I

<220>
<221> modified_base
<222> (12)...(12)
<223> I

220>
<221> modified_base
<222> (15)...(15)
<223> I

<220>
<221> modified_base
<222> (24)...(24)
<223> I

<400> 340
tcctcgatga acgatcaaca agtgattgat g

31

<210> 341
<211> 20

<212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 341
 tcctgaaaaa tggagcacgg 20

 <210> 342
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 342
 tcctgaagca agtgcattta cga 23

 <210> 343
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 343
 tcctgaccga cccattattc cctttatc 28

 <210> 344
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 344
 tcctgatgct caaagtgcct ttttagatcc ttt 33

 <210> 345
 <211> 34
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 345
 tcctgttacc cctgaagtag ttaatcaagt ttgt 34

 <210> 346
 <211> 35
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 346
 tcctgttatt cctgaagtag ttaatcaagt ttgtt 35

 <210> 347
 <211> 36
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 347
 tcctgttatt cctgaagtag ttaatcaagt ttgtta 36

 <210> 348
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 348
 tccttacttc gaactatgaa tcttttggaa g 31

 <210> 349
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 349
 tccttatagg gatggctatc agtaatgtt 29

 <210> 350
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 350
 tccttgaccg cctttccgat ac 22

 <210> 351
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 351
 tccttgcttt agttttaagt gcatgtaatt caa 33

 <210> 352
 <211> 32
 <212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 352

tcctttgata tattatgcga tggaaagggtg gt

32

<210> 353

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 353

tcctttgatg catgtaattg ctgcaaaagc

30

<210> 354

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 354

tcgaaagctt ttgcatatta tatcgagcca c

31

<210> 355

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 355

tcgaagtaca atacaagaca aaagaagg

28

<210> 356

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 356

tcgacaacac cattatctat ggtgtgaa

28

<210> 357

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 357

tcgacctttg gcaggaacta gac 23

<210> 358
 <211> 17
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 358
 tcgagcaggc gctgccg 17

<210> 359
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 359
 tcgagtatag ctaaaaaaat agtttatgac a 31

<210> 360
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 360
 tcgatctggt ttcagtctgt ttcagt 26

<210> 361
 <211> 28
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 361
 tcgatgaacg accaacaagt gattgatg 28

<210> 362
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 362
 tcgattagcg agcaacgaaa gcgc 24

<210> 363
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer
 <400> 363
 tcgcaaaaaa atccagctat tagc 24
 <210> 364
 <211> 26
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 364
 tcgccaatca aaactaaggg aatggc 26
 <210> 365
 <211> 19
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 365
 tcgcatcggt caccaaccg 19
 <210> 366
 <211> 16
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 366
 tcgcccgcga ggacgt 16
 <210> 367
 <211> 21
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 367
 tcgccgactt cgacggtgac c 21
 <210> 368
 <211> 22
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 368
 tcgccggcaa tgccattgga ta 22

<210> 369
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 369
 tcgcccgtgga aaaatcctac gct 23

<210> 370
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 370
 tcgcggtgca acaaaaacttt ctaaagtatg t 31

<210> 371
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 371
 tcgctacagg cccttttaga caag 24

<210> 372
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 372
 tcgctatctt atcggttgaga agggatt 27

<210> 373
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 373
 tcggaatctg atgttgacgt tggt 24

<210> 374
 <211> 22
 <212> DNA
 <213> Artificial Sequence

<220>
<223> Primer

<400> 374
tcggcgcac cttcatcgaa gt 22

<210> 375
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 375
tcggcgaaat ccgatttcct gaaaatga 28

<210> 376
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 376
tcggcgcccg cctgatga 18

<210> 377
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 377
tcgggtgatg atgcgcgtga agg 23

<210> 378
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 378
tcgggttagt aaaagaacgt attgctcaac c 31

<210> 379
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 379
tcgtacgtat tattaggtgc tggtcact 28

<210> 380
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 380
tcgtatggct caatggtgga g 21

<210> 381
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 381
tcgtattttt gattttttcc ctgataatgc 30

<210> 382
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 382
tcgtattttt gattttttcc ctgataatgc tc 32

<210> 383
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 383
tcgtgattat ggatggcaac gtgaa 25

<210> 384
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 384
tcgtgcccgc aatttgcata aagc 24

<210> 385
<211> 21
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 385

tcgtggcggc gtggttatcg a

21

<210> 386

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 386

tcgtgttgaa cggtgtcaaa tcaaagt

27

<210> 387

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 387

tcgttccttg aacacgatga cgc

23

<210> 388

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 388

tcgttttggtg gtggtagatg aaaaagg

27

<210> 389

<211> 11

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 389

tccaccctca a

11

<210> 390

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 390

tctaaaacac caggtcaccc agaag

25

<210> 391

<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 391
tctaaatggt cgtgcagttg cgtg 24

<210> 392
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 392
tctactgatt ttggtaatct tgcagcacag 30

<210> 393
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 393
tctagtaata ataggaccct cagc 24

<210> 394
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 394
tctcaagggtg atattggtgt aggtaactta a 31

<210> 395
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 395
totcattacg ttgcatogga aaca 24

<210> 396
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 396
tctcgtatgaa cgaccaacaa gtgattgatg 30

<210> 397
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 397
tctcgtggtg cacaagtaac ggatatta 28

<210> 398
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 398
tctgaaatga atagtgatag aactgtaggc ac 32

<210> 399
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 399
tctgaacatg ataatatctt tgaaatcggc tc 32

<210> 400
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 400
tctgaatgtc tatatggagg tacaacacta 30

<210> 401
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 401
tctgacacct gcccggtgc 19

<210> 402
<211> 20

<212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 402
 tctgcccggtg tcgttggtga 20

<210> 403
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 403
 tctggaggca caccaataa aaca 24

<210> 404
 <211> 21
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 404
 tctggataac ggctgcgcg g 21

<210> 405
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 405
 tctggcaggt atgcgtggtc tgatg 25

<210> 406
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 406
 tctggctaata actttggcaa cggt 24

<210> 407
 <211> 29
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 407 tctggtccaa caaaaggaac gattacagg	29
<210> 408	
<211> 25	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 408 tctgtcccta gtaacgagagg acogg	25
<210> 409	
<211> 23	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 409 tctgttccta gtaacgagagg acc	23
<210> 410	
<211> 26	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 410 tcttatgccca agaggacaga gtgagt	26
<210> 411	
<211> 29	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 411 tcttatgccca agaggacaga gtgagtact	29
<210> 412	
<211> 31	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 412 tcttattcca acttcaaacc gaactatgac g	31
<210> 413	
<211> 29	
<212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 413

tcttctcatc ctatggctat tatgcttgc

29

<210> 414

<211> 33

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 414

tcttgatact tgtaatgtgg gcgataaata tgt

33

<210> 415

<211> 32

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 415

tcttgcagca gtttatttga tgaacctaaa gt

32

<210> 416

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 416

tcttgctctt tcgtgagttc agtaaatg

28

<210> 417

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 417

tcttgtaact gtggctcaca cggctgttg g

31

<210> 418

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 418

tcttgtttat gctggtaaag cagatgg 27

<210> 419

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 419

tctttatggt ggagatgact gaaaccga 28

<210> 420

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 420

tctttcttga atgctggtgt acgtatcg 28

<210> 421

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 421

tctttgaaat cggctcagga aaagg 25

<210> 422

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 422

tctttgccat tgaagatgac ttaagc 26

<210> 423

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 423

tcttttaciaa aaggggaaaa agttgactt 29

<210> 424

<211> 34

<212> DNA

<213> Artificial Sequence

<220>
 <223> Primer

 <400> 424
 tgaaaaatgt ccaagaagca tagcaaaaaa agca 34

 <210> 425
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 425
 tgaaaagggt gaagtagcaa atggagatag 30

 <210> 426
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 426
 tgaaaagtat ggatttgaac aactcgtgaa ta 32

 <210> 427
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 427
 tgaaatctca ttacgttgca tcggaaa 27

 <210> 428
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 428
 tgaaattgct acaggccctt taggacaagg 30

 <210> 429
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 429
 tgaacgctgg tggcatgctt aacac 25

<210> 430
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 430
tgaacgtggt caaatcaaag ttggtgaaga 30

<210> 431
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 431
tgaacgtggt caaatcaaag ttggtgaaga a 31

<210> 432
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 432
tgaagcttgt tcttttagcag gacttca 27

<210> 433
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 433
tgaaggtgga cgtcacactc cattcttc 28

<210> 434
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 434
tgaagtagaa atgactgaac gtccga 26

<210> 435
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 435
tgaagtagaa ggtgcaaagc aagttaga 28

<210> 436
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 436
tgaagtgcgt gatgatatcg atgcacttga tgta 34

<210> 437
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 437
tgaatagtga tagaactgta ggcacaatcg t 31

<210> 438
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 438
tgaatgctta ttacctgca ctcccacaac t 31

<210> 439
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 439
tgaattagtt caatcatttg ttgaacgacg t 31

<210> 440
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 440
tgaattgctg ctatgaaagc tggctt 26

<210> 441
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 441
tgacagcgaa gaaggttaga cttgtcc 27

<210> 442
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 442
tgacatccgg ctcacgttat tatggt 26

<210> 443
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 443
tgacatccgg ctcacgttat tatggta 27

<210> 444
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 444
tgacatccgg ctcacgttat tatggtac 28

<210> 445
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 445
tgacatgata ataaccgatt gaccgaaga 29

<210> 446
<211> 22
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 446

tgacatgctt gtcggttcag gc

22

<210> 447

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 447

tgacatggac tccccctata taactcttga g

31

<210> 448

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 448

tgaccagggtg atggccatgt tcg

23

<210> 449

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 449

tgacctacag taagaggttc tgtaatgaac c

31

<210> 450

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 450

tgacgatctt cgcggtgact agt

23

<210> 451

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 451

tgacggccta tacggtgttg gtttct

26

<210> 452

<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 452
tgacgtcatc ggtaagtacc accc 24

<210> 453
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 453
tgagatggat ttaaacctgt tcaccg 27

<210> 454
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 454
tgagattgct gaacatttaa tgctgattga 30

<210> 455
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 455
tgagcaatgg ggctttgaaa gaatttttaa at 32

<210> 456
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 456
tgagctgcat caactgtatt ggatag 26

<210> 457
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 457
tgagctttta gttgactttt tcaacagc 28

<210> 458
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 458
tgaggaccgt gtcgcgctca 20

<210> 459
<211> 35
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 459
tgagggtttt atgcttaaag ttggttttat tggtt 35

<210> 460
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 460
tgagggtggtg gataactcaa ttgatgaagc 30

<210> 461
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 461
tgagtaacat ccatatttct gccatacgt 29

<210> 462
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 462
tgagtaagtt ccaccgcac gg 22

<210> 463
<211> 30

<212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 463
 tgagtcactt gaagttgata caaatcctct 30

 <210> 464
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 464
 tgagtgatga aggccttagg gttgtaaa 28

 <210> 465
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 465
 tgagtgccaa catatcagtg ctgaaga 27

 <210> 466
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 466
 tgagtttaac agttcaccat atgaaacagg 30

 <210> 467
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 467
 tgatacttca acgcctgctg ctttc 25

 <210> 468
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 468
tgatcactgg tgctgctcag atgga 25

<210> 469
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 469
tgatcatcgg tggataacg atttattagt 30

<210> 470
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 470
tgatcgttga gaagggattt gcgaaaaga 29

<210> 471
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 471
tgatctcaga atctaataat tgggacgaa 29

<210> 472
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 472
tgatcttaaa aatttcgcc aacttcattc 30

<210> 473
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 473
tgatgacttt ttagctaattg gtcaggcagc 30

<210> 474
<211> 29
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 474

tgatggcaag tggatagggt ataatacag

29

<210> 475

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 475

tgattaccat gagtggcaag caag

24

<210> 476

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 476

tgattattgt tatcctgtta tgccattga g

31

<210> 477

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 477

tgattccggt gcccgtaggt

19

<210> 478

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 478

tgattctggt gcccgtaggt

19

<210> 479

<211> 34

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 479

tgattttgct aaatttagag aaattgcgga tgaa 34

<210> 480
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 480
 tgcaaaatct gcaacgagct ttgg 24

<210> 481
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 481
 tgcaaggag gtactcagac cat 23

<210> 482
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 482
 tgcaagcaaa cgcacaatca gaagc 25

<210> 483
 <211> 20
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 483
 tgcaagcgcg accacatag 20

<210> 484
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 484
 tgcaagcttc tgggtctagc att 23

<210> 485
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 485
 tgcaagtggt acttcaacat gggg 24

 <210> 486
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 486
 tgcaagttaa gaaagctggt gcaggtttat 30

 <210> 487
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 487
 tgcaattgct ttagttttta gtgcatgtaa ttc 33

 <210> 488
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 488
 tgcacaatca gaagctaaga aagcgcaagc t 31

 <210> 489
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 489
 tgcacaagcc gttcttcaac aact 24

 <210> 490
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 490
 tgcacatcgt gtccaacgtc ac 22

<210> 491
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 491
 tgcacoggct attaagaatt actttgccaa ct 32

 <210> 492
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 492
 tgcacgatgc ggaatgggtc aca 23

 <210> 493
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 493
 tgcacgcoga ctatgttaag aacatgat 28

 <210> 494
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 494
 tgcacttata gcaaatggaa aattaagcaa 30

 <210> 495
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 495
 tgcagggaac agcttttaggc a 21

 <210> 496
 <211> 22
 <212> DNA
 <213> Artificial Sequence

<220>

<223> Primer

<400> 496

tgcatacaaa cagtcggagc ct

22

<210> 497

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 497

tgcataccgg taagttggca aca

23

<210> 498

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 498

tgcataattat atcgagccac agcatcg

27

<210> 499

<211> 33

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 499

tgcattattt ggaactattg caactgctaa tgc

33

<210> 500

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 500

tgccaagagg acagagtgg tactttga

28

<210> 501

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 501

tgccggacaa ttacgattca tcgagtatta a

31

<210> 502
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 502
tgccgtaatg ataggtgaag atatacaaag agt 33

<210> 503
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 503
tgccgtgttg aacgtggtca aat 23

<210> 504
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 504
tgccctagaag atcttaaaaa ttccgcctaa ctt 33

<210> 505
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 505
tgccctatctt ttgtctgata tagcacatat tgc 33

<210> 506
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 506
tgccctcgaag ctgaatataa ccaagtt 27

<210> 507
<211> 22
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 507
tgctgtagg gaatcctgct ga

22

<210> 508
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 508
tgctgtttct tagtacgaga ggacc

25

<210> 509
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 509
tgcgagctc ttggtatcga gtt

23

<210> 510
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 510
tgcgcggaag atgtaacggg

20

<210> 511
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 511
tgcggatcgt ttggtggttg tagatgaaaa

30

<210> 512
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 512
tgcgggtagg gagcttgagc

20

<210> 513

<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 513
tgcgtacaat acgctttatg aaattttaac a 31

<210> 514
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 514
tgcgtataaa aaacacagat ggcagca 27

<210> 515
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 515
tgcgtttacc gcaatgcgtg c 21

<210> 516
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 516
tgctacggta ggatctcctt atcctattg 29

<210> 517
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 517
tgctagtcaa tctatcattc cggttgatac 30

<210> 518
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 518
tgctagttat ggtacagagt ttgcgac 27

<210> 519
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 519
tgctatgggtg ttaccttccc tatgca 26

<210> 520
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 520
tgctcaaccc gatcctaaat tagacga 27

<210> 521
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 521
tgctcaatct aaacctaaag tcgaaga 27

<210> 522
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 522
tgctcgagtg attgactttg ctaaatttag aga 33

<210> 523
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 523
tgctcgtaag ggtctggcgg atac 24

<210> 524
<211> 29

<212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 524
 tgctcgtggt gcacaagtaa cggatatta 29

 <210> 525
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 525
 tgctgaggcc tggaccgatt atttac 26

 <210> 526
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 526
 tgctggtaac agagccttat aggcgca 27

 <210> 527
 <211> 19
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 527
 tgctggtacg ggtcgagga 19

 <210> 528
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 528
 tgctgggtgaa aataaccag atgtcgtcctt c 31

 <210> 529
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 529
 tgctgtagct tatcgcgaaa tgtctttgat tt 32

 <210> 530
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 530
 tgcttatatta cctgcactcc cacaactg 28

 <210> 531
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 531
 tgcttcagga atcaatgatg gagcag 26

 <210> 532
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 532
 tgcttcggat ccagcagcac ttcaata 27

 <210> 533
 <211> 19
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 533
 tgcttctggt gctagcatt 19

 <210> 534
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 534
 tgctttccta tggcttatcc aaatttagat cg 32

 <210> 535
 <211> 29
 <212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 535

tgcttttgat ggtgatgcag atcgtttgg 29

<210> 536

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 536

tggaagcca tgcgtctgac atct 24

<210> 537

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 537

tggaaggtg ttgcagctac tca 23

<210> 538

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 538

tggaatggc agctagaata gtagctaaaa t 31

<210> 539

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 539

tggaacaaaa tagtctctcg gattttgact 30

<210> 540

<211> 33

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 540

tggaacagga attaattctc atcctgatta tcc

33

<210> 541
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 541
tggaacgtta tcagggtgcc caaaaattcg

30

<210> 542
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 542
tggaactatt gcaactgcta atg

23

<210> 543
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 543
tggaacttga agctctcgct cttaaagatg

30

<210> 544
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 544
tggaagatct gggtcaggc

19

<210> 545
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 545
tggaagtcat caagcgcctg gc

22

<210> 546
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 546
tggataaaca aaacatgaag gaaaccactt 30

<210> 547
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 547
tggaatgatg ataagattt cgcagatagc ta 32

<210> 548
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 548
tggacaatag acaatcactt ggatttaca 29

<210> 549
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 549
tggacacata tcgtgagcaa tgaactga 28

<210> 550
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 550
tggacggcat cagattctc tac 23

<210> 551
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 551
tggactctc ggtggtgc 19

<210> 552
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 552
tggagcacgg cttctgac 19

<210> 553
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 553
tggagcttga agctatcgct cttaaagatg 30

<210> 554
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 554
tggagggtgc actccacacg aa 22

<210> 555
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 555
tggagggtgt tgtatgtatg gtggt 25

<210> 556
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 556
tggatatcca ccgaacacta gggttg 26

<210> 557
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 557
tggatggcat ggtgaaatgg atatgtc 27

<210> 558
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 558
tggatgggga ttagcgggta caatg 25

<210> 559
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 559
tggatgttaa ggggtatttt cccgaa 26

<210> 560
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 560
tggattagag accctggtag tcc 23

<210> 561
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 561
tggcacggcc atctccgtg 19

<210> 562
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 562
tggcactctt gcctttaata ttagtaaaact atca 34

<210> 563
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 563
 tggcagctag aatagtagct aaaatcccta c 31

 <210> 564
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 564
 tggcagttttt acaaggtgct gtttcato 28

 <210> 565
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 565
 tggcatttct tatgaagctt gttcttttagc a 31

 <210> 566
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 566
 tggccagcgc ttcggtgaaa tgga 24

 <210> 567
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 567
 tggcccgaaa gaagctgagc g 21

 <210> 568
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>

<223> Primer

<400> 568
 tggcctaatag ggcttaatat caatgaaaat tg 32

<210> 569
 <211> 22
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 569
 tggcgaaacct ggtgaacgaa gc 22

<210> 570
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 570
 tggcgagtggt atagggtata atacag 26

<210> 571
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 571
 tggcgtagta gagctattta cagacac 27

<210> 572
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 572
 tggcgaagtggt atagggtata atacag 26

<210> 573
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 573
 tggctccttg gtatgactct gcttc 25

<210> 574

<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 574
tggctgacat cctacatgac tgtga 25

<210> 575
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 575
tggcttatcc aaatttagat cgtgggtttta c 31

<210> 576
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 576
tgggacttga agctatcgct cttaaagatg 30

<210> 577
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 577
tgggatgaaa aagcggttctt ttatccatga 30

<210> 578
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 578
tgggattatt gttatcctgt tatgccattt gaga 34

<210> 579
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 579
tgggatttta aaaaacattg gtaacatcgc ag 32

<210> 580
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 580
tgggcaacag cagcggattg cgattgcgcg 30

<210> 581
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 581
tgggcagcgt ttcggcgaaa tgga 24

<210> 582
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 582
tgggcctaata gggcctaata tcaatgaaaa ttg 33

<210> 583
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 583
tgggcgatgc tgcgaaatgg ttaaaaga 28

<210> 584
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 584
tgggcgtgag caatgaactg attatac 27

<210> 585
<211> 18

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 585
tgggctgtgga acgtccac 18

<210> 586
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 586
tgggctcttt ctgccttaaa cacc 24

<210> 587
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 587
tgggctcttt ctgccttaaa cacct 25

<210> 588
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 588
tggggattca gccatcaaag cagctattga c 31

<210> 589
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 589
tggggattga tatcaccgat aagaagaa 28

<210> 590
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 590
tggggccttta aatattccaa ttgaagattt tca 33

<210> 591
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 591
tggggccttg ctttatagtt ttttacattt aag 33

<210> 592
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

220>
<221> modified_base
<222> (5)...(5)
<223> I

<220>
<221> modified_base
<222> (14)...(14)
<223> I

<400> 592
tgggtgatgc tgctaaatgg ttaaaaga 28

<210> 593
<211> 35
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 593
tgggtcgtgg ttttacagaa aatttcttat atatg 35

<210> 594
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 594
tgggtgacat tcatcaattt catcgttc 28

<210> 595
<211> 32
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 595

tggggtttaca catatcgtga gcaatgaact ga

32

<210> 596

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 596

tggttaatac caccacatg gtgac

25

<210> 597

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 597

tggtaacaga gccttatagg cgca

24

<210> 598

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 598

tggtaacaga gccttatagg cgcatatg

28

<210> 599

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 599

tggttaagagc gcaccggtaa gttggttaaca

30

<210> 600

<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 600

tggtacagag tttgcgac	18
<210> 601	
<211> 26	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 601	
tggtacatgt gccttcattg atgctg	26
<210> 602	
<211> 18	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 602	
tggtacagag tttgcgac	18
<210> 603	
<211> 24	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 603	
tggtactcac ttagcgggtt tccg	24
<210> 604	
<211> 24	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 604	
tggtatgata tgatgcctgc acca	24
<210> 605	
<211> 21	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 605	
tggtatgcgt ggtctgatgg c	21
<210> 606	
<211> 31	
<212> DNA	
<213> Artificial Sequence	

<220>
<223> Primer

<400> 606 31
tggtattcta ttttgcctgat aatgacctcg c

<210> 607
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 607 25
tggtcaaatac aaagttggtg aagaa

<210> 608
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 608 29
tggtcttatg ccaagaggac agagtgagt

<210> 609
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 609 24
tggtgactcg gcatgttatg aagc

<210> 610
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 610 30
tggtgacttc ataattgatg aagttgaagt

<210> 611
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 611 27
tggtgcgagt gcttatgctc gtattat

<210> 612
 <211> 13
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 612
 tgggtgctagc att 13

<210> 613
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 613
 tgggtgctttc tggcgcttaa acga 24

<210> 614
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 614
 tgggtggacat ttaacacatg gtgcaaa 27

<210> 615
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 615
 tgggtggtgaa atagatagga ctgctt 26

<210> 616
 <211> 13
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 616
 tgggtgctagc att 13

<210> 617
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 617
 tgggttatgc tcaggcgaac tccaac 26

 <210> 618
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 618
 tggttatgta ccaaatactt tgtctgaaga tgg 33

 <210> 619
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 619
 tggtttagat aattccttag gatctatgcy t 31

 <210> 620
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <220>
 <221> modified_base
 <222> (3)...(3)
 <223> I

 <220>
 <221> modified_base
 <222> (6)...(6)
 <223> I

 <220>
 <221> modified_base
 <222> (9)...(9)
 <223> I

 <220>
 <221> modified_base
 <222> (12)...(13)
 <223> I

 <220>
 <221> modified_base
 <222> (15)...(15)
 <223> I

<400> 620 tgtcctactg tttgtggttc tgtaatgaac c	31
<210> 621	
<211> 24	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 621 tgtaactatc acccgacgg tgat	24
<210> 622	
<211> 30	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 622 tgtaagctct acaaccaca aaaccttacg	30
<210> 623	
<211> 30	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 623 tgtaatgaac cctaatagacc atccacacgg	30
<210> 624	
<211> 31	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 624 tgtaccgct gaattaacga atttatacga c	31
<210> 625	
<211> 24	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 625 tgtaactcggg aagatcaccc cgca	24
<210> 626	
<211> 21	
<212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 626

tgtactgcta tccacctca a

21

<210> 627

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 627

tgtagccgct aagcactacc atcc

24

<210> 628

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 628

tgtagcttat cgcgaaatgt ctttgatttt

30

<210> 629

<211> 33

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 629

tgtatgggtg tgtaacgtta catgataata atc

33

<210> 630

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 630

tgtattaggg gcatacagtc ctcatac

27

<210> 631

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 631

tgtcaaagtg gcacgtttac tggc 24

<210> 632
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 632
tgtcatgggt aaatatcacc ctca 24

<210> 633
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 633
tgtccaagaa gcatagcaaa aaaagcaa 28

<210> 634
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 634
tgtcgatgca acgcaagaa cct 23

<210> 635
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 635
tgtcgggtaca cgatattctt cacga 25

<210> 636
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 636
tgtgaataaa tcacgattga ttgagca 27

<210> 637
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 637
tgtggagtaa cactgcatga aaacaa

26

<210> 638
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 638
tgtgggtcaaa tcaaagttgg tgaagaa

27

<210> 639
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 639
tgttcaagag ctatgatcttc aggca

25

<210> 640
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 640
tgttcaagag ctatgatcttc aggcaa

26

<210> 641
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 641
tgttcgctgt ttcacaaaca acattcca

28

<210> 642
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 642
tgttcttttag caggacttca caaacttgat aa

32

<210> 643
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 643
 tgttgaacgt ggtcaaatca aagttggtg 29

 <210> 644
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 644
 tgttgggagt attccttacc atttaagcac a 31

 <210> 645
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 645
 tgttgggtgct ttctggcgct taa 23

 <210> 646
 <211> 38
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 646
 tccttgttgt cctactgttt gtggttctgt aatgaacc 38

 <210> 647
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <220>
 <221> modified_base
 <222> (7)...(7)
 <223> I

 <220>
 <221> modified_base
 <222> (10)...(10)

<223> I

<220>
<221> modified_base
<222> (13)...(13)
<223> I

<220>
<221> modified_base
<222> (16)...(16)
<223> I

<220>
<221> modified_base
<222> (19)...(20)
<223> I

<220>
<221> modified_base
<222> (22)...(22)
<223> I

<400> 647
ttaaagtgtg ttttattggt tggcgcgga 29

<210> 648
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 648
ttaacatgaa ggaaccact ttgataatgg 30

<210> 649
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 649
ttaacgggta tcatggccca gatggg 26

<210> 650
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 650
ttaagtcccg caacgagcgc aa 22

<210> 651
<211> 22
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 651

ttaagtcccg caacgatcgc aa

22

<210> 652

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 652

ttaatttgcc aaaaatgcaa ccaggtag

28

<210> 653

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 653

ttacacatat cgtgagcaat gaactga

27

<210> 654

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 654

ttacaggaag tttaggtggt aatctaaaag g

31

<210> 655

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 655

ttactccatt attgcttggt tacactttcc

30

<210> 656

<211> 35

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 656

ttataactta ctgcaatcta ttcagttgct tgggtg 35

<210> 657
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 657
ttataccgga aacttccga aaggag 26

<210> 658
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 658
ttatcagcta gaccttttag gtaaagctaa gc 32

<210> 659
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 659
ttatcgctca ggcgaaactcc aac 23

<210> 660
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 660
ttatcgtttg tggagctagt gcttatgc 28

<210> 661
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 661
ttatgaagcg tggtcttttag caggacttca 30

<210> 662
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
 <223> Primer

 <400> 662
 ttatggatgg caacgtgaaa cgcggt 25

 <210> 663
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 663
 ttattgttat cctgttatgc c 21

 <210> 664
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 664
 ttattttacct gcactccac aactg 25

 <210> 665
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 665
 ttcaaaaact ccaggccatc ctgaaatttc aac 33

 <210> 666
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 666
 ttcaacaggt accaatgatt tgatctca 28

 <210> 667
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 667
 ttccccaccga tatcatgggt taccacgg 28

<210> 668
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 668
 ttccgtaagt cggctaaac agtcg 25

 <210> 669
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 669
 ttctctcttt tgaaagcgac gggt 24

 <210> 670
 <211> 17
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 670
 ttctctgggc gcctggc 17

 <210> 671
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 671
 ttctgacgc acccattatt ccctttatc 29

 <210> 672
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 672
 ttcgatgcaa cgcaagaac ct 22

 <210> 673
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
<223> Primer

<400> 673
ttcgccaatc aaaactaagg gaatggc 27

<210> 674
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 674
ttcgcggtc agcgcttcgg 20

<210> 675
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 675
ttctaaaaca ccaggtcacc cagaag 26

<210> 676
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 676
ttctatctcg ttggtttatt cggagtt 27

<210> 677
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 677
ttctgaatgt ctatatggag gtacaacact 30

<210> 678
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 678
ttgactgcc aggtcacgt g 21

<210> 679
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 679
 ttgactgagg cacaacacgg at 22

 <210> 680
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 680
 ttgagaagac atccggctca cggtattatg gta 33

 <210> 681
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 681
 ttgagggtat gcaccgtctt ttgtattctt t 31

 <210> 682
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 682
 ttgcaactgc tgatttagct caga 24

 <210> 683
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 683
 ttgcacaagc aaggcgctat tt 22

 <210> 684
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>

<223> Primer

<400> 684
 ttgccaatga tattcgttgg ttagcaag 28

<210> 685
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 685
 ttgccgcggtg tgcggaagta accgatatta c 31

<210> 686
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 686
 ttgcgaatag aacgatggct cgt 23

<210> 687
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 687
 ttgctcgtgg tgcacaagta acggatatta 30

<210> 688
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 688
 ttgctcgtgg tgcacaagta acggatatta c 31

<210> 689
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<220>
 <221> modified_base
 <222> (15)...(15)
 <223> I

<220>
<221> modified_base
<222> (29)...(29)
<223> I

<400> 689
ttgctctgttg tgcacaagta acggatatta c 31

<210> 690
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 690
ttgcttaaaag ttggttttat tgggtggcg 29

<210> 691
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 691
ttggtccttt ttatacgaaa gaagaagttg aa 32

<210> 692
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 692
ttgtaaatgc cgggtcctca gatcc 25

<210> 693
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 693
ttgtacacac cgcccgctcat ac 22

<210> 694
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 694
ttgttagcaca gcaaggcaaa tttcttgaaa c 31

<210> 695
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 695
ttgtatgtat ggtggtgtaa cgttacatga 30

<210> 696
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 696
ttgtatgtat ggtggtgtaa ctgagca 27

<210> 697
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 697
tttaagtccc gcaacgagcg caa 23

<210> 698
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 698
tttacacata tcgtgagcaa tgaactga 28

<210> 699
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 699
tttacactac ttttattcat tggcctaacg 30

<210> 700
<211> 19

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 700
tttacagctt tatgcaccg 19

<210> 701
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 701
tttcacacag cgtgtttata gttctacca 29

<210> 702
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 702
tttcacatgt aattttgata ttogcactga 30

<210> 703
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 703
tttcatttta tcaggagccc gaaatcga 28

<210> 704
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 704
tttctctcctt ttgaaagcga cgggtt 25

<210> 705
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 705 tttcgaaggg cctttcgacc tg	22
<210> 706	
<211> 23	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 706 tttogatgca acgcgaagaa cct	23
<210> 707	
<211> 27	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 707 tttgatttta cgcggtcctc caggctg	27
<210> 708	
<211> 32	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 708 tttgcggtg aagtaggtgc ctatcttttt gc	32
<210> 709	
<211> 32	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 709 ttttatgctt aaagttgggtt ttattgggtg gc	32
<210> 710	
<211> 24	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 710 ttttgaaggt gatccgtgcc aacg	24
<210> 711	
<211> 29	
<212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 711
aaactatattt tttagctata ctggaacac 29

<210> 712
<211> 19
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 712
aacatagcct tctcgtcc 19

<210> 713
<211> 21
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 713
aaattcgctt tcggtcatgt t 21

<210> 714
<211> 17
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 714
aaggaggtga tccagcc 17

<210> 715
<211> 29
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 715
aatcgacgac catcttgga agatttctc 29

<210> 716
<211> 21
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 716

acaaaaggca cgccatcacc c 21

<210> 717
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 717
acaaaaggta cgccgtcacc c 21

<210> 718
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (12)...(12)
<223> I

<400> 718
acaacacgag ctgacgac 18

<210> 719
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 719
acaacacgag ctgacgac 18

<210> 720
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (16)...(16)
<223> I

<400> 720
acaacacgag ctgacgac 18

<210> 721
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (14)...(14)
<223> I

<400> 721
acaacacgag ctgacgac

<210> 722
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (11)...(11)
<223> I

<400> 722
acaacacgag ctgacgac

<210> 723
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (10)...(10)
<223> I

<220>
<221> modified_base
<222> (13)...(13)
<223> I

<400> 723
acaacacgag ctgacgac

<210> 724
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (8)...(8)
<223> I

18

18

18


```

<220>
<221> modified_base
<222> (10)...(10)
<223> I

<220>
<221> modified_base
<222> (13)...(13)
<223> I

<400> 724
acaacacgag ctgacgac

<210> 725
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (8)...(8)
<223> I

<220>
<221> modified_base
<222> (10)...(10)
<223> I

<220>
<221> modified_base
<222> (13)...(13)
<223> I

<220>
<221> modified_base
<222> (16)...(16)
<223> I

<400> 725
acaacacgag ctgacgac

<210> 726
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 726
acaaccatgc accacatgtc

<210> 727
<211> 12
<212> DNA
<213> Artificial Sequence

```

18

18

20

<220>
 <223> Primer

 <400> 727
 acacgagctg ac 12

 <210> 728
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 728
 accactttta ataaggtttg tagctaac 28

 <210> 729
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 729
 acctgcaata tctaatagcac tcttacg 27

 <210> 730
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 730
 acctgcatcc ctaaactgac ttgc 24

 <210> 731
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 731
 accttggtac gacttcaccc ca 22

 <210> 732
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 732
 acgaactgga tgtcgccgtt 20

<210> 733
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 733
acgacacgag ctgacgac 18

<210> 734
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 734
acgacacgag ctgacgac 18

<210> 735
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 735
acgagctgac gacagccatg 20

<210> 736
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 736
acgccatcag gccacgcat 19

<210> 737
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 737
acgcgggcat gcagagatgc c 21

<210> 738
<211> 18
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 738
acggcacgag gtagtcgc 18

<210> 739
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 739
acggttacct tgttacgact 20

<210> 740
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 740
acgtccttca tcgcctctga 20

<210> 741
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 741
acgttttttcg ttttgaacga taatgct 27

<210> 742
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 742
actgctgcct cccgtag 17

<210> 743
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 743
acttagatgc tttcagcggc 20

<210> 744

<211> 20
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 744
 agacctcctg cgtgcaaagc 20

<210> 745
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 745
 agataaagaa tcacgaatat caattgttag c 31

<210> 746
 <211> 22
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 746
 agccgacatc gaggtgccaa ac 22

<210> 747
 <211> 28
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 747
 agctgctaga tgagcttctg ccatggcc 28

<210> 748
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 748
 aggatagatt tattttotgt tcg 23

<210> 749
 <211> 20
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 749
agtcacatccc ggtcctctcg 20

<210> 750
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 750
ataagccatg ttctgttcca tc 22

<210> 751
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 751
ataagccggg ttctgtcg 18

<210> 752
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 752
atatgattat cattgaactg cggccg 26

<210> 753
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 753
atcccctgct tctgctgcc 19

<210> 754
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 754
attcaagagc cattttctttt ggtaaaccac 30

<210> 755
<211> 22

<212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 755
 attgccacaga aatcaaatca tc 22

 <210> 756
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 756
 attgcttctt acttgcttag cataaatctt cca 33

 <210> 757
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 757
 attgtagcac gtgtgtagcc c 21

 <210> 758
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 758
 caagcggttt gcctcaaata gtca 24

 <210> 759
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 759
 caatctgctg acggatctga gc 22

 <210> 760
 <211> 15
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 760 caccgggcag gcgtc	15
<210> 761 <211> 20 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 761 cacggctacc ttgttagc	20
<210> 762 <211> 20 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 762 cagataaaga atcgctccag	20
<210> 763 <211> 25 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 763 catgacagcc aagacctcac ccacc	25
<210> 764 <211> 18 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 764 catgatggtc acaaccgg	18
<210> 765 <211> 20 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 765 ccaaacaccg ccgtcgatat	20
<210> 766 <211> 26 <212> DNA	

<213> Artificial Sequence
 <220>
 <223> Primer
 <400> 766
 ccaacgtttt ccacaacaga atcagc 26
 <210> 767
 <211> 27
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 767
 ccaagtgtg gtttaccoca tggagta 27
 <210> 768
 <211> 23
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 768
 ccacttttaa taaggtttgt agc 23
 <210> 769
 <211> 27
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 769
 ccagcagtta ctgtccctc atctttg 27
 <210> 770
 <211> 28
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 770
 ccataagggtc accgtcacca ttcaaagc 28
 <210> 771
 <211> 18
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 771

cctatgcagca cctgtctc 18

<210> 772
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 772
cccatttttt cagcatgct gaaaatctc 29

<210> 773
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 773
ccccogtcaa ttcotttgag t 21

<210> 774
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 774
ccctgtagta gaagaggtaa ccac 24

<210> 775
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 775
ccgacaagga atttcgctac c 21

<210> 776
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 776
ccgcggtcga attgcatgcc ttc 23

<210> 777
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 777
ccggtcctct cgtacta 17

<210> 778
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 778
ccgtgctcca tttttcag 18

<210> 779
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 779
cctacccaac gttcaccaag ggcag 25

<210> 780
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 780
cctcctgcgt gcaaagc 17

<210> 781
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 781
cctgtagtag aagaggtaac 20

<210> 782
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 782
ccttctccg aagttacg 18

<210> 783
 <211> 20
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 783
 ccttggttacg acttcacccc

20

<210> 784
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 784
 cgaacggcca gagtagtcaa cacg

24

<210> 785
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 785
 cgaacggcct gagtagtcaa cacg

24

<210> 786
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 786
 cgacttgacg gttaacatctt cctg

24

<210> 787
 <211> 21
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 787
 cgagttgcag actgcgatcc g

21

<210> 788
 <211> 21
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 788
 cgagttgcag actgcgatcc g 21

 <210> 789
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 789
 cgcaccatgc gtagagatga agtac 25

 <210> 790
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 790
 cgcaccgtgg gttgagatga agtac 25

 <210> 791
 <211> 18
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 791
 cgcatttcac cgctacac 18

 <210> 792
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 792
 cgcggtcggc tcggtgatga 20

 <210> 793
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 793
 cggctgctgg cacgaagtta g 21

<210> 794
<211> 15
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 794
cggcttcaag acccc 15

<210> 795
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 795
cggtaacgaac tggatgtcgc cggt 24

<210> 796
<211> 15
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 796
cgtaactccc aggcg 15

<210> 797
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 797
cgtataagct gcaccataag cttgtaatgc 30

<210> 798
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 798
cgtggactac cagggtatct a 21

<210> 799
<211> 21
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 799
ctatcgggtca gtcaggagta t 21

<210> 800
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 800
cttctacatt tttagccatc ac 22

<210> 801
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 801
ctttacgccc agtaattccg 20

<210> 802
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 802
ctttcgcttt ctcgaaactca accat 25

<210> 803
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 803
gaatatcaat ttgttagc 17

<210> 804
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 804
gaccccaacc tggccttttg tcgttga 27

<210> 805

<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 805
gaccgttata gttacggcc 19

<210> 806
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 806
gacgggcggt gtgtacaag 19

<210> 807
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 807
gacgggcggt gtgtacaag 19

<210> 808
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 808
gacgggcggt gtgtacaag 19

<210> 809
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 809
gacgtcatcc ccaccttcct c 21

<210> 810
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 810
gacgtcatcc ccacottcct cc 22
<210> 811
<211> 19
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 811
gagcatcagc gtgcgtgct 19
<210> 812
<211> 25
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 812
gagctgcgcc aacgaataaa tcgtc 25
<210> 813
<211> 29
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 813
gattggcgat aaagtgatat tttctaaaa 29
<210> 814
<211> 26
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 814
gcccaccaga aagactagca ggataa 26
<210> 815
<211> 21
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 815
gccgtccatc tgagcagcac c 21
<210> 816
<211> 21

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 816
gccgttcatt tgagcagcac c 21

<210> 817
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 817
gccttgcgac cgtactccc 19

<210> 818
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 818
gcgaccgtac tccccagg 18

<210> 819
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 819
gcgctccacg tcttcacgc 19

<210> 820
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 820
gcgtgacagg caggatttc 19

<210> 821
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 821
gcgtgacgac cttcttgaat tgtaatca 28
<210> 822
<211> 24
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 822
gcgttccaca gcttgttga gaag 24
<210> 823
<211> 18
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 823
gctgctggca cggagtta 18
<210> 824
<211> 26
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 824
gctgctttga tggtgaac cccttc 26
<210> 825
<211> 21
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 825
gctggattcg cctttgtac g 21
<210> 826
<211> 20
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 826
gcttacacac ccggcctatc 20
<210> 827
<211> 23
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 827

ggaatttacc agcgatagac acc

23

<210> 828

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 828

ggataattgg tcgtaacaag ggatagtgag

30

<210> 829

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 829

ggcatcacca ttctcttgct cttcg

25

<210> 830

<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 830

ggcgcgtactc cccagcgc

18

<210> 831

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 831

ggcgcgttgta cttaccgcac

20

<210> 832

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 832

gggtctacac ctgcacttgc ataac	25
<210> 833	
<211> 16	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 833	
gggtttccccc attcgg	16
<210> 834	
<211> 20	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 834	
ggtaaccctt gtctttgaat	20
<210> 835	
<211> 19	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 835	
ggtaagggtt ttgcggtg	19
<210> 836	
<211> 28	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 836	
ggtataacgc atgcagcaa aagattta	28
<210> 837	
<211> 28	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 837	
gtaacccttg tctttgaatt gtatttgc	28
<210> 838	
<211> 22	
<212> DNA	
<213> Artificial Sequence	

<220>
 <223> Primer

 <400> 838
 gtaagccatg ttttgttcca tc 22

 <210> 839
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 839
 gtatcctaac ctgtttgctc cc 22

 <210> 840
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 840
 gtccgacttg acggtcaaca ttctctg 27

 <210> 841
 <211> 18
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 841
 gtgcgccctt tctaactt 18

 <210> 842
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 842
 gtgctgggtt accccatgga gt 22

 <210> 843
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 843
 gttcaaatgc ctggataccc a 21

<210> 844
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 844
 gttgtcacca ggcattacca ttcc

24

<210> 845
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 845
 gttgtcgcca ggcataacca ttcc

24

<210> 846
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 846
 gtttcacgct tagatgcttt cagc

24

<210> 847
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 847
 gtttttcggt gcgtacgatg atgtc

25

<210> 848
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 848
 taaacgtccg ataccaatgg ttcgctc

27

<210> 849
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 849
 taaactatatt ttttagctat actcgaacac 30

 <210> 850
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 850
 taaagacacc gctgggttta aatgtgca 28

 <210> 851
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 851
 taaagagacg tttggtagtt catttgc 27

 <210> 852
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 852
 taaaggatag cggtaactaa atggtgagc cat 33

 <210> 853
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 853
 taaatgcact tgcttcaggg ccatat 26

 <210> 854
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 854
 taaattccgc aaagactttg gcattaggtg t 31

<210> 855
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 855
taacaaatcc cgtctgagtt cctcttgca 29

<210> 856
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 856
taacaaagtt accttcgcga tccactaa 28

<210> 857
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 857
taaccatttc gogtaagatt caa 23

<210> 858
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 858
taaccacccc aagatttate tttttgcc 29

<210> 859
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 859
taaccatttc gogtaagatt caa 23

<210> 860
<211> 39
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 860

taacccttgt ctttgaattg tatttgcaat taatcctgg

39

<210> 861

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 861

taaccgtttc caaaggctact gtattttgt

29

<210> 862

<211> 34

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 862

taaccgtttc caaaggctact gtattttggt tacc

34

<210> 863

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 863

taactcctct tccttcaaca ggtgga

26

<210> 864

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 864

taactgaccc aaagctgaaa gctttactg

29

<210> 865

<211> 33

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 865

taagacaagg ttttgtggat ttttagcgtt gtt

33

<210> 866

<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 866
taagagtgat gcgggctggt tcaaca 26

<210> 867
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 867
taagcaatac ctttacttgc accac 25

<210> 868
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 868
taagcaatac ctttacttgc accacct 27

<210> 869
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 869
taagcaatac ctttacttgc accacctg 28

<210> 870
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 870
taagcaccat ataagtctac tttttccct t 31

<210> 871
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 871
 taagccagca agagctgtat agttcca 27

 <210> 872
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 872
 taagctcccg tatcttgagt cgcctc 26

 <210> 873
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 873
 taagttacct tgcccgtaaa cca 23

 <210> 874
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 874
 taagttcctt cgctagtatg ttggctt 27

 <210> 875
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 875
 taatcgagca ccatcttgga aagatttctc 30

 <210> 876
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 876
 taatctggct gcggaagtga aat 23

 <210> 877
 <211> 25

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 877
taatctggct gcggaagtga aatcg 25

<210> 878
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 878
taatgccggg tagtgcaatc cattcttcta g 31

<210> 879
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 879
taatgcgata ctggcctgca agtc 24

<210> 880
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 880
tacaaccttc ggataatcag gatgagaatt aat 33

<210> 881
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 881
tacaacgtga taaacacgac cagaagc 27

<210> 882
<211> 35
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 882
tacaactggt tcaaaaacat taagctgtaa ttgtc 35

<210> 883
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 883
tacagcttta aagccagcaa aatgaattac ag 32

<210> 884
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 884
tacaggagca gcaggettca ag 22

<210> 885
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 885
tacatogttt cgccaagat caatca 26

<210> 886
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 886
tacatctcct tcgatagaaa tttoattgct atc 33

<210> 887
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 887
taccaaaagcg tgcacgatag ttgag 25

<210> 888
<211> 27
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 888
taccatctac ccaaacatta gcaccaa 27

<210> 889
<211> 22
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 889
taccacagtt cccctgacct tc 22

<210> 890
<211> 27
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 890
taccggaagc accagcgaca ttaatag 27

<210> 891
<211> 27
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 891
tacctgcatt aatcgcttgt tcatcaa 27

<210> 892
<211> 22
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 892
taccttacgg ccaaagctgt ct 22

<210> 893
<211> 24
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 893

taccttagga cggttatagt taag	24
<210> 894	
<211> 24	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 894	
taccttttcc acaacagaat cagc	24
<210> 895	
<211> 21	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 895	
tacgagctga cgacagccat g	21
<210> 896	
<211> 23	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 896	
tacgagctga cgacagccat gca	23
<210> 897	
<211> 22	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 897	
tacgcattac tcaccogtcc gc	22
<210> 898	
<211> 20	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 898	
tacgccaatca ggccacgcat	20
<210> 899	
<211> 27	
<212> DNA	
<213> Artificial Sequence	

<220>
<223> Primer

<400> 899
tacgctaagc cacgtccata ttatcat 27

<210> 900
<211> 35
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 900
tacgtatgta aattccgcaa agactttggc attag 35

<210> 901
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 901
tacgtgcct ttaacttggc tatattcagc 30

<210> 902
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 902
tacgttctac gatttcttca tcaggatcat c 31

<210> 903
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 903
tacgttttga ttttctgcag aacc 24

<210> 904
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 904
tacacctggc ttcgttttga tgatttgta 29

<210> 905
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 905
tactagacga cgggtcaggt aacc 24

<210> 906
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 906
tacttcagct tcgtccaata aaaaatcaca at 32

<210> 907
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 907
tactttaagg ggctatcttt accatgaacc t 31

<210> 908
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 908
tagagagtag ccatacttcac cgttgtc 27

<210> 909
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 909
tagcaccaat caccctttcc tgt 23

<210> 910
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
 <223> Primer

 <400> 910
 tagcagcaaa agttatcaca cctgcagt 28

 <210> 911
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 911
 tagcagctag ctcgtaacca gtgta 25

 <210> 912
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 912
 tagccatagc taccattgct tcataaatag a 31

 <210> 913
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 913
 tagcccagct gtttgagcaa ct 22

 <210> 914
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 914
 tagccgcggt cgaattgcat 20

 <210> 915
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 915
 tagccttggc aacatcagca aaact 25

<210> 916
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 916
tagccttttc tccggcgtag atct 24

<210> 917
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 917
tagcgatttc tactcctaga gttgaaattt cagg 34

<210> 918
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 918
tagctgctag atgagcttct gccatggcc 29

<210> 919
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 919
taggatgaaa gcattccgct ggc 23

<210> 920
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 920
taggatgagc attatcaggg aaagaatc 28

<210> 921
<211> 23
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer
 <400> 921
 taggattttt ccacggcggc atc 23
 <210> 922
 <211> 31
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 922
 taggcataac catttcagta ccttctggta a 31
 <210> 923
 <211> 28
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 923
 tagtatcacc acgtacacc ggatcagt 28
 <210> 924
 <211> 28
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <220>
 <221> modified_base
 <222> (17)...(17)
 <223> I
 <220>
 <221> modified_base
 <222> (20)...(20)
 <223> I
 <400> 924
 tagtatcacc acgtacacc ggatcagt 28
 <210> 925
 <211> 32
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 925
 tagtcctttc tgaattttac catcaaaggt ac 32
 <210> 926

<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 926
tagtcttttg gaacacgcgc tttaattaaa gt 32

<210> 927
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 927
tagtggtgta cctccatata gacattcaga 30

<210> 928
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 928
tagttgaagt tgcactatat actggttga 29

<210> 929
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 929
tataacgcac atcgtcaggg tga 23

<210> 930
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 930
tatagcacca tccatctgag cggcac 26

<210> 931
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 931
tatatgaaca ataccagttc cttctgagt 29
<210> 932
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 932
tatatgatta tcattgaact ggggocg 27
<210> 933
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 933
tatccattga accaaagtta ccttggcc 28
<210> 934
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 934
tatcccctgc ttctgctgcc 20
<210> 935
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 935
tatcgacaga tccaaagtta ccatgccc 28
<210> 936
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 936
tatggtctat ttcaatggca gttacga 27
<210> 937
<211> 24

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 937
tatgtgctca cgagtttgcg gcat 24

<210> 938
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 938
tatgtgtagt tgagcttact acatgagc 28

<210> 939
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 939
tattcttctgt tactcatgcc ataca 25

<210> 940
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 940
tattgccag aatcaaatc atc 23

<210> 941
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 941
tattgcggat caccatgatg atattcttgc 30

<210> 942
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 942 tattgctttt ttgtctatgc ttcttggaca t	31
<210> 943 <211> 24 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 943 tattggaaat accggcagca tctc	24
<210> 944 <211> 30 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 944 tatttgggtt tcattccact cagattcttg	30
<210> 945 <211> 32 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 945 tcaaaaacaa agaattcatt ttctgggtcca aa	32
<210> 946 <211> 34 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 946 tcaaaaacgca tttttacatc ttctgtaaag gcta	34
<210> 947 <211> 29 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 947 tcaaaaacttg ctctagacca tttaactcc	29
<210> 948 <211> 30 <212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 948
tcaaaatctt ttgattgat catacgagac 30

<210> 949
<211> 27
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 949
tcaaacgacg cgcacaccca tcaaaag 27

<210> 950
<211> 30
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 950
tcaaaagaacc agcacctaatt tcattcattta 30

<210> 951
<211> 30
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 951
tcaaaagaacc cgcacctaatt tcattcattta 30

<210> 952
<211> 26
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 952
tcaacaacac ctctctattc coactc 26

<210> 953
<211> 29
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 953

tcaacaatca gatagatgtc agacgcatg

29

<210> 954

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 954

tcaacaccag cgttacctaa agtacctt

28

<210> 955

<211> 35

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 955

tcaactgggt caaaaacatt aagttgtaat tgttc

35

<210> 956

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 956

tcaacttctg ccattaaaag taatgcc

28

<210> 957

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 957

tcaagcgatc taccgcatt acaa

24

<210> 958

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 958

tcaagcgcca tctctttcgg taatccacat

30

<210> 959

<211> 30

<212> DNA

<213> Artificial Sequence

<220>
<223> Primer

<400> 959
tcaagcgcca ttctcttttg taaaccacat 30

<210> 960
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 960
tcaagctata tgctacaact ggttcaaaaa c 31

<210> 961
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 961
tcaagctcta caccataaaa aaagctctca 30

<210> 962
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 962
tcaagggttct caccgtttac cttaggag 28

<210> 963
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 963
tcaagtgctt ttacttctat aggtttaage tc 32

<210> 964
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 964
tcaatacaga gtctacactt ggcttaggat 30

<210> 965
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 965
tcaatctcga ctttttgtgc cggta 25

<210> 966
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 966
tcacaaggac cattataatc aatgccaa 28

<210> 967
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 967
tcacaccaag tagtgcaagg atc 23

<210> 968
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 968
tcacacctgt aagtgagaaa aaggttgat 29

<210> 969
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 969
tcacaggttc tacttcatca ataattcca ttgc 34

<210> 970
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 970
tcaccagctt cagcgtagtc taataattta cgga 34

<210> 971
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 971
tcaccatgcg cccgttcaca ta 22

<210> 972
<211> 37
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 972
tcaccgataa ataaaatacc taaagttaat gccattg 37

<210> 973
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 973
tcacctacag ctttaaagcc agcaaaatg 29

<210> 974
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 974
tcacgatacc tgcacatca aattggtt 28

<210> 975
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 975
tcacgatcta aattggata agccatagga aa 32

<210> 976
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 976
tcacgcgcacg agtgccatcc attg

24

<210> 977
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 977
tcacgcgcacat catcaccagt ca

22

<210> 978
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 978
tcacgggcca gctcgtct

18

<210> 979
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 979
tcacgtcgtc cgacttcacg gtcagcat

28

<210> 980
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 980
tcagaatcga tgccaaatgc gtcac

26

<210> 981
<211> 30
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 981
tcagatataa atggaacaaa tggagccact 30

<210> 982

<211> 32

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 982
tcagcgtagt ctaataattt acggaacatt tc 32

<210> 983

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 983
tcagctgtta acggttcaa gaccc 25

<210> 984

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 984
tcaggtatga aacacgatta gtcctttct 29

<210> 985

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 985
tcagtttgca cttcaaaga aattgtgtt 29

<210> 986

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 986
tcataactag catttgtgct ttgaatgct 29

<210> 987

<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 987
tcataagggt tgcgttgtag attatcttta c 31

<210> 988
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 988
tcattctgggt taggattctgg ttgact 26

<210> 989
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 989
tcattctgggt tatggcgggt aagtt 25

<210> 990
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 990
tcattgacagc caagacctca cccacc 26

<210> 991
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 991
tcattgataga actacctggt tgcatttttg g 31

<210> 992
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 992
 tcatgtgcta atgttactgc tggatctg 28
 <210> 993
 <211> 31
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 993
 tcattaggta aaatgtctgg acatgatcca a 31
 <210> 994
 <211> 28
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 994
 tcatttatatt ctctgctttt ctcgctac 28
 <210> 995
 <211> 20
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 995
 tcattttgtgc tttgaatgct 20
 <210> 996
 <211> 28
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 996
 tccaaacgat ctgcatcacc atcaaaag 28
 <210> 997
 <211> 29
 <212> DNA
 <213> Artificial Sequence
 <220>
 <223> Primer
 <400> 997
 tccaaccag aaccacatac ttatttcac 29
 <210> 998
 <211> 27

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 998
tccaaccttt tccacaacag aatcagc 27

<210> 999
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 999
tccaagtgtt gggttaaccc atgg 24

<210> 1000
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1000
tccaagtgtt gggttaaccc atggag 26

<210> 1001
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1001
tccaagtttg acctaaacgt accatcgc 28

<210> 1002
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1002
tccacactgg attgtaattt acctgtttct tt 32

<210> 1003
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1003 tccaccacct caaagaccat gtggtg	26
<210> 1004	
<211> 23	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1004 tccagcaggt tctgacggaa acg	23
<210> 1005	
<211> 28	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1005 tccagcaggt actgtccct catctttg	28
<210> 1006	
<211> 30	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1006 tccaggcatt accatttcta ctcttctg	30
<210> 1007	
<211> 29	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1007 tccataaggt caccgtcacc attcaaagc	29
<210> 1008	
<211> 35	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<220>	
<221> modified_base	
<222> (21) ... (21)	
<223> I	
<220>	

<221> modified_base
<222> (30)...(30)
<223> I

<400> 1008
tccatacctt tatgcaactt tgtatcaact ggaat 35

<210> 1009
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1009
tccatattgt tgcataaaac ctgttggc 28

<210> 1010
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1010
tccatccata gaaccaaagt taccttg 27

<210> 1011
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1011
tccatcgag tcacgtttac tgttgg 26

<210> 1012
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1012
tccatcgcca gtttttgcac aatcgctaaa aa 32

<210> 1013
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1013

tccatctggt aaaccatcat ataccatgct atc 33

<210> 1014
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1014 26
tccatttcgg acacgtcggt gatcac

<210> 1015
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1015 27
tcccaatcta acttccacat accatct

<210> 1016
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1016 29
tcccaatott ttgattogat catacgaga

<210> 1017
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1017 27
tcccatacct atggcgataa ctgtcat

<210> 1018
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1018 30
tcccattttt tcacgcatgc tgaaaatatc

<210> 1019
<211> 15
<212> DNA
<213> Artificial Sequence

<220>
 <223> Primer

 <400> 1019
 tccccacatt cctcc 15

 <210> 1020
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1020
 tccccatttc cgcaaagaca ataaa 25

 <210> 1021
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1021
 tccccattta ataattccac ctactatcac act 33

 <210> 1022
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1022
 tccccctcatg tttaaattgat caggataaaa agc 33

 <210> 1023
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1023
 tcccccttta agcaccatta ctcattatag t 31

 <210> 1024
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1024
 tcccgaacaa tgagttgtat caactatttt tac 33

<210> 1025
 <211> 21
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1025
 tcccgctggc aaataaactc g 21

<210> 1026
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1026
 tcccgctag agattctgta tacga 25

<210> 1027
 <211> 28
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1027
 tcccgctctga gttctcttg catgatca 28

<210> 1028
 <211> 32
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1028
 tccctaataag tagaaataac tgcacagta gc 32

<210> 1029
 <211> 34
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1029
 tcccttattt ttctttctac taccttcgga taat 34

<210> 1030
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 1030
 tcccttcctt aatatgagaa ggaaaccact 30

 <210> 1031
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1031
 tccgaaactt gttttgtagc ttttaattga gc 32

 <210> 1032
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1032
 tccgaagttg ccttggccgt c 21

 <210> 1033
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1033
 tccgagacca gcgtagggtg aacg 24

 <210> 1034
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1034
 tccgataagc cggattctgt gc 22

 <210> 1035
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1035
 tccgcaaaga ctttggcatt aggtgtga 28

<210> 1036
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1036
 tccgccaaaaa aaactccctttt tcacagg 27

 <210> 1037
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1037
 tccgccttca aaatggtggc gagt 24

 <210> 1038
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1038
 tccggctaga gattctgtat acgaaaaatat c 31

 <210> 1039
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1039
 tccggctaga gattctgtat acgacaatat c 31

 <210> 1040
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1040
 tccggtaact gggtcagctc gaa 23

 <210> 1041
 <211> 34
 <212> DNA
 <213> Artificial Sequence

 <220>

<223> Primer

<400> 1041

tccgtagttt tgcataattt atggtctatt tcaa

34

<210> 1042

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1042

tccgtcatcg ctgacagaaa ctgagtt

27

<210> 1043

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1043

tccgtctatc cacaagttaa ttggtact

28

<210> 1044

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1044

tcctacccaa cgttcaccaa gggcag

26

<210> 1045

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1045

tcctccttgt gcctcaaaac gcattttta

29

<210> 1046

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1046

tcctctatgc aacttagtat caacaggaat

30

<210> 1047

<211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1047
 tcctcttggg ccacgcaag tttt 24

 <210> 1048
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1048
 tcctcttttc acaggtctta cttcattc 27

 <210> 1049
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1049
 tcctgaagat ctagtcttgg aatgggttact 30

 <210> 1050
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1050
 tcctgcaata tataatgcac tcttacg 27

 <210> 1051
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1051
 tcctgcagct ctacctgctc catta 25

 <210> 1052
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 1052
 tcctggccat cctgcaggat 20

 <210> 1053
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1053
 tcctgtttta tagccgcca gagtaag 27

 <210> 1054
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1054
 tccttcacgc gcatcatcac c 21

 <210> 1055
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1055
 tcctttgat gcctgatgga ccaggag 27

 <210> 1056
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1056
 tccttgcat acatcatgtc gtagca 26

 <210> 1057
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1057
 tccttggtct tcaaacgca tttttacatt ttc 33

 <210> 1058
 <211> 28

<212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1058
 tccttttaaaa taaccgctag tagctcct 28

 <210> 1059
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1059
 tccttttatgc aacttagtat caaccggaat 30

 <210> 1060
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1060
 tccttttatgc aacttggtat caacaggaat 30

 <210> 1061
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1061
 tccttttatgc aacttggtat caaccggaat 30

 <210> 1062
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1062
 tcctttcaat gttacagaaa actctacag 29

 <210> 1063
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 1063
tcgaaccgaa gttaccctga ccat 24
<210> 1064
<211> 30
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 1064
tcgaattcag ctaaatactt ttcagcatct 30
<210> 1065
<211> 26
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 1065
tcgacctgga ggacgacgta aatca 26
<210> 1066
<211> 25
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 1066
tcgacgacca tcttggaaag atttc 25
<210> 1067
<211> 25
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 1067
tcgagccgaa gttaccctgt ccgtc 25
<210> 1068
<211> 27
<212> DNA
<213> Artificial Sequence
<220>
<223> Primer
<400> 1068
tcgatccgca tcaccatcaa aagcaaa 27
<210> 1069
<211> 26
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1069

tcgatcgaac cgaagttacc ctgacc

26

<210> 1070

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1070

tcgatcgtga ctctctttat tttcagtt

28

<210> 1071

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1071

tcgatctcct tggcgccga

20

<210> 1072

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1072

tcgcaccgtg ggttgagatg aagtac

26

<210> 1073

<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1073

tcgcagcgtg cgtggcac

18

<210> 1074

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1074

tcgcaggctt acagaacgct ctctcta 26

<210> 1075
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1075
 tcgcagttca tcagcacgaa gcg 23

<210> 1076
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1076
 tcgccagcta gcacgatgto attttc 26

<210> 1077
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1077
 tcgccatagc taagttgttt attgtttcca t 31

<210> 1078
 <211> 22
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1078
 tcgcctggtg caggcatcat at 22

<210> 1079
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1079
 tcgcgctgta ttttctctcc gaga 24

<210> 1080
 <211> 18
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 1080
 tcgctacctt aggaccgt 18

 <210> 1081
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1081
 tcgctcagca ataattcact ataagcga 29

 <210> 1082
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1082
 tcgctctctc aagtgatcta aactggag 29

 <210> 1083
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1083
 tcgcttgagt gtagtcatga ttgcg 25

 <210> 1084
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1084
 tcggaaacaa agaattcatt ttctggcca aa 32

 <210> 1085
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1085
 tcggaaatat tctttcaata cctttatgca act 33

<210> 1086
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1086
tcggaaatat ttttcaata ctttatgca actt

34

<210> 1087
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (20)...(20)
<223> I

<220>
<221> modified_base
<222> (26)...(26)
<223> I

<400> 1087
tcggaaatat ttttcaata ctttatgca actt

34

<210> 1088
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1088
tcggactcgc ttctgctacg

20

<210> 1089
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1089
tcggataagc tgccacaagg

20

<210> 1090
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
 <223> Primer

 <400> 1090
 tgggcacac gccgtgctc 19

 <210> 1091
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1091
 tggcggaaca tggccatcac 20

 <210> 1092
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1092
 tcgggcgtag tttttagtaa ttaaatcaga agt 33

 <210> 1093
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1093
 tcgggtacgaa ctggatgtcg ccgtt 25

 <210> 1094
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1094
 tcgggtacgaa aaacggtagc ttgc 24

 <210> 1095
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1095
 tcgggtggtgg tagccgatct c 21

<210> 1096
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1096
tcggtttaag ctctacatga tcgtaaggat a 31

<210> 1097
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1097
tcggtttcag tcctctccac cataaagg 29

<210> 1098
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1098
tcgtatgacc agcttcggta ctacta 26

<210> 1099
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1099
tcgtcaacac taccattatt accatgcac tc 32

<210> 1100
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1100
tcgtccgact taacggtcag catttc 26

<210> 1101
<211> 29
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 1101
tcgtccgact taacggtcag catttctg 29

<210> 1102

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1102
tcgtccgact taacggtcag catttctgc a 31

<210> 1103

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1103
tcgtctcttc gaatctcga tataacc 26

<210> 1104

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1104
tcgtcgcgga cttcgaagcc 20

<210> 1105

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1105
tcgtcggact taacggtcag catttc 26

<210> 1106

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1106
tcgtcggact taacggtcag catttctg 29

<210> 1107

<211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1107
 tcgtcggact taacggtcag catttctctgc a 31

 <210> 1108
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1108
 tcgtcgtatt tatagtgacc agcaccta 28

 <210> 1109
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1109
 tcgtgcctaa caaatcccggt ctgagttc 28

 <210> 1110
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1110
 tcgtgggacta ccagggtatc ta 22

 <210> 1111
 <211> 17
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1111
 tcgtggggcct tgccgggt 17

 <210> 1112
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 1112
tcgttaatta atctggctgc ggaagtga 28

<210> 1113
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1113
tcgttgagat ggtttttacc ttogttg 27

<210> 1114
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1114
tcgtttaagc gccagaaagc accaa 25

<210> 1115
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1115
tcgtttcacc ctgtcatgcc g 21

<210> 1116
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1116
tccttcgtat aaaaaggacc aattgg 26

<210> 1117
<211> 35
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1117
tctacaacac ttgattgtaa ttgccttgt tcttt 35

<210> 1118
<211> 24

<212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1118
 tctagcggaa caacagttct gatg 24

 <210> 1119
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1119
 tctatagagt cggactttc ctggtga 27

 <210> 1120
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1120
 tctataggtg ctgtagtttg tttccgtct 30

 <210> 1121
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1121
 tctcacctac agctttaaag ccagcaa 27

 <210> 1122
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1122
 tctcacctac agctttaaag ccagcaaat g 31

 <210> 1123
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 1123 tctcatcccg atattacgc catga	25
<210> 1124 <211> 29 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1124 tctcatgaaa aaggctcagg agatacaag	29
<210> 1125 <211> 27 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1125 tctcttacc caccctttca cccttac	27
<210> 1126 <211> 31 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1126 tctctttcaa agcaccattg ctcattatag t	31
<210> 1127 <211> 24 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1127 tctgcatttt tgcgagcctg tcta	24
<210> 1128 <211> 28 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1128 tctgcctgag atgtcgaaaa aaacgttg	28
<210> 1129 <211> 26 <212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 1129

tctggccctt ccatacatgt atttag

26

<210> 1130

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1130

tctggctgcg gaagtgaat cgt

23

<210> 1131

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1131

tctgggtgac ctggtgtttt aga

23

<210> 1132

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1132

tctgtttcag ttgcaaattc

20

<210> 1133

<211> 34

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1133

tcttcacact tttagaatca accgttttat tgtc

34

<210> 1134

<211> 35

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1134

tcttcagcgt agtctaataa ttacggaac atttc 35

<210> 1135
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1135
tcttccaagg atagatttat ttctgttgcg 30

<210> 1136
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1136
tcttctgtaa aggggtggtt attattcacc cca 33

<210> 1137
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1137
tcttcttctt tcttataaaa aggaccaatt gg 32

<210> 1138
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1138
tcttcttgaa aaattgtgtt ccggaac 28

<210> 1139
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1139
tcttcttctg tataaaaagg accaattggt t 31

<210> 1140
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1140
tcttgacagc atccggttg 18

<210> 1141
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1141
tcttgagcat tgggttttac ttgttttgca ta 32

<210> 1142
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1142
tcttgagcca tacgtaccat tgc 23

<210> 1143
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1143
tcttggccta ggatgaaaat atagtgggtg ta 32

<210> 1144
<211> 39
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1144
tctttaagtt cttccaagga tagatttatt tcttgttcg 39

<210> 1145
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1145
tcttttcttt gettaatttt ccatttgcga t 31

<210> 1146
 <211> 14
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1146
 tgttactgct gcat 14

<210> 1147
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1147
 tgaacatttg cgacggtata cccat 25

<210> 1148
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1148
 tgaatatgta atgcaaacca gtctttgtca t 31

<210> 1149
 <211> 21
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1149
 tgaatcttga aacaccatac g 21

<210> 1150
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1150
 tgaatcttga aacaccatac gtaacg 26

<210> 1151
 <211> 34
 <212> DNA
 <213> Artificial Sequence

<220>
<223> Primer

<400> 1151
tgaattatgc aagaagtgat caattttctc acga 34

<210> 1152
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1152
tgaattctttt caaagcacca ttgctcatta tagt 34

<210> 1153
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1153
tgacaggaca caatctgcat gaagtctgag 30

<210> 1154
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1154
tgacccaaag ctgaaagctt tactg 25

<210> 1155
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1155
tgaccccaac ctggcctttt gtggtga 28

<210> 1156
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1156
tgaccgttat agttacggcc 20

<210> 1157
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1157
tgacggcatc gataccaccg tc 22

<210> 1158
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1158
tgacgtcatc cccaccttcc 20

<210> 1159
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1159
tgacgtcatc cccaccttcc tc 22

<210> 1160
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1160
tgacgtcatg cccaccttcc 20

<210> 1161
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1161
tgacgtcatg gccaccttcc 20

<210> 1162
<211> 32
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 1162
tgacttaaac gtaccatgc ttcatatata ga 32

<210> 1163

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1163
tgactttcct ccccttatc agtctcc 27

<210> 1164

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1164
tgagatgtcg aaaaaaacgt tggcaaaata c 31

<210> 1165

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1165
tgagatgttg atgatttacc agttccgatt g 31

<210> 1166

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1166
tgagcatcag cgtgcgtgct 20

<210> 1167

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1167
tgagcatttt tatatccatc tccaccat 28

<210> 1168

<211> 32
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1168 32
 tgagccatac gaacaatggg ttcataaaca gc

<210> 1169
 <211> 32
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1169 32
 tgagccatga gtaccatggc ttcataacat gc

<210> 1170
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1170 26
 tgagcgtgtg gaaaaggact tggatg

<210> 1171
 <211> 29
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1171 29
 tgagctgggt ctatatgaac aataccagt

<210> 1172
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1172 30
 tgagtcacc tccacaatgt atagttcaga

<210> 1173
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1173
tgagtcgggt tcactttacc tggca 25

<210> 1174
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1174
tgagtcotaca cttgggttag gatgaaa 27

<210> 1175
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1175
tgagttaaaa tgcgattgat ttcagtttcc aa 32

<210> 1176
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1176
tgagtttgaa ccatttcaga gcgaatatct ac 32

<210> 1177
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1177
tgagtttgca cttcaaaaga aattgtgt 28

<210> 1178
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1178
tgataaaaag cactaagcga tgaacagc 29

<210> 1179
<211> 27

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1179
tgataatgaa gggaaacctt ttccaag 27

<210> 1180
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1180
tgatatgaa ctggtgtacc ataatagttg cc 32

<210> 1181
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1181
tgatcctgaa tgtttatatc tttaacgcct 30

<210> 1182
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1182
tgatctccat ggcgcggatc tt 22

<210> 1183
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1183
tgatgcgggc tggttcaac 19

<210> 1184
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1184
tgatgcgggc tggttcaaca agag 24

<210> 1185
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1185
tgatgggtcta tttcaatggc agttacgaaa 30

<210> 1186
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1186
tgattatcag cggaagtag 19

<210> 1187
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1187
tgattcaaat gcagaaccat caaactcg 28

<210> 1188
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1188
tgattcgatc atacgagaca ttaaaactga g 31

<210> 1189
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1189
tgattggcga taaagtgata ttttctaaaa 30

<210> 1190
<211> 22
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1190

tgattgtttt gcagctgatt gt

22

<210> 1191

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1191

tgattgtttt gcagctgatt gt

22

<210> 1192

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1192

tgcaaaagta acggttacat ctgctccaat

30

<210> 1193

<211> 33

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1193

tgcaacaatt aatgctccga caattaaagg att

33

<210> 1194

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1194

tgcaactcat ctggtttagg atct

24

<210> 1195

<211> 32

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1195

tgcaactgaa tagattgcag taagttataa gc 32

<210> 1196
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1196
 tgcaagagca accctagtgt tcg 23

<210> 1197
 <211> 29
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1197
 tgcaagggaa acctagaatt acaaacctt 29

<210> 1198
 <211> 29
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1198
 tgcaatgtgt gctatgtcag caaaaagat 29

<210> 1199
 <211> 18
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1199
 tgcacctgcg gtcgagcg 18

<210> 1200
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1200
 tgcacgcaaa cgctttactt cagc 24

<210> 1201
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 1201
 tgcacgtctg tttcagttgc aaattc 26

 <210> 1202
 <211> 13
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1202
 tgcagctgat tgt 13

 <210> 1203
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1203
 tgcataggga aggtaacacc atagtt 26

 <210> 1204
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1204
 tgcatacaca tttccttgtc cttcg 25

 <210> 1205
 <211> 35
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1205
 tgcatagaagc ataaaaactg tatcaagtc tttta 35

 <210> 1206
 <211> 36
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1206
 tgcatagtta ctcaaatcat cataaacaat taaagc 36

<210> 1207
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1207
 tgcatgtac cgaagtagtt cacattgtt 29

 <210> 1208
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1208
 tgccaagtgc tggtttacct catgg 25

 <210> 1209
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1209
 tgccactttg acaactcctg ttgctg 26

 <210> 1210
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1210
 tgccagcgac agaccatcgt a 21

 <210> 1211
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1211
 tgccagctta gtcatacgga cttc 24

 <210> 1212
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 1212
 tgccagtttc cacatttcac gttcgtg 27

 <210> 1213
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1213
 tgccatacgt accatcggtt cataaacagc 30

 <210> 1214
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1214
 tgccatagca aagcctacag catt 24

 <210> 1215
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1215
 tgccatccat aatcacgcc a tactgacg 28

 <210> 1216
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1216
 tgccatttcc atgtactctt ctctaacatt 30

 <210> 1217
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1217
 tgcccaccag aaagactagc aggataa 27

<210> 1218
 <211> 20
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1218
 tgcccaggta caacctgcat 20

<210> 1219
 <211> 27
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1219
 tgccccattg ctcatgatag tagctac 27

<210> 1220
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1220
 tgccctttct aaaagtcttg agtgaagata 30

<210> 1221
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1221
 tgcccttttg taaaagcagg gctat 25

<210> 1222
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1222
 tgccgataag ccgattctg tgc 23

<210> 1223
 <211> 28
 <212> DNA
 <213> Artificial Sequence

<220>

<223> Primer

<400> 1223
tgccgtaaca tagaagttac cgttgatt

28

<210> 1224
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1224
tgccgtaact aacataagag aattatgcaa gaa

33

<210> 1225
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1225
tgccgtatatac gaaaatatct tatcatttag cgt

33

<210> 1226
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1226
tgccctaaca atcccgtctg agttc

25

<210> 1227
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1227
tgccctcgcg aacctaccg

20

<210> 1228
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1228
tgccctcgctg aacccaccg

20

<210> 1229

<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1229
tgcgaggaac ttcacgtcct gc 22

<210> 1230
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1230
tgcgatggta ggtatcttag caatcattct 30

<210> 1231
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1231
tgcgcgagct tttatttggg ttct 24

<210> 1232
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1232
tgcgctaatt cttcaacttc ttctttcgt 29

<210> 1233
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1233
tgcgctatca acgattttga caatatatgt ga 32

<210> 1234
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1234
tgcggcagca ctatcaccat cca 23

<210> 1235
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1235
tgcgggctgg ttcaacaaga g 21

<210> 1236
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1236
tgcgggtgat acttaccgag tac 23

<210> 1237
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1237
tgcggtctgg cgcataatagg ta 22

<210> 1238
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1238
tgcgtagtct aataatttac ggaacatttc 30

<210> 1239
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1239
tgcgtgacga cttctctgaa ttgtaatca 29

<210> 1240
<211> 23

<212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1240 23
 tgcgtggact accagggtat cta

<210> 1241
 <211> 13
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1241 13
 tgcagctgat tgt

<210> 1242
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1242 30
 tgctaaagtc ttgagccata cgaacaatgg

<210> 1243
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1243 24
 tgctagacct ttaogtgac cgtg

<210> 1244
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1244 23
 tgctaggcca tcaggccacg cat

<210> 1245
 <211> 34
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1245 tgctatatgc tacaactggg tcaaaaacat taag	34
<210> 1246 <211> 29 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1246 tgctcacctg ctacaacaag tccagcaat	29
<210> 1247 <211> 29 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1247 tgctcttacc tcacggttc acccttacc	29
<210> 1248 <211> 29 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1248 tgctgcttcc gcattggttaa ttgcttcaa	29
<210> 1249 <211> 27 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1249 tgctgctttg atggctgaat ccccttc	27
<210> 1250 <211> 22 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1250 tgctggattc gcctttgcta cg	22
<210> 1251 <211> 21 <212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 1251

tgctgtaggg aaatcagggc c

21

<210> 1252

<211> 18

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1252

tgcttagatg ctttcagc

18

<210> 1253

<211> 34

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1253

tgcttcaaaa cgcatttttta cattttcggt aaag

34

<210> 1254

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1254

tgcttcagca cggccaccaa cttctag

27

<210> 1255

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1255

tgcttcagcg tagtctaata atttacggaa c

31

<210> 1256

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1256

tgcttctctt cegggtcggc

20

<210> 1257
 <211> 32
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1257
 tgcttgctca aatcatcata aacaattaaa gc

32

<210> 1258
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1258
 tgcttgctct ttcaagcagt ctggaatgaa g

31

<210> 1259
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1259
 tgcttggtgg cttcttcgtc gaa

23

<210> 1260
 <211> 32
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1260
 tgctttgtaa tctagtctct gaatagtaac ca

32

<210> 1261
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1261
 tggaaaactc atgaaattaa agtgaaagga

30

<210> 1262
 <211> 29
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 1262 29
 tggaaaccgg ctaagtgagt accaccatc

 <210> 1263
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1263 30
 tggaaacaccg tctttaatta aagtatctcc

 <210> 1264
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1264 24
 tggaaatttac cagcgataga cacc

 <210> 1265
 <211> 20
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1265 20
 tggaccacgc cgaagaacgg

 <210> 1266
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1266 30
 tggacgatat tcacggttta ccacttata

 <210> 1267
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1267 31
 tggactaata acaatgagct cattgtactg a

<210> 1268
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1268
 tggataattg gtcgtaacaa gggatagtga g 31

 <210> 1269
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1269
 tggatagacg tcatatgaag gtgtgct 27

 <210> 1270
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1270
 tggatcactg cttacgaact cagcttc 27

 <210> 1271
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1271
 tggatgtgct cacgagtctg tggcat 26

 <210> 1272
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1272
 tggcaacagc tcaacacctt tgg 23

 <210> 1273
 <211> 26
 <212> DNA
 <213> Artificial Sequence

<220>
<223> Primer

<400> 1273
tggcaccgtg ggttgagatg aagtac 26

<210> 1274
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1274
tggcaccgagc ctgacctgt 19

<210> 1275
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1275
tggcagcaat agtttgacgt acaaatgcac acat 34

<210> 1276
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1276
tggcatcacc atttccttgt ccttcg 26

<210> 1277
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1277
tggccacttt tatcagcaac cttacagtc 29

<210> 1278
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1278
tggccgtact ccccgagcg 19

<210> 1279
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1279
tggcgatgca ctggcttgag

20

<210> 1280
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1280
tggctcataa gacgcgcttg taga

24

<210> 1281
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1281
tggctgcgga agtgaaatcg ta

22

<210> 1282
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1282
tggctgcttc taagccaac

19

<210> 1283
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1283
tggcttgaga atttaggatac cggcac

26

<210> 1284
<211> 30
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 1284

tgggacgtaa tcgtataaat tcatacatttc

30

<210> 1285

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1285

tgggataaca tgggttgga tataagcaga aacatc

36

<210> 1286

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1286

tgggatggag gtgtagaagg tggtatcatc

30

<210> 1287

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1287

tgggatggag gtgtagaagg tggtatcatc

30

<210> 1288

<211> 32

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1288

tgggcacccat ttatccacaa attgattggt at

32

<210> 1289

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1289

tggggacttc cttaccactt ttagtatcta a

31

<210> 1290

<211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1290
 tggggatatg gaggtgtaga aggtgttatc atc 33

 <210> 1291
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1291
 tggggtaaga cgcggttagc atgtatt 27

 <210> 1292
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1292
 tgggtacgaa ctggatgtcg cggtt 25

 <210> 1293
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1293
 tgggtagggt tttatctgtg aogcctt 27

 <210> 1294
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1294
 tgggtctaca cctgcacttg cataac 26

 <210> 1295
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 1295
tgggtgctgg ttaccocat ggag 24
<210> 1296
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1296
tgggttgctg tgcagattat ctttaccaa 29
<210> 1297
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1297
tgggtttctgc gcttagatgc tttca 25
<210> 1298
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1298
tggtaacct tgtctttg 18
<210> 1299
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1299
tggtaacct tgtctttgaa ttgtatttgc a 31
<210> 1300
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1300
tggtacaaca tegttagctt taccacttgc acg 33
<210> 1301
<211> 32

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1301 32
tggtagacct ggttctgttt tgatgatttg ta

<210> 1302
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1302 30
tggtagctca acttcatcca ttatgaagtc

<210> 1303
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1303 31
tggtagatctc gtttaattaat ctggctgagg a

<210> 1304
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1304 21
tggtagtagt acctcctttg c

<210> 1305
<211> 31
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1305 31
tggtagggat ctttagcaatc attctaatag c

<210> 1306
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1306 tggtgttcta gatatagattg aggtagtggg ga	32
<210> 1307	
<211> 24	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1307 tggttagaag tcgtaacgtg gacc	24
<210> 1308	
<211> 25	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1308 tggttcaaca agagtgcgg ttgca	25
<210> 1309	
<211> 31	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1309 tggttcttac ttgctttgca taaactttcc a	31
<210> 1310	
<211> 31	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1310 tggttgtagt tcctgtagtt gttgcattaa c	31
<210> 1311	
<211> 30	
<212> DNA	
<213> Artificial Sequence	
<220>	
<223> Primer	
<400> 1311 tggtttgtca gaatcacgtt ctggagttgg	30
<210> 1312	
<211> 28	
<212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 1312

tgtaaaagca gggctataat aaggactc

28

<210> 1313

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1313

tgtaaatgcc gcaaagactt tggcattag

29

<210> 1314

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1314

tgtaaccocc gtctttgaat tgtatttgc

29

<210> 1315

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1315

tgtaattaac cgaaggttct gtagaagtat g

31

<210> 1316

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1316

tgtacaagga ccattataat caatgcca

28

<210> 1317

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1317

tgtacaataa ggagtcacct tatgtcctt a 31

<210> 1318

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1318

tgtacacccat ttatccacaa attgattggt 30

<210> 1319

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1319

tgtaggcaag tgcataagaa attgataca 29

<210> 1320

<211> 27

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1320

tgtcaatatg aagggtctct gtggata 27

<210> 1321

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1321

tgtcaccagc ttcagcgtag tctaataa 28

<210> 1322

<211> 19

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1322

tgtcaactccc gacacgcca 19

<210> 1323

<211> 29

<212> DNA

<213> Artificial Sequence

<220>
<223> Primer

<400> 1323
tgtcagctaa gctaataacg ttgtagag 29

<210> 1324
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1324
tgtcatcaag caccctaaaa tgaact 26

<210> 1325
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1325
tgtccgactt gacggtoaac atttcctg 28

<210> 1326
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1326
tgtccgactt gacggtcagc atttcctg 28

<210> 1327
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1327
tgtccgactt gacggttagc atttcctg 28

<210> 1328
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1328
tgtccagca tctgttcctg c 21

<210> 1329
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1329
tgtctattgt cgattgttac ctgtacagt 29

<210> 1330
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1330
tgtgaacatt tgcgacggta tacccat 27

<210> 1331
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1331
tgtgaagaac tttcaaatct gtgaatcca 29

<210> 1332
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1332
tgtgatatgg aggtgtagaa ggtg 24

<210> 1333
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1333
tgtgatatgg aggtgtagaa ggtgtta 27

<210> 1334
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
 <223> Primer

 <400> 1334
 tgtgcaggca tcatgtcata ccaa 24

 <210> 1335
 <211> 31
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1335
 tgtgctgctt tcgcatgggt aattgcttca a 31

 <210> 1336
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1336
 tgtgctgggt taccocatgg ag 22

 <210> 1337
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1337
 tgtgctgggt taccocatgg agt 23

 <210> 1338
 <211> 15
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1338
 tgtgctttga atgct 15

 <210> 1339
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1339
 tgtgcttttt ttgctgccat agcaaagc 28

<210> 1340
<211> 26
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1340
tgtggcgcgat ttcaccacct gctoct 26

<210> 1341
<211> 18
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1341
tgtgtgtgtcg ccgcgcag 18

<210> 1342
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1342
tgttaacggc ttcaagaccc 20

<210> 1343
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1343
tgttaagtgt gbtgcggctg tctttatt 28

<210> 1344
<211> 36
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1344
tgttaatggt aacccttgct tttgaattgt atttgc 36

<210> 1345
<211> 22
<212> DNA
<213> Artificial Sequence

<220>

<223> Primer

<400> 1345

tggttactcac ccgtctgccca ct

22

<210> 1346

<211> 14

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1346

tggttactgct ggat

14

<210> 1347

<211> 34

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1347

tgttcatggtt taaatgatca ggataaaaag cact

34

<210> 1348

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1348

tgttccaata gcagttccgc ccaaattgat

30

<210> 1349

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1349

tgttctggat tgattgcaca atcaccaaag

30

<210> 1350

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1350

tgttcttgat acacctgggtt tcgttttgat

30

<210> 1351

<211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1351
 tggtgaagct gtacttgacc tgattttacg 30

 <210> 1352
 <211> 19
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1352
 tgttgaccat gcttcttag 19

 <210> 1353
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1353
 tgttgtgccc cagtcaaata tctaaata 28

 <210> 1354
 <211> 25
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1354
 tgtttgtgat gcatttgcg agcta 25

 <210> 1355
 <211> 33
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1355
 tgttttatgt gtagttgagc ttactacatg agc 33

 <210> 1356
 <211> 29
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 1356
tgttttgtat ccaagtgtg gtttaccoc 29

<210> 1357
<211> 11
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1357
tactcatgcc a 11

<210> 1358
<211> 11
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1358
tattcttcgt t 11

<210> 1359
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1359
ttacttctaa cccactc 17

<210> 1360
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1360
ttaatctggc tgcggaagtg aaatcgt 27

<210> 1361
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1361
ttaccatctt caaataccg aacagtaa 28

<210> 1362
<211> 26

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1362
ttaccgagca ggttctgacg gaaacg

26

<210> 1363
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1363
ttacgccatc aggccacgca

20

<210> 1364
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1364
ttactcaccg gtccgcc

17

<210> 1365
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1365
ttactcaccg gtccgccgct

20

<210> 1366
<211> 34
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1366
ttacttcctt accactttta gtatctaaag cata

34

<210> 1367
<211> 17
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1367 ttacttcttaa cccactc	17
<210> 1368 <211> 21 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1368 ttagaagtcg taacgtggac c	21
<210> 1369 <211> 22 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1369 ttagatgctt tcagcactta tc	22
<210> 1370 <211> 26 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1370 ttcaaaaacct tgctctcgcc aaacaa	26
<210> 1371 <211> 24 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1371 ttcaaaaagtt gctcgagacc attg	24
<210> 1372 <211> 23 <212> DNA <213> Artificial Sequence	
<220> <223> Primer	
<400> 1372 ttcaaaaatgc ggaggcgtat gtg	23
<210> 1373 <211> 22 <212> DNA	

<213> Artificial Sequence

<220>

<223> Primer

<400> 1373

ttcaacaaga gttgccgttg ca

22

<210> 1374

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1374

ttcaacacto tcacctacag ctttaaag

28

<210> 1375

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1375

ttcaagtgtc tgctcaccat tgtc

24

<210> 1376

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1376

ttcaggtaca gcaggtggtt caggat

26

<210> 1377

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1377

ttcaggtcca tcgggttcac gcc

23

<210> 1378

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1378

ttcataagca atacctttac ttgcaccac 29

<210> 1379
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1379
ttcattttct ggtccaaagt aagcagtatc 30

<210> 1380
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1380
ttccaagtgc tggtttacc catgg 25

<210> 1381
<211> 33
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1381
ttccaccttg gatacctgga aaaatagctg aat 33

<210> 1382
<211> 36
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1382
ttccatttca actaattcta ataattcttc atcgtc 36

<210> 1383
<211> 28
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1383
ttccctgac cttcgattaa aggatagc 28

<210> 1384
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
 <223> Primer

 <400> 1384
 ttcgcgcatc caggagaagt acatggt 27

 <210> 1385
 <211> 15
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1385
 ttcgctcgcc gctac 15

 <210> 1386
 <211> 18
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1386
 ttcgctctcg gcttggcc 16

 <210> 1387
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1387
 ttcggtataa cgcacgcag ca 22

 <210> 1388
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1388
 ttcggtgctgg attttgtcct tgtcct 26

 <210> 1389
 <211> 21
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1389
 ttcggtgctta gatgctttca g 21

<210> 1390
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1390
 ttctgagcta aatcagcagt tgca

24

<210> 1391
 <211> 23
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1391
 ttctgcgaat caatgcacg ctg

23

<210> 1392
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1392
 ttctgcttga ggaatagtgc gtgg

24

<210> 1393
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1393
 ttctgggtga cctgggtgtt taga

24

<210> 1394
 <211> 31
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1394
 ttcttccaag gatagattta tttcttggtc g

31

<210> 1395
 <211> 24
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

 <400> 1395
 ttcttgaacg cgaggtttcg attg 24

 <210> 1396
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1396
 ttgacatcgt ccctcttcac ag 22

 <210> 1397
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1397
 ttgacatttg catgcttcaa agcctg 26

 <210> 1398
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1398
 ttgacgtcat cccaccttc ctc 23

 <210> 1399
 <211> 24
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1399
 ttgacgttgc atgttcgagc ccat 24

 <210> 1400
 <211> 30
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1400
 ttgcaatcga catatccatt tcaccatgcc 30

<210> 1401
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1401
ttgcacgtct gtttcagttg caaatc 27

<210> 1402
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1402
ttgcacgtct gtttcagttg caaatc 27

<210> 1403
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1403
ttgcacggg ttggttaagtc 20

<210> 1404
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1404
ttgccacttt gacaactcct gttgctg 27

<210> 1405
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1405
ttgccatagc aaagcctaca gcatt 25

<210> 1406
<211> 30
<212> DNA
<213> Artificial Sequence

<220>

<400> FLAME1

<400> 1406
 ttgccattca tggattttaa gtgtagcaga

30

<210> 1407
 <211> 21
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1407
 ttgcgccata cgtaccatcg t

21

<210> 1408
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1408
 ttgcgttgca gattatcttt accaa

25

<210> 1409
 <211> 25
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1409
 ttgctgccat agcaaagcct acagc

25

<210> 1410
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1410
 ttgctgcttt cgcctggtta atcgcttcaa

30

<210> 1411
 <211> 30
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Primer

<400> 1411
 ttgctgcttt cgcctggtta attgcttcaa

30

<210> 1412

<211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1412
 ttggacctgt aatcacgtga atactgg 27

 <210> 1413
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1413
 ttggccatca gaccacgcat ac 22

 <210> 1414
 <211> 22
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1414
 ttggccatca ggccacgcat ac 22

 <210> 1415
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1415
 ttggcgacgg tatacccata gctttata 28

 <210> 1416
 <211> 17
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1416
 ttggtgcgct tggecgta 17

 <210> 1417
 <211> 32
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

<400> 1417
 ttggtttctta cttgttttgc ataaactttc ca 32

 <210> 1418
 <211> 35
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1418
 ttgtacattt gaaacaatat gcatgacatg tgaat 35

 <210> 1419
 <211> 23
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1419
 ttgtcagact catcgcgaaac atc 23

 <210> 1420
 <211> 28
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1420
 ttgtgatatg gaggtgtaga aggtgtta 28

 <210> 1421
 <211> 26
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1421
 ttgtgattgt ttgcagctg attgtg 26

 <210> 1422
 <211> 27
 <212> DNA
 <213> Artificial Sequence

 <220>
 <223> Primer

 <400> 1422
 ttgtggccga ttccaccacc tgetcct 27

 <210> 1423
 <211> 21

<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1423
ttgttaacgg cttcaagacc c 21

<210> 1424
<211> 32
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1424
ttgtttattg ttccatatg ctacacactt tc 32

<210> 1425
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1425
tttaagcgcc agaaagcacc aac 23

<210> 1426
<211> 25
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1426
tttacctcgc ctttcacccc ttacc 25

<210> 1427
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1427
tttagctact attotagctg ccatttcca 29

<210> 1428
<211> 27
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1428
tttatgacca gcttcggtac tactaaa 27

<210> 1429
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1429
tttatggtct atttcaatgg cagttacgaa 30

<210> 1430
<211> 37
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<220>
<221> modified_base
<222> (23)...(23)
<223> I

<220>
<221> modified_base
<222> (32)...(32)
<223> I

<400> 1430
tttcaatacc tttatgcaac ttgtatcaa ctggaat 37

<210> 1431
<211> 30
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1431
tttcacagca tgcacgtctg tttagttgc 30

<210> 1432
<211> 35
<212> DNA
<213> Artificial Sequence

<220>
<223> Primer

<400> 1432
tttccccgat cttaaatttg ataagccata ggaaa 35

<210> 1433
<211> 27
<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1433

tttcgatgc aacgaatga gatttca

27

<210> 1434

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1434

tttcgtgctt agatgctttc ag

22

<210> 1435

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1435

tttcttgaag agtatgagct gctcgttaag

30

<210> 1436

<211> 22

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1436

tttgcacctt acgcctaaag ct

22

<210> 1437

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1437

tttgctcatg atctgcatga agcataaa

28

<210> 1438

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1438

tttgctctcc gccaaagttt ccac

24

<210> 1439

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1439

tttggacctg taatcagctg aatactgg

28

<210> 1440

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1440

tttgtgaaac agcgaacatt ttcttggtta

29

<210> 1441

<211> 20

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1441

ttttccagcc atgcagcgac

20

<210> 1442

<211> 34

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<220>

<221> modified_base

<222> (29)...(29)

<223> I

<400> 1442

ttttcccttt atgcaactta gtatcaactg gaat

34

<210> 1443

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 1443

ttttgctcat gatctgcatg aagcataaa

29

<210> 1444
 <211> 2118
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Concatenation of *A. baumannii* genes

<220>
 <221> misc_feature
 <222> 447-486, 778-817, 1162-1201, 1495-1534, 1928-1967, 2115-2118
 <223> n = A,T,C or G

```

<400> 1444
cgcgcggttaa aactaaagaa gaagatatag cattagaaaa agatttgctg tctgatgaaa 60
aagagattgc tgaacattta atgctgattg atcttgggcg aaacgatgta gggcggtgat 120
cgaaaatagg taaagtccaa gtccgcggtc aaatgggtgat cgagcggtat tcacatgtca 180
tgcatattgt ttcaaatgta caaggtgaag tgcgtgatga tatcgatgca cttgatgtat 240
ttaaagccac cttccagcca ggaacgttat caggtgcccc aaaaattcgt gcaatggaaa 300
ttattgatga agtagaacct gtgaaaaggg caatcgctac ttgtgttttc ggttatttgg 360
tagggcatgg tgaaatggat atgtcgattg caatcgctac ttgtgttttc cgtgataaaa 420
aggtgtatgt acaggctggt gcaggggnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn 480
nnnnngggaa tctggcggtt tagtttcaga tgaactcatt atcggtttag taaaagaacg 540
tatgtctcaa cctgactcgg tgaatggttg tattttcgac ggcttccacc gcaactatcc 600
tcaagcagaa gctttggaaa aagaagggat cagcattgat catgtaattg aaattgatgt 660
accctgatga gaaatcgtaa aacgtcttct tggctcgtgt cagcatcctg cttctggctg 720
tgtttatcac gttgtatata atccacctaa agtggaaagt aaagatgatg tccacggnnn 780
nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnncg tcaaccgctg aaaattacgt 840
aaccttaaaa ctgtgaagt tttagaaaaa acttttaaat ctggtgatac tttagaagct 900
gctgacatcg tagaagtaga aatgaactac ctatacaacg atggcgaagt gtggcacttc 960
atggaccagc aagagcttcga acaaattgca gctgacaaaa ctgcaatggg tgatgctgtc 1020
aaatgggttaa aagacgactc aaatgaaaca tgtacaatca tgttatcaa cggcggtcct 1080
ttaacgttaa atgcacctaa cttcgttgta ttgaaagtg ttgaaactga tccggcgcta 1140
cgtggtgata cttctgtgtg tnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn 1200
ntcgtgcccg yaatttgcatt aaagctgcgc gccttgtagc acagcaaggc aaatttctctg 1260
aaaactctaga agaattgatt gcaactaccg gcattggctg ctogaccgca ggtgcactca 1320
tgtctttagg tttaactcag tatggcgtga ttatggatgg caacgtgaaa cgcgcttttag 1380
ccgctttctt tgccattgaa gatgacttaa gcaaaaccaca gcacgaacgt gaaatgtgga 1440
aactggctga agagctttgt cccaccacaac gcaatcatga ctacactcaa gcaannnnnn 1500
nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnttaaaa acactagcgg taagcttaaa 1560
caagattgcc aatgatattc gttgggttagc aagtgttcca cgttgcggct tcggcgaaat 1620
cgtatttcct gaaaatgaac ctgggttcaag tatcatgccca ggttaagtga acccgactca 1680
aagtgaagcc atgaccatgg ttgtgtgtca agtacttggc aacgatacca ctattaatgt 1740
cgctgggtgct tctggttaact tcgagctcaa tgatttatgt ccagtgatgt cttataactg 1800
actgcaatct attcagttgc ttggtgatgc atgtaaatgt tttaaatgat atgtgcaggt 1860
agggatcgag ccaaatcgtg agaaaatgta tcaatttctg cataattctc ttatgttagt 1920
tacggcannn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn ggttatgtac 1980
caaaactatt gtctgaagat ggtgacccat tagacgtact tgttgtaact ccacatcctg 2040
ttgctgcggc ttctgttaatt cgttgcggcc cagtggcgaa attaaactag gaagacgacg 2100
gtggatctga tgcnnnnn

```

<210> 1445
 <211> 102
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Polynucleotide

<400> 1445
 ttttaagtcgcc gcaacgagcg caacccttga tcttagttgt ttagttgggc actctaaggt 60
 gactgcgcgt gacaaaccgg aggaaggtgg ggatgaogtc aa 102

<210> 1446
 <211> 94
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Polynucleotide

<400> 1446
 tagaacacgg atggcggaagg cgactttctg gtctgtaact gacactgaga aagcgtgggg 60
 agcaaacagg attagatacc ctggtagtcc acga 94

<210> 1447
 <211> 108
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Polynucleotide

<400> 1447
 tggattagag accttggttag tccacgcggt aaacgatgag tgctaagtgt tagaggcctt 60
 tagtgtctgaa gttaacgcat taagcactcc gcctgggggag tacggcca 108

<210> 1448
 <211> 108
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Polynucleotide

<400> 1448
 ttctgatgca acgcgaagaa ccttaccagg tcttgacatc ctctgacaac cctagcttct 60
 ccttcggggag cagagtgaca ggtggtgcat ggctgtcgtc agctcgta 108

<210> 1449
 <211> 95
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Polynucleotide

<400> 1449
 tctgcacact gcccggtgct ggaaggttaa ggagaggggt tagcgtact ctgaactgaa 60
 gccccagtaa acggcgggcg taactataac ggtca 95

<210> 1450
 <211> 117
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Polynucleotide

<400> 1450

tctgttttta gtacgagagg accgggatgg acgcaccggg accagttggt ctgccaaggg 60
catagctggg tagctatgtg cggaagggat aagtgcgtgaa agcatctaag caccgaaa 117

<210> 1451
<211> 100
<212> DNA
<213> Artificial Sequence

<220>
<223> Calibration Polynucleotide

<400> 1451
tgattattgt tatocgtgta tgccatttga gatttttgag tggatttga gttattgttc 60
caggattaat tgcaaatata attcaaagac aagggttaca 100

<210> 1452
<211> 112
<212> DNA
<213> Artificial Sequence

<220>
<223> Calibration Polynucleotide

<400> 1452
tcgaagtaca atacaagaca aaagaaggta aaattactgt tttaggggaa aaattcaaga 60
aatatagaag tgatggtctaa aaatgtagaa ggggtcttga agccgttaac aa 112

<210> 1453
<211> 100
<212> DNA
<213> Artificial Sequence

<220>
<223> Calibration Polynucleotide

<400> 1453
ttgctcgtgg tgcacaagta acggatatta caatcattgt tgttcagct gatgacggcg 60
taataaacag ttgaagcaat taaccatgcg aaagcagcaa 100

<210> 1454
<211> 114
<212> DNA
<213> Artificial Sequence

<220>
<223> Calibration Polynucleotide

<400> 1454
tagcttttgc atattatata gagccacagc atcgtgatgt ttacagctt tatgcaccgg 60
aagcttttaa tggataaatt taacgaacaa gaaataaata tacccttgga agaa 114

<210> 1455
<211> 116
<212> DNA
<213> Artificial Sequence

<220>
<223> Calibration Polynucleotide

<400> 1455
tgacctacag taagagggtc tgtaatgaac ctaaatgacc atccacacgg tgggtggtgaa 60

ggtagatctc ctatcggaag gtccacgtac tccatggggg aaaccagcac ttggaa 116

<210> 1456

<211> 70

<212> DNA

<213> Artificial Sequence

<220>

<223> Calibration Polynucleotide

<400> 1456

tccacacggt ggtggtgaag gtatgctcc tctcggaag tccacgtact ccatggggta 60
aaccagcaca 70

<210> 1457

<211> 82

<212> DNA

<213> Artificial Sequence

<220>

<223> Calibration Polynucleotide

<400> 1457

ttatcgctca ggcgaaactcc aacctggatg atgaaggccg ctttttagaa ggtgacttgt 60
cgtagcaaaag gcgaatccag ca 82

<210> 1458

<211> 87

<212> DNA

<213> Artificial Sequence

<220>

<223> Calibration Polynucleotide

<400> 1458

tgggcagcgt ttccgcaaaa tgggaagtggc tcgaagcgta tggcgcttcg tacgtgctgc 60
aggaaatggt gaccgtcaag tcggaca 87

<210> 1459

<211> 97

<212> DNA

<213> Artificial Sequence

<220>

<223> Calibration Polynucleotide

<400> 1459

tcaggagtcg ttcaactcga tctacatgat ggccgaccgc cgggggttcg gcggtgcaga 60
ttcgtcagct ggccgcatg cgtggcctga tggcgta 97

<210> 1460

<211> 117

<212> DNA

<213> Artificial Sequence

<220>

<223> Calibration Polynucleotide

<400> 1460

tctggcaggt atgcgtggtc tgatggccaa tccatctggt cgtatcatcg aacttccaat 60
caagtttccg tgaagggtta acagtacttg agtacttcat ctcaaccac ggtgcga 117

<210> 1461
 <211> 98
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Polynucleotide

<400> 1461
 tcaagcaaac gcacaatcag aagctaagaa agcgcaagct tctggaaagc acaaatgcta 60
 gttatggtag agaatttgca actgaaacag acgtgcaa 98

<210> 1462
 <211> 99
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Polynucleotide

<400> 1462
 tccacacggc gtttttcaac aactaccgtg ttctacttcc gtacgacgga cgtgacgggg 60
 tcgatcgagc tgcogaagga caaggaaatg gtgatgcaa 99

<210> 1463
 <211> 111
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Calibration Sequence

<400> 1463
 tcgtggcgcc gtggttatcg aacctatgct gaccgatcaa tggtagctgc acaccgcccc 60
 ccaaagtgcg gattgaagcc gtagagaacg ggcacatcca gttcgtaccg a 111

<210> 1464
 <211> 2100
 <212> DNA
 <213> Artificial Sequence

<220>
 <223> Combination Calibration Polynucleotide

<400> 1464
 gaagtagaga tatggaggaa caccagtggc gaaggcgact ttctggctctg taactgacac 60
 tgagaaagcg tggggagcaa acaggattag ataccctggg agtccacgcc gtaaacgatg 120
 agtgctaagt gttagaggcc tttagtgctg aagttaacgc attaacgact ccgcctgggg 180
 agtacggccg caaggctgaa actcaaaagga attgacgggg cacaagcggg ggagcatgtg 240
 gtttaattcg aagcaacgag aagaacctta ccaggtcttg acatcctctg acaaccctag 300
 ctctcccttc gggagcagag tgacaggtgg tgcattggtg tcgtcagctc gtgtcgttag 360
 atgttgggtt aagtcccgca acgagcgcaa ccttgatctc tagttgttta gttgggcact 420
 ctaaggtgac tgccgggtgac aaaccggagg aagtgggga tgacgtcaaa tcatcatgcc 480
 ccagtagcgt gagggaagg tgaaaagcac ccgggaaggg gagtgaaga gatcctgaaa 540
 ccgtgtgcca tagtcaagagc ccgttaacgg gtgatggcgt gccttttgta gaatgaaccg 600
 gcgagttata agatccgtag tcaaaaggga aacagcccag accgccagct aaggtcccaa 660
 agtctgtatt gaaaaggatg tggagtgtct tagacaacta ggatgtggc ttagaagcag 720
 ccaccattta aagagtata ggggtgacac ctgcccgggt ctggaaggtt aaggagaggg 780
 gtttagcgtaa ctctgaactg aagccccagt aaacggcggc cgtaactata accgtccata 840

ggtagcgaaa	gaaatttgag	aggagctgtc	cttagtacga	gaggaccggg	atggacgcac	900
cggtaccagt	tgttctgccca	agggcatagc	tgggtagcta	tgtgcggaag	ggataaagtgc	960
tgaagcagtc	taagcatgaa	gccccctca	agatgagagc	agtaaaacaa	gcaaaccgcac	1020
aatcagaagc	taagaaagcg	caagcttctg	gaaagcacaa	atgctagtta	tggtacagaa	1080
tttgcaactg	aaacagacgt	gcatgctgtg	aaatttgcga	aagcttttgc	atattatattc	1140
gagccacagc	atcgtgatgt	tttacagctt	tatgcaccgg	aagcttttaa	tggataaatt	1200
taacgaacaa	gaaataaatc	tatccttgga	agaacttaaa	gatcaacgga	tgctgggcaag	1260
atataaaaa	taagataaaa	cagcactatc	aacactggag	cgattcttta	tctgaagaag	1320
gaagagcgat	gaaaacaaac	aagtacaata	caagacaaaa	gaaggtaaaa	tactctgtttt	1380
aggggaaaaa	ttcaagaaat	atagaagtga	tggctaaaaa	tgtagaaggg	gtcttgaaac	1440
cgttaacagc	tgttatggcg	acogtggcgg	cgtggttatc	gaaccatgc	tgaccgatca	1500
atggtagctg	cacacogccc	cccaaagtgc	cgattgaagc	cgtagagAAC	ggcgagatcc	1560
agttcgtccc	taaacagtac	ggcaacttcg	ttatcgctca	ggcgaaactcc	aaactggatg	1620
atgaaggccg	cttcttagaa	ggtgacttgt	cgtagcaaaag	gcgaatcaag	ccctgtctagc	1680
cacaactatg	cgtgctcgtg	gtgcacaagt	aacggatatt	acaatcattg	ttgttgccagc	1740
tgatgacggc	gtaataaaca	gttgaagoga	ttaaccatgc	gaaagcagca	ggagtaccaa	1800
ctttactcag	cttgctggta	tgctgtggtct	gatggccaat	ccatctggtc	gtatcatcga	1860
acttccaatc	aagtcttcgt	gaaggtttaa	cagtacttga	gtacttcatc	tctacgcattg	1920
gtgcgcgtaa	aggttcattgg	agtaagacct	acagtaagag	gttctgtaat	gaaccctaat	1980
gaccattccac	acogtgggtg	tgaaggtaga	tctctatcgc	gaaagtcacc	gtactccattg	2040
gggtaaacca	gcacttggat	acaaaacaag	cgcagttcgg	cgccacgcgc	ttcgggtgaaa	2100